



**Off-Base Drinking Water Sample Results,  
Level 2 Laboratory Report, Level 4 Laboratory Report,  
Electronic Data Deliverable, Data Validation Report,  
and the Sample Location Figure, SDG J24224-1**

*Naval Air Station Whidbey Island  
Oak Harbor, Washington*

June 2019

# TestAmerica

THE LEADER IN ENVIRONMENTAL TESTING


## ANALYTICAL REPORT

TestAmerica Laboratories, Inc.  
TestAmerica Sacramento  
880 Riverside Parkway  
West Sacramento, CA 95605  
Tel: (916)373-5600

TestAmerica Job ID: 320-24224-1  
Client Project/Site: Whidbey Island

For:  
CH2M Hill Constructors, Inc.  
1100 NE Circle Blvd  
Corvallis, Oregon 97330

Attn: Tiffany Hill

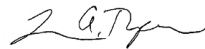

---

Authorized for release by:  
12/20/2016 4:36:23 PM

Laura Turpen, Project Manager I  
(916)374-4414  
[laura.turpen@testamericainc.com](mailto:laura.turpen@testamericainc.com)

### LINKS

Review your project  
results through  
**TotalAccess**

Have a Question?



Visit us at:  
[www.testamericainc.com](http://www.testamericainc.com)

*The test results in this report meet all 2003 NELAC and 2009 TNI requirements for accredited parameters, exceptions are noted in this report. This report may not be reproduced except in full, and with written approval from the laboratory. For questions please contact the Project Manager at the e-mail address or telephone number listed on this page.*

*This report has been electronically signed and authorized by the signatory. Electronic signature is intended to be the legally binding equivalent of a traditionally handwritten signature.*

*Results relate only to the items tested and the sample(s) as received by the laboratory.*

- 1
- 2
- 3
- 4
- 5
- 6
- 7
- 8
- 9
- 10
- 11
- 12
- 13
- 14
- 15



# Table of Contents

Cover Page . . . . .	1
Table of Contents . . . . .	2
Definitions/Glossary . . . . .	3
Case Narrative . . . . .	4
Detection Summary . . . . .	6
Client Sample Results . . . . .	7
Surrogate Summary . . . . .	8
QC Sample Results . . . . .	9
QC Association Summary . . . . .	11
Lab Chronicle . . . . .	12
Certification Summary . . . . .	13
Method Summary . . . . .	14
Sample Summary . . . . .	15
Chain of Custody . . . . .	16
Receipt Checklists . . . . .	17

# Definitions/Glossary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

## Qualifiers

### LCMS

Qualifier	Qualifier Description
U	Undetected at the Limit of Detection.
J	Estimated: The analyte was positively identified; the quantitation is an estimation

## Glossary

Abbreviation	These commonly used abbreviations may or may not be present in this report.
α	Listed under the "D" column to designate that the result is reported on a dry weight basis
%R	Percent Recovery
CFL	Contains Free Liquid
CNF	Contains no Free Liquid
DER	Duplicate error ratio (normalized absolute difference)
Dil Fac	Dilution Factor
DL, RA, RE, IN	Indicates a Dilution, Re-analysis, Re-extraction, or additional Initial metals/anion analysis of the sample
DLC	Decision level concentration
MDA	Minimum detectable activity
EDL	Estimated Detection Limit
MDC	Minimum detectable concentration
MDL	Method Detection Limit
ML	Minimum Level (Dioxin)
NC	Not Calculated
ND	Not detected at the reporting limit (or MDL or EDL if shown)
PQL	Practical Quantitation Limit
QC	Quality Control
RER	Relative error ratio
RL	Reporting Limit or Requested Limit (Radiochemistry)
RPD	Relative Percent Difference, a measure of the relative difference between two points
TEF	Toxicity Equivalent Factor (Dioxin)
TEQ	Toxicity Equivalent Quotient (Dioxin)

# Case Narrative

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

**Job ID: 320-24224-1**

**Laboratory: TestAmerica Sacramento**

**Narrative**

## CASE NARRATIVE

**Client: CH2M Hill Constructors, Inc.**

**Project: Whidbey Island**

**Report Number: 320-24224-1**

With the exceptions noted as flags or footnotes, standard analytical protocols were followed in the analysis of the samples and no problems were encountered or anomalies observed. In addition all laboratory quality control samples were within established control limits, with any exceptions noted below. Each sample was analyzed to achieve the lowest possible reporting limit within the constraints of the method. In some cases, due to interference or analytes present at high concentrations, samples were diluted. For diluted samples, the reporting limits are adjusted relative to the dilution required.

TestAmerica Sacramento attests to the validity of the laboratory data generated by TestAmerica facilities reported herein. All analyses performed by TestAmerica facilities were done using established laboratory SOPs that incorporate QA/QC procedures described in the applicable methods. TestAmerica's operations groups have reviewed the data for compliance with the laboratory QA/QC plan, and data have been found to be compliant with laboratory protocols unless otherwise noted below.

TestAmerica utilizes USEPA approved methods and DOD QSM, where applicable, in all analytical work. The samples presented in this report were analyzed for the parameter(s) listed on the analytical methods summary page in accordance with the method(s) indicated. A summary of QC data for these analyses is included at the back of the report.

All parameters for which TestAmerica Sacramento has certification were evaluated to the QSM specified reporting convention or to the client specified format if different from QSM. Parameters not certified under QSM, if any, were evaluated to the detection limit (DL) and include qualified results where applicable.

The sample(s) that contain constituents flagged with U are undetected. The result associated with this flag is the limit of detection (LOD).

Calculations are performed before rounding to avoid round-off errors in calculated results.

All holding times were met and proper preservation noted for the methods performed on these samples, unless otherwise detailed in the individual sections below.

All solid sample results are reported on an "as received" basis unless otherwise indicated by the presence of a % solids value in the method header.

This laboratory report is confidential and is intended for the sole use of TestAmerica and its client.

### **RECEIPT**

The samples were received on 12/09/2016; the samples arrived in good condition, properly preserved and on ice. The temperature of the coolers at receipt was 0.9 C.

The container label and the COC do not match for sample WI-AF-3FB30-1216 (320-24224-2). The container label lists WI-AF-FB30-1216, while the COC lists WI-AF-3FB30-1216. The sample was logged in according to the COC.

An extended TAT was requested by the client via an email on December 8. Samples received on December 8 and 9 were requested to have a due date of December 19, 2016.

### **PFOA/PFOS**

Samples WI-AF-3RW30-1216 (320-24224-1) and WI-AF-3FB30-1216 (320-24224-2) were analyzed for PFOA/PFOS in accordance with

# Case Narrative

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

---

## Job ID: 320-24224-1 (Continued)

---

### Laboratory: TestAmerica Sacramento (Continued)

537. The samples were prepared on 12/12/2016 and analyzed on 12/18/2016.

No analytical or quality issues were noted, other than those described above or in the Definitions/Glossary page.

1

2

3

4

5

6

7

8

9

10

11

12

13

14

15

# Detection Summary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

**Client Sample ID: WI-AF-3RW30-1216**

**Lab Sample ID: 320-24224-1**

No Detections.

**Client Sample ID: WI-AF-3FB30-1216**

**Lab Sample ID: 320-24224-2**

No Detections.

1

2

3

4

5

6

7

8

9

10

11

12

13

14

15

This Detection Summary does not include radiochemical test results.

TestAmerica Sacramento

# Client Sample Results

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

**Client Sample ID: WI-AF-3RW30-1216**

**Lab Sample ID: 320-24224-1**

**Date Collected: 12/07/16 14:45**

**Matrix: Water**

**Date Received: 12/09/16 09:50**

**Method: 537 - Perfluorinated Alkyl Acids (LC/MS)**

Analyte	Result	Qualifier	LOQ	DL	Unit	D	Prepared	Analyzed	Dil Fac
Perfluorooctanesulfonic acid (PFOS)	0.048	U	0.060	0.015	ug/L		12/12/16 10:03	12/18/16 09:01	1
Perfluorooctanoic acid (PFOA)	0.024	U	0.030	0.0094	ug/L		12/12/16 10:03	12/18/16 09:01	1
Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.14	0.047	ug/L		12/12/16 10:03	12/18/16 09:01	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C2 PFHxA	100		70 - 130				12/12/16 10:03	12/18/16 09:01	1
13C2 PFDA	117		70 - 130				12/12/16 10:03	12/18/16 09:01	1

**Client Sample ID: WI-AF-3FB30-1216**

**Lab Sample ID: 320-24224-2**

**Date Collected: 12/07/16 14:46**

**Matrix: Water**

**Date Received: 12/09/16 09:50**

**Method: 537 - Perfluorinated Alkyl Acids (LC/MS)**

Analyte	Result	Qualifier	LOQ	DL	Unit	D	Prepared	Analyzed	Dil Fac
Perfluorooctanesulfonic acid (PFOS)	0.047	U	0.058	0.015	ug/L		12/12/16 10:03	12/18/16 10:30	1
Perfluorooctanoic acid (PFOA)	0.023	U	0.029	0.0091	ug/L		12/12/16 10:03	12/18/16 10:30	1
Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.14	0.046	ug/L		12/12/16 10:03	12/18/16 10:30	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C2 PFHxA	121		70 - 130				12/12/16 10:03	12/18/16 10:30	1
13C2 PFDA	126		70 - 130				12/12/16 10:03	12/18/16 10:30	1



# Surrogate Summary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

## Method: 537 - Perfluorinated Alkyl Acids (LC/MS)

Matrix: Water

Prep Type: Total/NA

### Percent Surrogate Recovery (Acceptance Limits)

Lab Sample ID	Client Sample ID	Percent Surrogate Recovery (Acceptance Limits)	
		3C2 PFHx (70-130)	3C2 PFDA (70-130)
320-24224-1	WI-AF-3RW30-1216	100	117
320-24224-1 MS	WI-AF-3RW30-1216	99	122
320-24224-1 MSD	WI-AF-3RW30-1216	96	113
320-24224-2	WI-AF-3FB30-1216	121	126
LCS 320-141642/2-A	Lab Control Sample	118	123
MB 320-141642/1-A	Method Blank	110	111

#### Surrogate Legend

13C2 PFHxA = 13C2 PFHxA

13C2 PFDA = 13C2 PFDA

# QC Sample Results

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

## Method: 537 - Perfluorinated Alkyl Acids (LC/MS)

**Lab Sample ID: MB 320-141642/1-A**  
**Matrix: Water**  
**Analysis Batch: 142886**

**Client Sample ID: Method Blank**  
**Prep Type: Total/NA**  
**Prep Batch: 141642**

Analyte	MB Result	MB Qualifier	LOQ	DL	Unit	D	Prepared	Analyzed	Dil Fac
Perfluorooctanesulfonic acid (PFOS)	0.048	U	0.060	0.016	ug/L		12/12/16 10:03	12/19/16 13:19	1
Perfluorooctanoic acid (PFOA)	0.024	U	0.030	0.0094	ug/L		12/12/16 10:03	12/19/16 13:19	1
Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.14	0.048	ug/L		12/12/16 10:03	12/19/16 13:19	1

Surrogate	MB %Recovery	MB Qualifier	Limits	Prepared	Analyzed	Dil Fac
13C2 PFHxA	110		70 - 130	12/12/16 10:03	12/19/16 13:19	1
13C2 PFDA	111		70 - 130	12/12/16 10:03	12/19/16 13:19	1

**Lab Sample ID: LCS 320-141642/2-A**  
**Matrix: Water**  
**Analysis Batch: 142886**

**Client Sample ID: Lab Control Sample**  
**Prep Type: Total/NA**  
**Prep Batch: 141642**

Analyte	Spike Added	LCS Result	LCS Qualifier	Unit	D	%Rec	Limits
Perfluorooctanesulfonic acid (PFOS)	0.160	0.150		ug/L		94	70 - 130
Perfluorooctanoic acid (PFOA)	0.0811	0.0704		ug/L		87	70 - 130
Perfluorobutanesulfonic acid (PFBS)	0.359	0.334		ug/L		93	70 - 130

Surrogate	LCS %Recovery	LCS Qualifier	Limits
13C2 PFHxA	118		70 - 130
13C2 PFDA	123		70 - 130

**Lab Sample ID: 320-24224-1 MS**  
**Matrix: Water**  
**Analysis Batch: 142809**

**Client Sample ID: WI-AF-3RW30-1216**  
**Prep Type: Total/NA**  
**Prep Batch: 141642**

Analyte	Sample Result	Sample Qualifier	Spike Added	MS Result	MS Qualifier	Unit	D	%Rec	Limits
Perfluorooctanesulfonic acid (PFOS)	0.048	U	0.0388	0.0334	J	ug/L		86	70 - 130
Perfluorooctanoic acid (PFOA)	0.024	U	0.0192	0.0189	J	ug/L		98	70 - 130
Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.0870	0.0700	J	ug/L		80	70 - 130

Surrogate	MS %Recovery	MS Qualifier	Limits
13C2 PFHxA	99		70 - 130
13C2 PFDA	122		70 - 130

**Lab Sample ID: 320-24224-1 MSD**  
**Matrix: Water**  
**Analysis Batch: 142809**

**Client Sample ID: WI-AF-3RW30-1216**  
**Prep Type: Total/NA**  
**Prep Batch: 141642**

Analyte	Sample Result	Sample Qualifier	Spike Added	MSD Result	MSD Qualifier	Unit	D	%Rec	Limits	RPD	RPD Limit
Perfluorooctanesulfonic acid (PFOS)	0.048	U	0.0403	0.0337	J	ug/L		84	70 - 130	1	30
Perfluorooctanoic acid (PFOA)	0.024	U	0.0200	0.0165	J	ug/L		83	70 - 130	14	30
Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.0903	0.0706	J	ug/L		78	70 - 130	1	30

TestAmerica Sacramento

# QC Sample Results

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

## Method: 537 - Perfluorinated Alkyl Acids (LC/MS) (Continued)

Lab Sample ID: 320-24224-1 MSD

Matrix: Water

Analysis Batch: 142809

Client Sample ID: WI-AF-3RW30-1216

Prep Type: Total/NA

Prep Batch: 141642

Surrogate	MSD		Limits
	%Recovery	Qualifier	
13C2 PFHxA	96		70 - 130
13C2 PFDA	113		70 - 130

- 1
- 2
- 3
- 4
- 5
- 6
- 7
- 8
- 9
- 10
- 11
- 12
- 13
- 14
- 15

# QC Association Summary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

## LCMS

### Prep Batch: 141642

Lab Sample ID	Client Sample ID	Prep Type	Matrix	Method	Prep Batch
320-24224-1	WI-AF-3RW30-1216	Total/NA	Water	537	
320-24224-2	WI-AF-3FB30-1216	Total/NA	Water	537	
MB 320-141642/1-A	Method Blank	Total/NA	Water	537	
LCS 320-141642/2-A	Lab Control Sample	Total/NA	Water	537	
320-24224-1 MS	WI-AF-3RW30-1216	Total/NA	Water	537	
320-24224-1 MSD	WI-AF-3RW30-1216	Total/NA	Water	537	

### Analysis Batch: 142809

Lab Sample ID	Client Sample ID	Prep Type	Matrix	Method	Prep Batch
320-24224-1	WI-AF-3RW30-1216	Total/NA	Water	537	141642
320-24224-2	WI-AF-3FB30-1216	Total/NA	Water	537	141642
320-24224-1 MS	WI-AF-3RW30-1216	Total/NA	Water	537	141642
320-24224-1 MSD	WI-AF-3RW30-1216	Total/NA	Water	537	141642

### Analysis Batch: 142886

Lab Sample ID	Client Sample ID	Prep Type	Matrix	Method	Prep Batch
MB 320-141642/1-A	Method Blank	Total/NA	Water	537	141642
LCS 320-141642/2-A	Lab Control Sample	Total/NA	Water	537	141642

# Lab Chronicle

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

**Client Sample ID: WI-AF-3RW30-1216**

**Lab Sample ID: 320-24224-1**

**Date Collected: 12/07/16 14:45**

**Matrix: Water**

**Date Received: 12/09/16 09:50**

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	537			251.6 mL	1.0 mL	141642	12/12/16 10:03	NS1	TAL SAC
Total/NA	Analysis	537		1			142809	12/18/16 09:01	JRB	TAL SAC

**Client Sample ID: WI-AF-3FB30-1216**

**Lab Sample ID: 320-24224-2**

**Date Collected: 12/07/16 14:46**

**Matrix: Water**

**Date Received: 12/09/16 09:50**

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	537			257.5 mL	1.0 mL	141642	12/12/16 10:03	NS1	TAL SAC
Total/NA	Analysis	537		1			142809	12/18/16 10:30	JRB	TAL SAC

## Laboratory References:

TAL SAC = TestAmerica Sacramento, 880 Riverside Parkway, West Sacramento, CA 95605, TEL (916)373-5600

# Certification Summary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

## Laboratory: TestAmerica Sacramento

Unless otherwise noted, all analytes for this laboratory were covered under each certification below.

Authority	Program	EPA Region	Certification ID	Expiration Date
A2LA	DoD ELAP		2928-01	01-31-17

Analysis Method	Prep Method	Matrix	Analyte
-----------------	-------------	--------	---------

- 1
- 2
- 3
- 4
- 5
- 6
- 7
- 8
- 9
- 10
- 11
- 12
- 13
- 14
- 15

# Method Summary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

---

Method	Method Description	Protocol	Laboratory
537	Perfluorinated Alkyl Acids (LC/MS)	EPA	TAL SAC

---

**Protocol References:**

EPA = US Environmental Protection Agency

**Laboratory References:**

TAL SAC = TestAmerica Sacramento, 880 Riverside Parkway, West Sacramento, CA 95605, TEL (916)373-5600



# Sample Summary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

Lab Sample ID	Client Sample ID	Matrix	Collected	Received
320-24224-1	WI-AF-3RW30-1216	Water	12/07/16 14:45	12/09/16 09:50
320-24224-2	WI-AF-3FB30-1216	Water	12/07/16 14:46	12/09/16 09:50

1

2

3

4

5

6

7

8

9

10

11

12

13

14


15



Regulatory Program:  DW  NPDES  RCRA  Other:  Site Contact: Eric Epple  Date: 12/8/2016  
 Project Manager: Katie Tipping  Lab Contact: Laura Turpen  Carrier: FedEx  
 Tel/Fax: (757) 671-6258 Analysis Turnaround Time:  CALENDAR DAYS  WORKING DAYS  
 TAT if different from Below:  7-Day  2 weeks  1 week  2 days  1 day

Client Contact: Tiffany Hill  
 Project Chemist: 1100 NE Circle Blvd Ste 300 Corvallis, OR 97330  
 (541) 768-3109  
 (541) 908-3794  
 Project Name: CTO-08  
 Site: NAS Whidbey Island  
 P.O.#: 100067106050 - 679580.06.F.I.F.S

Sample Identification	Sample Date	Sample Time	Sample Type (C=Comp, G=Grab)	Matrix	# of Cont.	Filtered Sample (Y/N)	Perform MS/MSD (Y/N)	USEPA Method 537 (FOA, FOS, and FBS)	Date	Carrier	COC No: 7 of 1 COCs	Sampler:	For Lab Use Only: Walk-in Client: Lab Sampling: Job / SDG No.:	Sample Specific Notes:
WI-AF-3RW30-1216	12/07/16	1445	G	DW	2	N	X							
WI-AF-3RW30-1216-MS	12/07/16	1445	G	DW	2	N	X							
WI-AF-3RW30-1216-SD	12/07/16	1445	G	DW	2	N	X							
* WI-AF-3FB30-1216	12/07/16	1446	G	DW	2	N	X							
* container ID WI-AF-FB30-1216 006/12/16														

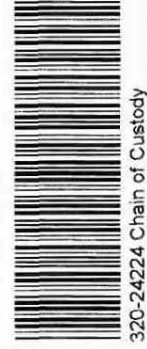

Preservation Used: 1= Ice, 2= HCl; 3= H2SO4; 4= HNO3; 5= NaOH; 6= Other - Trizma  
 Possible Hazard Identification: Please List any EPA Waste Codes for the sample in the Comments Section if the lab is to dispose of the sample.  
 Non-Hazard  Flammable  Skin Irritant  Poison B  Unknown  
 Return to Client  Disposal by Lab  Archive for \_\_\_\_\_ Months

Special Instructions/QC Requirements & Comments:  
 Sample Disposal (A fee may be assessed if samples are retained longer than 1 month)  
 Cooler Temp. (°C): Obs'd: \_\_\_\_\_ Cor'd: \_\_\_\_\_ Therm ID No.: \_\_\_\_\_

Relinquished by:	Date/Time:	Received by:	Date/Time:	Company:
Eric Epple	12-8-16/1600	Laura Turpen	12-9-16 09:50	JA SAE

## Login Sample Receipt Checklist

Client: CH2M Hill Constructors, Inc.

Job Number: 320-24224-1

**Login Number: 24224**

**List Source: TestAmerica Sacramento**

**List Number: 1**

**Creator: Turpen, Troy**

Question	Answer	Comment
Radioactivity wasn't checked or is <math>\leq</math> background as measured by a survey meter.	True	
The cooler's custody seal, if present, is intact.	True	(E88)
Sample custody seals, if present, are intact.	N/A	
The cooler or samples do not appear to have been compromised or tampered with.	True	
Samples were received on ice.	True	
Cooler Temperature is acceptable.	True	
Cooler Temperature is recorded.	True	
COC is present.	True	
COC is filled out in ink and legible.	True	
COC is filled out with all pertinent information.	True	
Is the Field Sampler's name present on COC?	N/A	
There are no discrepancies between the containers received and the COC.	False	Refer to Job Narrative for details.
Samples are received within Holding Time (excluding tests with immediate HTs)	True	
Sample containers have legible labels.	True	
Containers are not broken or leaking.	True	
Sample collection date/times are provided.	True	
Appropriate sample containers are used.	True	
Sample bottles are completely filled.	True	
Sample Preservation Verified.	N/A	
There is sufficient vol. for all requested analyses, incl. any requested MS/MSDs	True	
Containers requiring zero headspace have no headspace or bubble is <math><6\text{mm}</math> (1/4").	True	
Multiphasic samples are not present.	True	
Samples do not require splitting or compositing.	True	
Residual Chlorine Checked.	N/A	

## ANALYTICAL REPORT

Job Number: 320-24224-1  
Job Description: Whidbey Island

For:  
CH2M Hill Constructors, Inc.  
1100 NE Circle Blvd  
Corvallis, OR 97330  
Attention: Tiffany Hill



Approved for release.  
Laura Turpen  
Project Manager I  
12/20/2016 4:37 PM

---

Laura Turpen, Project Manager I  
880 Riverside Parkway, West Sacramento, CA, 95605  
(916)374-4414  
laura.turpen@testamericainc.com  
12/20/2016

# Table of Contents

Cover Title Page . . . . .	1
Data Summaries . . . . .	4
Definitions . . . . .	4
Case Narrative . . . . .	5
Detection Summary . . . . .	6
Client Sample Results . . . . .	7
Default Detection Limits . . . . .	8
Surrogate Summary . . . . .	9
QC Sample Results . . . . .	10
QC Association . . . . .	12
Chronicle . . . . .	13
Certification Summary . . . . .	14
Method Summary . . . . .	15
Sample Summary . . . . .	16
Manual Integration Summary . . . . .	17
Reagent Traceability . . . . .	19
COAs . . . . .	28
Organic Sample Data . . . . .	78
LCMS . . . . .	78
Method 537 DOD . . . . .	78
Method 537 DOD QC Summary . . . . .	79
Method 537 DOD Sample Data . . . . .	89
Standards Data . . . . .	97
Method 537 DOD ICAL Data . . . . .	97
Method 537 DOD CCAL Data . . . . .	120
Raw QC Data . . . . .	150

# Table of Contents

Method 537 DOD Blank Data .....	150
Method 537 DOD LCS/LCSD Data .....	154
Method 537 DOD MS/MSD Data .....	158
Method 537 DOD Run Logs .....	166
Method 537 DOD Prep Data .....	171
Shipping and Receiving Documents .....	185
Client Chain of Custody .....	186
Sample Receipt Checklist .....	187

# Definitions/Glossary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

---

## Qualifiers

---

### LCMS

Qualifier	Qualifier Description
U	Undetected at the Limit of Detection.
J	Estimated: The analyte was positively identified; the quantitation is an estimation

---

## Glossary

---

Abbreviation	These commonly used abbreviations may or may not be present in this report.
α	Listed under the "D" column to designate that the result is reported on a dry weight basis
%R	Percent Recovery
CFL	Contains Free Liquid
CNF	Contains no Free Liquid
DER	Duplicate error ratio (normalized absolute difference)
Dil Fac	Dilution Factor
DL, RA, RE, IN	Indicates a Dilution, Re-analysis, Re-extraction, or additional Initial metals/anion analysis of the sample
DLC	Decision level concentration
MDA	Minimum detectable activity
EDL	Estimated Detection Limit
MDC	Minimum detectable concentration
MDL	Method Detection Limit
ML	Minimum Level (Dioxin)
NC	Not Calculated
ND	Not detected at the reporting limit (or MDL or EDL if shown)
PQL	Practical Quantitation Limit
QC	Quality Control
RER	Relative error ratio
RL	Reporting Limit or Requested Limit (Radiochemistry)
RPD	Relative Percent Difference, a measure of the relative difference between two points
TEF	Toxicity Equivalent Factor (Dioxin)
TEQ	Toxicity Equivalent Quotient (Dioxin)

## CASE NARRATIVE

Client: CH2M Hill Constructors, Inc.

Project: Whidbey Island

Report Number: 320-24224-1

With the exceptions noted as flags or footnotes, standard analytical protocols were followed in the analysis of the samples and no problems were encountered or anomalies observed. In addition all laboratory quality control samples were within established control limits, with any exceptions noted below. Each sample was analyzed to achieve the lowest possible reporting limit within the constraints of the method. In some cases, due to interference or analytes present at high concentrations, samples were diluted. For diluted samples, the reporting limits are adjusted relative to the dilution required.

TestAmerica Sacramento attests to the validity of the laboratory data generated by TestAmerica facilities reported herein. All analyses performed by TestAmerica facilities were done using established laboratory SOPs that incorporate QA/QC procedures described in the applicable methods. TestAmerica's operations groups have reviewed the data for compliance with the laboratory QA/QC plan, and data have been found to be compliant with laboratory protocols unless otherwise noted below.

TestAmerica utilizes USEPA approved methods and DOD QSM, where applicable, in all analytical work. The samples presented in this report were analyzed for the parameter(s) listed on the analytical methods summary page in accordance with the method(s) indicated. A summary of QC data for these analyses is included at the back of the report.

All parameters for which TestAmerica Sacramento has certification were evaluated to the QSM specified reporting convention or to the client specified format if different from QSM. Parameters not certified under QSM, if any, were evaluated to the detection limit (DL) and include qualified results where applicable.

The sample(s) that contain constituents flagged with U are undetected. The result associated with this flag is the limit of detection (LOD).

Calculations are performed before rounding to avoid round-off errors in calculated results.

All holding times were met and proper preservation noted for the methods performed on these samples, unless otherwise detailed in the individual sections below.

All solid sample results are reported on an "as received" basis unless otherwise indicated by the presence of a % solids value in the method header.

This laboratory report is confidential and is intended for the sole use of TestAmerica and its client.

### **RECEIPT**

The samples were received on 12/09/2016; the samples arrived in good condition, properly preserved and on ice. The temperature of the coolers at receipt was 0.9 C.

The container label and the COC do not match for sample WI-AF-3FB30-1216 (320-24224-2). The container label lists WI-AF-FB30-1216, while the COC lists WI-AF-~~3~~FB30-1216. The sample was logged in according to the COC.

An extended TAT was requested by the client via an email on December 8. Samples received on December 8 and 9 were requested to have a due date of December 19, 2016.

### **PFOA/PFOS**

Samples WI-AF-3RW30-1216 (320-24224-1) and WI-AF-3FB30-1216 (320-24224-2) were analyzed for PFOA/PFOS in accordance with 537. The samples were prepared on 12/12/2016 and analyzed on 12/18/2016.

No analytical or quality issues were noted, other than those described above or in the Definitions/Glossary page.

## Detection Summary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

---

**Client Sample ID: WI-AF-3RW30-1216**

**Lab Sample ID: 320-24224-1**

No Detections.

---

**Client Sample ID: WI-AF-3FB30-1216**

**Lab Sample ID: 320-24224-2**

No Detections.

This Detection Summary does not include radiochemical test results.

TestAmerica Sacramento



# Client Sample Results

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

## Client Sample ID: WI-AF-3RW30-1216

Date Collected: 12/07/16 14:45

Date Received: 12/09/16 09:50

## Lab Sample ID: 320-24224-1

Matrix: Water

### Method: 537 - Perfluorinated Alkyl Acids (LC/MS)

Analyte	Result	Qualifier	LOQ	DL	Unit	D	Prepared	Analyzed	Dil Fac
Perfluorooctanesulfonic acid (PFOS)	0.048	U	0.060	0.015	ug/L		12/12/16 10:03	12/18/16 09:01	1
Perfluorooctanoic acid (PFOA)	0.024	U	0.030	0.0094	ug/L		12/12/16 10:03	12/18/16 09:01	1
Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.14	0.047	ug/L		12/12/16 10:03	12/18/16 09:01	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C2 PFHxA	100		70 - 130				12/12/16 10:03	12/18/16 09:01	1
13C2 PFDA	117		70 - 130				12/12/16 10:03	12/18/16 09:01	1

## Client Sample ID: WI-AF-3FB30-1216

Date Collected: 12/07/16 14:46

Date Received: 12/09/16 09:50

## Lab Sample ID: 320-24224-2

Matrix: Water

### Method: 537 - Perfluorinated Alkyl Acids (LC/MS)

Analyte	Result	Qualifier	LOQ	DL	Unit	D	Prepared	Analyzed	Dil Fac
Perfluorooctanesulfonic acid (PFOS)	0.047	U	0.058	0.015	ug/L		12/12/16 10:03	12/18/16 10:30	1
Perfluorooctanoic acid (PFOA)	0.023	U	0.029	0.0091	ug/L		12/12/16 10:03	12/18/16 10:30	1
Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.14	0.046	ug/L		12/12/16 10:03	12/18/16 10:30	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C2 PFHxA	121		70 - 130				12/12/16 10:03	12/18/16 10:30	1
13C2 PFDA	126		70 - 130				12/12/16 10:03	12/18/16 10:30	1

# Default Detection Limits

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

## Method: 537 - Perfluorinated Alkyl Acids (LC/MS)

Prep: 537

Analyte	LOQ	DL	Units	Method
Perfluorobutanesulfonic acid (PFBS)	0.14	0.048	ug/L	537
Perfluorooctanesulfonic acid (PFOS)	0.060	0.016	ug/L	537
Perfluorooctanoic acid (PFOA)	0.030	0.0094	ug/L	537

# Surrogate Summary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

## Method: 537 - Perfluorinated Alkyl Acids (LC/MS)

Matrix: Water

Prep Type: Total/NA

Lab Sample ID	Client Sample ID	Percent Surrogate Recovery (Acceptance Limits)	
		3C2 PFHx (70-130)	3C2 PFD (70-130)
320-24224-1	WI-AF-3RW30-1216	100	117
320-24224-1 MS	WI-AF-3RW30-1216	99	122
320-24224-1 MSD	WI-AF-3RW30-1216	96	113
320-24224-2	WI-AF-3FB30-1216	121	126
LCS 320-141642/2-A	Lab Control Sample	118	123
MB 320-141642/1-A	Method Blank	110	111

### Surrogate Legend

13C2 PFHxA = 13C2 PFHxA

13C2 PFDA = 13C2 PFDA

# QC Sample Results

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

## Method: 537 - Perfluorinated Alkyl Acids (LC/MS)

**Lab Sample ID: MB 320-141642/1-A**  
**Matrix: Water**  
**Analysis Batch: 142886**

**Client Sample ID: Method Blank**  
**Prep Type: Total/NA**  
**Prep Batch: 141642**

Analyte	MB MB		LOQ	DL	Unit	D	Prepared	Analyzed	Dil Fac
	Result	Qualifier							
Perfluorooctanesulfonic acid (PFOS)	0.048	U	0.060	0.016	ug/L		12/12/16 10:03	12/19/16 13:19	1
Perfluorooctanoic acid (PFOA)	0.024	U	0.030	0.0094	ug/L		12/12/16 10:03	12/19/16 13:19	1
Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.14	0.048	ug/L		12/12/16 10:03	12/19/16 13:19	1

Surrogate	MB MB		Limits	Prepared	Analyzed	Dil Fac
	%Recovery	Qualifier				
13C2 PFHxA	110		70 - 130	12/12/16 10:03	12/19/16 13:19	1
13C2 PFDA	111		70 - 130	12/12/16 10:03	12/19/16 13:19	1

**Lab Sample ID: LCS 320-141642/2-A**  
**Matrix: Water**  
**Analysis Batch: 142886**

**Client Sample ID: Lab Control Sample**  
**Prep Type: Total/NA**  
**Prep Batch: 141642**

Analyte	Spike Added	LCS Result	LCS Qualifier	Unit	D	%Rec	Limits
Perfluorooctanoic acid (PFOA)	0.0811	0.0704		ug/L		87	70 - 130
Perfluorobutanesulfonic acid (PFBS)	0.359	0.334		ug/L		93	70 - 130

Surrogate	LCS LCS		Limits
	%Recovery	Qualifier	
13C2 PFHxA	118		70 - 130
13C2 PFDA	123		70 - 130

**Lab Sample ID: 320-24224-1 MS**  
**Matrix: Water**  
**Analysis Batch: 142809**

**Client Sample ID: WI-AF-3RW30-1216**  
**Prep Type: Total/NA**  
**Prep Batch: 141642**

Analyte	Sample Result	Sample Qualifier	Spike Added	MS Result	MS Qualifier	Unit	D	%Rec	Limits
Perfluorooctanoic acid (PFOA)	0.024	U	0.0192	0.0189	J	ug/L		98	70 - 130
Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.0870	0.0700	J	ug/L		80	70 - 130

Surrogate	MS MS		Limits
	%Recovery	Qualifier	
13C2 PFHxA	99		70 - 130
13C2 PFDA	122		70 - 130

**Lab Sample ID: 320-24224-1 MSD**  
**Matrix: Water**  
**Analysis Batch: 142809**

**Client Sample ID: WI-AF-3RW30-1216**  
**Prep Type: Total/NA**  
**Prep Batch: 141642**

Analyte	Sample Result	Sample Qualifier	Spike Added	MSD Result	MSD Qualifier	Unit	D	%Rec	Limits	RPD	
										RPD	Limit
Perfluorooctanesulfonic acid (PFOS)	0.048	U	0.0403	0.0337	J	ug/L		84	70 - 130	1	30
Perfluorooctanoic acid (PFOA)	0.024	U	0.0200	0.0165	J	ug/L		83	70 - 130	14	30
Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.0903	0.0706	J	ug/L		78	70 - 130	1	30

TestAmerica Sacramento

# QC Sample Results

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

## Method: 537 - Perfluorinated Alkyl Acids (LC/MS) (Continued)

Lab Sample ID: 320-24224-1 MSD

Matrix: Water

Analysis Batch: 142809

Client Sample ID: WI-AF-3RW30-1216

Prep Type: Total/NA

Prep Batch: 141642

<i>Surrogate</i>	<i>MSD</i> <i>%Recovery</i>	<i>MSD</i> <i>Qualifier</i>	<i>Limits</i>
13C2 PFHxA	96		70 - 130
13C2 PFDA	113		70 - 130

# QC Association Summary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

## LCMS

### Prep Batch: 141642

Lab Sample ID	Client Sample ID	Prep Type	Matrix	Method	Prep Batch
320-24224-1	WI-AF-3RW30-1216	Total/NA	Water	537	
320-24224-2	WI-AF-3FB30-1216	Total/NA	Water	537	
MB 320-141642/1-A	Method Blank	Total/NA	Water	537	
LCS 320-141642/2-A	Lab Control Sample	Total/NA	Water	537	
320-24224-1 MS	WI-AF-3RW30-1216	Total/NA	Water	537	
320-24224-1 MSD	WI-AF-3RW30-1216	Total/NA	Water	537	

### Analysis Batch: 142809

Lab Sample ID	Client Sample ID	Prep Type	Matrix	Method	Prep Batch
320-24224-1	WI-AF-3RW30-1216	Total/NA	Water	537	141642
320-24224-2	WI-AF-3FB30-1216	Total/NA	Water	537	141642
320-24224-1 MS	WI-AF-3RW30-1216	Total/NA	Water	537	141642
320-24224-1 MSD	WI-AF-3RW30-1216	Total/NA	Water	537	141642

### Analysis Batch: 142886

Lab Sample ID	Client Sample ID	Prep Type	Matrix	Method	Prep Batch
MB 320-141642/1-A	Method Blank	Total/NA	Water	537	141642
LCS 320-141642/2-A	Lab Control Sample	Total/NA	Water	537	141642

# Lab Chronicle

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

**Client Sample ID: WI-AF-3RW30-1216**

**Date Collected: 12/07/16 14:45**

**Date Received: 12/09/16 09:50**

**Lab Sample ID: 320-24224-1**

**Matrix: Water**

Prep Type	Batch Type	Batch Method	Run	Dilution Factor	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	537			141642	12/12/16 10:03	NS1	TAL SAC
Total/NA	Analysis	537		1	142809	12/18/16 09:01	JRB	TAL SAC

**Client Sample ID: WI-AF-3FB30-1216**

**Date Collected: 12/07/16 14:46**

**Date Received: 12/09/16 09:50**

**Lab Sample ID: 320-24224-2**

**Matrix: Water**

Prep Type	Batch Type	Batch Method	Run	Dilution Factor	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	537			141642	12/12/16 10:03	NS1	TAL SAC
Total/NA	Analysis	537		1	142809	12/18/16 10:30	JRB	TAL SAC

**Laboratory References:**

TAL SAC = TestAmerica Sacramento, 880 Riverside Parkway, West Sacramento, CA 95605, TEL (916)373-5600

# Certification Summary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

## Laboratory: TestAmerica Sacramento

Unless otherwise noted, all analytes for this laboratory were covered under each certification below.

<u>Authority</u>	<u>Program</u>	<u>EPA Region</u>	<u>Certification ID</u>	<u>Expiration Date</u>
A2LA	DoD ELAP		2928-01	01-31-17

<u>Analysis Method</u>	<u>Prep Method</u>	<u>Matrix</u>	<u>Analyte</u>
------------------------	--------------------	---------------	----------------



# Method Summary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

---

---

<b>Method</b>	<b>Method Description</b>	<b>Protocol</b>	<b>Laboratory</b>
537	Perfluorinated Alkyl Acids (LC/MS)	EPA	TAL SAC

**Protocol References:**

EPA = US Environmental Protection Agency

**Laboratory References:**

TAL SAC = TestAmerica Sacramento, 880 Riverside Parkway, West Sacramento, CA 95605, TEL (916)373-5600

# Sample Summary

Client: CH2M Hill Constructors, Inc.  
Project/Site: Whidbey Island

TestAmerica Job ID: 320-24224-1

---

---

<b>Lab Sample ID</b>	<b>Client Sample ID</b>	<b>Matrix</b>	<b>Collected</b>	<b>Received</b>
320-24224-1	WI-AF-3RW30-1216	Water	12/07/16 14:45	12/09/16 09:50
320-24224-2	WI-AF-3FB30-1216	Water	12/07/16 14:46	12/09/16 09:50

LCMS MANUAL INTEGRATION SUMMARY

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Instrument ID: A6 Analysis Batch Number: 140688

Lab Sample ID: STD 320-140688/2 IC Client Sample ID: \_\_\_\_\_

Date Analyzed: 12/05/16 17:26 Lab File ID: 05DEC2016A6A\_004.d GC Column: Acquity ID: 2.1(mm)

COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION		
		REASON	ANALYST	DATE
Perfluoroheptanoic acid	19.37	Split Peak	barnettj	12/06/16 10:00
Perfluorooctanoic acid (PFOA)	20.05	Split Peak	barnettj	12/06/16 10:00

Lab Sample ID: STD 320-140688/3 IC Client Sample ID: \_\_\_\_\_

Date Analyzed: 12/05/16 17:55 Lab File ID: 05DEC2016A6A\_005.d GC Column: Acquity ID: 2.1(mm)

COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION		
		REASON	ANALYST	DATE
Perfluoroheptanoic acid	19.38	Split Peak	barnettj	12/06/16 10:03
Perfluorooctanoic acid (PFOA)	20.05	Split Peak	barnettj	12/06/16 10:03

Lab Sample ID: CCV 320-140688/9 CCVL Client Sample ID: \_\_\_\_\_

Date Analyzed: 12/05/16 20:53 Lab File ID: 05DEC2016A6A\_011.d GC Column: Acquity ID: 2.1(mm)

COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION		
		REASON	ANALYST	DATE
Perfluoroheptanoic acid	19.38	Split Peak	barnettj	12/06/16 10:08
Perfluorooctanoic acid (PFOA)	20.05	Split Peak	barnettj	12/06/16 10:08

LCMS MANUAL INTEGRATION SUMMARY

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Instrument ID: A6 Analysis Batch Number: 142223

Lab Sample ID: CCV 320-142223/2 CCVL Client Sample ID: \_\_\_\_\_

Date Analyzed: 12/15/16 08:03 Lab File ID: 15DEC2016A6A\_002.d GC Column: Acquity ID: 2.1(mm)

COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION		
		REASON	ANALYST	DATE
Perfluoroheptanoic acid	19.38	Split Peak	barnettj	12/15/16 15:05
Perfluorooctanoic acid (PFOA)	20.05	Split Peak	barnettj	12/15/16 15:05

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Sacramento

Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
<b>LC537-ICV_00017</b>	01/13/17	08/09/16	MeOH/H2O, Lot 067374	10 mL	LC537-IS_00018	200 uL	13C2-PFOA	10 ng/mL
.LC537-IS_00018	01/13/17	07/13/16	Methanol, Lot 090285	10000 uL	LCM2PFOA_00004	100 uL	13C2-PFOA	0.5 ug/mL
..LCM2PFOA_00004	03/19/17	Wellington Laboratories, Lot M2PFOA0312			LCMPFOS_00013	300 uL	13C4 PFOS	1.434 ug/mL
..LCMPFOS_00013	01/22/21	Wellington Laboratories, Lot MPFOS0116			(Purchased Reagent)		13C2-PFOA	50 ug/mL
<b>LC537-ICV_00017</b>	01/13/17	08/09/16	MeOH/H2O, Lot 067374	10 mL	LC537-SU_00017	500 uL	13C2 PFDA	10 ng/mL
					LC537ICIM_00013	25 uL	13C2 PFHxA	10 ng/mL
							Perfluorobutanesulfonic acid (PFBS)	114.77 ng/mL
							Perfluorooctanoic acid (PFOA)	25.0965 ng/mL
							Perfluorooctanesulfonic acid (PFOS)	27.2389 ng/mL
.LC537-SU_00017	01/19/17	07/19/16	Methanol, Lot 104453	25000 uL	LCMPFDA_00008	100 uL	13C2 PFDA	0.2 ug/mL
..LCMPFDA_00008	08/19/20	Wellington Laboratories, Lot MPFDA0815			LCMPFHxA_00009	100 uL	13C2 PFHxA	0.2 ug/mL
..LCMPFHxA_00009	04/09/20	Wellington Laboratories, Lot MPFHxA0415			(Purchased Reagent)		13C2 PFDA	50 ug/mL
..LC537ICIM_00013	02/05/17	08/09/16	Methanol, Lot 090285	25 mL	LC537-PFBS2_00005	0.5 mL	Perfluorobutanesulfonic acid (PFBS)	45.908 ug/mL
					LC537-PFOA2_00007	0.13 mL	Perfluorooctanoic acid (PFOA)	10.0386 ug/mL
					LC537-PFOS2_00005	0.22 mL	Perfluorooctanesulfonic acid (PFOS)	10.8956 ug/mL
..LC537-PFBS2_00005	03/01/17	02/29/16	Methanol, Lot 090285	10 mL	LC537_PFBS2_00001	0.023 g	Perfluorobutanesulfonic acid (PFBS)	2295.4 ug/mL
...LC537_PFBS2_00001	08/09/17	Santa Cruz Biotechnology, Lot H0112			(Purchased Reagent)		Perfluorobutanesulfonic acid (PFBS)	0.998 g/g
..LC537-PFOA2_00007	07/25/17	08/05/16	Methanol, Lot 090285	10 mL	LC537_PFOA2_00001	0.0195 g	Perfluorooctanoic acid (PFOA)	1930.5 ug/mL
...LC537_PFOA2_00001	07/25/17	Afla Aesar, Lot D24Y026			(Purchased Reagent)		Perfluorooctanoic acid (PFOA)	0.99 g/g
..LC537-PFOS2_00005	03/01/17	02/29/16	Methanol, Lot 090285	10 mL	LC537_PFOS2_00001	0.0159 g	Perfluorooctanesulfonic acid (PFOS)	1238.13 ug/mL
...LC537_PFOS2_00001	07/26/17	Sigma, Lot BCBF5116V			(Purchased Reagent)		Perfluorooctanesulfonic acid (PFOS)	0.7787 g/g
<b>LC537-IS_00026</b>	03/19/17	12/05/16	Methanol, Lot 090285	10000 uL	LCM2PFOA_00003	100 uL	13C2-PFOA	0.5 ug/mL
.LCM2PFOA_00003	03/19/17	Wellington Laboratories, Lot M2PFOA0312			LCMPFOS_00018	300 uL	13C4 PFOS	1.434 ug/mL
..LCMPFOS_00018	08/03/21	Wellington Laboratories, Lot MPFOS0816			(Purchased Reagent)		13C2-PFOA	50 ug/mL
					(Purchased Reagent)		13C4 PFOS	47.8 ug/mL
<b>LC537-L1_00015</b>	01/13/17	07/28/16	MeOH/H2O, Lot 090285	5 mL	LC537-IS_00018	100 uL	13C2-PFOA	10 ng/mL
					LC537-MSP_00012	24.4 uL	13C4 PFOS	28.68 ng/mL
							Perfluorobutanesulfonic acid (PFBS)	8.76058 ng/mL
							Perfluoroheptanoic acid	0.993847 ng/mL
							Perfluorohexanesulfonic acid	2.9532 ng/mL
							Perfluorononanoic acid	1.91737 ng/mL
							Perfluorooctanoic acid (PFOA)	1.9793 ng/mL
							Perfluorooctanesulfonic acid (PFOS)	3.91048 ng/mL
					LC537-SU_00017	250 uL	13C2 PFDA	10 ng/mL

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Sacramento

Job No.: 320-24224-1

SDG No.:

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
.LC537-IS_00018	01/13/17	07/13/16	Methanol, Lot 090285	10000 uL	LCM2PFOA_00004	100 uL	13C2 PFHxA	10 ng/mL
..LCM2PFOA_00004	03/19/17	Wellington Laboratories, Lot M2PFOA0312			(Purchased Reagent)		13C2-PFOA	0.5 ug/mL
..LCMPFOS_00013	01/22/21	Wellington Laboratories, Lot MPFOS0116			(Purchased Reagent)		13C4 PFOS	1.434 ug/mL
.LC537-MSP_00012	01/28/17	07/28/16	Methanol, Lot 090285	10000 uL	LC537SPIM_00012	200 uL	13C2 PFOS	47.8 ug/mL
							Perfluorobutanesulfonic acid (PFBS)	1795.2 ng/mL
							Perfluoroheptanoic acid	203.657 ng/mL
							Perfluorohexanesulfonic acid	605.164 ng/mL
							Perfluorononanoic acid	392.904 ng/mL
							Perfluorooctanoic acid (PFOA)	405.594 ng/mL
							Perfluorooctanesulfonic acid (PFOS)	801.328 ng/mL
..LC537SPIM_00012	01/28/17	07/28/16	Methanol, Lot 104453	10 mL	LC537-PFBS_00006	0.44 mL	Perfluorobutanesulfonic acid (PFBS)	89.76 ug/mL
					LC537-PFHpA_00010	0.1 mL	Perfluoroheptanoic acid	10.1829 ug/mL
					LC537-PFHxS_00008	0.3 mL	Perfluorohexanesulfonic acid	30.2582 ug/mL
					LC537-PFNA_00008	0.2 mL	Perfluorononanoic acid	19.6452 ug/mL
					LC537-PFOA_00009	0.098 mL	Perfluorooctanoic acid (PFOA)	20.2797 ug/mL
					LC537-PFOS_00006	0.4 mL	Perfluorooctanesulfonic acid (PFOS)	40.0664 ug/mL
...LC537-PFBS_00006	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537_PFBS_00002	0.0102 g	Perfluorobutanesulfonic acid (PFBS)	2040 ug/mL
....LC537_PFBS_00002	04/01/18	Sigma, Lot MKBP8842V			(Purchased Reagent)		Perfluorobutanesulfonic acid (PFBS)	1 g/g
...LC537-PFHpA_00010	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537_PFHpA_00002	0.0072 g	Perfluoroheptanoic acid	1018.29 ug/mL
...LC537_PFHpA_00002	04/01/18	Aldrich, Lot BCM2579V			(Purchased Reagent)		Perfluoroheptanoic acid	0.99 g/g
...LC537-PFHxS_00008	07/28/17	07/28/16	Methanol, Lot 090285	5.5 mL	LC537_PFHxS_00002	0.0061 g	Perfluorohexanesulfonic acid	1008.61 ug/mL
...LC537_PFHxS_00002	04/01/18	Sigma, Lot BCBL3545V			(Purchased Reagent)		Perfluorohexanesulfonic acid	0.9094 g/g
...LC537-PFNA_00008	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537_PFNA_00002	0.0051 g	Perfluorononanoic acid	982.26 ug/mL
...LC537_PFNA_00002	04/01/18	TCI America, Lot QN44F			(Purchased Reagent)		Perfluorononanoic acid	0.963 g/g
...LC537-PFOA_00009	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537_PFOA_00002	0.0145 g	Perfluorooctanoic acid (PFOA)	2069.36 ug/mL
...LC537_PFOA_00002	11/04/18	Fluka, Lot SZBD308XV			(Purchased Reagent)		Perfluorooctanoic acid (PFOA)	0.999 g/g
...LC537-PFOS_00006	07/28/17	07/28/16	Methanol, Lot 090285	6 mL	LC537_PFOS_00002	0.0066 g	Perfluorooctanesulfonic acid (PFOS)	1001.66 ug/mL
....LC537_PFOS_00002	08/09/17	Fluka, Lot SZBC222XV			(Purchased Reagent)		Perfluorooctanesulfonic acid (PFOS)	0.9106 g/g
.LC537-SU_00017	01/19/17	07/19/16	Methanol, Lot 104453	25000 uL	LCMPFDA_00008	100 uL	13C2 PFDA	0.2 ug/mL
..LCMPFDA_00008	08/19/20	Wellington Laboratories, Lot MPFDA0815			(Purchased Reagent)		13C2 PFHxA	0.2 ug/mL
..LCMPFHxA_00009	04/09/20	Wellington Laboratories, Lot MPFHxA0415			(Purchased Reagent)		13C2 PFDA	50 ug/mL
							13C2 PFHxA	50 ug/mL
<b>LC537-L2_00014</b>	01/13/17	07/28/16	MeOH/H2O, Lot 090285	5 mL	LC537-HSP_00010	34 uL	Perfluorobutanesulfonic acid (PFBS)	22.8888 ng/mL
							Perfluoroheptanoic acid	2.59663 ng/mL
							Perfluorohexanesulfonic acid	7.71585 ng/mL
							Perfluorononanoic acid	5.00953 ng/mL
							Perfluorooctanoic acid (PFOA)	5.17132 ng/mL
							Perfluorooctanesulfonic acid (PFOS)	10.2169 ng/mL

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Sacramento

Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
					LC537-IS_00018	100 uL	13C2-PFOA	10 ng/mL
							13C4 PFOS	28.68 ng/mL
					LC537-SU_00017	250 uL	13C2 PFDA	10 ng/mL
							13C2 PFHxA	10 ng/mL
.LC537-HSP_00010	01/28/17	07/28/16	Methanol, Lot 090285	10000 uL	LC537SPIM_00012	375 uL	Perfluorobutanesulfonic acid (PFBS)	3366 ng/mL
							Perfluoroheptanoic acid	381.857 ng/mL
							Perfluorohexanesulfonic acid	1134.68 ng/mL
							Perfluorononanoic acid	736.695 ng/mL
							Perfluorooctanoic acid (PFOA)	760.489 ng/mL
							Perfluorooctanesulfonic acid (PFOS)	1502.49 ng/mL
..LC537SPIM_00012	01/28/17	07/28/16	Methanol, Lot 104453	10 mL	LC537-PFBS_00006	0.44 mL	Perfluorobutanesulfonic acid (PFBS)	89.76 ug/mL
					LC537-PFHpA_00010	0.1 mL	Perfluoroheptanoic acid	10.1829 ug/mL
					LC537-PFHxS_00008	0.3 mL	Perfluorohexanesulfonic acid	30.2582 ug/mL
					LC537-PFNA_00008	0.2 mL	Perfluorononanoic acid	19.6452 ug/mL
					LC537-PFOA_00009	0.098 mL	Perfluorooctanoic acid (PFOA)	20.2797 ug/mL
					LC537-PFOS_00006	0.4 mL	Perfluorooctanesulfonic acid (PFOS)	40.0664 ug/mL
...LC537-PFBS_00006	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537_PFBS_00002	0.0102 g	Perfluorobutanesulfonic acid (PFBS)	2040 ug/mL
....LC537_PFBS_00002	04/01/18		Sigma, Lot MKBP8842V		(Purchased Reagent)		Perfluorobutanesulfonic acid (PFBS)	1 g/g
...LC537-PFHpA_00010	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537_PFHpA_00002	0.0072 g	Perfluoroheptanoic acid	1018.29 ug/mL
....LC537_PFHpA_00002	04/01/18		Aldrich, Lot BCM2579V		(Purchased Reagent)		Perfluoroheptanoic acid	0.99 g/g
...LC537-PFHxS_00008	07/28/17	07/28/16	Methanol, Lot 090285	5.5 mL	LC537_PFHxS_00002	0.0061 g	Perfluorohexanesulfonic acid	1008.61 ug/mL
....LC537_PFHxS_00002	04/01/18		Sigma, Lot BCBL3545V		(Purchased Reagent)		Perfluorohexanesulfonic acid	0.9094 g/g
...LC537-PFNA_00008	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537_PFNA_00002	0.0051 g	Perfluorononanoic acid	982.26 ug/mL
....LC537_PFNA_00002	04/01/18		TCI America, Lot QN44F		(Purchased Reagent)		Perfluorononanoic acid	0.963 g/g
...LC537-PFOA_00009	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537_PFOA_00002	0.0145 g	Perfluorooctanoic acid (PFOA)	2069.36 ug/mL
....LC537_PFOA_00002	11/04/18		Fluka, Lot SZBD308XV		(Purchased Reagent)		Perfluorooctanoic acid (PFOA)	0.999 g/g
...LC537-PFOS_00006	07/28/17	07/28/16	Methanol, Lot 090285	6 mL	LC537_PFOS_00002	0.0066 g	Perfluorooctanesulfonic acid (PFOS)	1001.66 ug/mL
....LC537_PFOS_00002	08/09/17		Fluka, Lot SZBC222XV		(Purchased Reagent)		Perfluorooctanesulfonic acid (PFOS)	0.9106 g/g
.LC537-IS_00018	01/13/17	07/13/16	Methanol, Lot 090285	10000 uL	LCM2PFOA_00004	100 uL	13C2-PFOA	0.5 ug/mL
..LCM2PFOA_00004	03/19/17		Wellington Laboratories, Lot M2PFOA0312		LCMPFOS_00013	300 uL	13C4 PFOS	1.434 ug/mL
..LCMPFOS_00013	01/22/21		Wellington Laboratories, Lot MPFOS0116		(Purchased Reagent)		13C2-PFOA	50 ug/mL
.LC537-SU_00017	01/19/17	07/19/16	Methanol, Lot 104453	25000 uL	LCMPFDA_00008	100 uL	13C4 PFOS	47.8 ug/mL
..LCMPFDA_00008	08/19/20		Wellington Laboratories, Lot MPFDA0815		LCMPFHxA_00009	100 uL	13C2 PFDA	0.2 ug/mL
..LCMPFHxA_00009	04/09/20		Wellington Laboratories, Lot MPFHxA0415		(Purchased Reagent)		13C2 PFHxA	0.2 ug/mL
..LCMPFHxA_00009	04/09/20		Wellington Laboratories, Lot MPFHxA0415		(Purchased Reagent)		13C2 PFDA	50 ug/mL
..LCMPFHxA_00009	04/09/20		Wellington Laboratories, Lot MPFHxA0415		(Purchased Reagent)		13C2 PFHxA	50 ug/mL
<b>LC537-L3_00016</b>	01/28/17	11/07/16	MeOH/H2O, Lot 090285	5 mL	LC537-HSP_00010	67 uL	Perfluorobutanesulfonic acid (PFBS)	45.1044 ng/mL
							Perfluoroheptanoic acid	5.11689 ng/mL
							Perfluorohexanesulfonic acid	15.2048 ng/mL
							Perfluorononanoic acid	9.87171 ng/mL

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Sacramento

Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
							Perfluorooctanoic acid (PFOA)	10.1905 ng/mL
							Perfluorooctanesulfonic acid (PFOS)	20.1334 ng/mL
					LC537-IS_00024	100 uL	13C2-PFOA	10 ng/mL
							13C4 PFOS	28.68 ng/mL
					LC537-SU_00020	250 uL	13C2 PFDA	10 ng/mL
							13C2 PFHxA	10 ng/mL
.LC537-HSP_00010	01/28/17	07/28/16	Methanol, Lot 090285	10000 uL	LC537SPIM_00012	375 uL	Perfluorobutanesulfonic acid (PFBS)	3366 ng/mL
							Perfluoroheptanoic acid	381.857 ng/mL
							Perfluorohexanesulfonic acid	1134.68 ng/mL
							Perfluorononanoic acid	736.695 ng/mL
							Perfluorooctanoic acid (PFOA)	760.489 ng/mL
							Perfluorooctanesulfonic acid (PFOS)	1502.49 ng/mL
..LC537SPIM_00012	01/28/17	07/28/16	Methanol, Lot 104453	10 mL	LC537-PFBS_00006	0.44 mL	Perfluorobutanesulfonic acid (PFBS)	89.76 ug/mL
					LC537-PFHpA_00010	0.1 mL	Perfluoroheptanoic acid	10.1829 ug/mL
					LC537-PFHxS_00008	0.3 mL	Perfluorohexanesulfonic acid	30.2582 ug/mL
					LC537-PFNA_00008	0.2 mL	Perfluorononanoic acid	19.6452 ug/mL
					LC537-PFOA_00009	0.098 mL	Perfluorooctanoic acid (PFOA)	20.2797 ug/mL
					LC537-PFOS_00006	0.4 mL	Perfluorooctanesulfonic acid (PFOS)	40.0664 ug/mL
...LC537-PFBS_00006	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537_PFBS_00002	0.0102 g	Perfluorobutanesulfonic acid (PFBS)	2040 ug/mL
....LC537_PFBS_00002	04/01/18		Sigma, Lot MKBP8842V			(Purchased Reagent)	Perfluorobutanesulfonic acid (PFBS)	1 g/g
...LC537-PFHpA_00010	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537 PFHpA_00002	0.0072 g	Perfluoroheptanoic acid	1018.29 ug/mL
....LC537 PFHpA_00002	04/01/18		Aldrich, Lot BCBM2579V			(Purchased Reagent)	Perfluoroheptanoic acid	0.99 g/g
...LC537-PFHxS_00008	07/28/17	07/28/16	Methanol, Lot 090285	5.5 mL	LC537 PFHxS_00002	0.0061 g	Perfluorohexanesulfonic acid	1008.61 ug/mL
....LC537 PFHxS_00002	04/01/18		Sigma, Lot BCBL3545V			(Purchased Reagent)	Perfluorohexanesulfonic acid	0.9094 g/g
...LC537-PFNA_00008	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537 PFNA_00002	0.0051 g	Perfluorononanoic acid	982.26 ug/mL
....LC537 PFNA_00002	04/01/18		TCI America, Lot QN44F			(Purchased Reagent)	Perfluorononanoic acid	0.963 g/g
...LC537-PFOA_00009	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537 PFOA_00002	0.0145 g	Perfluorooctanoic acid (PFOA)	2069.36 ug/mL
....LC537 PFOA_00002	11/04/18		Fluka, Lot SZBD308XV			(Purchased Reagent)	Perfluorooctanoic acid (PFOA)	0.999 g/g
...LC537-PFOS_00006	07/28/17	07/28/16	Methanol, Lot 090285	6 mL	LC537_PFOS_00002	0.0066 g	Perfluorooctanesulfonic acid (PFOS)	1001.66 ug/mL
....LC537_PFOS_00002	08/09/17		Fluka, Lot SZBC222XV			(Purchased Reagent)	Perfluorooctanesulfonic acid (PFOS)	0.9106 g/g
.LC537-IS_00024	03/19/17	11/03/16	Methanol, Lot 090285	10000 uL	LCM2PFOA_00003	100 uL	13C2-PFOA	0.5 ug/mL
					LCMPFOS_00018	300 uL	13C4 PFOS	1.434 ug/mL
..LCM2PFOA_00003	03/19/17		Wellington Laboratories, Lot M2PFOA0312			(Purchased Reagent)	13C2-PFOA	50 ug/mL
..LCMPFOS_00018	08/03/21		Wellington Laboratories, Lot MPFOS0816			(Purchased Reagent)	13C4 PFOS	47.8 ug/mL
.LC537-SU_00020	04/07/17	10/07/16	Methanol, Lot 104453	25000 uL	LCMPFDA_00008	100 uL	13C2 PFDA	0.2 ug/mL
					LCMPFHxA_00009	100 uL	13C2 PFHxA	0.2 ug/mL
..LCMPFDA_00008	08/19/20		Wellington Laboratories, Lot MPFDA0815			(Purchased Reagent)	13C2 PFDA	50 ug/mL
..LCMPFHxA_00009	04/09/20		Wellington Laboratories, Lot MPFHxA0415			(Purchased Reagent)	13C2 PFHxA	50 ug/mL
<b>LC537-L4_00015</b>	01/13/17	07/28/16	MeOH/H2O, Lot 090285	5 mL	LC537-HSP_00010	135 uL	Perfluorobutanesulfonic acid (PFBS)	90.882 ng/mL



REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Sacramento

Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
							Perfluoroheptanoic acid	10.3101 ng/mL
							Perfluorohexanesulfonic acid	30.6364 ng/mL
							Perfluorononanoic acid	19.8908 ng/mL
							Perfluorooctanoic acid (PFOA)	20.5332 ng/mL
							Perfluorooctanesulfonic acid (PFOS)	40.5672 ng/mL
					LC537-IS_00018	100 uL	13C2-PFOA	10 ng/mL
							13C4 PFOS	28.68 ng/mL
					LC537-SU_00017	250 uL	13C2 PFDA	10 ng/mL
							13C2 PFHxA	10 ng/mL
.LC537-HSP_00010	01/28/17	07/28/16	Methanol, Lot 090285	10000 uL	LC537SPIM_00012	375 uL	Perfluorobutanesulfonic acid (PFBS)	3366 ng/mL
							Perfluoroheptanoic acid	381.857 ng/mL
							Perfluorohexanesulfonic acid	1134.68 ng/mL
							Perfluorononanoic acid	736.695 ng/mL
							Perfluorooctanoic acid (PFOA)	760.489 ng/mL
							Perfluorooctanesulfonic acid (PFOS)	1502.49 ng/mL
..LC537SPIM_00012	01/28/17	07/28/16	Methanol, Lot 104453	10 mL	LC537-PFBS_00006	0.44 mL	Perfluorobutanesulfonic acid (PFBS)	89.76 ug/mL
					LC537-PFHpA_00010	0.1 mL	Perfluoroheptanoic acid	10.1829 ug/mL
					LC537-PFHxS_00008	0.3 mL	Perfluorohexanesulfonic acid	30.2582 ug/mL
					LC537-PFNA_00008	0.2 mL	Perfluorononanoic acid	19.6452 ug/mL
					LC537-PFOA_00009	0.098 mL	Perfluorooctanoic acid (PFOA)	20.2797 ug/mL
					LC537-PFOS_00006	0.4 mL	Perfluorooctanesulfonic acid (PFOS)	40.0664 ug/mL
...LC537-PFBS_00006	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537_PFBS_00002	0.0102 g	Perfluorobutanesulfonic acid (PFBS)	2040 ug/mL
....LC537_PFBS_00002	04/01/18		Sigma, Lot MKBP8842V		(Purchased Reagent)		Perfluorobutanesulfonic acid (PFBS)	1 g/g
...LC537-PFHpA_00010	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537 PFHpA_00002	0.0072 g	Perfluoroheptanoic acid	1018.29 ug/mL
....LC537 PFHpA_00002	04/01/18		Aldrich, Lot BCBM2579V		(Purchased Reagent)		Perfluoroheptanoic acid	0.99 g/g
...LC537-PFHxS_00008	07/28/17	07/28/16	Methanol, Lot 090285	5.5 mL	LC537 PFHxS_00002	0.0061 g	Perfluorohexanesulfonic acid	1008.61 ug/mL
....LC537 PFHxS_00002	04/01/18		Sigma, Lot BCBL3545V		(Purchased Reagent)		Perfluorohexanesulfonic acid	0.9094 g/g
...LC537-PFNA_00008	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537 PFNA_00002	0.0051 g	Perfluorononanoic acid	982.26 ug/mL
....LC537 PFNA_00002	04/01/18		TCI America, Lot QN44F		(Purchased Reagent)		Perfluorononanoic acid	0.963 g/g
...LC537-PFOA_00009	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537 PFOA_00002	0.0145 g	Perfluorooctanoic acid (PFOA)	2069.36 ug/mL
....LC537 PFOA_00002	11/04/18		Fluka, Lot SZBD308XV		(Purchased Reagent)		Perfluorooctanoic acid (PFOA)	0.999 g/g
...LC537-PFOS_00006	07/28/17	07/28/16	Methanol, Lot 090285	6 mL	LC537_PFOS_00002	0.0066 g	Perfluorooctanesulfonic acid (PFOS)	1001.66 ug/mL
....LC537_PFOS_00002	08/09/17		Fluka, Lot SZBC222XV		(Purchased Reagent)		Perfluorooctanesulfonic acid (PFOS)	0.9106 g/g
.LC537-IS_00018	01/13/17	07/13/16	Methanol, Lot 090285	10000 uL	LCM2PFOA_00004	100 uL	13C2-PFOA	0.5 ug/mL
					LCMPFOS_00013	300 uL	13C4 PFOS	1.434 ug/mL
..LCM2PFOA_00004	03/19/17		Wellington Laboratories, Lot M2PFOA0312		(Purchased Reagent)		13C2-PFOA	50 ug/mL
..LCMPFOS_00013	01/22/21		Wellington Laboratories, Lot MPFOS0116		(Purchased Reagent)		13C4 PFOS	47.8 ug/mL
.LC537-SU_00017	01/19/17	07/19/16	Methanol, Lot 104453	25000 uL	LCMPFDA_00008	100 uL	13C2 PFDA	0.2 ug/mL
					LCMPFHxA_00009	100 uL	13C2 PFHxA	0.2 ug/mL
..LCMPFDA_00008	08/19/20		Wellington Laboratories, Lot MPFDA0815		(Purchased Reagent)		13C2 PFDA	50 ug/mL

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Sacramento

Job No.: 320-24224-1

SDG No.:

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration		
					Reagent ID	Volume Added				
..LCMPFHxA_00009	04/09/20	Wellington Laboratories, Lot MPFHxA0415			(Purchased Reagent)		13C2 PFHxA	50 ug/mL		
LC537-L5_00017	01/28/17	11/07/16	MeOH/H2O, Lot 090285	5 mL	LC537-HSP_00010	200 uL	Perfluorobutanesulfonic acid (PFBS)	134.64 ng/mL		
							Perfluoroheptanoic acid	15.2743 ng/mL		
							Perfluorohexanesulfonic acid	45.3873 ng/mL		
							Perfluorononanoic acid	29.4678 ng/mL		
							Perfluorooctanoic acid (PFOA)	30.4196 ng/mL		
					Perfluorooctanesulfonic acid (PFOS)	60.0996 ng/mL				
					LC537-IS_00024	100 uL	13C2-PFOA	10 ng/mL		
							13C4 PFOS	28.68 ng/mL		
					LC537-SU_00020	250 uL	13C2 PFDA	10 ng/mL		
							13C2 PFHxA	10 ng/mL		
.LC537-HSP_00010	01/28/17	07/28/16	Methanol, Lot 090285	10000 uL	LC537SPIM_00012	375 uL	Perfluorobutanesulfonic acid (PFBS)	3366 ng/mL		
							Perfluoroheptanoic acid	381.857 ng/mL		
							Perfluorohexanesulfonic acid	1134.68 ng/mL		
							Perfluorononanoic acid	736.695 ng/mL		
							Perfluorooctanoic acid (PFOA)	760.489 ng/mL		
Perfluorooctanesulfonic acid (PFOS)	1502.49 ng/mL									
..LC537SPIM_00012	01/28/17	07/28/16	Methanol, Lot 104453	10 mL	LC537-PFBS_00006	0.44 mL	Perfluorobutanesulfonic acid (PFBS)	89.76 ug/mL		
							LC537-PFHpA_00010	0.1 mL	Perfluoroheptanoic acid	10.1829 ug/mL
							LC537-PFHxS_00008	0.3 mL	Perfluorohexanesulfonic acid	30.2582 ug/mL
							LC537-PFNA_00008	0.2 mL	Perfluorononanoic acid	19.6452 ug/mL
							LC537-PFOA_00009	0.098 mL	Perfluorooctanoic acid (PFOA)	20.2797 ug/mL
LC537-PFOS_00006	0.4 mL	Perfluorooctanesulfonic acid (PFOS)	40.0664 ug/mL							
...LC537-PFBS_00006	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537_PFBS_00002	0.0102 g	Perfluorobutanesulfonic acid (PFBS)	2040 ug/mL		
....LC537_PFBS_00002	04/01/18	Sigma, Lot MKBP8842V			(Purchased Reagent)		Perfluorobutanesulfonic acid (PFBS)	1 g/g		
...LC537-PFHpA_00010	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537 PFHpA_00002	0.0072 g	Perfluoroheptanoic acid	1018.29 ug/mL		
...LC537 PFHpA_00002	04/01/18	Aldrich, Lot BCM2579V			(Purchased Reagent)		Perfluoroheptanoic acid	0.99 g/g		
...LC537-PFHxS_00008	07/28/17	07/28/16	Methanol, Lot 090285	5.5 mL	LC537 PFHxS_00002	0.0061 g	Perfluorohexanesulfonic acid	1008.61 ug/mL		
...LC537 PFHxS_00002	04/01/18	Sigma, Lot BCBL3545V			(Purchased Reagent)		Perfluorohexanesulfonic acid	0.9094 g/g		
...LC537-PFNA_00008	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537 PFNA_00002	0.0051 g	Perfluorononanoic acid	982.26 ug/mL		
...LC537 PFNA_00002	04/01/18	TCI America, Lot QN44F			(Purchased Reagent)		Perfluorononanoic acid	0.963 g/g		
...LC537-PFOA_00009	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537 PFOA_00002	0.0145 g	Perfluorooctanoic acid (PFOA)	2069.36 ug/mL		
...LC537 PFOA_00002	11/04/18	Fluka, Lot SZBD308XV			(Purchased Reagent)		Perfluorooctanoic acid (PFOA)	0.999 g/g		
...LC537-PFOS_00006	07/28/17	07/28/16	Methanol, Lot 090285	6 mL	LC537_PFOS_00002	0.0066 g	Perfluorooctanesulfonic acid (PFOS)	1001.66 ug/mL		
....LC537_PFOS_00002	08/09/17	Fluka, Lot SZBC222XV			(Purchased Reagent)		Perfluorooctanesulfonic acid (PFOS)	0.9106 g/g		
.LC537-IS_00024	03/19/17	11/03/16	Methanol, Lot 090285	10000 uL	LCM2PFOA_00003	100 uL	13C2-PFOA	0.5 ug/mL		
					LCMPFOS_00018	300 uL	13C4 PFOS	1.434 ug/mL		
..LCM2PFOA_00003	03/19/17	Wellington Laboratories, Lot M2PFOA0312			(Purchased Reagent)		13C2-PFOA	50 ug/mL		
..LCMPFOS_00018	08/03/21	Wellington Laboratories, Lot MPFOS0816			(Purchased Reagent)		13C4 PFOS	47.8 ug/mL		

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Sacramento

Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
.LC537-SU_00020	04/07/17	10/07/16	Methanol, Lot 104453	25000 uL	LCMPFDA_00008	100 uL	13C2 PFDA	0.2 ug/mL
					LCMPFHxA_00009	100 uL	13C2 PFHxA	0.2 ug/mL
..LCMPFDA_00008	08/19/20	Wellington Laboratories, Lot MPFDA0815			(Purchased Reagent)		13C2 PFDA	50 ug/mL
..LCMPFHxA_00009	04/09/20	Wellington Laboratories, Lot MPFHxA0415			(Purchased Reagent)		13C2 PFHxA	50 ug/mL
<b>LC537-L6_00014</b>	01/13/17	07/28/16	MeOH/H2O, Lot 090285	5 mL	LC537-HSP_00010	265 uL	Perfluorobutanesulfonic acid (PFBS)	178.398 ng/mL
							Perfluoroheptanoic acid	20.2384 ng/mL
							Perfluorohexanesulfonic acid	60.1382 ng/mL
							Perfluorononanoic acid	39.0448 ng/mL
							Perfluorooctanoic acid (PFOA)	40.3059 ng/mL
							Perfluorooctanesulfonic acid (PFOS)	79.632 ng/mL
					LC537-IS_00018	100 uL	13C2-PFOA	10 ng/mL
							13C4 PFOS	28.68 ng/mL
					LC537-SU_00017	250 uL	13C2 PFDA	10 ng/mL
							13C2 PFHxA	10 ng/mL
.LC537-HSP_00010	01/28/17	07/28/16	Methanol, Lot 090285	10000 uL	LC537SPIM_00012	375 uL	Perfluorobutanesulfonic acid (PFBS)	3366 ng/mL
							Perfluoroheptanoic acid	381.857 ng/mL
							Perfluorohexanesulfonic acid	1134.68 ng/mL
							Perfluorononanoic acid	736.695 ng/mL
							Perfluorooctanoic acid (PFOA)	760.489 ng/mL
							Perfluorooctanesulfonic acid (PFOS)	1502.49 ng/mL
..LC537SPIM_00012	01/28/17	07/28/16	Methanol, Lot 104453	10 mL	LC537-PFBS_00006	0.44 mL	Perfluorobutanesulfonic acid (PFBS)	89.76 ug/mL
					LC537-PFHpA_00010	0.1 mL	Perfluoroheptanoic acid	10.1829 ug/mL
					LC537-PFHxS_00008	0.3 mL	Perfluorohexanesulfonic acid	30.2582 ug/mL
					LC537-PFNA_00008	0.2 mL	Perfluorononanoic acid	19.6452 ug/mL
					LC537-PFOA_00009	0.098 mL	Perfluorooctanoic acid (PFOA)	20.2797 ug/mL
					LC537-PFOS_00006	0.4 mL	Perfluorooctanesulfonic acid (PFOS)	40.0664 ug/mL
...LC537-PFBS_00006	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537_PFBS_00002	0.0102 g	Perfluorobutanesulfonic acid (PFBS)	2040 ug/mL
....LC537_PFBS_00002	04/01/18	Sigma, Lot MKBP8842V			(Purchased Reagent)		Perfluorobutanesulfonic acid (PFBS)	1 g/g
...LC537-PFHpA_00010	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537 PFHpA_00002	0.0072 g	Perfluoroheptanoic acid	1018.29 ug/mL
....LC537 PFHpA_00002	04/01/18	Aldrich, Lot BCM2579V			(Purchased Reagent)		Perfluoroheptanoic acid	0.99 g/g
...LC537-PFHxS_00008	07/28/17	07/28/16	Methanol, Lot 090285	5.5 mL	LC537 PFHxS_00002	0.0061 g	Perfluorohexanesulfonic acid	1008.61 ug/mL
....LC537 PFHxS_00002	04/01/18	Sigma, Lot BCBL3545V			(Purchased Reagent)		Perfluorohexanesulfonic acid	0.9094 g/g
...LC537-PFNA_00008	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537 PFNA_00002	0.0051 g	Perfluorononanoic acid	982.26 ug/mL
....LC537 PFNA_00002	04/01/18	TCI America, Lot QN44F			(Purchased Reagent)		Perfluorononanoic acid	0.963 g/g
...LC537-PFOA_00009	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537 PFOA_00002	0.0145 g	Perfluorooctanoic acid (PFOA)	2069.36 ug/mL
....LC537 PFOA_00002	11/04/18	Fluka, Lot SZBD308XV			(Purchased Reagent)		Perfluorooctanoic acid (PFOA)	0.999 g/g
...LC537-PFOS_00006	07/28/17	07/28/16	Methanol, Lot 090285	6 mL	LC537_PFOS_00002	0.0066 g	Perfluorooctanesulfonic acid (PFOS)	1001.66 ug/mL
....LC537_PFOS_00002	08/09/17	Fluka, Lot SZBC222XV			(Purchased Reagent)		Perfluorooctanesulfonic acid (PFOS)	0.9106 g/g
.LC537-IS_00018	01/13/17	07/13/16	Methanol, Lot 090285	10000 uL	LCM2PFOA_00004	100 uL	13C2-PFOA	0.5 ug/mL

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Sacramento

Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
..LCM2PFOA_00004	03/19/17		Wellington Laboratories, Lot M2PFOA0312		LCMPFOS_00013	300 uL	13C4 PFOS	1.434 ug/mL
..LCMPFOS_00013	01/22/21		Wellington Laboratories, Lot MPFOS0116		(Purchased Reagent)		13C2-PFOA	50 ug/mL
..LC537-SU_00017	01/19/17	07/19/16	Methanol, Lot 104453	25000 uL	LCMPFDA_00008	100 uL	13C2 PFDA	0.2 ug/mL
..LCMPFDA_00008	08/19/20		Wellington Laboratories, Lot MPFDA0815		LCMPFHxA_00009	100 uL	13C2 PFHxA	0.2 ug/mL
..LCMPFHxA_00009	04/09/20		Wellington Laboratories, Lot MPFHxA0415		(Purchased Reagent)		13C2 PFDA	50 ug/mL
..LCMPFHxA_00009	04/09/20		Wellington Laboratories, Lot MPFHxA0415		(Purchased Reagent)		13C2 PFHxA	50 ug/mL
<b>LC537-LSP_00016</b>	05/04/17	11/04/16	Methanol, Lot 090285	10000 uL	LC537SPIM_00015	50 uL	Perfluorobutane Sulfonate	448.8 ng/mL
							Perfluorobutanesulfonic acid (PFBS)	448.8 ng/mL
							Perfluoroheptanoic acid	53.7429 ng/mL
							Perfluorohexanesulfonic acid	151.291 ng/mL
							Perfluorononanoic acid	101.553 ng/mL
							Perfluorooctanoic acid (PFOA)	99.234 ng/mL
							Perfluorooctanesulfonic acid (PFOS)	200.332 ng/mL
.LC537SPIM_00015	05/04/17	11/04/16	Methanol, Lot 104453	10000 uL	LC537-PFBS_00006	440 uL	Perfluorobutane Sulfonate	89.76 ug/mL
							Perfluorobutanesulfonic acid (PFBS)	89.76 ug/mL
					LC537-PFHpA_00011	100 uL	Perfluoroheptanoic acid	10.7486 ug/mL
					LC537-PFHxS_00008	300 uL	Perfluorohexanesulfonic acid	30.2582 ug/mL
					LC537-PFNA_00009	200 uL	Perfluorononanoic acid	20.3105 ug/mL
					LC537-PFOA_00010	100 uL	Perfluorooctanoic acid (PFOA)	19.8468 ug/mL
					LC537-PFOS_00006	400 uL	Perfluorooctanesulfonic acid (PFOS)	40.0664 ug/mL
..LC537-PFBS_00006	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537_PFBS_00002	0.0102 g	Perfluorobutane Sulfonate	2040 ug/mL
							Perfluorobutanesulfonic acid (PFBS)	2040 ug/mL
...LC537_PFBS_00002	04/01/18		Sigma, Lot MKBP8842V		(Purchased Reagent)		Perfluorobutane Sulfonate	1 g/g
							Perfluorobutanesulfonic acid (PFBS)	1 g/g
..LC537-PFHpA_00011	11/04/17	11/04/16	Methanol, Lot 090285	7 mL	LC537 PFHpA_00002	0.0076 g	Perfluoroheptanoic acid	1074.86 ug/mL
...LC537 PFHpA_00002	04/01/18		Aldrich, Lot BCM2579V		(Purchased Reagent)		Perfluoroheptanoic acid	0.99 g/g
..LC537-PFHxS_00008	07/28/17	07/28/16	Methanol, Lot 090285	5.5 mL	LC537 PFHxS_00002	0.0061 g	Perfluorohexanesulfonic acid	1008.61 ug/mL
...LC537 PFHxS_00002	04/01/18		Sigma, Lot BCBL3545V		(Purchased Reagent)		Perfluorohexanesulfonic acid	0.9094 g/g
..LC537-PFNA_00009	11/04/17	11/04/16	Methanol, Lot 090285	5.5 mL	LC537 PFNA_00002	0.0058 g	Perfluorononanoic acid	1015.53 ug/mL
...LC537 PFNA_00002	04/01/18		TCI America, Lot QN44F		(Purchased Reagent)		Perfluorononanoic acid	0.963 g/g
..LC537-PFOA_00010	11/04/17	11/04/16	Methanol, Lot 090285	7.5 mL	LC537 PFOA_00002	0.0149 g	Perfluorooctanoic acid (PFOA)	1984.68 ug/mL
...LC537 PFOA_00002	11/04/18		Fluka, Lot SZBD308XV		(Purchased Reagent)		Perfluorooctanoic acid (PFOA)	0.999 g/g
..LC537-PFOS_00006	07/28/17	07/28/16	Methanol, Lot 090285	6 mL	LC537_PFOS_00002	0.0066 g	Perfluorooctanesulfonic acid (PFOS)	1001.66 ug/mL
...LC537_PFOS_00002	08/09/17		Fluka, Lot SZBC222XV		(Purchased Reagent)		Perfluorooctanesulfonic acid (PFOS)	0.9106 g/g
<b>LC537-MSP_00014</b>	03/14/17	09/14/16	Methanol, Lot 090285	10000 uL	LC537SPIM_00013	200 uL	Perfluorobutane Sulfonate	1795.2 ng/mL
							Perfluorobutanesulfonic acid (PFBS)	1795.2 ng/mL
							Perfluoroheptanoic acid	203.657 ng/mL
							Perfluorohexanesulfonic acid	605.164 ng/mL

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Sacramento

Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
							Perfluorononanoic acid	392.904 ng/mL
							Perfluorooctanoic acid (PFOA)	405.594 ng/mL
							Perfluorooctanesulfonic acid (PFOS)	801.328 ng/mL
.LC537SPIM_00013	03/14/17	09/14/16	Methanol, Lot 104453	10000 uL	LC537-PFBS_00006	440 uL	Perfluorobutane Sulfonate	89760 ng/mL
							Perfluorobutanesulfonic acid (PFBS)	89760 ng/mL
					LC537-PFHpA_00010	100 uL	Perfluoroheptanoic acid	10182.9 ng/mL
					LC537-PFHxS_00008	300 uL	Perfluorohexanesulfonic acid	30258.2 ng/mL
					LC537-PFNA_00008	200 uL	Perfluorononanoic acid	19645.2 ng/mL
					LC537-PFOA_00009	98 uL	Perfluorooctanoic acid (PFOA)	20279.7 ng/mL
					LC537-PFOS_00006	400 uL	Perfluorooctanesulfonic acid (PFOS)	40066.4 ng/mL
..LC537-PFBS_00006	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537_PFBS_00002	0.0102 g	Perfluorobutane Sulfonate	2040 ug/mL
							Perfluorobutanesulfonic acid (PFBS)	2040 ug/mL
...LC537_PFBS_00002	04/01/18		Sigma, Lot MKBP8842V			(Purchased Reagent)	Perfluorobutane Sulfonate	1 g/g
							Perfluorobutanesulfonic acid (PFBS)	1 g/g
..LC537-PFHpA_00010	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537 PFHpA_00002	0.0072 g	Perfluoroheptanoic acid	1018.29 ug/mL
..LC537 PFHpA_00002	04/01/18		Aldrich, Lot BCM2579V			(Purchased Reagent)	Perfluoroheptanoic acid	0.99 g/g
..LC537-PFHxS_00008	07/28/17	07/28/16	Methanol, Lot 090285	5.5 mL	LC537 PFHxS_00002	0.0061 g	Perfluorohexanesulfonic acid	1008.61 ug/mL
..LC537 PFHxS_00002	04/01/18		Sigma, Lot BCBL3545V			(Purchased Reagent)	Perfluorohexanesulfonic acid	0.9094 g/g
..LC537-PFNA_00008	07/28/17	07/28/16	Methanol, Lot 090285	5 mL	LC537 PFNA_00002	0.0051 g	Perfluorononanoic acid	982.26 ug/mL
..LC537 PFNA_00002	04/01/18		TCI America, Lot QN44F			(Purchased Reagent)	Perfluorononanoic acid	0.963 g/g
..LC537-PFOA_00009	07/28/17	07/28/16	Methanol, Lot 090285	7 mL	LC537 PFOA_00002	0.0145 g	Perfluorooctanoic acid (PFOA)	2069.36 ug/mL
..LC537 PFOA_00002	11/04/18		Fluka, Lot SZBD308XV			(Purchased Reagent)	Perfluorooctanoic acid (PFOA)	0.999 g/g
..LC537-PFOS_00006	07/28/17	07/28/16	Methanol, Lot 090285	6 mL	LC537_PFOS_00002	0.0066 g	Perfluorooctanesulfonic acid (PFOS)	1001.66 ug/mL
...LC537_PFOS_00002	08/09/17		Fluka, Lot SZBC222XV			(Purchased Reagent)	Perfluorooctanesulfonic acid (PFOS)	0.9106 g/g
<b>LC537-SU_00024</b>	06/05/17	12/05/16	Methanol, Lot 104453	20000 uL	LCMPFDA_00008	80 uL	13C2 PFDA	0.2 ug/mL
					LCMPFHxA_00009	80 uL	13C2 PFHxA	0.2 ug/mL
.LCMPFDA_00008	08/19/20		Wellington Laboratories, Lot MPFDA0815			(Purchased Reagent)	13C2 PFDA	50 ug/mL
.LCMPFHxA_00009	04/09/20		Wellington Laboratories, Lot MPFHxA0415			(Purchased Reagent)	13C2 PFHxA	50 ug/mL

Reagent

---

**LC537\_PFB\_00002**

7: 4/1/15 SPV

3050 Spruce Street, Saint Louis, MO 63103, USA

Website: www.sigmaaldrich.com

Email USA: techserv@sial.com


Outside USA: eurtechserv@sial.com

## Certificate of Analysis

Product Name:

Nonafluorobutane-1-sulfonic acid - 97%

**Product Number:** 562629  
**Batch Number:** MKBP8842V  
**Brand:** ALDRICH  
**CAS Number:** 375-73-5  
**MDL Number:** MFCD01320794  
**Formula:** C4HF9O3S  
**Formula Weight:** 300.10 g/mol  
**Storage Temperature:** Store at 2 - 8 °C  
**Quality Release Date:** 11 OCT 2013

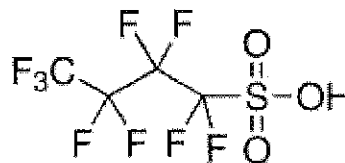

PFBS

Test	Specification	Result
Appearance (Color)	Colorless	Colorless
Appearance (Form)	Liquid	Liquid
Infrared Spectrum	Conforms to Structure	Conforms
Fluorine NMR Spectrum	Conforms to Structure	Conforms
Purity (Titration by NaOH)	96.5 - 103.5 %	101.6 %

*Jamie Gleason*

Jamie Gleason, Manager  
 Quality Control  
 Milwaukee, Wisconsin US


Sigma-Aldrich warrants, that at the time of the quality release or subsequent retest date this product conformed to the information contained in this publication. The current Specification sheet may be available at Sigma-Aldrich.com. For further inquiries, please contact Technical Service. Purchaser must determine the suitability of the product for its particular use. See reverse side of invoice or packing slip for additional terms and conditions of sale.

Reagent

---

**LC537\_PFB2\_00001**



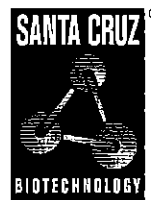

*The Power to Question*

# CERTIFICATE OF ANALYSIS

Catalog Number: sc-236187  
Product Name: Nonafluorobutane-1-sulfonic acid  
CAS Number: 375-73-5  
Molecular Formula:  $C_4HF_9O_3S$   
Molecular Weight: 300.10  
Lot Number: H0112

Test		Result
Refractive Index	1.3200 to 1.3290	1.3219
Purity (Titration)	min. 98.0%	99.8%

Test Conditions: Refractive Index: n<sub>20/D</sub>

Reagent

---

**LC537\_PFHpA\_00002**

R: 4/1/15 4V

### Certificate of Analysis

**Product Name:** PERFLUOROHEPTANOIC ACID  
 99 %

**Product Number:** 342041

**Batch Number:** BCBM2579V

**Brand:** Aldrich

**CAS Number:** 375-85-9

**Formula:**  $CF_3(CF_2)_5CO_2H$

**Formula Weight:** 364.06

**Quality Release Date:** 06 DEC 2013

**Recommended Retest Date:** OCT 2018

PFHpA

TEST	SPECIFICATION	RESULT
APPEARANCE (COLOR)	COLORLESS OR WHITE	WHITE
APPEARANCE (FORM)	LIQUID OR SOLID	SOLID
TITRATION	98.5 - 101.5 %	99.8 %
TITRATION (METHOD)	-	BACK TITRATION
PURITY (GC AREA %)	≥ 98.5 %	99.5 %
INFRARED SPECTRUM	CONFORMS TO STRUCTURE	CONFORMS

Dr. Claudia Geitner  
Manager Quality Control  
Buchs, Switzerland

Sigma-Aldrich warrants that at the time of the quality release or subsequent retest date this product conformed to the information contained in this publication. The current specification sheet may be available at Sigma-Aldrich.com. For further inquiries, please contact Technical Service. Purchaser must determine the suitability of the product for its particular use. See reverse side of invoice or packing slip for additional terms and conditions of sale.

Reagent

---

**LC537\_PFHxS\_00002**

r: 4/1/15 stw

### Certificate of Analysis

**Product Name:** TRIDECAFLUOROHEXANE-1-SULFONIC ACID POTASSIUM SALT  
 >= 98.0 % T

**Product Number:** 50929

**Batch Number:** BCBL3545V

**Brand:** Aldrich

**CAS Number:** 3871-99-6

**Formula:** C<sub>6</sub>F<sub>13</sub>KO<sub>3</sub>S

**Formula Weight:** 438.20

**Quality Release Date:** 20 JUN 2013

PFH<sub>13</sub>S-K

TEST	SPECIFICATION	RESULT
APPEARANCE (COLOR)	WHITE TO FAINT BEIGE	WHITE
APPEARANCE (FORM)	POWDER OR CRYSTALS	POWDER
TITRATION (ION EXCHANGE)	≥ 98.0 %	99.5 %
INFRARED SPECTRUM	CONFORMS TO STRUCTURE	CONFORMS

Dr. Claudia Geitner  
Manager Quality Control  
Buchs, Switzerland

$$MW_{corr} = \frac{(k_{form}) - (K) + (H)}{438.20 (k_{form})} = \frac{(438.20 - 39.10 + 1.01)}{438.20 (k_{form})} = 0.91307 \text{ (anion form)}$$

$$Purity = 90.94 \% \text{ w/m.w correction}$$

stw 4/1/15

Sigma-Aldrich warrants that at the time of the quality release or subsequent retest date this product conformed to the information contained in this publication. The current specification sheet may be available at Sigma-Aldrich.com. For further inquiries, please contact Technical Service. Purchaser must determine the suitability of the product for its particular use. See reverse side of invoice or packing slip for additional terms and conditions of sale.

Reagent

---

**LC537\_PENA\_00002**

R: 4/1/15 SKV



### Certificate of Analysis

Apr 2, 2015 (JST)

TOKYO CHEMICAL INDUSTRY CO.,LTD.  
4-10-1 Nihonbashi-Honcho, Chuo-ku, Tokyo 103-0023 Japan

Chemical Name: Heptadecafluorononanoic Acid		
Product Number: H0843 CAS: 375-95-1	Lot: QN44F	

Tests	Results	Specifications
Purity(GC)	96.3 %	min. 95.0 %
Purity(Neutralization titration)	98.1 %	min. 95.0 %
Melting point	63.3 deg-C	62.0 to 67.0 deg-C

TCI Lot numbers are 4-5 characters in length.  
Characters listed after the first 4-5 characters are control numbers for internal purpose only.

**Customer service:**

TCI AMERICA  
Tel: +1-800-423-8616 / +1-503-283-1681  
Fax: +1-888-520-1075 / +1-503-283-1987  
E-mail: Sales-US@TCIchemicals.com

PFNA

Reagent

---

**LC537\_PFOA\_00002**



3/21/15

# SIGMA-ALDRICH

## CERTIFICATE OF ANALYSIS

Sigma-Aldrich Laborchemikalien GmbH D-30918 Seelze  
Telefon: +49 5137 8238-150

Seelze, 13.11.2013/505378/13/24029
Order-No.:
Customer-No.:
Order-Code:
Quantity:
Production Date: 04.Nov.2013
Expiry Date: 04.Nov.2018

Article/Product: 33824	Batch : SZBD308XV	PFOA
Pentadecafluorooctanoic acid OEKANAL®		

### Reference Material (RM)

#### 1. General Information

Formula: C<sub>8</sub>HF<sub>15</sub>O<sub>2</sub>  
CAS-No.: [335-67-1]  
Usage : PFOA

Molar mass: 414.07 g/Mole  
Recomm. storage temp.: roomtemp.

The estimated uncertainty of a single measurement of the assay can be expected to be 0.5 % relative (confidence level = 95%, n= 6) whereby the assay measurements are calculated by 100% minus found impurities.

#### 2. Batch Analysis

identity (GC-MS)  
Assay (GCMS)  
Date of Analysis

complying  
99.4 %  
13.Nov.2013

#### 3. Advice and Remarks

- The expiry date is based on the current knowledge and holds only for proper storage conditions in the originally closed flasks/ packages.
- Whenever the container is opened for removal of aliquot portions of the substance, the person handling the substance must assure, that the integrity of the substance is maintained and proper records of all its handlings are kept. Special care has to be taken to avoid any contamination or adulteration of the substance.
- We herewith confirm that the delivery is effected according to the technical delivery conditions agreed.
- Particular properties of the products or the suitability for a particular area of application are not assured.
- We guarantee a proper quality within our General Conditions of Sales.

Sigma-Aldrich Laborchemikalien GmbH  
Quality Management SA-LC

This document was produced electronically and is valid without a signature

**GC/MS-Method**

**Analytical Department**

**Article:** Pentadecafluorooctanoic acid OEKANAL

**Article-No.:** 33824

**Batch:** SZBD308XV

**Column:** XTI-5 (Restek); 30 m; fs cap.; I.D.:0.25 mm; 1 µm df

**Injector:** Split mode

**Injection:** approx. 1 µl of reaction mixture with MSTFA (approx. 10 mg + 200 µl MSTFA)

**Inj.-temp.:** 280°C

**Oven-temp.:** 40°C (for 2 min) to 320°C (6°C/min) hold for 2 min

**Split:** 1:100

**Flow:** 1 ml He/min (Constant flow mode)

**Detector:** MSD


**Mass range:** 10-600 amu (Scan mode)

**Evaluation:** Purity: Total Ion Chromatogram  
(MSTFA and reaction products blinded out in report)

Identity: Mass spectrum complies


**Operator:** Ahrens / 2013-11-13

**Total Ion Chromatogram:**

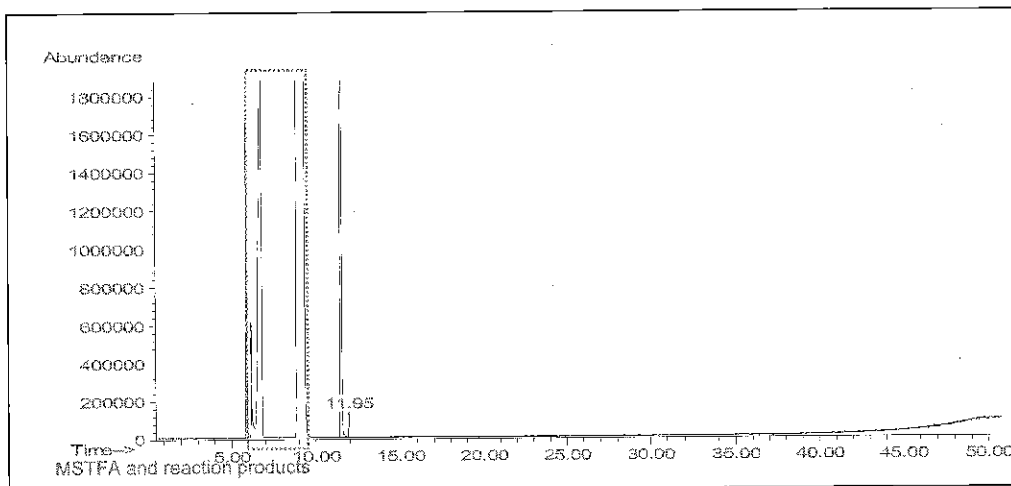

Ret. time	Area	Area-%	Com
11.54	565.1670	99.4	Pentadecafluorooctanoic acid (as TMS-ester)
11.95	3.6792	0.64	

**Mass spectrum (rt = 11.54 min):**

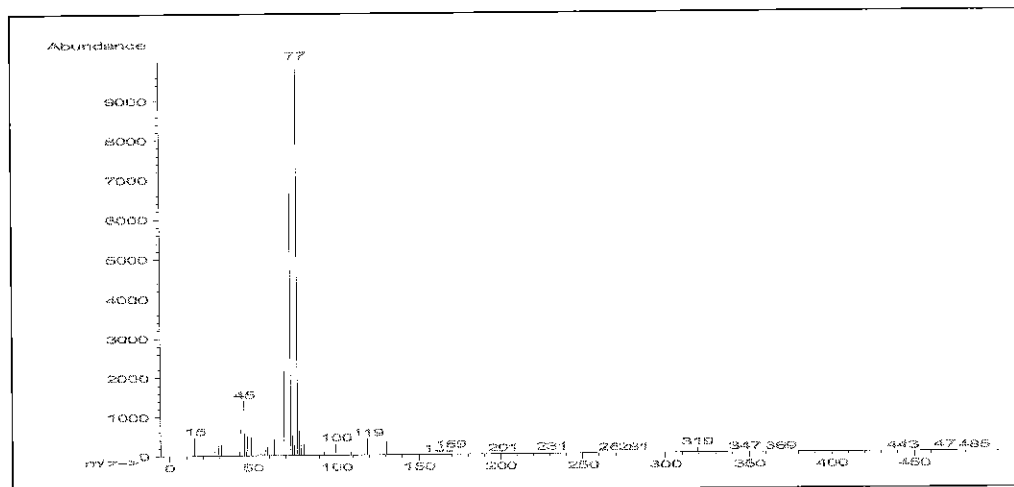

Reagent

---

**LC537\_PFOA2\_00001**

# Certificate of Analysis

**Alfa Aesar**  
A Johnson Matthey Company

Product No.: L08862  
Product: Perfluorooctanoic acid, 95%  
Lot No.: D24Y026

PFOA

Appearance White solid  
Melting point 58 - 60°C  
Assay 99 %  
Identity Matches reference

This document has been electronically generated and does not require a signature.

[www.alfa.com](http://www.alfa.com)

**NORTH AMERICA**  
Tel: +1-800-343-0660 or  
+1-978-521-6300  
Fax: +1-800-322-4757  
Email: [info@alfa.com](mailto:info@alfa.com)

**GERMANY**  
Tel: 00800 4566 4566 or  
+49 721 84007 280  
Fax: 00800 4577 4577 or  
+49 721 84007 300  
Email: [Eurosales@alfa.com](mailto:Eurosales@alfa.com)

**UNITED KINGDOM**  
Tel: 0800-801812 or  
+44 (0)1524-850506  
Fax: +44 (0)1524-850608  
Email: [UKsales@alfa.com](mailto:UKsales@alfa.com)

**FRANCE**  
Tel: 0800 03 51 47 or  
+33 (0)3 8862 2690  
Fax: 0800 10 20 67 or  
+33 (0)3 8862 6864  
Email: [frventes@alfa.com](mailto:frventes@alfa.com)

**INDIA**  
Tel: +91 8008 812424 or  
+91 8008 812525 or  
+91 8008 812626  
Fax: +91 8418 260060  
Email: [India@alfa.com](mailto:India@alfa.com)

**CHINA**  
Tel: +86 (010) 8567-8600  
Fax: +86 (010) 8567-8601  
Email: [saleschina@alfa-asia.com](mailto:saleschina@alfa-asia.com)

**KOREA**  
Tel: +82-2-3140-6000  
Fax: +82-2-3140-6002  
Email: [saleskorea@alfa-asia.com](mailto:saleskorea@alfa-asia.com)

Reagent

---

**LC537\_PFOs\_00002**

**SIGMA-ALDRICH®****CERTIFICATE OF ANALYSIS**

Sigma-Aldrich Laborchemikalien GmbH D-30918 Seelze  
 Telefon: +49 5137 8238-150

Seelze, 13.08.2012/419060/12/17583

Order-No.:

Customer-No.:

Order-Code:

Quantity:

Production Date: 09.Aug.2012

Expiry Date: 09.Aug.2017 - *err date*

Article/Product: 33829

Batch : SZBC222XV

Heptadecafluorooctanesulfonic acid potassium salt OEKANAL®

*PFOS-K<sup>+</sup>*

**Reference Material (RM)****1. General Information**

Formula: C8F17KO3S  
 CAS-No.: [2795-39-3]  
 Usage : PFOS

Molar mass: 538.22 g/Mole  
 Recomm. storage temp.: roomtemp.

The estimated uncertainty of a single measurement of the assay can be expected to be 0.5 % relative (confidence level = 95%, n= 6) whereby the assay measurements are calculated by 100% minus found impurities.

**2. Batch Analysis**

Identity  
 Assay (LC-MS)  
 Date of Analysis

complying  
 98 %  
 10.Aug.2012

$$\text{PW-Correction: } \frac{538.22 - 39.10 + 1.01}{538.22} = \frac{500.13}{538.22} = 0.92923$$

*Purity = 91.06%*

**3. Advice and Remarks**

- The minimum shelf life is based on the current knowledge and holds only for proper storage conditions in the originally closed flasks/ packages.
- Whenever the container is opened for removal of aliquot portions of the substance, the person handling the substance must assure, that the integrity of the substance is maintained and proper records of all its handlings are kept. Special care has to be taken to avoid any contamination or adulteration of the substance.
- We herewith confirm that the delivery is effected according to the technical delivery conditions agreed.
- Particular properties of the products or the suitability for a particular area of application are not assured.
- We guarantee a proper quality within our General Conditions of Sales.

Sigma-Aldrich Laborchemikalien GmbH  
 Quality Management SA-LC

Reagent

---

**LC537\_PFOs2\_00001**

Certificate of Analysis

Inv 820  
12LCMS 0579

Product Name: HEPTADEC AFLUORO OCTANESULFONIC ACID TETRAETHYLAMMONIUM SALT  
98 %  
Product Number: 365289  
Product Brand: Aldrich  
Molecular Formula: C<sub>16</sub>H<sub>20</sub>F<sub>17</sub>NO<sub>3</sub>S  
Molecular Mass: 629.37  
CAS Number: 56773-42-3

TEST	SPECIFICATION	LOT BCBF5116V RESULTS
APPEARANCE (COLOR)	OFF-WHITE TO WHITE	WHITE
APPEARANCE (FORM)	POWDER, LUMPS OR CHUNKS	POWDER WITH LUMPS
CARBON CONTENT	29.77 % - 31.29 %	30.52
INFRARED SPECTRUM	CONFORMS TO STRUCTURE	CONFORMS

QC RELEASE DATE 13/APR/11



$$\text{Mw correction} = \frac{500.125}{629.37} = 0.7946$$

~~79.46%~~ Oct 7-26-12

*E. Schwarzler*

Purity + Mw Correction = 77.87%

Edeltraud Schwarzler, Manager  
Quality Control  
Buchs, Switzerland

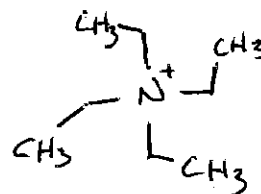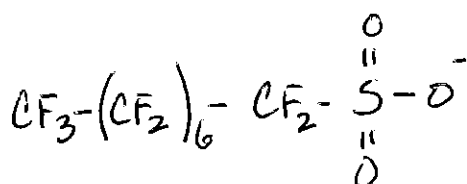

	<u>C<sub>8</sub>F<sub>17</sub>SO<sub>3</sub>H</u>	<u>C<sub>8</sub>H<sub>20</sub>N</u>
C = 12.011	96.088	96.088
F = 18.998	322.966	-
S = 32.066	32.066	-
O = 15.999	47.997	-
H = 1.008	1.008	20.160
N = 14.007	-	14.007
	<u>500.125</u>	<u>130.255</u> →

Sigma-Aldrich warrants, that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser must determine the suitability of the product for its particular use. See reverse side of invoice for additional terms and conditions of sale. The values given on the 'Certificate of Analysis' are the results determined at the time of analysis.



## Certificate of Origin

**Product Name:** Heptadecafluorooctanesulfonic acid tetraethylammonium salt  
98 %  
**Product Number:** 365289  
**Product Brand:** Aldrich  
**Lot:** BCBF5116V  
**Molecular Formula:**  $C_{16}H_{20}F_{17}NO_3S$   
**Molecular Mass:** 629.37  
**CAS Number:** 56773-42-3  
**Date of Issue:** 30-MAR-11

---

**Country of Origin** China

---

<b>product is of synthetic origin</b>	yes
<b>only synthetic materials used in the manufacturing process</b>	yes
<b>compounds of animal origin used</b>	no
<b>genetically modified organisms used</b>	no
<b>allergenic materials used</b>	no
<b>procedures in place to avoid cross contamination with residue of animal, human, GMO or allergenes in manufacturing process</b>	yes

Sigma-Aldrich has quality systems and procedures in place for monitoring the production process, traceability and batch consistency.

Document issued by Sigma-Aldrich Corporation "Sigma-Aldrich". This document is valid without signature and has been produced digitally.

This information is to be used for the purpose of determining animal or other biological origin only and not to be confused with "Country of Origin" for import/export purposes. Data provided on this document are property of Sigma-Aldrich.

This information is considered accurate and reliable as of the date appearing on the document and is presented in good faith.

Sigma-Aldrich shall not be held liable for any damage resulting from handling or from processing the above product(s). This document does not make any warranty, express or implied, of fitness for any particular use of the product(s). Purchaser must determine the suitability of the product(s) for its use under the applicable law and regulations.


For further questions please contact your local Sigma-Aldrich representative.

*We are committed to the success of our Customers, Employees and Shareholders through leadership in Life Science, High Technology and Service.*

Reagent

---


**LCM2PFOA\_00003**

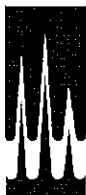

# WELLINGTON LABORATORIES

## CERTIFICATE OF ANALYSIS DOCUMENTATION

**PRODUCT CODE:** M2PFOA **LOT NUMBER:** M2PFOA0312  
**COMPOUND:** Perfluoro-n-[1,2-<sup>13</sup>C<sub>2</sub>]octanoic acid  
**STRUCTURE:** **CAS #:** Not available

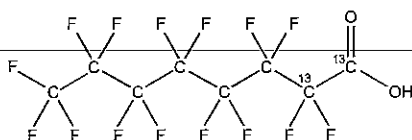

**MOLECULAR FORMULA:** <sup>13</sup>C<sub>2</sub><sup>12</sup>C<sub>6</sub>HF<sub>15</sub>O<sub>2</sub> **MOLECULAR WEIGHT:** 416.05  
**CONCENTRATION:** 50 ± 2.5 µg/ml **SOLVENT(S):** Methanol  
Water (<1%)  
**CHEMICAL PURITY:** >98% **ISOTOPIC PURITY:** ≥99%<sup>13</sup>C  
**LAST TESTED:** (mm/dd/yyyy) 03/19/2012 (1,2-<sup>13</sup>C<sub>2</sub>)  
**EXPIRY DATE:** (mm/dd/yyyy) 03/19/2017  
**RECOMMENDED STORAGE:** Store ampoule in a cool, dark place

### DOCUMENTATION/ DATA ATTACHED:

Figure 1: LC/MS Data (TIC and Mass Spectrum)  
Figure 2: LC/MS/MS Data (Selected MRM Transitions)

### ADDITIONAL INFORMATION:

- See page 2 for further details.
- Contains 4 mole eq. of NaOH to prevent conversion of the carboxylic acid to the methyl ester.

**FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE**

Certified By: \_\_\_\_\_

B.G. Chittim

Date: 01/09/2013  
(mm/dd/yyyy)

Wellington Laboratories Inc., 345 Southgate Dr. Guelph ON N1G 3M5 CANADA  
519-822-2436 • Fax: 519-822-2849 • info@well-labs.com

#### **INTENDED USE:**

The products prepared by Wellington Laboratories Inc. are for laboratory use only. They are designed to be used as reference standards for the identification and/or quantification of specific chemical compound(s).

#### **HAZARDS:**

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Material Safety Data Sheets (MSDSs) are available upon request.

#### **SYNTHESIS / CHARACTERIZATION:**

Where possible, all of our products are synthesized using single-product, unambiguous routes. They are then characterized, and their structures and purities confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, x-ray crystallography and melting point. Isotopic purities of mass-labelled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

#### **HOMOGENEITY:**

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given diluent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS and/or LC/MS/MS. The relative response factors of the analyte of interest in each solution are required to be <5% RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers.

#### **UNCERTAINTY:**

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty,  $u_c(y)$ , of a value  $y$  and the uncertainty of the independent parameters  $x_1, x_2, \dots, x_n$  on which it depends is:

$$u_c(y(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{i=1}^n u(y, x_i)^2}$$

where  $x$  is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of  $\pm 5\%$  (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all our products.

#### **TRACEABILITY:**

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external, ISO/IEC 17025:2005 accredited calibration company. In addition, their calibration is verified prior to each weighing using NIST and/or NRC traceable external weights. All volumetric glassware used is of Class A tolerance and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to international interlaboratory studies has also been established.

#### **EXPIRY DATE / PERIOD OF VALIDITY:**

Ongoing stability studies of this product have demonstrated stability in its composition and concentration for the period of time specified by the expiry date in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

#### **LIMITED WARRANTY:**

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.



#### **QUALITY MANAGEMENT:**

This product was produced using a Quality Management System registered to ISO 9001:2008 by SAI Global, ISO/IEC 17025:2005 by the Canadian Association for Laboratory Accreditation Inc. (CALA; A 1226), and ISO GUIDE 34:2009 by ACLASS (certificate number AR-1523).



\*\*For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at [www.well-labs.com](http://www.well-labs.com) or contact us directly at [info@well-labs.com](mailto:info@well-labs.com)\*\*

**Figure 1: M2PFOA; LC/MS Data (TIC and Mass Spectrum)**

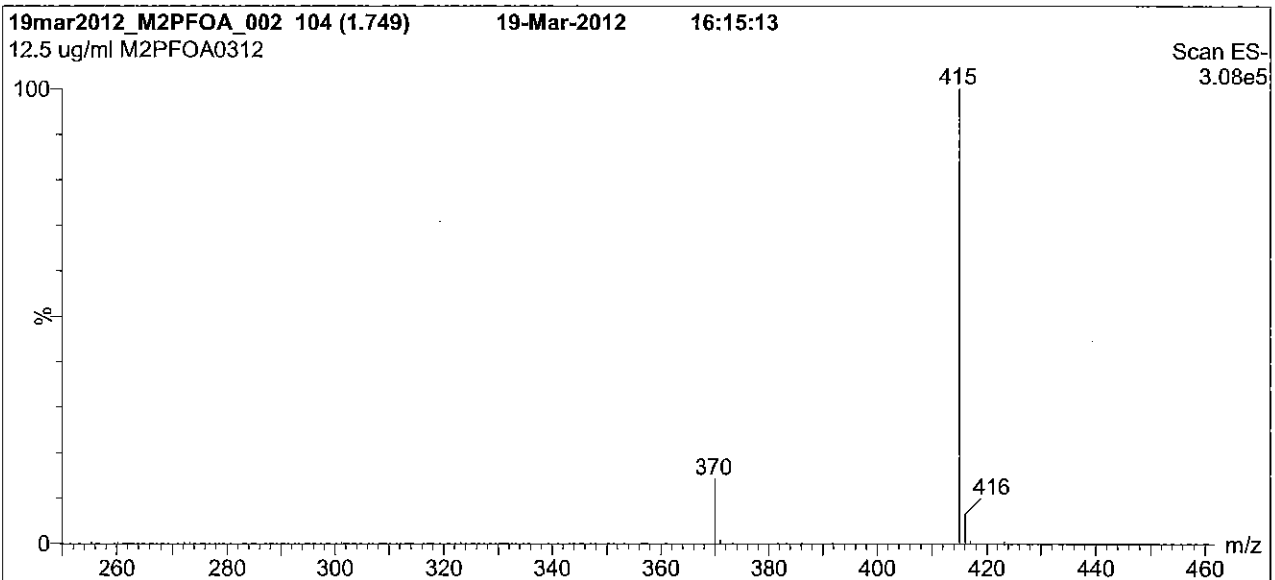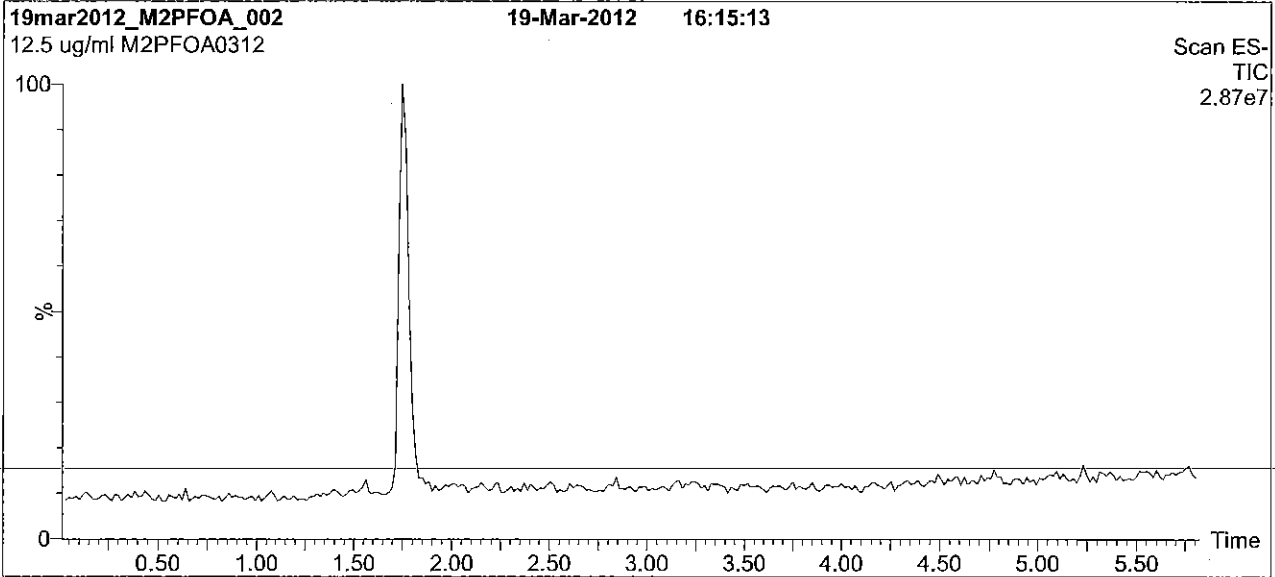

**Conditions for Figure 1:**

**LC:** Waters Acquity Ultra Performance LC  
**MS:** Micromass Quattro *micro* API MS

**Chromatographic Conditions**

**Column:** Acquity UPLC BEH Shield RP<sub>18</sub>  
1.7  $\mu$ m, 2.1 x 100 mm

**Mobile phase:** Gradient  
Start: 60% (80:20 MeOH:ACN) / 40% H<sub>2</sub>O  
(both with 10 mM NH<sub>4</sub>OAc buffer)  
Ramp to 90% organic over 6.5 min and hold for 2 min  
before returning to initial conditions in 0.5 min.  
Time: 10 min


**Flow:** 300  $\mu$ l/min

**MS Parameters**

**Experiment:** Full Scan (250 - 850 amu)

**Source:** Electrospray (negative)  
Capillary Voltage (kV) = 2.00  
Cone Voltage (V) = 15.00  
Cone Gas Flow (l/hr) = 100  
Desolvation Gas Flow (l/hr) = 750

**Figure 2: M2PFOA; LC/MS/MS Data (Selected MRM Transitions)**

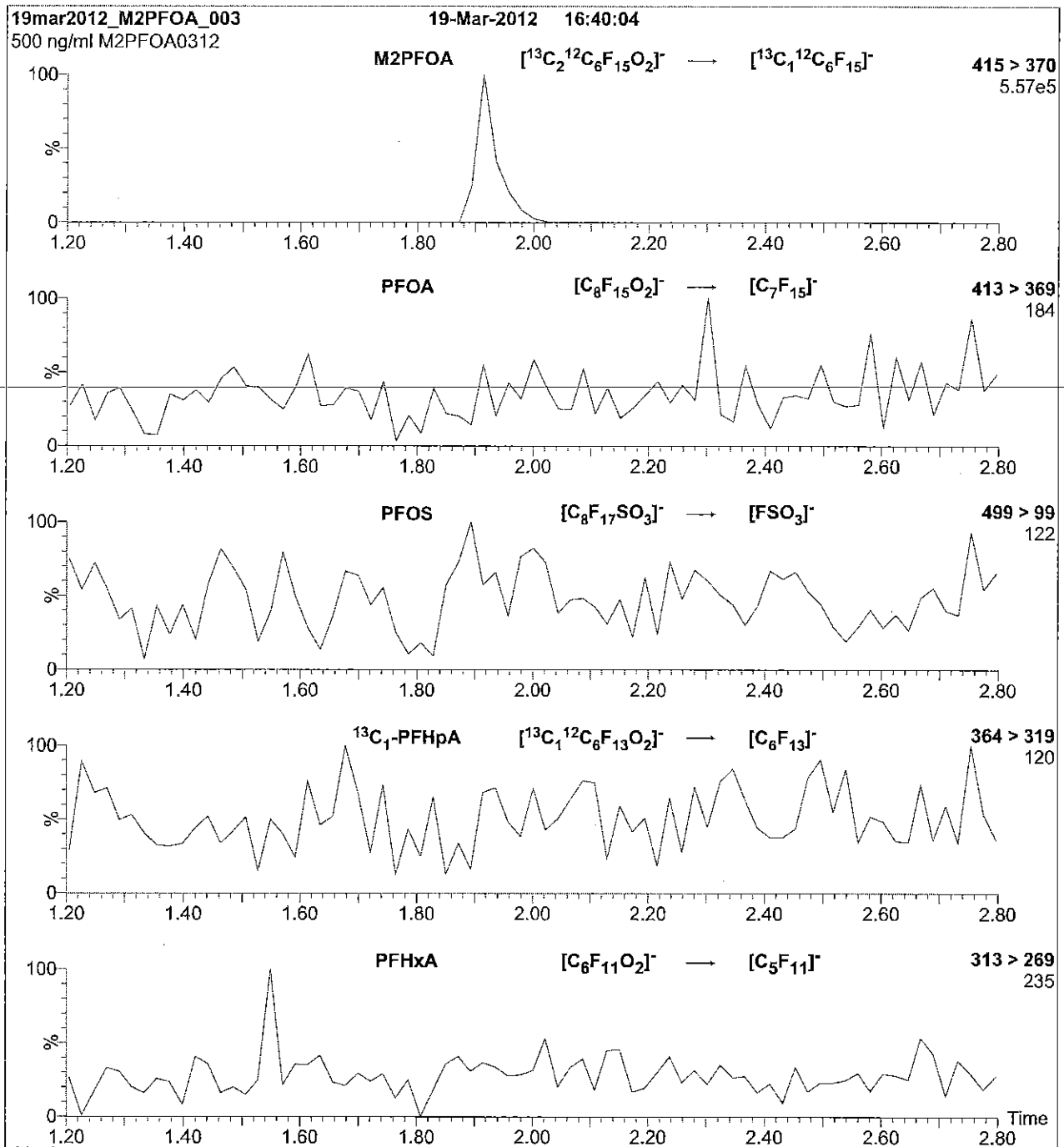

**Conditions for Figure 2:**

Injection: Direct loop injection  
10  $\mu\text{l}$  (500 ng/ml M2PFOA)

Mobile phase: Isocratic 70% (80:20 MeOH:ACN) / 30%  $\text{H}_2\text{O}$   
(both with 10 mM  $\text{NH}_4\text{OAc}$  buffer)

Flow: 300  $\mu\text{l}/\text{min}$


**MS Parameters**

Collision Gas (mbar) =  $3.35\text{e-}3$   
Collision Energy (eV) = 11

Reagent

---


**LCM2PFOA\_00004**

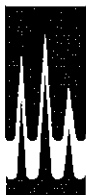

# WELLINGTON LABORATORIES

## CERTIFICATE OF ANALYSIS DOCUMENTATION

**PRODUCT CODE:** M2PFOA **LOT NUMBER:** M2PFOA0312  
**COMPOUND:** Perfluoro-n-[1,2-<sup>13</sup>C<sub>2</sub>]octanoic acid  
**STRUCTURE:** **CAS #:** Not available

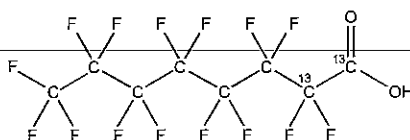

**MOLECULAR FORMULA:** <sup>13</sup>C<sub>2</sub><sup>12</sup>C<sub>6</sub>HF<sub>15</sub>O<sub>2</sub> **MOLECULAR WEIGHT:** 416.05  
**CONCENTRATION:** 50 ± 2.5 µg/ml **SOLVENT(S):** Methanol  
Water (<1%)  
**CHEMICAL PURITY:** >98% **ISOTOPIC PURITY:** ≥99%<sup>13</sup>C  
(1,2-<sup>13</sup>C<sub>2</sub>)  
**LAST TESTED:** (mm/dd/yyyy) 03/19/2012  
**EXPIRY DATE:** (mm/dd/yyyy) 03/19/2017  
**RECOMMENDED STORAGE:** Store ampoule in a cool, dark place

### DOCUMENTATION/ DATA ATTACHED:

Figure 1: LC/MS Data (TIC and Mass Spectrum)  
Figure 2: LC/MS/MS Data (Selected MRM Transitions)

### ADDITIONAL INFORMATION:

- See page 2 for further details.
- Contains 4 mole eq. of NaOH to prevent conversion of the carboxylic acid to the methyl ester.

**FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE**

Certified By:

B.G. Chittim

Date: 01/09/2013  
(mm/dd/yyyy)

Wellington Laboratories Inc., 345 Southgate Dr. Guelph ON N1G 3M5 CANADA  
519-822-2436 • Fax: 519-822-2849 • info@well-labs.com



#### **INTENDED USE:**

The products prepared by Wellington Laboratories Inc. are for laboratory use only. They are designed to be used as reference standards for the identification and/or quantification of specific chemical compound(s).

#### **HAZARDS:**

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Material Safety Data Sheets (MSDSs) are available upon request.

#### **SYNTHESIS / CHARACTERIZATION:**

Where possible, all of our products are synthesized using single-product, unambiguous routes. They are then characterized, and their structures and purities confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, x-ray crystallography and melting point. Isotopic purities of mass-labelled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

#### **HOMOGENEITY:**

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given diluent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS and/or LC/MS/MS. The relative response factors of the analyte of interest in each solution are required to be <5% RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers.

#### **UNCERTAINTY:**

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty,  $u_c(y)$ , of a value  $y$  and the uncertainty of the independent parameters  $x_1, x_2, \dots, x_n$  on which it depends is:

$$u_c(y(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{j=1}^n u(y, x_j)^2}$$

where  $x$  is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of  $\pm 5\%$  (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all our products.

#### **TRACEABILITY:**

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external, ISO/IEC 17025:2005 accredited calibration company. In addition, their calibration is verified prior to each weighing using NIST and/or NRC traceable external weights. All volumetric glassware used is of Class A tolerance and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to international interlaboratory studies has also been established.

#### **EXPIRY DATE / PERIOD OF VALIDITY:**

Ongoing stability studies of this product have demonstrated stability in its composition and concentration for the period of time specified by the expiry date in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

#### **LIMITED WARRANTY:**

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.



#### **QUALITY MANAGEMENT:**

This product was produced using a Quality Management System registered to ISO 9001:2008 by SAI Global, ISO/IEC 17025:2005 by the Canadian Association for Laboratory Accreditation Inc. (CALA; A 1226), and ISO GUIDE 34:2009 by ACLASS (certificate number AR-1523).



\*\*For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at [www.well-labs.com](http://www.well-labs.com) or contact us directly at [info@well-labs.com](mailto:info@well-labs.com)\*\*

**Figure 1: M2PFOA; LC/MS Data (TIC and Mass Spectrum)**

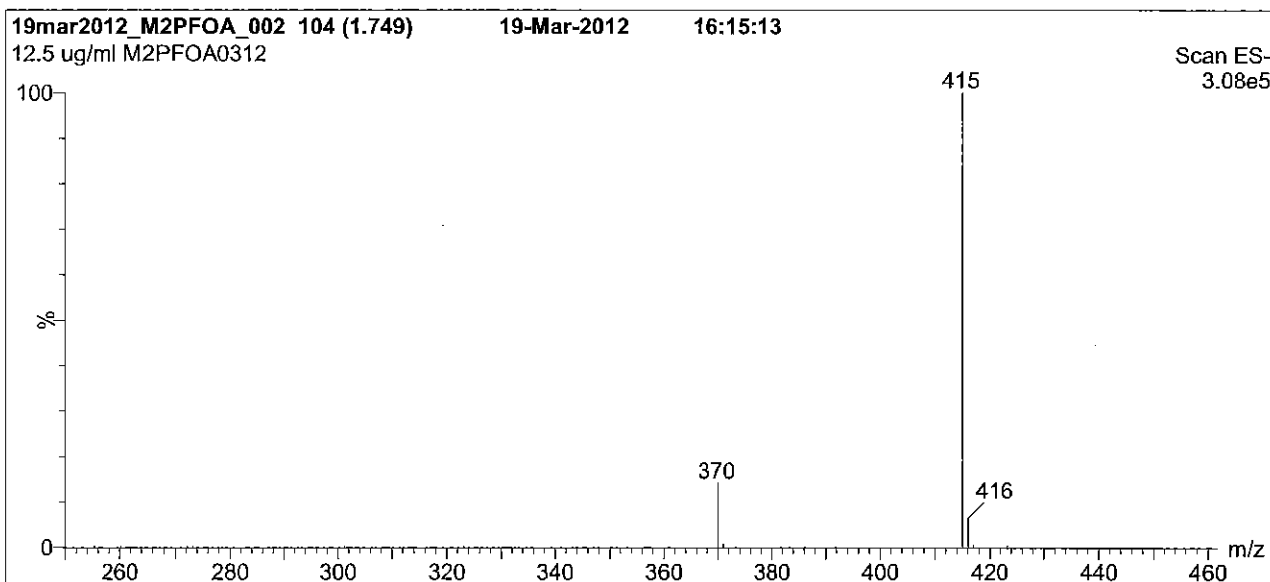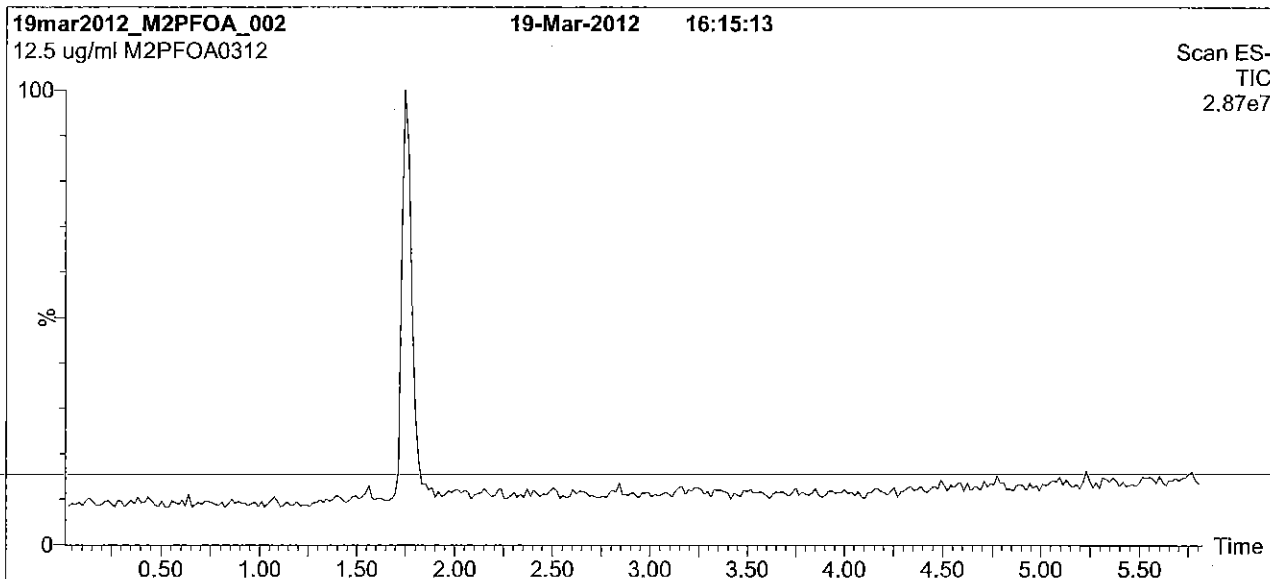

**Conditions for Figure 1:**

**LC:** Waters Acquity Ultra Performance LC  
**MS:** Micromass Quattro *micro* API MS

**Chromatographic Conditions**

**Column:** Acquity UPLC BEH Shield RP<sub>18</sub>  
1.7  $\mu$ m, 2.1 x 100 mm

**Mobile phase:** Gradient  
Start: 60% (80:20 MeOH:ACN) / 40% H<sub>2</sub>O  
(both with 10 mM NH<sub>4</sub>OAc buffer)  
Ramp to 90% organic over 6.5 min and hold for 2 min  
before returning to initial conditions in 0.5 min.  
Time: 10 min


**Flow:** 300  $\mu$ l/min

**MS Parameters**

**Experiment:** Full Scan (250 - 850 amu)

**Source:** Electrospray (negative)  
Capillary Voltage (kV) = 2.00  
Cone Voltage (V) = 15.00  
Cone Gas Flow (l/hr) = 100  
Desolvation Gas Flow (l/hr) = 750

**Figure 2: M2PFOA; LC/MS/MS Data (Selected MRM Transitions)**

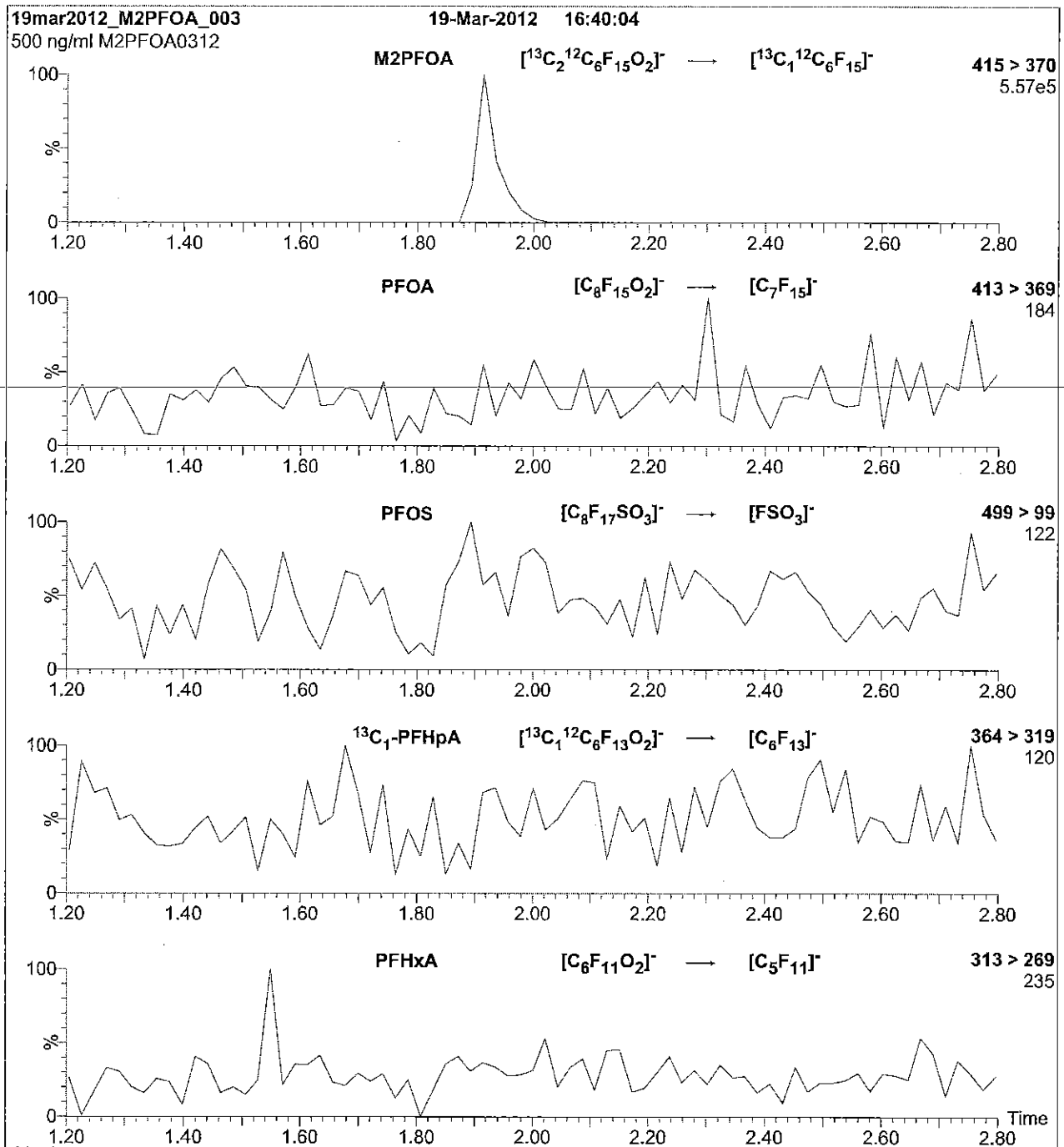

**Conditions for Figure 2:**

Injection: Direct loop injection  
10  $\mu\text{l}$  (500 ng/ml M2PFOA)

Mobile phase: Isocratic 70% (80:20 MeOH:ACN) / 30%  $\text{H}_2\text{O}$   
(both with 10 mM  $\text{NH}_4\text{OAc}$  buffer)

Flow: 300  $\mu\text{l}/\text{min}$

**MS Parameters**

Collision Gas (mbar) =  $3.35\text{e-}3$   
Collision Energy (eV) = 11

Reagent

---

**LCMPFDA\_00008**



605243

ID: LCMPFDA\_00008

Exp: 08/19/20 Pptd: CBW

13C2-Perfluorodecanoic acid

Rec. 3/29/16 JEB ✓



**WELLINGTON**  
LABORATORIES

**CERTIFICATE OF ANALYSIS**  
DOCUMENTATION

**PRODUCT CODE:**


MPFDA

**LOT NUMBER:**

MPFDA0815

**COMPOUND:**Perfluoro-n-[1,2-<sup>13</sup>C<sub>2</sub>]decanoic acid**STRUCTURE:****CAS #:**

Not available

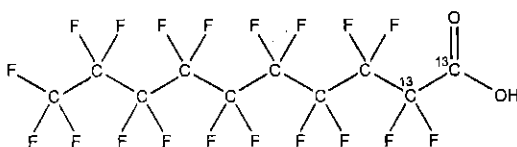**MOLECULAR FORMULA:**<sup>13</sup>C<sub>2</sub><sup>12</sup>C<sub>8</sub>HF<sub>19</sub>O<sub>2</sub>**MOLECULAR WEIGHT:**

516.07

**CONCENTRATION:**

50 ± 2.5 µg/ml

**SOLVENT(S):**

Methanol

Water (&lt;1%)

**CHEMICAL PURITY:**

&gt;98%

**ISOTOPIC PURITY:**≥99% <sup>13</sup>C**LAST TESTED:** (mm/dd/yyyy)

08/19/2015

(1,2-<sup>13</sup>C<sub>2</sub>)**EXPIRY DATE:** (mm/dd/yyyy)

08/19/2020

**RECOMMENDED STORAGE:**

Store ampoule in a cool, dark place

**DOCUMENTATION/ DATA ATTACHED:**

Figure 1: LC/MS Data (TIC and Mass Spectrum)

Figure 2: LC/MS/MS Data (Selected MRM Transitions)

**ADDITIONAL INFORMATION:**

- See page 2 for further details.
- Contains 4 mole eq. of NaOH to prevent conversion of the carboxylic acid to the methyl ester.
- Contains < 0.1% of <sup>13</sup>C<sub>1</sub>-PFNA.

**FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE**

Certified By:

B.G. Chittim

Date: 08/21/2015

(mm/dd/yyyy)

Wellington Laboratories Inc., 345 Southgate Dr. Guelph ON N1G 3M5 CANADA  
519-822-2436 • Fax: 519-822-2849 • info@well-labs.com

### **INTENDED USE:**

The products prepared by Wellington Laboratories Inc. are for laboratory use only. This certified reference material (CRM) was designed to be used as a standard for the identification and/or quantification of the specific chemical compound it contains.

### **HAZARDS:**

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection, and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Safety Data Sheets (SDSs) are available upon request.

### **SYNTHESIS / CHARACTERIZATION:**

Where possible, all of our products are synthesized using single-product unambiguous routes. They are then characterized, and their structures and purities confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, SFC/UV/MS/MS, x-ray crystallography, and melting point. Isotopic purities of mass-labelled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

### **HOMOGENEITY:**

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given diluent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS, LC/MS/MS and/or SFC/UV/MS/MS. The relative response factors of the analyte of interest in each solution are required to be <5% RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers.

### **UNCERTAINTY:**

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty,  $u_c(y)$ , of a value  $y$  and the uncertainty of the independent parameters  $x_1, x_2, \dots, x_n$  on which it depends is:

$$u_c(y(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{i=1}^n u(y, x_i)^2}$$

where  $x$  is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of  $\pm 5\%$  (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all of our products.

### **TRACEABILITY:**

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external ISO/IEC 17025 accredited calibration company. In addition, their calibration is verified prior to each weighing using NIST and/or NRC traceable external weights. All volumetric glassware used is of Class A tolerance and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to international interlaboratory studies has also been established.

### **EXPIRY DATE / PERIOD OF VALIDITY:**


Ongoing stability studies of this product have demonstrated stability in its composition and concentration, until the specified expiry date, in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

### **LIMITED WARRANTY:**

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.



### **QUALITY MANAGEMENT:**

This product was produced using a Quality Management System registered to the latest versions of ISO 9001 by SAI Global, ISO/IEC 17025 by the Canadian Association for Laboratory Accreditation Inc. (CALA; A 1226), and ISO GUIDE 34 by ANSI-ASQ National Accreditation Board (ANAB; AR-1523).

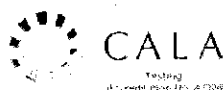

\*\*For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at [www.well-labs.com](http://www.well-labs.com) or contact us directly at [info@well-labs.com](mailto:info@well-labs.com)\*\*

**Figure 1: MPFDA; LC/MS Data (TIC and Mass Spectrum)**

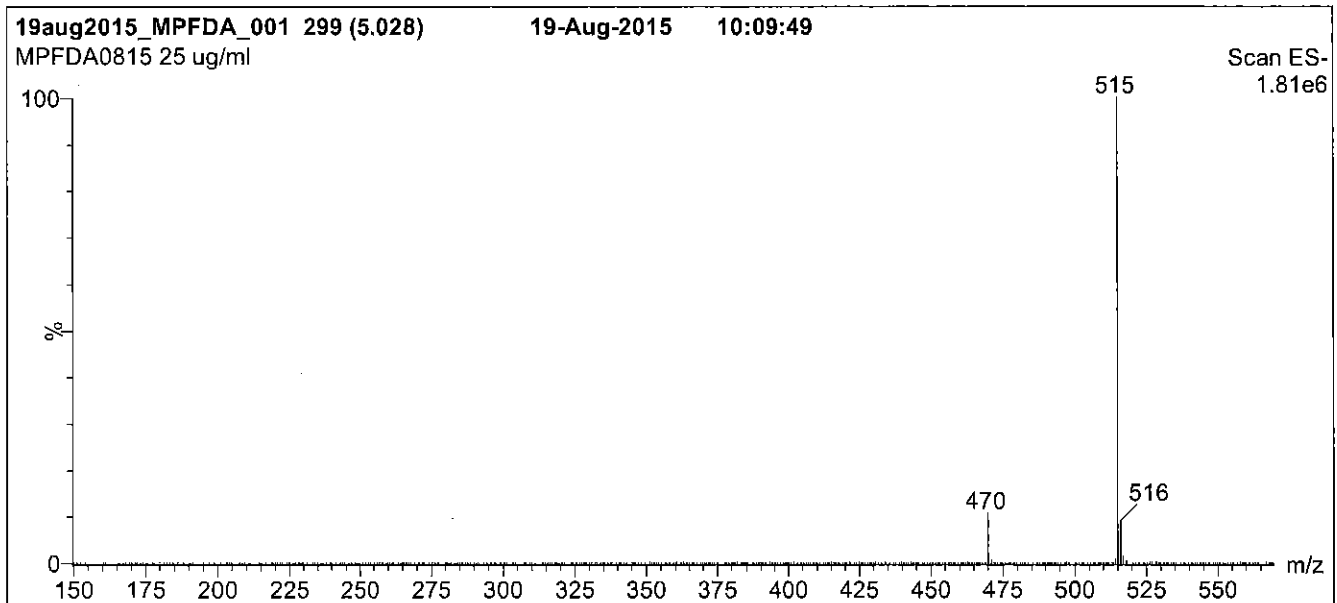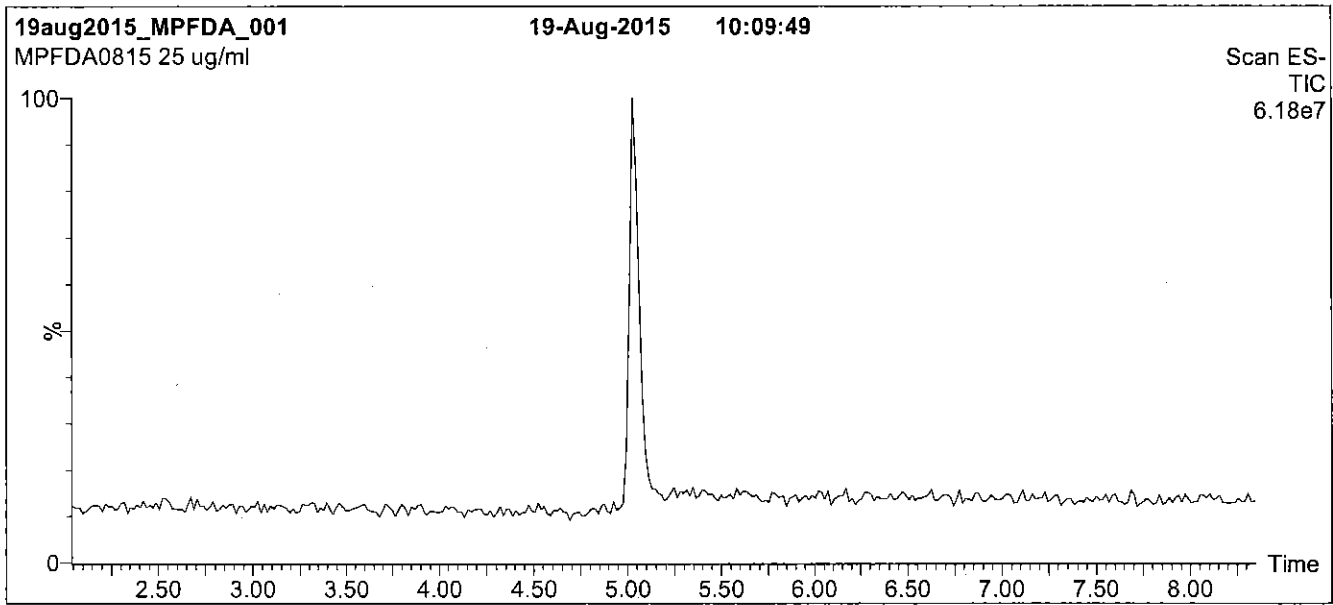

**Conditions for Figure 1:**

**LC:** Waters Acquity Ultra Performance LC  
**MS:** Micromass Quattro *micro* API MS

**Chromatographic Conditions**

Column: Acquity UPLC BEH Shield RP<sub>18</sub>  
1.7  $\mu$ m, 2.1 x 100 mm

Mobile phase: Gradient  
Start: 50% (80:20 MeOH:ACN) / 50% H<sub>2</sub>O  
(both with 10 mM NH<sub>4</sub>OAc buffer)  
Ramp to 90% organic over 7 min and hold for 2 min  
before returning to initial conditions in 0.5 min.  
Time: 10 min


Flow: 300  $\mu$ l/min

**MS Parameters**

Experiment: Full Scan (150 - 850 amu)

Source: Electrospray (negative)  
Capillary Voltage (kV) = 2.00  
Cone Voltage (V) = 15.00  
Cone Gas Flow (l/hr) = 50  
Desolvation Gas Flow (l/hr) = 750

**Figure 2: MPFDA; LC/MS/MS Data (Selected MRM Transitions)**

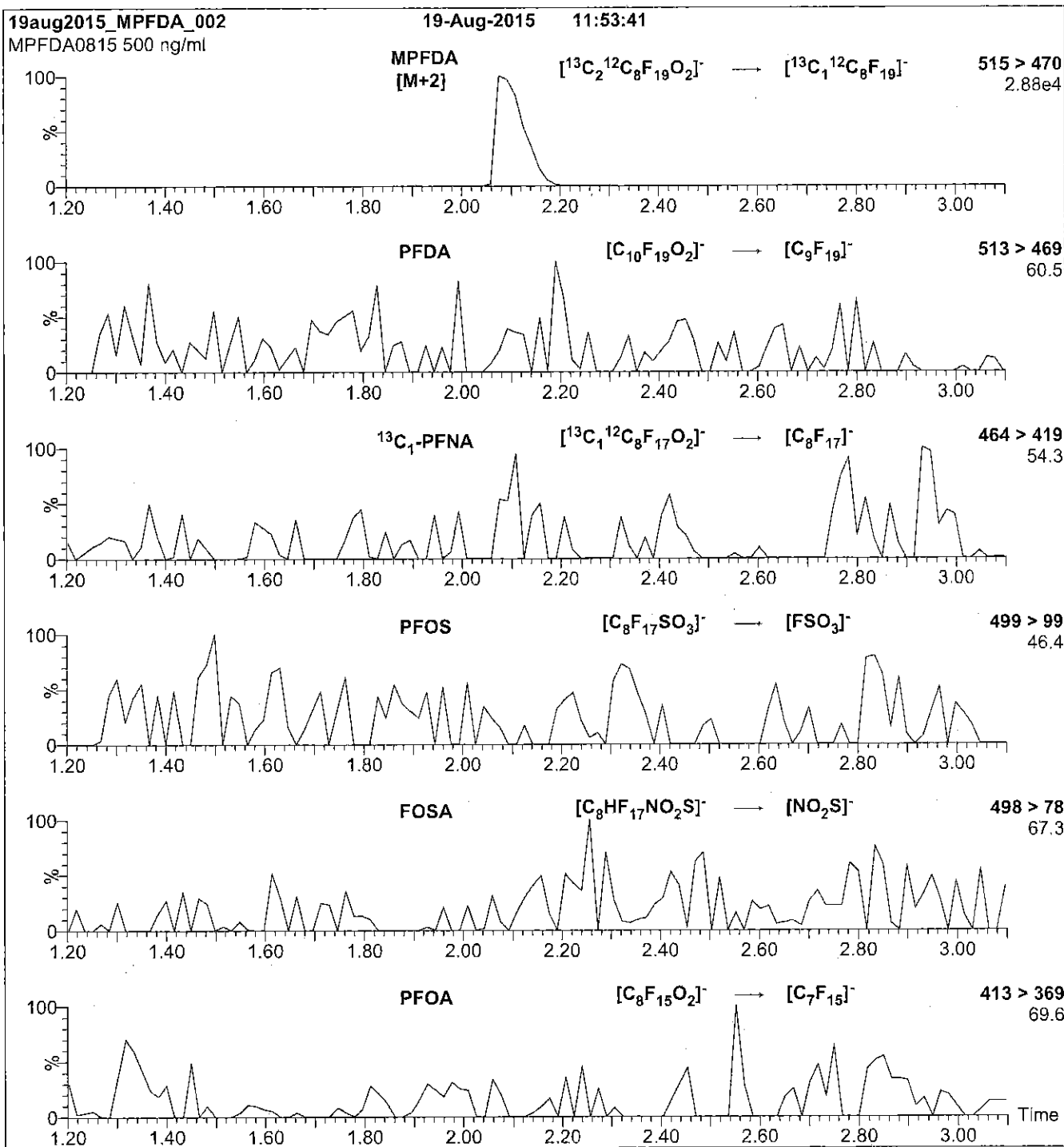

**Conditions for Figure 2:**

Injection: Direct loop injection  
 10  $\mu$ l (500 ng/ml MPFDA)

Mobile phase: Isocratic 80% (80:20 MeOH:ACN) / 20% H<sub>2</sub>O  
 (both with 10 mM NH<sub>4</sub>OAc buffer)

Flow: 300  $\mu$ l/min

**MS Parameters**

Collision Gas (mbar) = 3.35e-3  
 Collision Energy (eV) = 13



Reagent

---

**LCMPFHxA\_00009**



605244  
 ID: LCMPFHxA\_00009  
 Exp: 04/09/20 Prpd: CBW  
 13C2-Perfluorohexanoic ac

Rec. 3/29/16 JRB ✓



# WELLINGTON LABORATORIES


## CERTIFICATE OF ANALYSIS DOCUMENTATION

**PRODUCT CODE:** MPFHxA  
**COMPOUND:** Perfluoro-n-[1,2-<sup>13</sup>C<sub>2</sub>]hexanoic acid

**LOT NUMBER:** MPFHxA0415

**STRUCTURE:**

**CAS #:** Not available

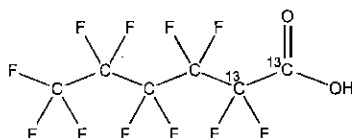

**MOLECULAR FORMULA:** <sup>13</sup>C<sub>2</sub><sup>12</sup>C<sub>4</sub>HF<sub>11</sub>O<sub>2</sub>  
**CONCENTRATION:** 50 ± 2.5 µg/ml

**MOLECULAR WEIGHT:** 316.04  
**SOLVENT(S):** Methanol  
 Water (<1%)

**CHEMICAL PURITY:** >98%  
**LAST TESTED:** (mm/dd/yyyy) 04/09/2015  
**EXPIRY DATE:** (mm/dd/yyyy) 04/09/2020

**ISOTOPIC PURITY:** ≥99%<sup>13</sup>C  
 (1,2-<sup>13</sup>C<sub>2</sub>)

**RECOMMENDED STORAGE:** Store ampoule in a cool, dark place

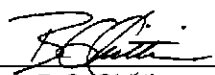
**DOCUMENTATION/ DATA ATTACHED:**

Figure 1: LC/MS Data (TIC and Mass Spectrum)  
 Figure 2: LC/MS/MS Data (Selected MRM Transitions)

**ADDITIONAL INFORMATION:**

- See page 2 for further details.
- Contains 4 mole eq. of NaOH to prevent conversion of the carboxylic acid to the methyl ester.
- Contains < 0.1% of perfluoro-n-hexanoic acid and ~ 0.3% of perfluoro-n-octanoic acid.

**FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE**

Certified By:   
 B.G. Chittim

Date: 04/14/2015  
 (mm/dd/yyyy)

Wellington Laboratories Inc., 345 Southgate Dr. Guelph ON N1G 3M5 CANADA  
 519-822-2436 • Fax: 519-822-2849 • info@well-labs.com

### **INTENDED USE:**

The products prepared by Wellington Laboratories Inc. are for laboratory use only. This certified reference material (CRM) was designed to be used as a standard for the identification and/or quantification of the specific chemical compound it contains.

### **HAZARDS:**

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection, and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Safety Data Sheets (SDSs) are available upon request.

### **SYNTHESIS / CHARACTERIZATION:**

Where possible, all of our products are synthesized using single-product unambiguous routes. They are then characterized, and their structures and purities confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, SFC/UV/MS/MS, x-ray crystallography, and melting point. Isotopic purities of mass-labelled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

### **HOMOGENEITY:**

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given diluent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS, LC/MS/MS and/or SFC/UV/MS/MS. The relative response factors of the analyte of interest in each solution are required to be <5% RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers.

### **UNCERTAINTY:**

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty,  $u_c(y)$ , of a value  $y$  and the uncertainty of the independent parameters  $x_1, x_2, \dots, x_n$  on which it depends is:

$$u_c(y(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{i=1}^n u(y, x_i)^2}$$

where  $x$  is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of  $\pm 5\%$  (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all of our products.

### **TRACEABILITY:**

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external ISO/IEC 17025 accredited calibration company. In addition, their calibration is verified prior to each weighing using NIST and/or NRC traceable external weights. All volumetric glassware used is of Class A tolerance and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to international interlaboratory studies has also been established.

### **EXPIRY DATE / PERIOD OF VALIDITY:**

Ongoing stability studies of this product have demonstrated stability in its composition and concentration, until the specified expiry date, in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

### **LIMITED WARRANTY:**

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.



### **QUALITY MANAGEMENT:**

This product was produced using a Quality Management System registered to the latest versions of ISO 9001 by SAI Global, ISO/IEC 17025 by the Canadian Association for Laboratory Accreditation Inc. (CALA; A 1226), and ISO GUIDE 34 by ANSI-ASQ National Accreditation Board (ANAB; AR-1523).



\*\*For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at [www.well-labs.com](http://www.well-labs.com) or contact us directly at [info@well-labs.com](mailto:info@well-labs.com)\*\*

**Figure 1: MPFHxA; LC/MS Data (TIC and Mass Spectrum)**

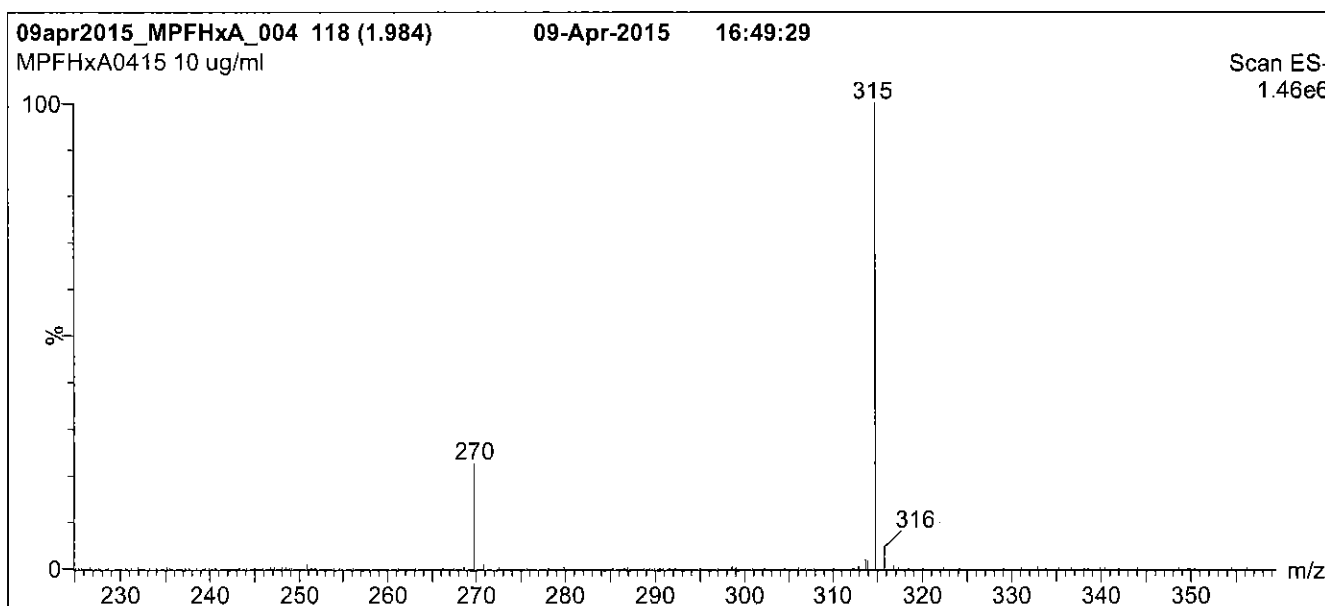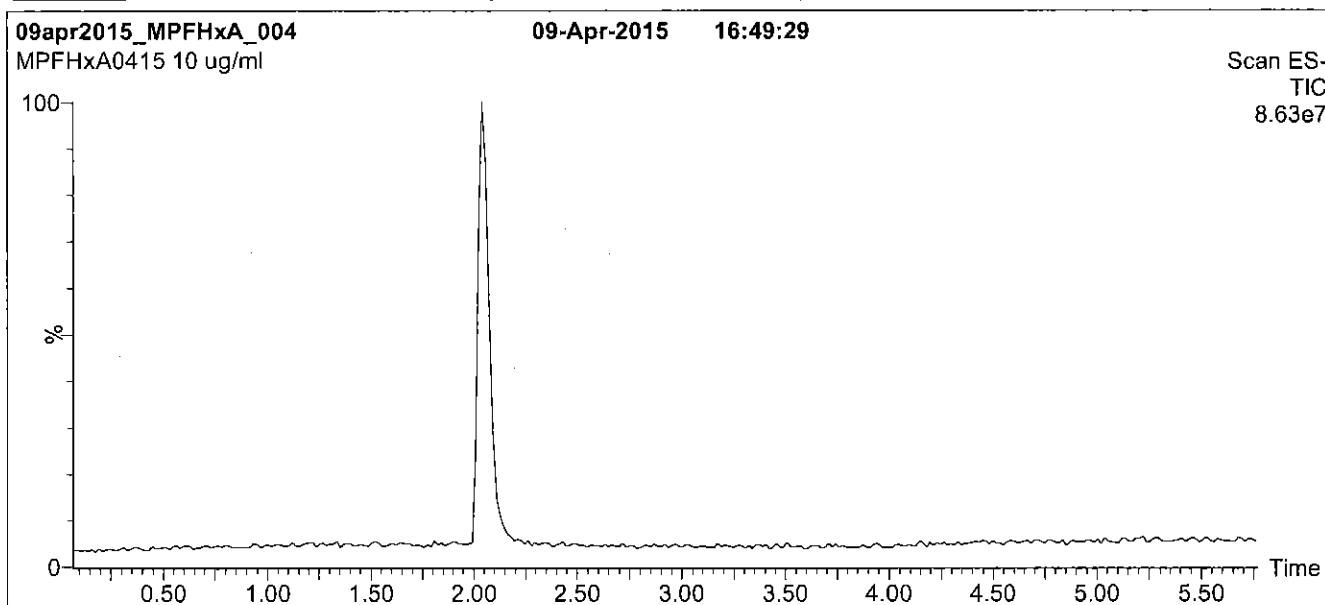

**Conditions for Figure 1:**

**LC:** Waters Acquity Ultra Performance LC  
**MS:** Micromass Quattro *micro* API MS

**Chromatographic Conditions**

Column: Acquity UPLC BEH Shield RP<sub>18</sub>  
1.7  $\mu$ m, 2.1 x 100 mm


Mobile phase: Gradient  
Start: 50% (80:20 MeOH:ACN) / 50% H<sub>2</sub>O  
(both with 10 mM NH<sub>4</sub>OAc buffer)  
Ramp to 90% organic over 7 min and hold for 2 min  
before returning to initial conditions over 0.5 min.  
Time: 10 min

Flow: 300  $\mu$ l/min

**MS Parameters**

Experiment: Full Scan (225 - 850 amu)  
Source: Electrospray (negative)  
Capillary Voltage (kV) = 2.00  
Cone Voltage (V) = 15.00  
Cone Gas Flow (l/hr) = 100  
Desolvation Gas Flow (l/hr) = 750

**Figure 2: MPFHxA; LC/MS/MS Data (Selected MRM Transitions)**

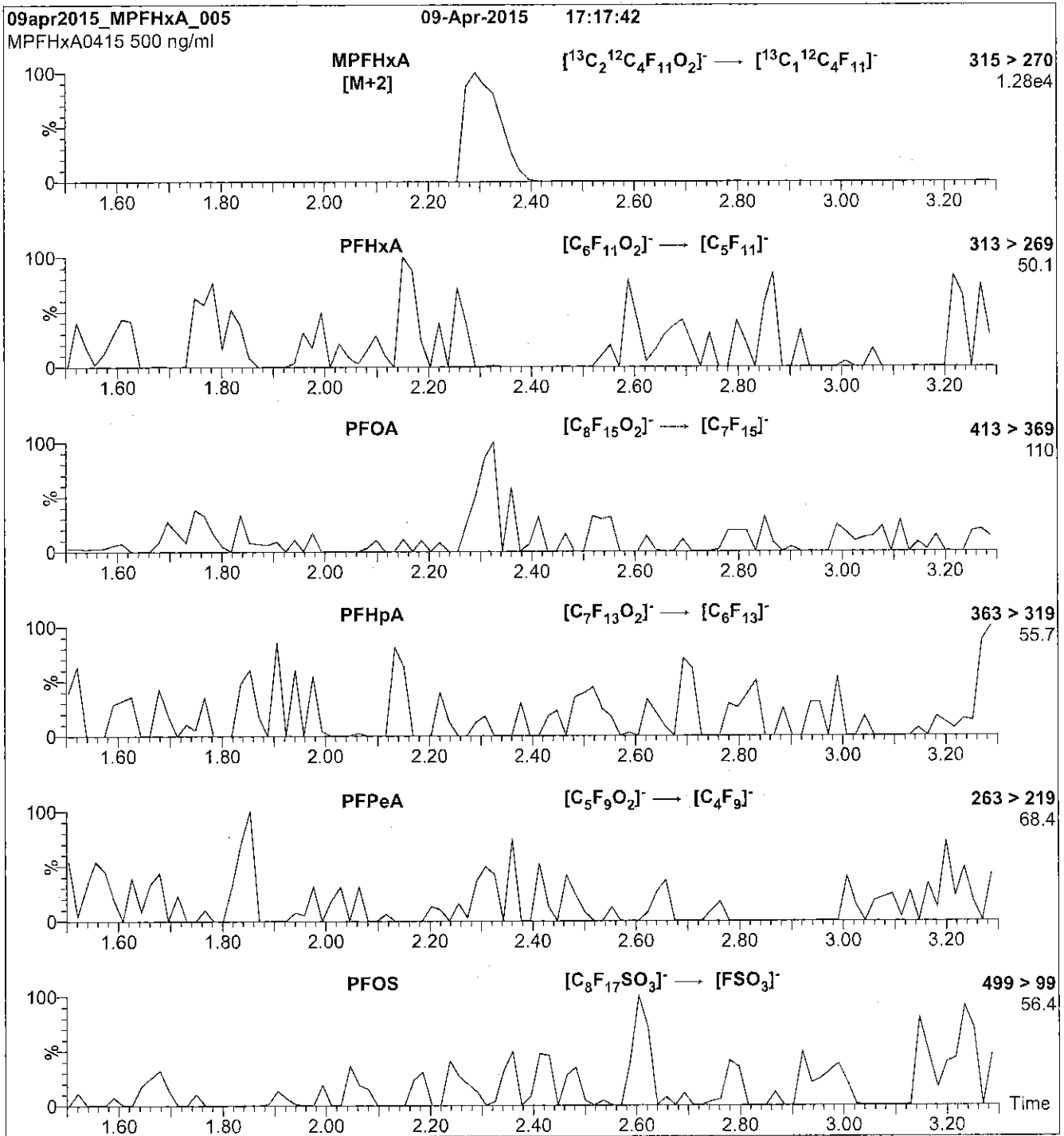

**Conditions for Figure 2:**

**Injection:** Direct loop injection  
10  $\mu$ l (500 ng/ml MPFHxA)

**Mobile phase:** Isocratic 80% (80:20 MeOH:ACN) / 20% H<sub>2</sub>O  
(both with 10 mM NH<sub>4</sub>OAc buffer)

**Flow:** 300  $\mu$ l/min

**MS Parameters**

Collision Gas (mbar) = 3.20e-3  
Collision Energy (eV) = 10

Reagent

---

**LCMPFOS\_00013**

605227  
ID: LCMFOS\_00012  
Exp: 01/22/21 Prpd: CBW  
13C4-Perfluorooctanesulfo

Rec 3/29/16 JRB ✓

606228  
ID: LCMFOS\_00013  
Exp: 01/22/21 Prpd: CBW  
13C4-Perfluorooctanesulfo




# WELLINGTON LABORATORIES

## CERTIFICATE OF ANALYSIS DOCUMENTATION

**PRODUCT CODE:** MPFOS **LOT NUMBER:** MPFOS0116  
**COMPOUND:** Sodium perfluoro-1-[1,2,3,4-<sup>13</sup>C<sub>4</sub>]octanesulfonate

**STRUCTURE:** **CAS #:** Not available

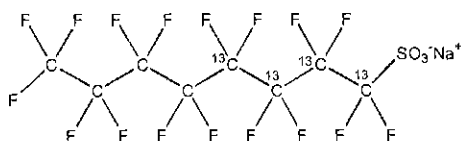

**MOLECULAR FORMULA:** <sup>13</sup>C<sub>4</sub><sup>12</sup>C<sub>4</sub>F<sub>17</sub>SO<sub>3</sub>Na **MOLECULAR WEIGHT:** 526.08  
**CONCENTRATION:** 50.0 ± 2.5 µg/ml (Na salt) **SOLVENT(S):** Methanol  
47.8 ± 2.4 µg/ml (MPFOS anion)  
**CHEMICAL PURITY:** >98% **ISOTOPIC PURITY:** ≥99% <sup>13</sup>C  
(1,2,3,4-<sup>13</sup>C<sub>4</sub>)  
**LAST TESTED:** (mm/dd/yyyy) 01/22/2016  
**EXPIRY DATE:** (mm/dd/yyyy) 01/22/2021  
**RECOMMENDED STORAGE:** Store ampoule in a cool, dark place

### DOCUMENTATION/ DATA ATTACHED:

Figure 1: LC/MS Data (TIC and Mass Spectrum)  
Figure 2: LC/MS/MS Data (Selected MRM Transitions)

### ADDITIONAL INFORMATION:

- See page 2 for further details.
- Contains ~ 0.8% Sodium perfluoro-1-[1,2,3-<sup>13</sup>C<sub>3</sub>]heptanesulfonate.

**FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE**

Certified By:

B.G. Chittim

Date: 02/01/2016  
(mm/dd/yyyy)

Wellington Laboratories Inc., 345 Southgate Dr. Guelph ON N1G 3M5 CANADA  
519-822-2436 • Fax: 519-822-2849 • info@well-labs.com

### **INTENDED USE:**

The products prepared by Wellington Laboratories Inc. are for laboratory use only. This certified reference material (CRM) was designed to be used as a standard for the identification and/or quantification of the specific chemical compound it contains.

### **HAZARDS:**

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection, and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Safety Data Sheets (SDSs) are available upon request.

### **SYNTHESIS / CHARACTERIZATION:**

Where possible, all of our products are synthesized using single-product unambiguous routes. They are then characterized, and their structures and purities confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, SFC/UV/MS/MS, x-ray crystallography, and melting point. Isotopic purities of mass-labelled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

### **HOMOGENEITY:**

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given diluent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS, LC/MS/MS and/or SFC/UV/MS/MS. The relative response factors of the analyte of interest in each solution are required to be <5% RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers.

### **UNCERTAINTY:**

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty,  $u_c(y)$ , of a value  $y$  and the uncertainty of the independent parameters  $x_1, x_2, \dots, x_n$  on which it depends is:

$$u_c(y(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{i=1}^n u(y, x_i)^2}$$

where  $x$  is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of  $\pm 5\%$  (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all of our products.

### **TRACEABILITY:**

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external ISO/IEC 17025 accredited calibration company. In addition, their calibration is verified prior to each weighing using NIST and/or NRC traceable external weights. All volumetric glassware used is of Class A tolerance and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to international interlaboratory studies has also been established.

### **EXPIRY DATE / PERIOD OF VALIDITY:**

Ongoing stability studies of this product have demonstrated stability in its composition and concentration, until the specified expiry date, in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

### **LIMITED WARRANTY:**

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.

### **QUALITY MANAGEMENT:**



This product was produced using a Quality Management System registered to the latest versions of ISO 9001 by SAI Global, ISO/IEC 17025 by the Canadian Association for Laboratory Accreditation Inc. (CALA; A 1226), and ISO GUIDE 34 by ANSI-ASQ National Accreditation Board (ANAB; AR-1523).



\*\*For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at [www.well-labs.com](http://www.well-labs.com) or contact us directly at [info@well-labs.com](mailto:info@well-labs.com)\*\*



**Figure 1: MPFOS; LC/MS Data (TIC and Mass Spectrum)**

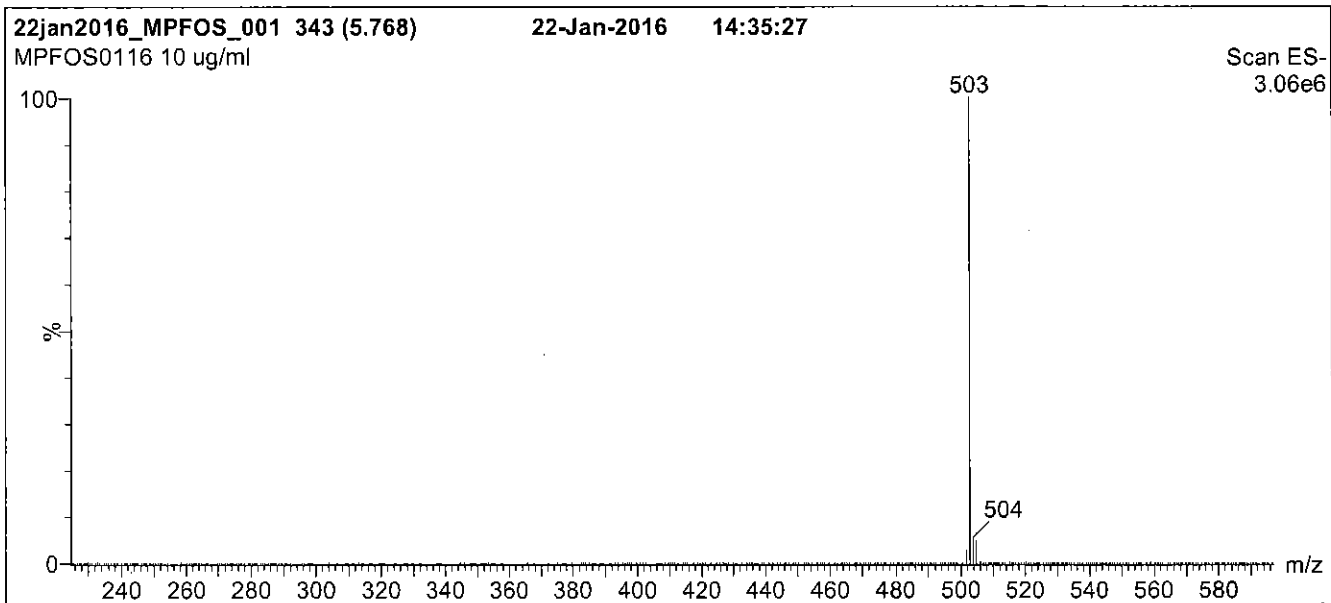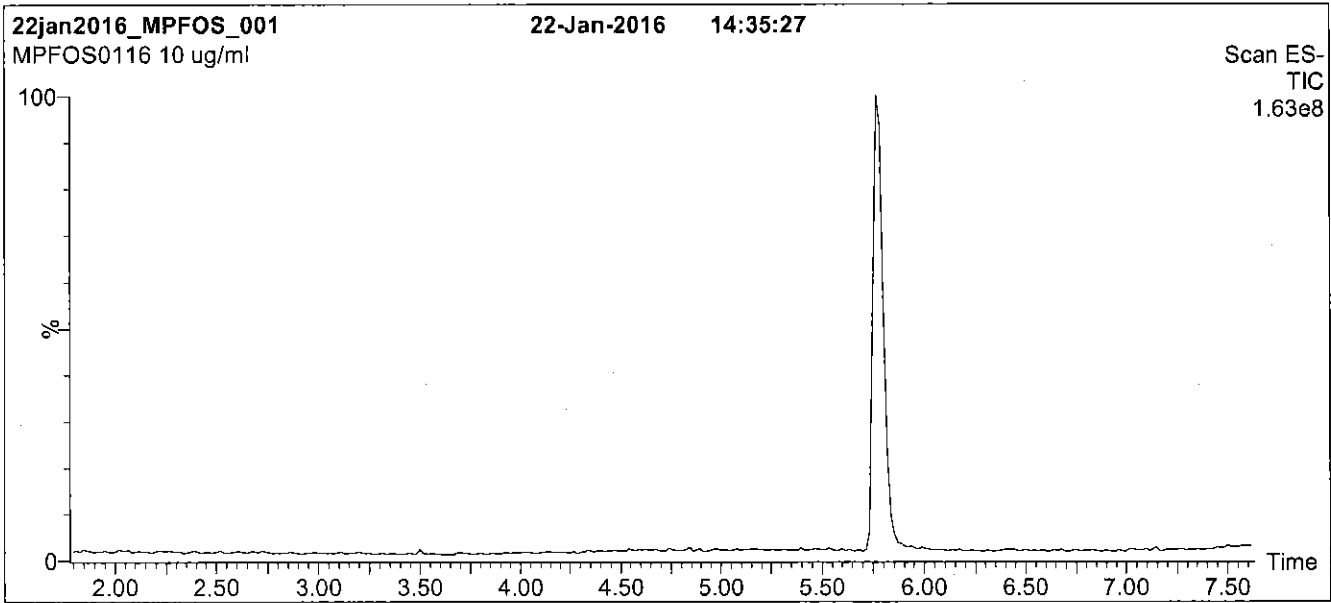

**Conditions for Figure 1:**

**LC:** Waters Acquity Ultra Performance LC  
**MS:** Micromass Quattro *micro* API MS

**Chromatographic Conditions**

**Column:** Acquity UPLC BEH Shield RP<sub>18</sub>  
1.7  $\mu$ m, 2.1 x 100 mm

**Mobile phase:** Gradient  
Start: 55% (80:20 MeOH:ACN) / 45% H<sub>2</sub>O  
(both with 10 mM NH<sub>4</sub>OAc buffer)  
Ramp to 90% organic over 7 min and hold for 2 min  
before returning to initial conditions in 0.5 min.  
Time: 10 min


**Flow:** 300  $\mu$ l/min

**MS Parameters**

Experiment: Full Scan (225 - 850 amu)

Source: Electrospray (negative)  
Capillary Voltage (kV) = 2.00  
Cone Voltage (V) = 60.00  
Cone Gas Flow (l/hr) = 50  
Desolvation Gas Flow (l/hr) = 750

**Figure 2: MPFOS; LC/MS/MS Data (Selected MRM Transitions)**

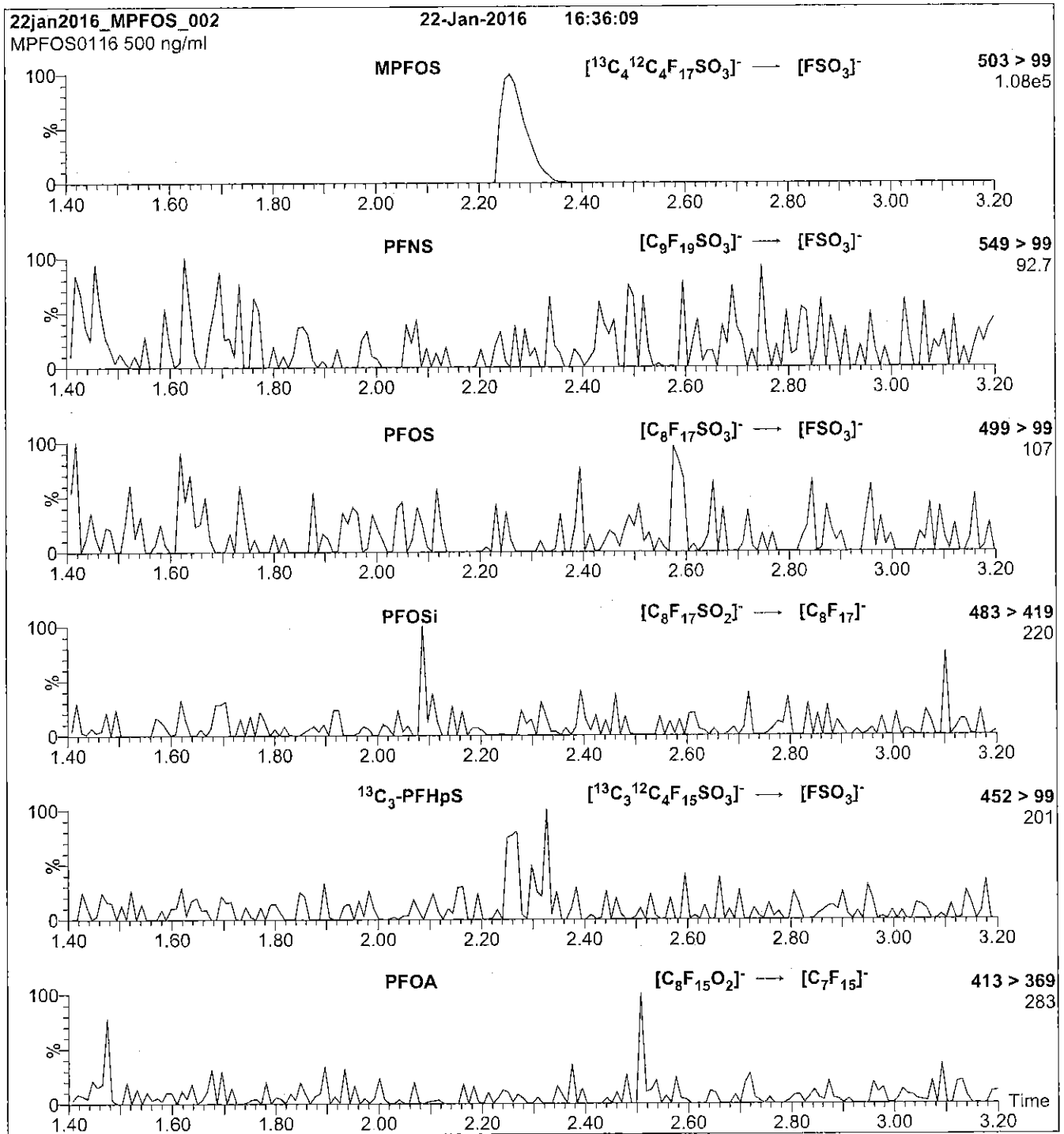

**Conditions for Figure 2:**

Injection: Direct loop injection  
10  $\mu\text{l}$  (500 ng/ml MPFOS)

Mobile phase: Isocratic 80% (80:20 MeOH:ACN) / 20%  $\text{H}_2\text{O}$   
(both with 10 mM  $\text{NH}_4\text{OAc}$  buffer)

Flow: 300  $\mu\text{l}/\text{min}$

**MS Parameters**

Collision Gas (mbar) = 3.70e-3  
Collision Energy (eV) = 40

Reagent

---

**LCMPFOS\_00018**

R: SBC 9/22/16



738686  
ID: LCMFOS\_00018  
Exp: 08/03/21 Papi: SBC  
13C4-Perfluorooctanesulfo




# WELLINGTON LABORATORIES

## CERTIFICATE OF ANALYSIS DOCUMENTATION

**PRODUCT CODE:** MPFOS **LOT NUMBER:** MPFOS0816  
**COMPOUND:** Sodium perfluoro-1-[1,2,3,4-<sup>13</sup>C]<sub>4</sub>octanesulfonate

**STRUCTURE:** **CAS #:** Not available

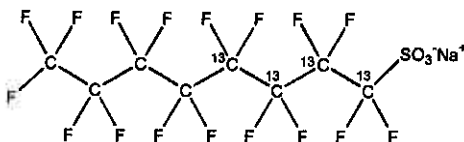

**MOLECULAR FORMULA:** <sup>13</sup>C<sub>4</sub><sup>12</sup>C<sub>4</sub>F<sub>17</sub>SO<sub>3</sub>Na **MOLECULAR WEIGHT:** 526.08  
**CONCENTRATION:** 50.0 ± 2.5 µg/ml (Na salt) **SOLVENT(S):** Methanol  
47.8 ± 2.4 µg/ml (MPFOS anion)  
**CHEMICAL PURITY:** >98% **ISOTOPIC PURITY:** ≥99% <sup>13</sup>C  
**LAST TESTED:** (mm/dd/yyyy) 08/03/2016 (1,2,3,4-<sup>13</sup>C<sub>4</sub>)  
**EXPIRY DATE:** (mm/dd/yyyy) 08/03/2021  
**RECOMMENDED STORAGE:** Store ampoule in a cool, dark place


### DOCUMENTATION/ DATA ATTACHED:

- Figure 1: LC/MS Data (TIC and Mass Spectrum)
- Figure 2: LC/MS/MS Data (Selected MRM Transitions)

### ADDITIONAL INFORMATION:

- See page 2 for further details.
- Contains ~ 0.8% Sodium perfluoro-1-[1,2,3-<sup>13</sup>C]<sub>3</sub>heptanesulfonate.

**FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE**

**Certified By:**   
B.G. Chittim **Date:** 08/05/2016  
(mm/dd/yyyy)

Wellington Laboratories Inc., 345 Southgate Dr. Guelph ON N1G 3M5 CANADA  
519-822-2436 • Fax: 519-822-2849 • info@well-labs.com

### **INTENDED USE:**

The products prepared by Wellington Laboratories Inc. are for laboratory use only. This certified reference material (CRM) was designed to be used as a standard for the identification and/or quantification of the specific chemical compound it contains.

### **HAZARDS:**

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection, and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Safety Data Sheets (SDSs) are available upon request.

### **SYNTHESIS / CHARACTERIZATION:**

Where possible, all of our products are synthesized using single-product unambiguous routes. They are then characterized, and their structures and purities confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, SFC/UV/MS/MS, x-ray crystallography, and melting point. Isotopic purities of mass-labelled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

### **HOMOGENEITY:**

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given diluent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS, LC/MS/MS and/or SFC/UV/MS/MS. The relative response factors of the analyte of interest in each solution are required to be <5% RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers.

### **UNCERTAINTY:**

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty,  $u_c(y)$ , of a value  $y$  and the uncertainty of the independent parameters  $x_1, x_2, \dots, x_n$  on which it depends is:

$$u_c(y(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{i=1}^n u(y, x_i)^2}$$

where  $x$  is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of  $\pm 5\%$  (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all of our products.

### **TRACEABILITY:**

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external ISO/IEC 17025 accredited calibration company. In addition, their calibration is verified prior to each weighing using NIST and/or NRC traceable external weights. All volumetric glassware used is of Class A tolerance and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to international interlaboratory studies has also been established.

### **EXPIRY DATE / PERIOD OF VALIDITY:**


Ongoing stability studies of this product have demonstrated stability in its composition and concentration, until the specified expiry date, in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

### **LIMITED WARRANTY:**

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.



### **QUALITY MANAGEMENT:**

This product was produced using a Quality Management System registered to the latest versions of ISO 9001 by SAI Global, ISO/IEC 17025 by the Canadian Association for Laboratory Accreditation Inc. (CALA; A 1226), and ISO GUIDE 34 by ANSI-ASQ National Accreditation Board (ANAB; AR-1523).

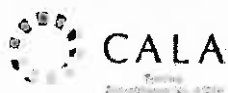

\*\*For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at [www.well-labs.com](http://www.well-labs.com) or contact us directly at [info@well-labs.com](mailto:info@well-labs.com)\*\*

**Figure 1: MPFOS; LC/MS Data (TIC and Mass Spectrum)**

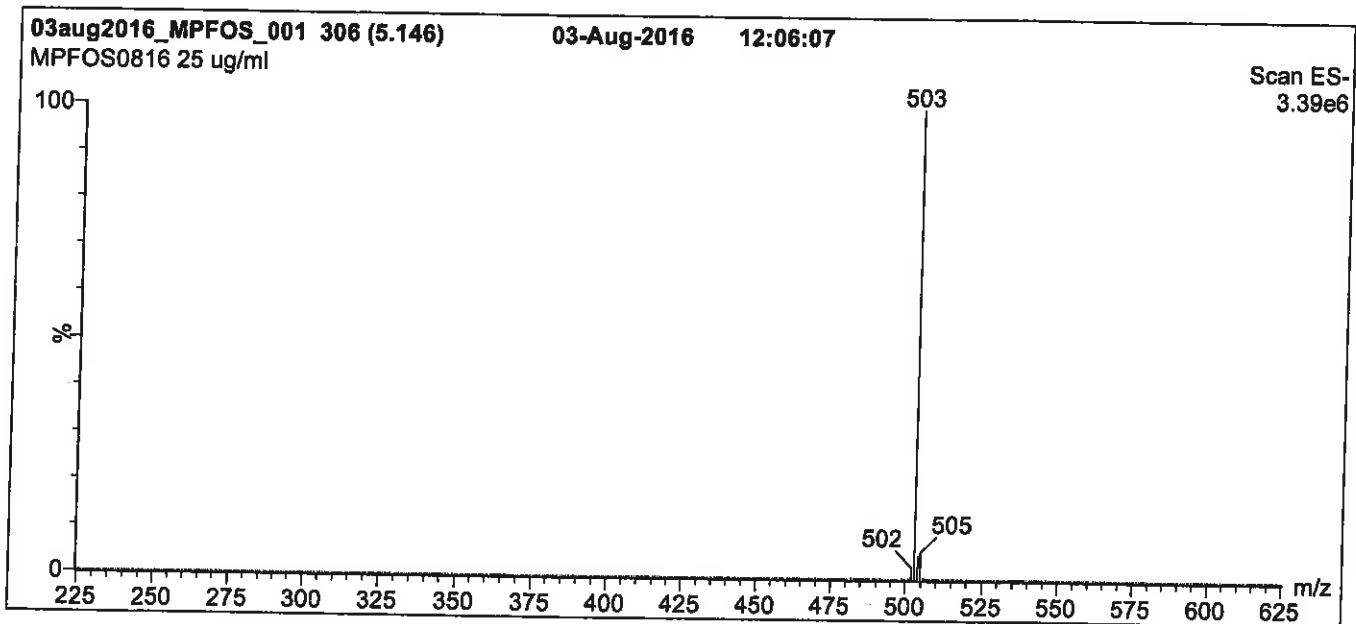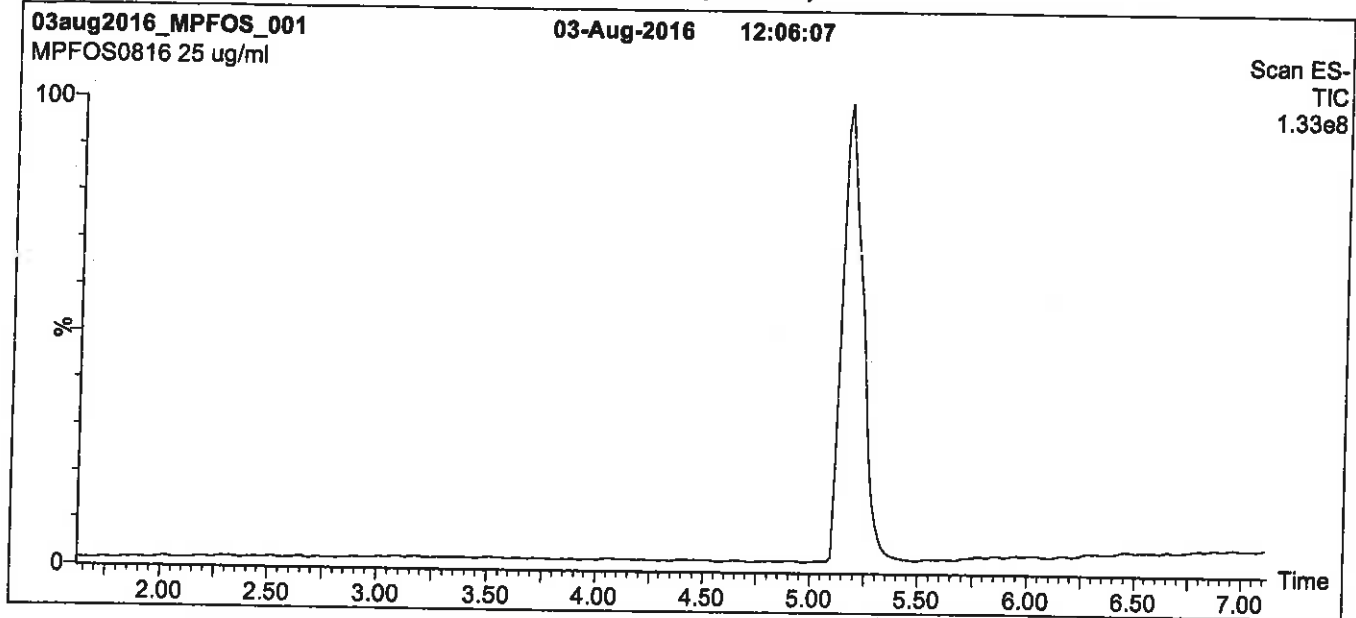

**Conditions for Figure 1:**

**LC:** Waters Acquity Ultra Performance LC  
**MS:** Micromass Quattro *micro* API MS

**Chromatographic Conditions**

**Column:** Acquity UPLC BEH Shield RP<sub>18</sub>  
 1.7  $\mu$ m, 2.1 x 100 mm

**Mobile phase:** Gradient  
 Start: 45% (80:20 MeOH:ACN) / 55% H<sub>2</sub>O  
 (both with 10 mM NH<sub>4</sub>OAc buffer)  
 Ramp to 90% organic over 7 min and hold for 1.5 min  
 before returning to initial conditions in 0.5 min.  
 Time: 10 min


**Flow:** 300  $\mu$ l/min

**MS Parameters**

**Experiment:** Full Scan (225 - 850 amu)

**Source:** Electrospray (negative)  
 Capillary Voltage (kV) = 3.00  
 Cone Voltage (V) = 60.00  
 Cone Gas Flow (l/hr) = 50  
 Desolvation Gas Flow (l/hr) = 750

**Figure 2: MPFOS; LC/MS/MS Data (Selected MRM Transitions)**

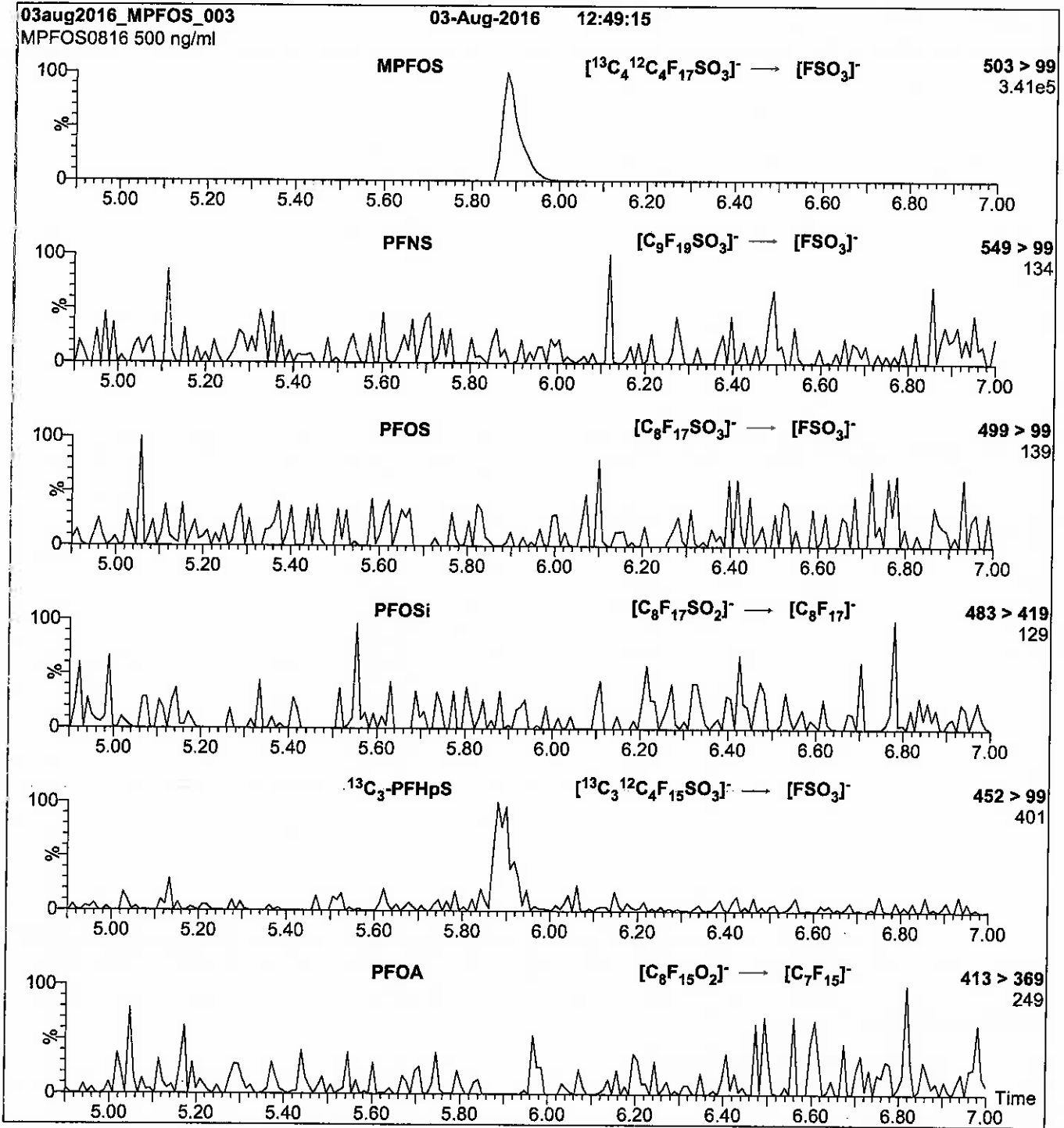

**Conditions for Figure 2:**

Injection: Direct loop injection  
10  $\mu\text{l}$  (500 ng/ml MPFOS)

Mobile phase: Isocratic 80% (80:20 MeOH:ACN) / 20%  $\text{H}_2\text{O}$   
(both with 10 mM  $\text{NH}_4\text{OAc}$  buffer)

Flow: 300  $\mu\text{l}/\text{min}$

**MS Parameters**

Collision Gas (mbar) = 3.46e-3  
Collision Energy (eV) = 40

# Method 537 DOD

---

Perfluorinated Alkyl Acids (LC/MS)  
by Method 537 DOD



FORM II  
LCMS SURROGATE RECOVERY

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Matrix: Water Level: Low

GC Column (1): Acquity ID: 2.1 (mm)

Client Sample ID	Lab Sample ID	PFHxA #	PFDA #
WI-AF-3RW30-1216	320-24224-1	100	117
WI-AF-3FB30-1216	320-24224-2	121	126
	MB 320-141642/1-A	110	111
	LCS 320-141642/2-A	118	123
WI-AF-3RW30-1216 MS	320-24224-1 MS	99	122
WI-AF-3RW30-1216 MSD	320-24224-1 MSD	96	113

PFHxA = 13C2 PFHxA  
PFDA = 13C2 PFDA

QC LIMITS  
70-130  
70-130

# Column to be used to flag recovery values

FORM III  
LCMS LAB CONTROL SAMPLE RECOVERY

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Matrix: Water Level: Low Lab File ID: 19DEC2016A6A\_011.d  
 Lab ID: LCS 320-141642/2-A Client ID: \_\_\_\_\_

COMPOUND	SPIKE ADDED (ug/L)	LCS CONCENTRATION (ug/L)	LCS % REC	QC LIMITS REC	#
Perfluorooctanesulfonic acid (PFOS)	0.160	0.150	94	70-130	
Perfluorooctanoic acid (PFOA)	0.0811	0.0704	87	70-130	
Perfluorobutanesulfonic acid (PFBS)	0.359	0.334	93	70-130	

# Column to be used to flag recovery and RPD values

FORM III  
LCMS MATRIX SPIKE RECOVERY

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Matrix: Water Level: Low Lab File ID: 15DEC2016A6A\_150.d  
 Lab ID: 320-24224-1 MS Client ID: WI-AF-3RW30-1216 MS

COMPOUND	SPIKE ADDED (ug/L)	SAMPLE CONCENTRATION (ug/L)	MS CONCENTRATION (ug/L)	MS % REC	QC LIMITS REC	#
Perfluorooctanesulfonic acid (PFOS)	0.0388	0.048 U	0.0334 J	86	70-130	
Perfluorooctanoic acid (PFOA)	0.0192	0.024 U	0.0189 J	98	70-130	
Perfluorobutanesulfonic acid (PFBS)	0.0870	0.11 U	0.0700 J	80	70-130	

# Column to be used to flag recovery and RPD values

FORM III  
LCMS MATRIX SPIKE DUPLICATE RECOVERY

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Matrix: Water Level: Low Lab File ID: 15DEC2016A6A\_151.d  
 Lab ID: 320-24224-1 MSD Client ID: WI-AF-3RW30-1216 MSD

COMPOUND	SPIKE ADDED (ug/L)	MSD CONCENTRATION (ug/L)	MSD % REC	% RPD	QC LIMITS		#
					RPD	REC	
Perfluorooctanesulfonic acid (PFOS)	0.0403	0.0337 J	84	1	30	70-130	
Perfluorooctanoic acid (PFOA)	0.0200	0.0165 J	83	14	30	70-130	
Perfluorobutanesulfonic acid (PFBS)	0.0903	0.0706 J	78	1	30	70-130	

# Column to be used to flag recovery and RPD values

FORM IV  
LCMS METHOD BLANK SUMMARY

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Lab File ID: 19DEC2016A6A\_010.d Lab Sample ID: MB 320-141642/1-A  
 Matrix: Water Date Extracted: 12/12/2016 10:03  
 Instrument ID: A6 Date Analyzed: 12/19/2016 13:19  
 Level: (Low/Med) Low

THIS METHOD BLANK APPLIES TO THE FOLLOWING SAMPLES:

CLIENT SAMPLE ID	LAB SAMPLE ID	LAB FILE ID	DATE ANALYZED
WI-AF-3RW30-1216	320-24224-1	15DEC2016A6 A 149.d	12/18/2016 09:01
WI-AF-3RW30-1216 MS	320-24224-1 MS	15DEC2016A6 A 150.d	12/18/2016 09:30
WI-AF-3RW30-1216 MSD	320-24224-1 MSD	15DEC2016A6 A 151.d	12/18/2016 10:00
WI-AF-3FB30-1216	320-24224-2	15DEC2016A6 A 152.d	12/18/2016 10:30
	LCS 320-141642/2-A	19DEC2016A6 A 011.d	12/19/2016 13:48

FORM VIII  
LCMS INTERNAL STANDARD AREA AND RETENTION TIME SUMMARY

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Instrument ID: A6 Calibration Start Date: 12/05/2016 17:26  
 GC Column: Acquity ID: 2.1(mm) Calibration End Date: 12/05/2016 19:54  
 Calibration ID: 26888

	13PFOA		PFOS		AREA #	RT #
	AREA #	RT #	AREA #	RT #		
INITIAL CALIBRATION MEAN AREA AND MEAN RT	965911	20.05	2046916	20.67		
UPPER LIMIT	1448867	20.55	3070374	21.17		
LOWER LIMIT	482956	19.55	1023458	20.17		
LAB SAMPLE ID	CLIENT SAMPLE ID					
CCV 320-140688/9 CCVL		1025187	20.05	2358079	20.67	
ICV 320-140688/11		877210	20.05	2015178	20.67	
CCV 320-142223/2 CCVL		615516	20.04	1606471	20.67	
CCV 320-142809/35 CCVIS		756954	20.00	1618248	20.62	
320-24224-1	WI-AF-3RW30-1216	591533	20.01	1579694	20.62	
320-24224-1 MS	WI-AF-3RW30-1216 MS	572717	19.99	1554853	20.60	
320-24224-1 MSD	WI-AF-3RW30-1216 MSD	580380	19.97	1580820	20.60	
320-24224-2	WI-AF-3FB30-1216	572204	19.97	1516952	20.60	
CCV 320-142809/47 CCVIS		974301	20.00	1988143	20.62	
CCV 320-142886/7 CCVIS		783659	20.01	1682080	20.62	
MB 320-141642/1-A		747991	20.01	2054845	20.63	
LCS 320-141642/2-A		715286	20.02	1813273	20.63	
CCV 320-142886/20 CCVIS		802443	20.01	1835426	20.62	

13PFOA = 13C2-PFOA

PFOS = 13C4 PFOS

Area Limit = 50%-150% of internal standard area  
 RT Limit = ± 0.5 minutes of internal standard RT

# Column used to flag values outside QC limits

FORM VIII  
LCMS INTERNAL STANDARD AREA AND RETENTION TIME SUMMARY

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Sample No.: CCV 320-142809/35 Date Analyzed: 12/18/2016 08:02  
 Instrument ID: A6 GC Column: Acquity ID: 2.1 (mm)  
 Lab File ID (Standard): 15DEC2016A6A\_147.d Heated Purge: (Y/N) N  
 Calibration ID: 26888

	13PFOA		PFOS		AREA #	RT #
	AREA #	RT #	AREA #	RT #		
12/24 HOUR STD	756954	20.00	1618248	20.62		
UPPER LIMIT	1059736	20.50	2265547	21.12		
LOWER LIMIT	529868	19.50	1132774	20.12		
LAB SAMPLE ID	CLIENT SAMPLE ID					
320-24224-1	WI-AF-3RW30-1216	591533	20.01	1579694	20.62	
320-24224-1 MS	WI-AF-3RW30-1216 MS	572717	19.99	1554853	20.60	
320-24224-1 MSD	WI-AF-3RW30-1216 MSD	580380	19.97	1580820	20.60	
320-24224-2	WI-AF-3FB30-1216	572204	19.97	1516952	20.60	

13PFOA = 13C2-PFOA  
 PFOS = 13C4 PFOS

Area Limit = 70%-140% of internal standard area  
 RT Limit = ± 0.5 minutes of internal standard RT

# Column used to flag values outside QC limits

FORM VIII  
LCMS INTERNAL STANDARD AREA AND RETENTION TIME SUMMARY

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Sample No.: CCV 320-142809/47 Date Analyzed: 12/18/2016 13:57  
 Instrument ID: A6 GC Column: Acquity ID: 2.1 (mm)  
 Lab File ID (Standard): 15DEC2016A6A\_159.d Heated Purge: (Y/N) N  
 Calibration ID: 26888

	13PFOA		PFOS		AREA #	RT #
	AREA #	RT #	AREA #	RT #		
12/24 HOUR STD	974301	20.00	1988143	20.62		
UPPER LIMIT	1364021	20.50	2783400	21.12		
LOWER LIMIT	682011	19.50	1391700	20.12		
LAB SAMPLE ID	CLIENT SAMPLE ID					
320-24224-1	WI-AF-3RW30-1216	591533	20.01	1579694	20.62	
320-24224-1 MS	WI-AF-3RW30-1216 MS	572717	19.99	1554853	20.60	
320-24224-1 MSD	WI-AF-3RW30-1216 MSD	580380	19.97	1580820	20.60	
320-24224-2	WI-AF-3FB30-1216	572204	19.97	1516952	20.60	

13PFOA = 13C2-PFOA  
 PFOS = 13C4 PFOS

Area Limit = 70%-140% of internal standard area  
 RT Limit = ± 0.5 minutes of internal standard RT

# Column used to flag values outside QC limits



FORM VIII  
LCMS INTERNAL STANDARD AREA AND RETENTION TIME SUMMARY

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Sample No.: CCV 320-142886/7 Date Analyzed: 12/19/2016 11:51  
 Instrument ID: A6 GC Column: Acquity ID: 2.1 (mm)  
 Lab File ID (Standard): 19DEC2016A6A\_007.d Heated Purge: (Y/N) N  
 Calibration ID: 26888

	13PFOA		PFOS		AREA #	RT #
	AREA #	RT #	AREA #	RT #		
12/24 HOUR STD	783659	20.01	1682080	20.62		
UPPER LIMIT	1097123	20.51	2354912	21.12		
LOWER LIMIT	548561	19.51	1177456	20.12		
LAB SAMPLE ID	CLIENT SAMPLE ID					
MB 320-141642/1-A	747991	20.01	2054845	20.63		
LCS 320-141642/2-A	715286	20.02	1813273	20.63		

13PFOA = 13C2-PFOA  
 PFOS = 13C4 PFOS

Area Limit = 70%-140% of internal standard area  
 RT Limit = ± 0.5 minutes of internal standard RT

# Column used to flag values outside QC limits

FORM VIII  
LCMS INTERNAL STANDARD AREA AND RETENTION TIME SUMMARY

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Sample No.: CCV 320-142886/20 Date Analyzed: 12/19/2016 18:32  
 Instrument ID: A6 GC Column: Acquity ID: 2.1 (mm)  
 Lab File ID (Standard): 19DEC2016A6A\_020.d Heated Purge: (Y/N) N  
 Calibration ID: 26888

	13PFOA		PFOS		AREA #	RT #
	AREA #	RT #	AREA #	RT #		
12/24 HOUR STD	802443	20.01	1835426	20.62		
UPPER LIMIT	1123420	20.51	2569596	21.12		
LOWER LIMIT	561710	19.51	1284798	20.12		
LAB SAMPLE ID	CLIENT SAMPLE ID					
MB 320-141642/1-A	747991	20.01	2054845	20.63		
LCS 320-141642/2-A	715286	20.02	1813273	20.63		

13PFOA = 13C2-PFOA  
 PFOS = 13C4 PFOS

Area Limit = 70%-140% of internal standard area  
 RT Limit = ± 0.5 minutes of internal standard RT

# Column used to flag values outside QC limits

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Client Sample ID: WI-AF-3RW30-1216 Lab Sample ID: 320-24224-1  
 Matrix: Water Lab File ID: 15DEC2016A6A\_149.d  
 Analysis Method: 537 Date Collected: 12/07/2016 14:45  
 Extraction Method: 537 Date Extracted: 12/12/2016 10:03  
 Sample wt/vol: 251.6(mL) Date Analyzed: 12/18/2016 09:01  
 Con. Extract Vol.: 1.0(mL) Dilution Factor: 1  
 Injection Volume: 10(uL) GC Column: Acquity ID: 2.1(mm)  
 % Moisture: \_\_\_\_\_ GPC Cleanup: (Y/N) N  
 Analysis Batch No.: 142809 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	LOQ	LOD	DL
1763-23-1	Perfluorooctanesulfonic acid (PFOS)	0.048	U	0.060	0.048	0.015
335-67-1	Perfluorooctanoic acid (PFOA)	0.024	U	0.030	0.024	0.0094
375-73-5	Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.14	0.11	0.047

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL00993	13C2 PFHxA	100		70-130
STL00996	13C2 PFDA	117		70-130

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_149.d  
 Lims ID: 320-24224-A-1-A  
 Client ID: WI-AF-3RW30-1216  
 Sample Type: Client  
 Inject. Date: 18-Dec-2016 09:01:16 ALS Bottle#: 44 Worklist Smp#: 37  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: 320-24224-a-1-a  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 19-Dec-2016 15:25:13 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK025

First Level Reviewer: westendorfc Date: 19-Dec-2016 09:11:36


Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
--------	----	--------	--------	--------	----------	--------------	-----	-------

\$ 2 13C2 PFHxA	315.0 > 270.0	18.557	18.548	0.009	1.000	690996	10.0	10132
4 Perfluoroheptanoic acid	363.0 > 319.0	19.344	19.344	0.0	1.000	868	0.0121	0.4
* 5 13C2-PFOA	415.0 > 370.0	20.011	19.999	0.012		591533	10.0	15144
* 8 13C4 PFOS	503.0 > 80.0	20.619	20.619	0.0		1579694	28.7	41517
\$ 10 13C2 PFDA	515.0 > 470.0	21.427	21.418	0.009	1.000	603897	11.7	19072


TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_149.d  
Injection Date: 18-Dec-2016 09:01:16 Instrument ID: A6  
Lims ID: 320-24224-A-1-A Lab Sample ID: 320-24224-1  
Client ID: WI-AF-3RW30-1216  
Operator ID: CBW ALS Bottle#: 44 Worklist Smp#: 37  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Method: 537\_A6 Limit Group: LC 537 ICAL


1 Perfluorobutanesulfonic acid (ND) \$ 2 13C2 PFHxA 3 Perfluorohexanesulfonic acid (ND)

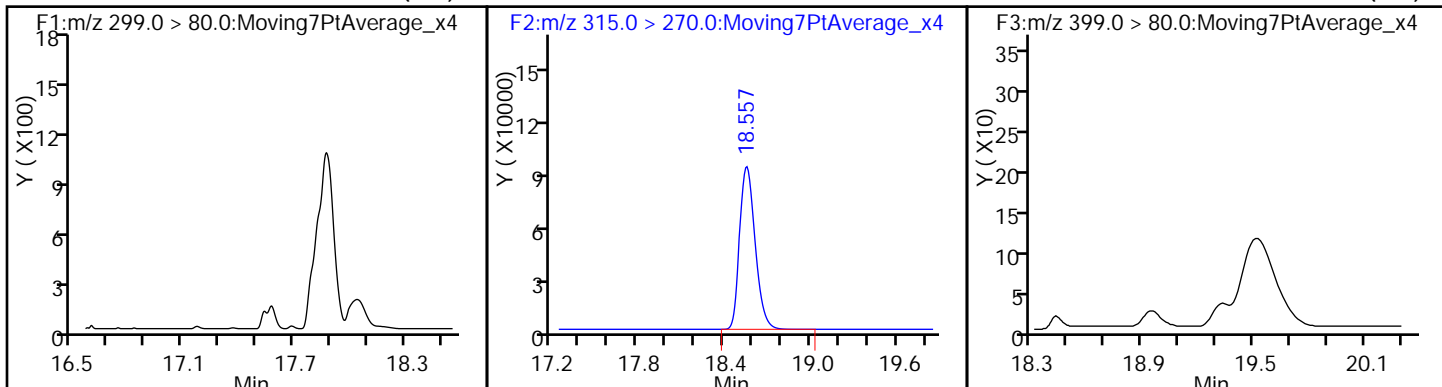


4 Perfluoroheptanoic acid \* 5 13C2-PFOA 6 Perfluorooctanoic acid (ND)

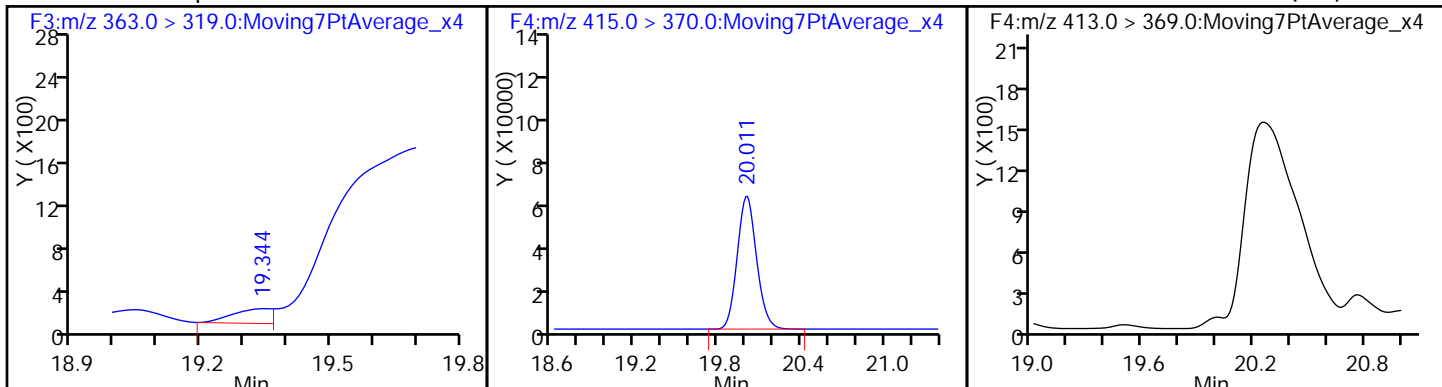

7 Perfluorooctane sulfonic acid (ND) \* 8 13C4 PFOS 9 Perfluorononanoic acid (ND)

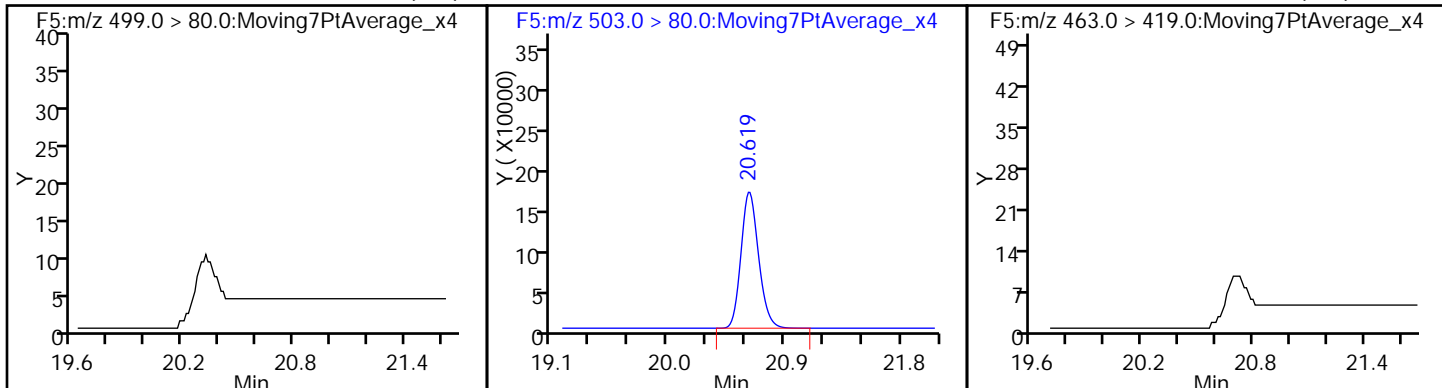

\$ 10 13C2 PFDA

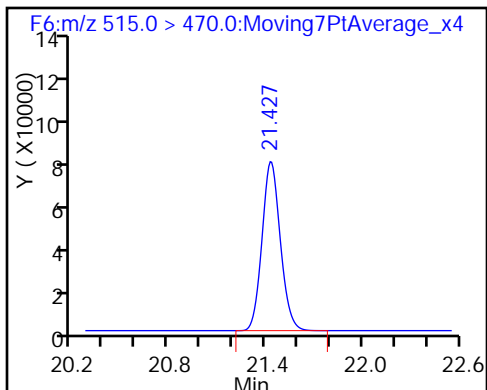

TestAmerica Sacramento  
Recovery Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_149.d  
 Lims ID: 320-24224-A-1-A  
 Client ID: WI-AF-3RW30-1216  
 Sample Type: Client  
 Inject. Date: 18-Dec-2016 09:01:16 ALS Bottle#: 44 Worklist Smp#: 37  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: 320-24224-a-1-a  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35°C  
 Operator ID: CBW Instrument ID: A6  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 19-Dec-2016 15:25:13 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK025

First Level Reviewer: westendorfc Date: 19-Dec-2016 09:11:36

Compound	Amount Added	Amount Recovered	% Rec.
\$ 2 13C2 PFHxA	10.0	10.0	100.14
\$ 10 13C2 PFDA	10.0	11.7	116.50

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Client Sample ID: WI-AF-3FB30-1216 Lab Sample ID: 320-24224-2  
 Matrix: Water Lab File ID: 15DEC2016A6A\_152.d  
 Analysis Method: 537 Date Collected: 12/07/2016 14:46  
 Extraction Method: 537 Date Extracted: 12/12/2016 10:03  
 Sample wt/vol: 257.5 (mL) Date Analyzed: 12/18/2016 10:30  
 Con. Extract Vol.: 1.0 (mL) Dilution Factor: 1  
 Injection Volume: 10 (uL) GC Column: Acquity ID: 2.1 (mm)  
 % Moisture: \_\_\_\_\_ GPC Cleanup: (Y/N) N  
 Analysis Batch No.: 142809 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	LOQ	LOD	DL
1763-23-1	Perfluorooctanesulfonic acid (PFOS)	0.047	U	0.058	0.047	0.015
335-67-1	Perfluorooctanoic acid (PFOA)	0.023	U	0.029	0.023	0.0091
375-73-5	Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.14	0.11	0.046

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL00993	13C2 PFHxA	121		70-130
STL00996	13C2 PFDA	126		70-130

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_152.d  
 Lims ID: 320-24224-A-2-A  
 Client ID: WI-AF-3FB30-1216  
 Sample Type: Client  
 Inject. Date: 18-Dec-2016 10:30:05 ALS Bottle#: 47 Worklist Smp#: 40  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: 320-24224-a-2-a  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 19-Dec-2016 15:25:13 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK025

First Level Reviewer: westendorfc Date: 19-Dec-2016 09:12:22

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
--------	----	--------	--------	--------	----------	--------------	-----	-------

\$ 2 13C2 PFHxA	315.0 > 270.0	18.530	18.548	-0.018	1.000	804794	12.1	26981
* 5 13C2-PFOA	415.0 > 370.0	19.972	19.999	-0.027		572204	10.0	14366
* 8 13C4 PFOS	503.0 > 80.0	20.596	20.619	-0.023		1516952	28.7	39684
\$ 10 13C2 PFDA	515.0 > 470.0	21.400	21.418	-0.018	1.000	631176	12.6	19893






TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_152.d  
Injection Date: 18-Dec-2016 10:30:05 Instrument ID: A6  
Lims ID: 320-24224-A-2-A Lab Sample ID: 320-24224-2  
Client ID: WI-AF-3FB30-1216  
Operator ID: CBW ALS Bottle#: 47 Worklist Smp#: 40  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Method: 537\_A6 Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid (ND) \$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid (ND)

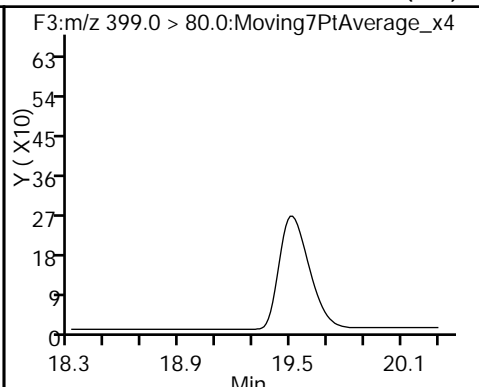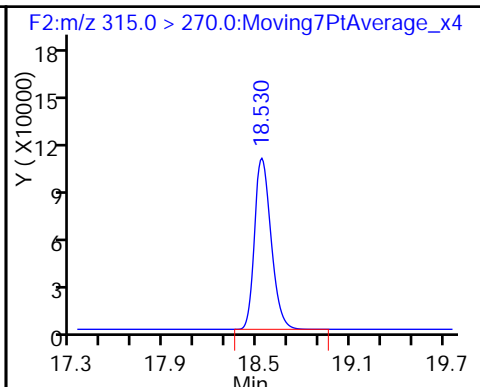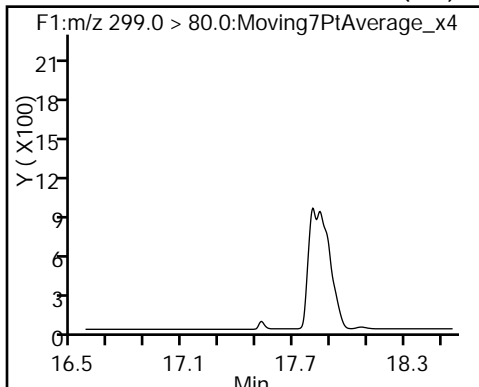

4 Perfluoroheptanoic acid (ND)

\* 5 13C2-PFOA




6 Perfluorooctanoic acid (ND)

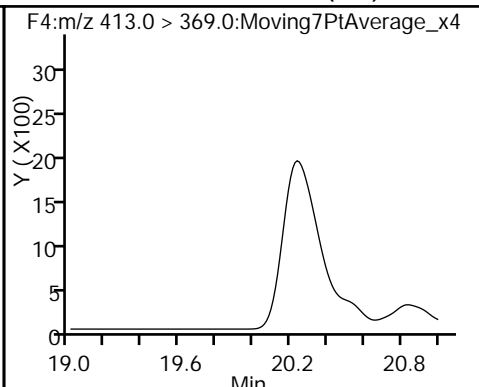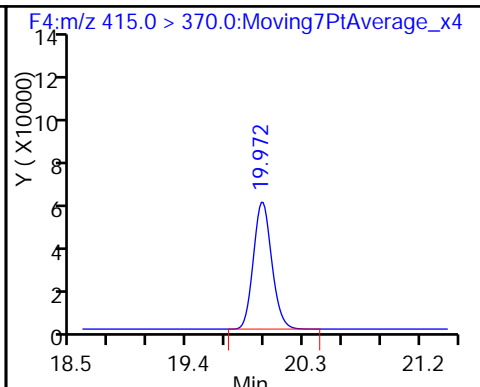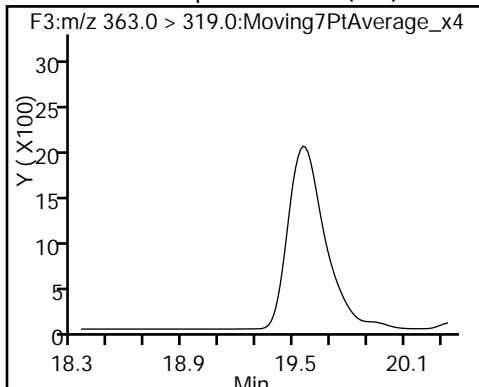


7 Perfluorooctane sulfonic acid (ND)

\* 8 13C4 PFOS

9 Perfluorononanoic acid (ND)

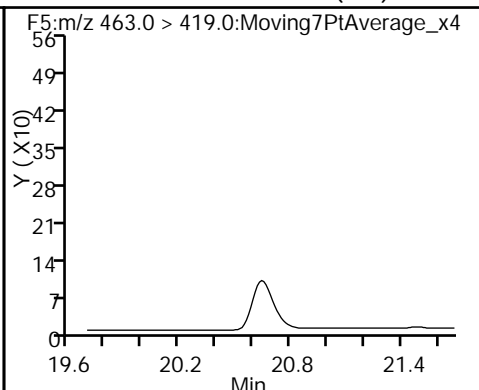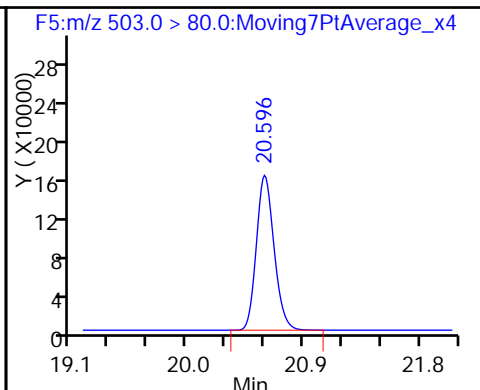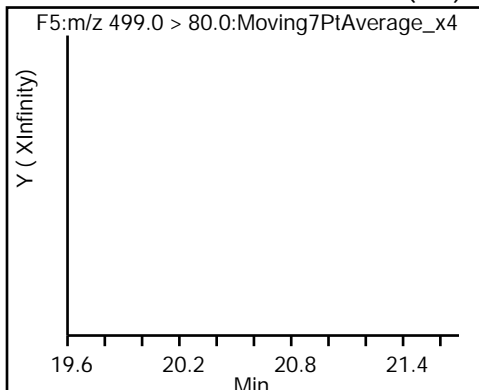

\$ 10 13C2 PFDA

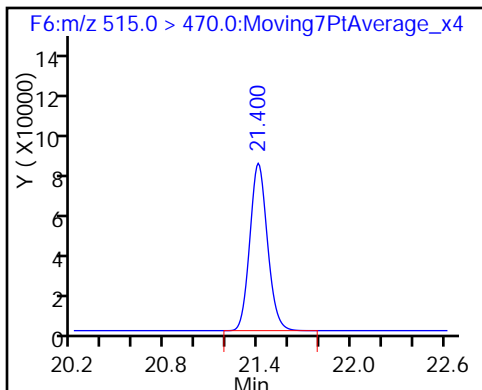

TestAmerica Sacramento  
Recovery Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_152.d  
 Lims ID: 320-24224-A-2-A  
 Client ID: WI-AF-3FB30-1216  
 Sample Type: Client  
 Inject. Date: 18-Dec-2016 10:30:05 ALS Bottle#: 47 Worklist Smp#: 40  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: 320-24224-a-2-a  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 19-Dec-2016 15:25:13 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK025

First Level Reviewer: westendorfc Date: 19-Dec-2016 09:12:22

Compound	Amount Added	Amount Recovered	% Rec.
\$ 2 13C2 PFHxA	10.0	12.1	120.57
\$ 10 13C2 PFDA	10.0	12.6	125.88

FORM VI  
LCMS BY INTERNAL STANDARD - INITIAL CALIBRATION DATA  
CURVE EVALUATION

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1 Analy Batch No.: 140688

SDG No.: \_\_\_\_\_

Instrument ID: A6 GC Column: Acquity ID: 2.1(mm) Heated Purge: (Y/N) N

Calibration Start Date: 12/05/2016 17:26 Calibration End Date: 12/05/2016 19:54 Calibration ID: 26888

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD 320-140688/2	05DEC2016A6A_004.d
Level 2	STD 320-140688/3	05DEC2016A6A_005.d
Level 3	STD 320-140688/4	05DEC2016A6A_006.d
Level 4	STD 320-140688/5	05DEC2016A6A_007.d
Level 5	STD 320-140688/6	05DEC2016A6A_008.d
Level 6	STD 320-140688/7	05DEC2016A6A_009.d

ANALYTE	RRF					CURVE TYPE	COEFFICIENT			#	MIN RRF	%RSD	#	MAX %RSD	R <sup>2</sup> OR COD	#	MIN R <sup>2</sup> OR COD
	LVL 1 LVL 6	LVL 2	LVL 3	LVL 4	LVL 5		B	M1	M2								
Perfluorobutanesulfonic acid (PFBS)	0.7247 0.6563	0.6525	0.7178	0.7256	0.7321	Ave		0.7015			5.2		30.0				
Perfluorohexanesulfonic acid	0.8344 0.8930	0.7757	0.9290	0.9478	1.0082	Ave		0.8980			9.3		30.0				
Perfluoroheptanoic acid	1.4137 1.1078	1.1891	1.2161	1.1975	1.1665	Ave		1.2151			8.6		30.0				
Perfluorooctanoic acid (PFOA)	0.9720 1.0610	0.9049	1.0674	1.1235	1.1136	Ave		1.0404			8.2		30.0				
Perfluorooctanesulfonic acid (PFOS)	0.8855 1.0951	0.9020	1.0711	1.0966	1.2136	Ave		1.0440			12.1		30.0				
Perfluorononanoic acid	0.9735 1.1655	0.9961	1.1929	1.2321	1.2453	Ave		1.1342			10.5		30.0				
13C2 PFHxA	1.0366 1.2091	1.0515	1.1929	1.2298	1.2791	Ave		1.1665			8.5		30.0				
13C2 PFDA	0.8084 0.9456	0.7439	0.8674	0.9054	0.9868	Ave		0.8763			10.2		30.0				

Note: The m1 coefficient is the same as Ave RRF for an Ave curve type.

FORM VI  
 LCMS BY INTERNAL STANDARD - INITIAL CALIBRATION DATA  
 RESPONSE AND CONCENTRATION

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1 Analy Batch No.: 140688

SDG No.: \_\_\_\_\_

Instrument ID: A6 GC Column: Acquity ID: 2.1(mm) Heated Purge: (Y/N) N

Calibration Start Date: 12/05/2016 17:26 Calibration End Date: 12/05/2016 19:54 Calibration ID: 26888

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD 320-140688/2	05DEC2016A6A_004.d
Level 2	STD 320-140688/3	05DEC2016A6A_005.d
Level 3	STD 320-140688/4	05DEC2016A6A_006.d
Level 4	STD 320-140688/5	05DEC2016A6A_007.d
Level 5	STD 320-140688/6	05DEC2016A6A_008.d
Level 6	STD 320-140688/7	05DEC2016A6A_009.d

ANALYTE	IS REF	CURVE TYPE	RESPONSE					CONCENTRATION (NG/ML)				
			LVL 1 LVL 6	LVL 2	LVL 3	LVL 4	LVL 5	LVL 1 LVL 6	LVL 2	LVL 3	LVL 4	LVL 5
Perfluorobutanesulfonic acid (PFBS)	PFOS	Ave	437563 7753569	1227165	2489398	4401661	6630132	8.76 178	22.9	45.1	90.9	135
Perfluorohexanesulfonic acid	PFOS	Ave	169827 3556638	491809	1086082	1938237	3077974	2.95 60.1	7.72	15.2	30.6	45.4
Perfluoroheptanoic acid	13PF OA	Ave	126557 2032288	324913	658044	1121930	1727957	0.994 20.2	2.60	5.12	10.3	15.3
Perfluorooctanoic acid (PFOA)	13PF OA	Ave	173304 3876381	492431	1150281	2096404	3285195	1.98 40.3	5.17	10.2	20.5	30.4
Perfluorooctanesulfonic acid (PFOS)	PFOS	Ave	238662 5775285	757269	1658139	2969550	4906017	3.91 79.6	10.2	20.1	40.6	60.1
Perfluorononanoic acid	13PF OA	Ave	168128 4124664	525061	1245341	2227031	3558831	1.92 39.0	5.01	9.87	19.9	29.5
13C2 PFHxA	13PF OA	Ave	933751 1095977	1106485	1261522	1117585	1240474	10.0 10.0	10.0	10.0	10.0	10.0
13C2 PFDA	13PF OA	Ave	728204 857144	782778	917302	822787	957025	10.0 10.0	10.0	10.0	10.0	10.0

Curve Type Legend:

Ave = Average ISTD

FORM VI  
 LCMS BY INTERNAL STANDARD - INITIAL CALIBRATION DATA  
 READBACK PERCENT ERROR

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1 Analy Batch No.: 140688

SDG No.: \_\_\_\_\_

Instrument ID: A6 GC Column: Acquity ID: 2.1(mm) Heated Purge: (Y/N) N

Calibration Start Date: 12/05/2016 17:26 Calibration End Date: 12/05/2016 19:54 Calibration ID: 26888

Calibration Files:


LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD 320-140688/2	05DEC2016A6A_004.d
Level 2	STD 320-140688/3	05DEC2016A6A_005.d
Level 3	STD 320-140688/4	05DEC2016A6A_006.d
Level 4	STD 320-140688/5	05DEC2016A6A_007.d
Level 5	STD 320-140688/6	05DEC2016A6A_008.d
Level 6	STD 320-140688/7	05DEC2016A6A_009.d

ANALYTE	PERCENT ERROR						PERCENT ERROR LIMIT					
	LVL 1 #	LVL 2 #	LVL 3 #	LVL 4 #	LVL 5 #	LVL 6 #	LVL 1	LVL 2	LVL 3	LVL 4	LVL 5	LVL 6
Perfluorobutanesulfonic acid (PFBS)	3.3	-7.0	2.3	3.4	4.4	-6.4	50	50	50	50	50	50
Perfluorohexanesulfonic acid	-7.1	-13.6	3.4	5.5	12.3	-0.6	50	50	50	50	50	50
Perfluoroheptanoic acid	16.3	-2.1	0.1	-1.5	-4.0	-8.8	50	50	50	50	50	50
Perfluorooctanoic acid (PFOA)	-6.6	-13.0	2.6	8.0	7.0	2.0	50	50	50	50	50	50
Perfluorooctanesulfonic acid (PFOS)	-15.2	-13.6	2.6	5.0	16.2	4.9	50	50	50	50	50	50
Perfluorononanoic acid	-14.2	-12.2	5.2	8.6	9.8	2.8	50	50	50	50	50	50
13C2 PFHxA	-11.1	-9.9	2.3	5.4	9.7	3.7	30	30	30	30	30	30
13C2 PFDA	-7.7	-15.1	-1.0	3.3	12.6	7.9	30	30	30	30	30	30

L4

05DEC2016A6A\_007 Sm (Mn, 2x3)


1: MRM of 1 Channel ES-  
299 > 80 (PFBS)  
7.74e5  
Area

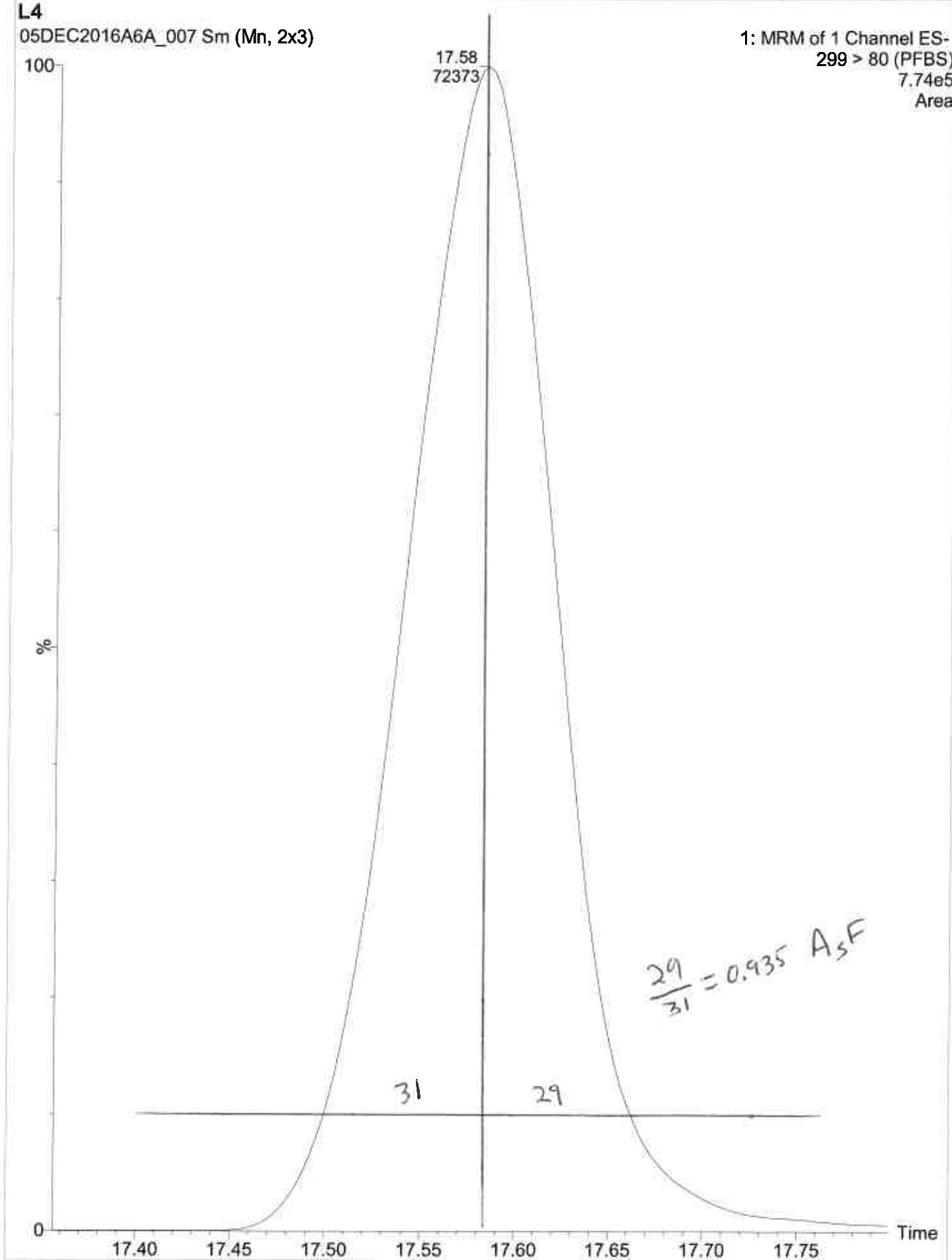

L4

05DEC2016A6A\_007 Sm (Mn, 2x3)

2: MRM of 1 Channel ES-  
315 > 270 (13C2-PFHxA)  
1.80e5  
Area

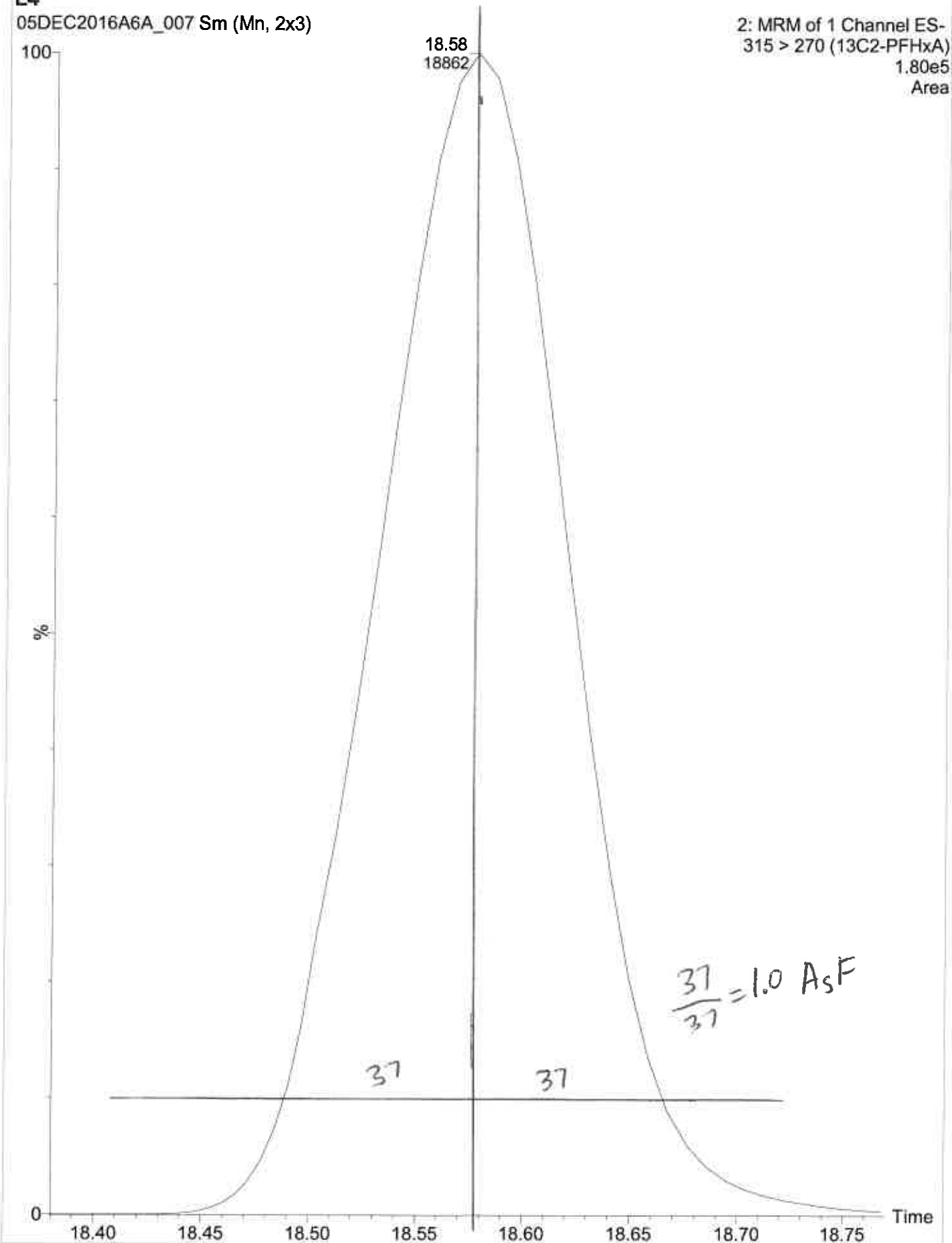

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_004.d  
 Lims ID: STD L1  
 Client ID:  
 Sample Type: IC Calib Level: 1  
 Inject. Date: 05-Dec-2016 17:26:03 ALS Bottle#: 1 Worklist Smp#: 2  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: STD L1 L1  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist: chrom-537\_\_A6\*sub3  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 06-Dec-2016 16:35:34 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK024

First Level Reviewer: barnettj Date: 06-Dec-2016 10:00:02

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.576	17.581	-0.005	1.000	437563	9.05	466
\$ 2 13C2 PFHxA	315.0 > 270.0	18.567	18.567	0.0	1.000	933751	8.89	30467
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.332	19.342	-0.010	1.000	169827	2.74	4140
4 Perfluoroheptanoic acid	363.0 > 319.0	19.368	19.378	-0.010	1.000	126557	1.16	45.1 M
* 5 13C2-PFOA	415.0 > 370.0	20.047	20.047	0.0		900764	10.0	23392
6 Perfluorooctanoic acid	413.0 > 369.0	20.047	20.047	0.0	1.000	173304	1.85	35.0 M
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.619	20.619	0.0	1.000	238662	3.32	2941
* 8 13C4 PFOS	503.0 > 80.0	20.667	20.669	-0.002		1976615	28.7	40886
9 Perfluorononanoic acid	463.0 > 419.0	20.738	20.748	-0.010	1.000	168128	1.65	6043
\$ 10 13C2 PFDA	515.0 > 470.0	21.471	21.474	-0.003	1.000	728204	9.23	22953



**QC Flag Legend**

Review Flags

M - Manually Integrated

**Reagents:**

LC537-L1\_00015

Amount Added: 1.00

Units: mL

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_004.d

Injection Date: 05-Dec-2016 17:26:03

Instrument ID: A6

Lims ID: STD L1

Client ID:

Operator ID: CBW

ALS Bottle#: 1

Worklist Smp#: 2

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

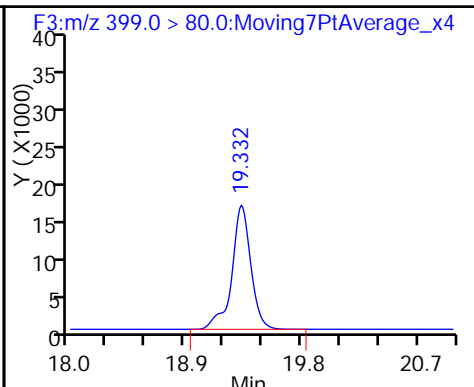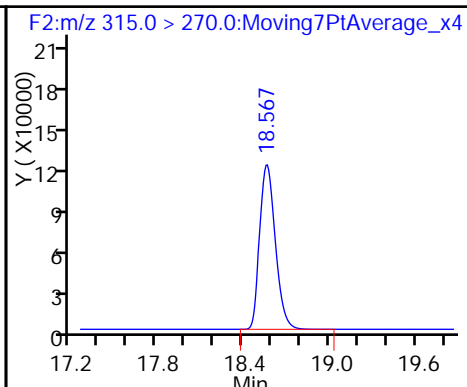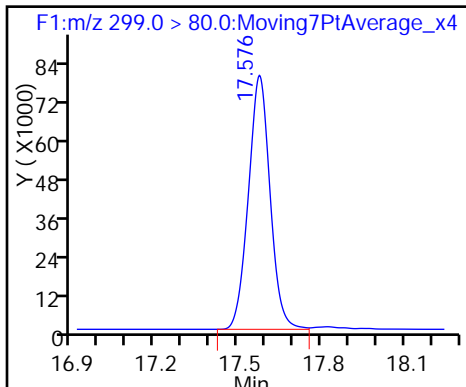

4 Perfluoroheptanoic acid (M)

\* 5 13C2-PFOA




6 Perfluorooctanoic acid (M)

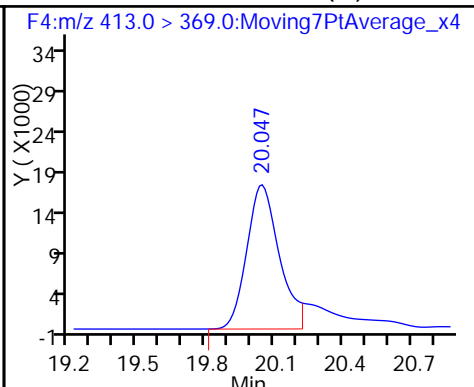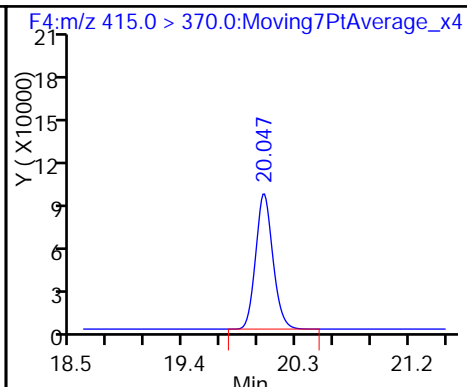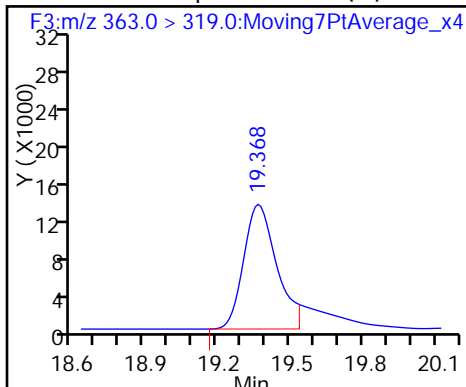


7 Perfluorooctane sulfonic acid

\* 8 13C4 PFOS

9 Perfluorononanoic acid

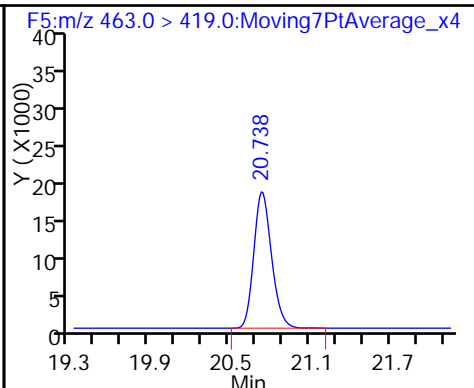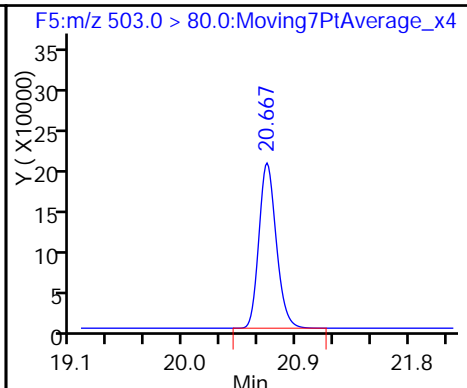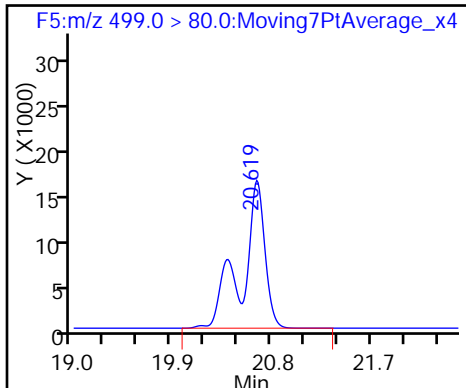

\$ 10 13C2 PFDA

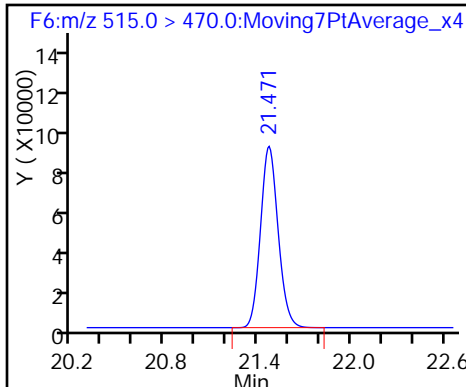

TestAmerica Sacramento


Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_004.d  
Injection Date: 05-Dec-2016 17:26:03 Instrument ID: A6  
Lims ID: STD L1  
Client ID:  
Operator ID: CBW ALS Bottle#: 1 Worklist Smp#: 2  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Method: 537\_\_A6 Limit Group: LC 537 ICAL  
Column: Acquity BEH C18 ( 2.10 mm) Detector F3:MRM

4 Perfluoroheptanoic acid, CAS: 375-85-9

Signal: 1


RT: 19.37  
Area: 155591  
Amount: 1.476072  
Amount Units: ng/ml

Processing Integration Results

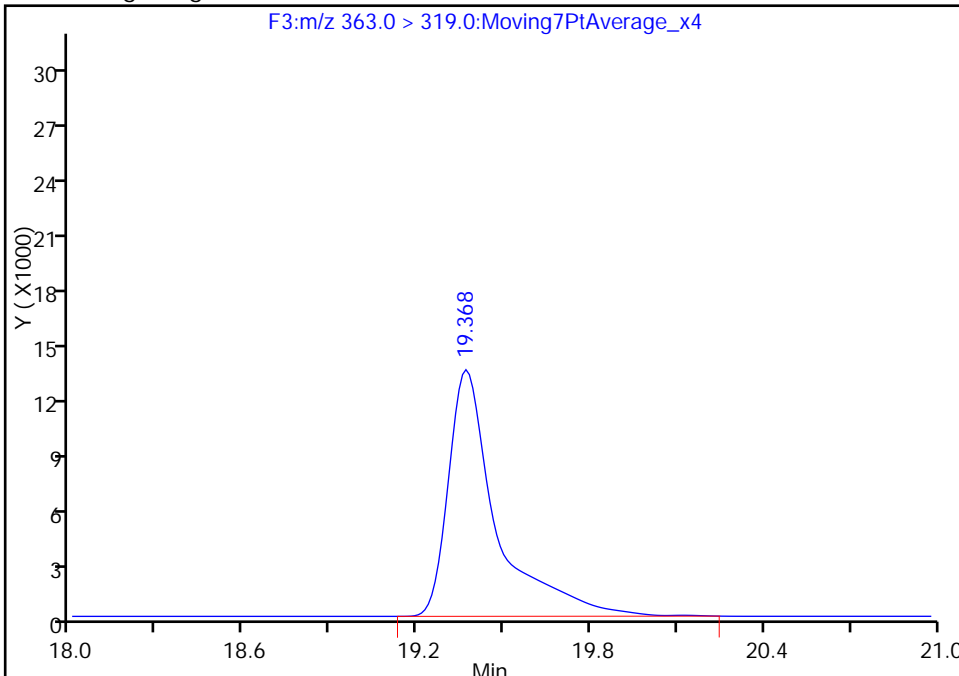

RT: 19.37  
Area: 126557  
Amount: 1.156251  
Amount Units: ng/ml

Manual Integration Results

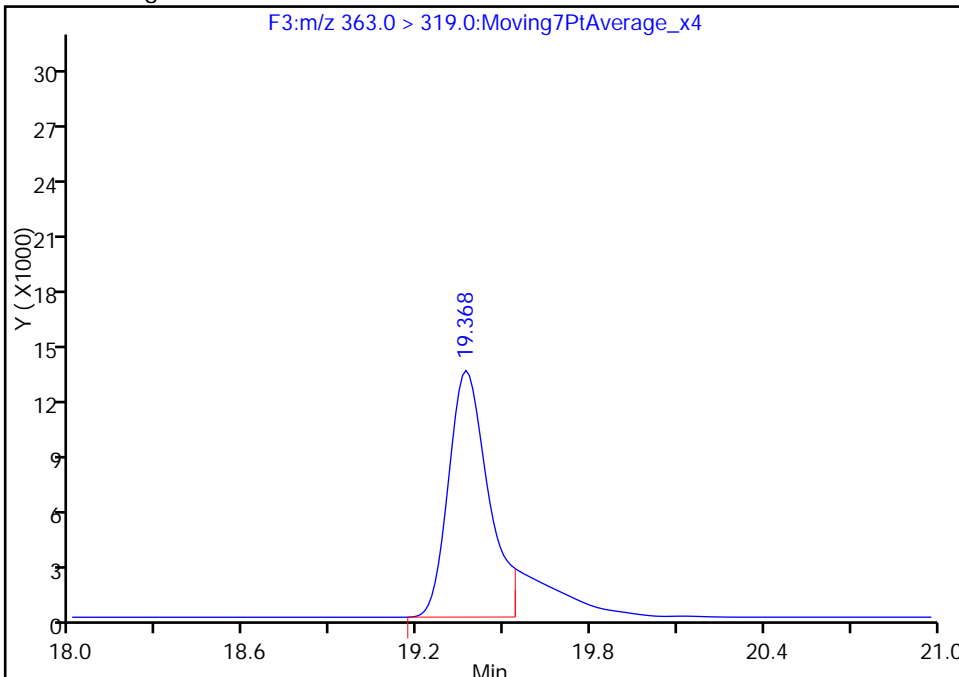

Reviewer: barnettj, 06-Dec-2016 10:00:02  
Audit Action: Manually Integrated

Audit Reason: Split Peak

TestAmerica Sacramento


Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_004.d  
Injection Date: 05-Dec-2016 17:26:03 Instrument ID: A6  
Lims ID: STD L1  
Client ID:  
Operator ID: CBW ALS Bottle#: 1 Worklist Smp#: 2  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Method: 537\_\_A6 Limit Group: LC 537 ICAL  
Column: Acquity BEH C18 ( 2.10 mm) Detector F4:MRM

6 Perfluorooctanoic acid, CAS: 335-67-1

Signal: 1


RT: 20.05  
Area: 186490  
Amount: 1.959453  
Amount Units: ng/ml

Processing Integration Results

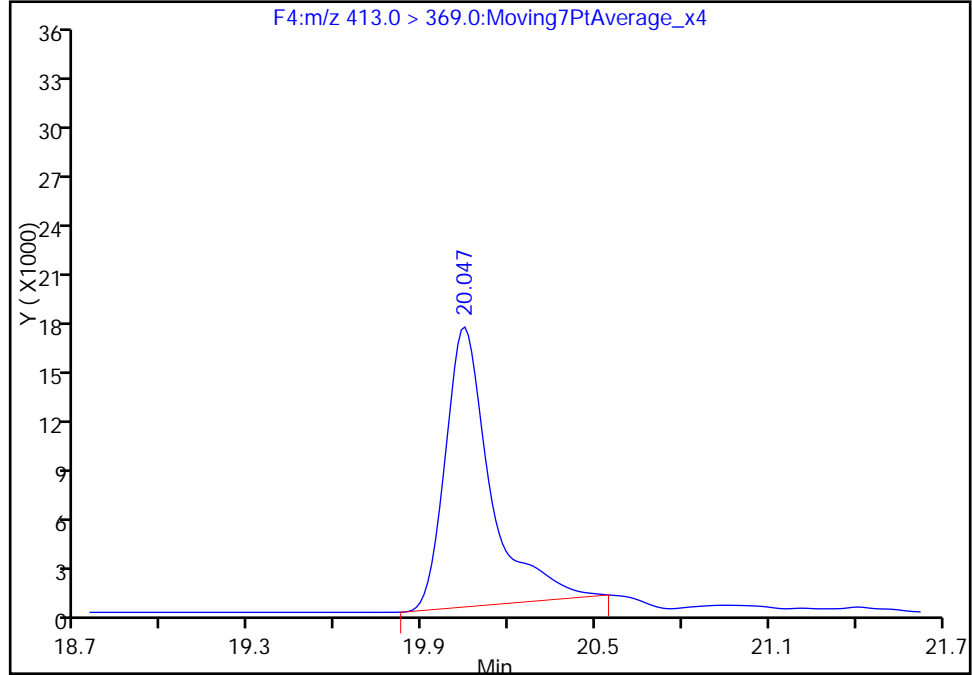

RT: 20.05  
Area: 173304  
Amount: 1.849212  
Amount Units: ng/ml

Manual Integration Results

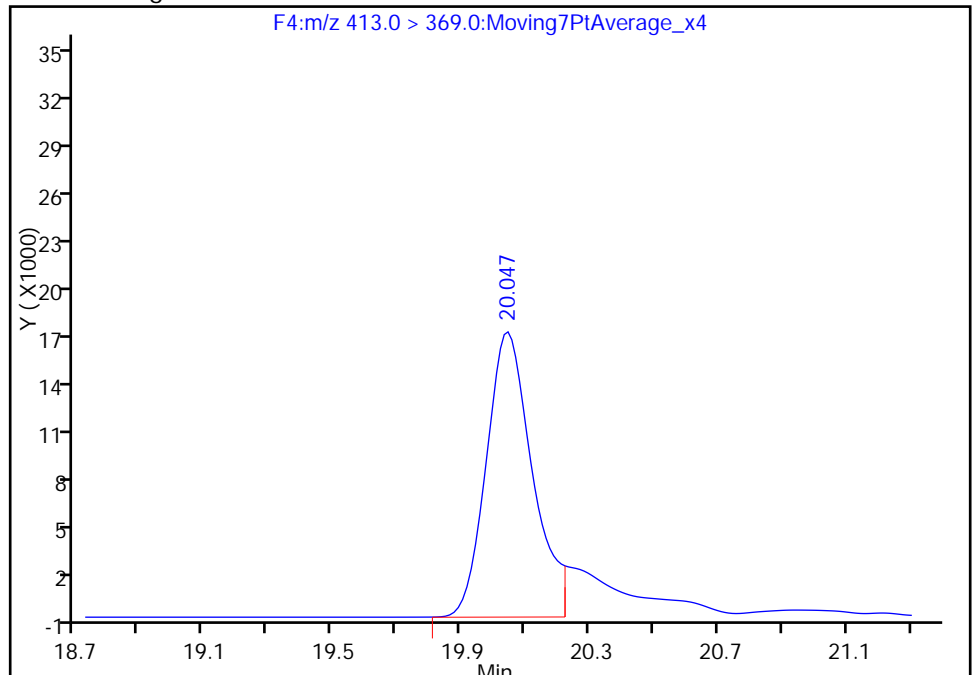

Reviewer: barnettj, 06-Dec-2016 10:00:02  
Audit Action: Manually Integrated

Audit Reason: Split Peak

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_005.d  
 Lims ID: STD L2  
 Client ID:  
 Sample Type: IC Calib Level: 2  
 Inject. Date: 05-Dec-2016 17:55:38 ALS Bottle#: 2 Worklist Smp#: 3  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: STD L2 L2  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist: chrom-537\_\_A6\*sub3  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 06-Dec-2016 16:35:35 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK024

First Level Reviewer: barnettj Date: 06-Dec-2016 09:58:24

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.582	17.581	0.001	1.000	1227165	21.3	5055
\$ 2 13C2 PFHxA	315.0 > 270.0	18.567	18.567	0.0	1.000	1106485	9.01	35678
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.344	19.342	0.002	1.000	491809	6.67	11495
4 Perfluoroheptanoic acid	363.0 > 319.0	19.380	19.378	0.002	1.000	324913	2.54	155 M
* 5 13C2-PFOA	415.0 > 370.0	20.047	20.047	0.0		1052273	10.0	27645
6 Perfluorooctanoic acid	413.0 > 369.0	20.047	20.047	0.0	1.000	492431	4.50	100 M
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.619	20.619	0.0	1.000	757269	8.83	8449
* 8 13C4 PFOS	503.0 > 80.0	20.667	20.669	-0.002		2356620	28.7	30757
9 Perfluorononanoic acid	463.0 > 419.0	20.750	20.748	0.002	1.000	525061	4.40	13911
\$ 10 13C2 PFDA	515.0 > 470.0	21.480	21.474	0.006	1.000	782778	8.49	24678

### QC Flag Legend

Review Flags

M - Manually Integrated

### Reagents:

LC537-L2\_00014

Amount Added: 1.00

Units: mL

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_005.d

Injection Date: 05-Dec-2016 17:55:38

Instrument ID: A6

Lims ID: STD L2

Client ID:

Operator ID: CBW

ALS Bottle#: 2

Worklist Smp#: 3

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

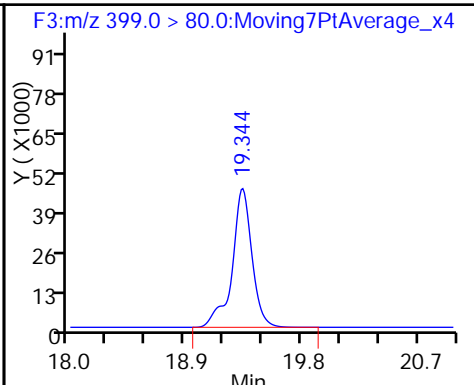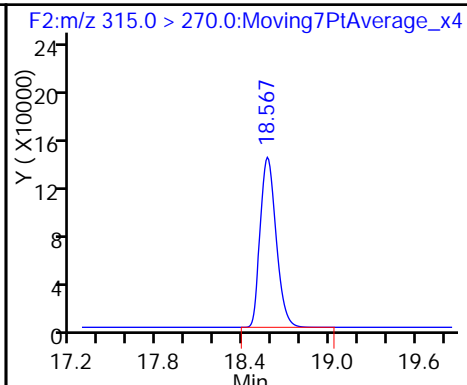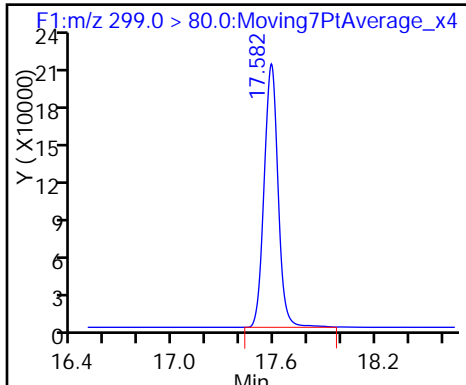

4 Perfluoroheptanoic acid (M)

\* 5 13C2-PFOA




6 Perfluorooctanoic acid (M)

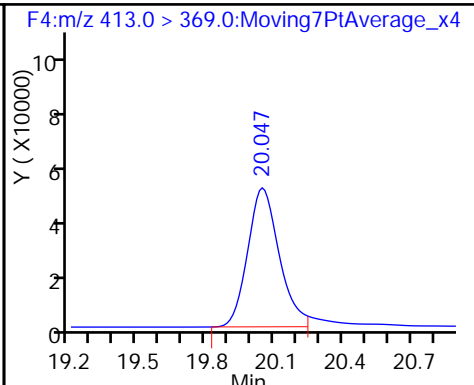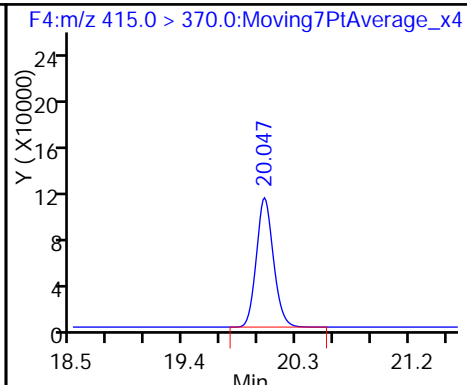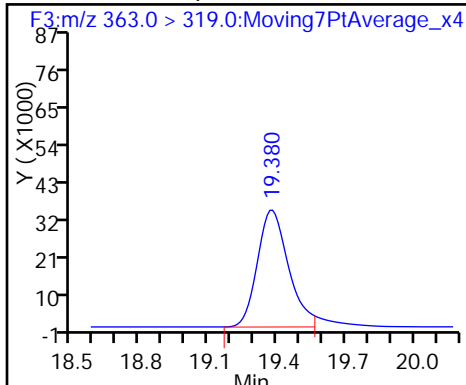


7 Perfluorooctane sulfonic acid

\* 8 13C4 PFOS

9 Perfluorononanoic acid

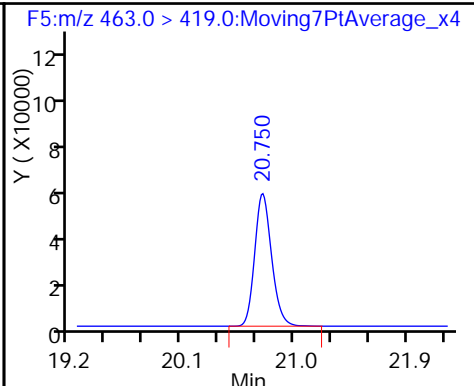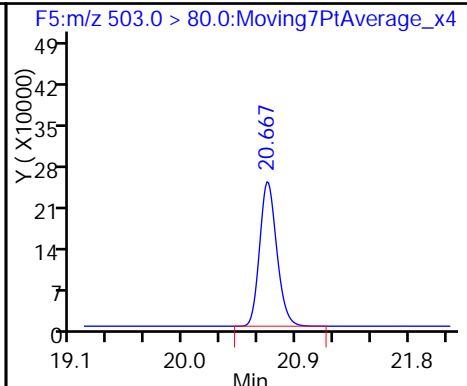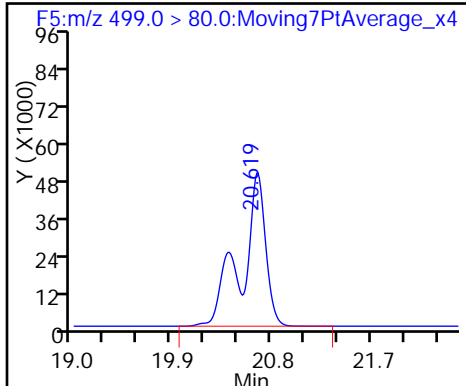

\$ 10 13C2 PFDA

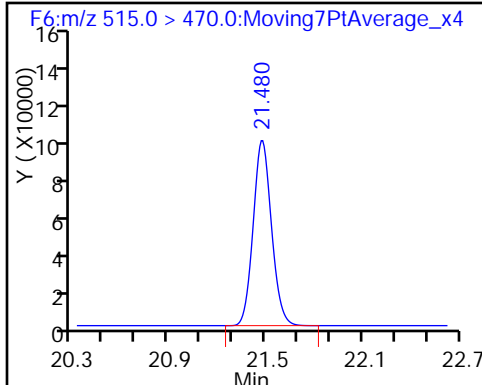

TestAmerica Sacramento


Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_005.d  
Injection Date: 05-Dec-2016 17:55:38 Instrument ID: A6  
Lims ID: STD L2  
Client ID:  
Operator ID: CBW ALS Bottle#: 2 Worklist Smp#: 3  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Method: 537\_\_A6 Limit Group: LC 537 ICAL  
Column: Acquity BEH C18 ( 2.10 mm) Detector F3:M/RM

4 Perfluoroheptanoic acid, CAS: 375-85-9

Signal: 1


RT: 19.38  
Area: 344811  
Amount: 2.670013  
Amount Units: ng/ml

Processing Integration Results

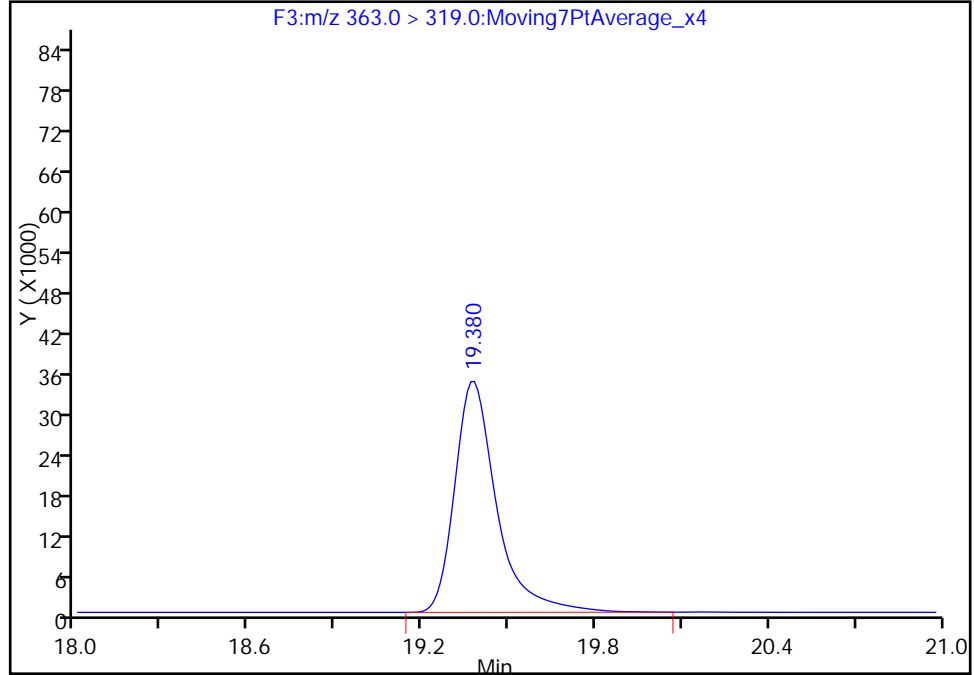

RT: 19.38  
Area: 324913  
Amount: 2.541065  
Amount Units: ng/ml

Manual Integration Results

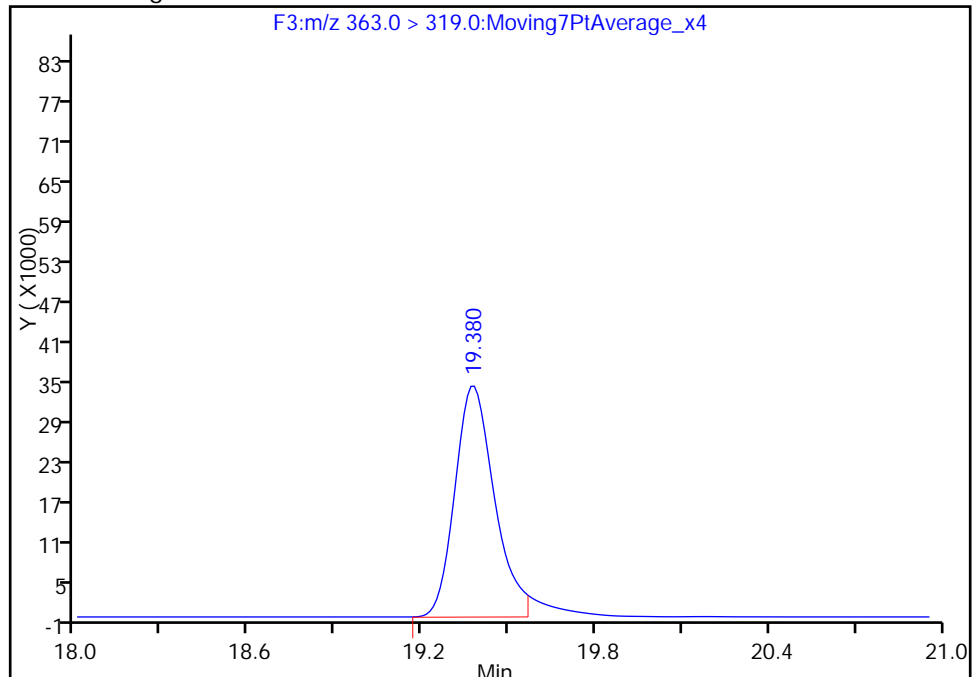

Reviewer: barnettj, 06-Dec-2016 10:03:30  
Audit Action: Manually Integrated

Audit Reason: Split Peak



TestAmerica Sacramento


Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_005.d  
Injection Date: 05-Dec-2016 17:55:38 Instrument ID: A6  
Lims ID: STD L2  
Client ID:  
Operator ID: CBW ALS Bottle#: 2 Worklist Smp#: 3  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Method: 537\_\_A6 Limit Group: LC 537 ICAL  
Column: Acquity BEH C18 ( 2.10 mm) Detector F4:MRM

6 Perfluorooctanoic acid, CAS: 335-67-1

Signal: 1


RT: 20.05  
Area: 504990  
Amount: 4.595586  
Amount Units: ng/ml

Processing Integration Results

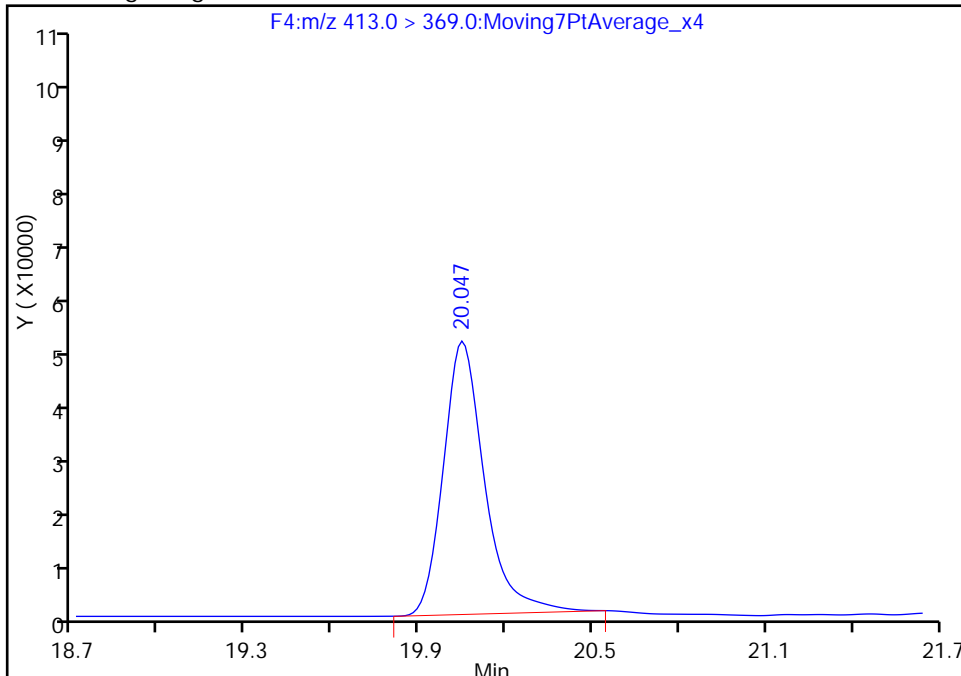

RT: 20.05  
Area: 492431  
Amount: 4.497863  
Amount Units: ng/ml

Manual Integration Results

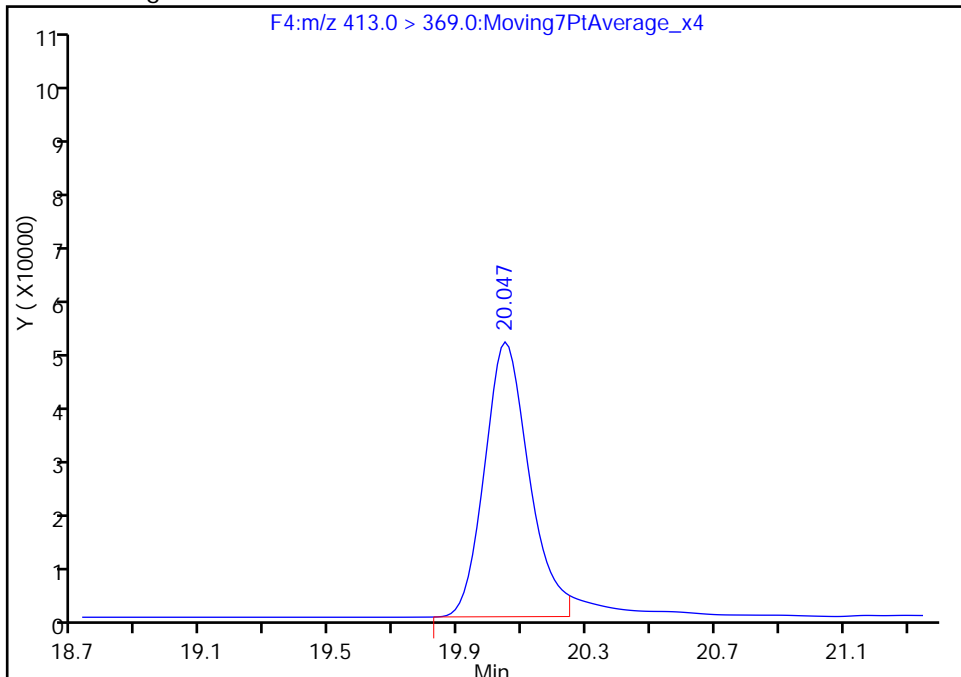

Reviewer: barnettj, 06-Dec-2016 10:03:30  
Audit Action: Manually Integrated

Audit Reason: Split Peak

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_006.d  
 Lims ID: STD L3  
 Client ID:  
 Sample Type: IC Calib Level: 3  
 Inject. Date: 05-Dec-2016 18:25:13 ALS Bottle#: 3 Worklist Smp#: 4  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: STD L3 L3  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist: chrom-537\_\_A6\*sub3  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 06-Dec-2016 16:35:36 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK024

First Level Reviewer: barnettj Date: 06-Dec-2016 09:58:05

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.582	17.581	0.001	1.000	2489398	46.2	1804
\$ 2 13C2 PFHxA	315.0 > 270.0	18.567	18.567	0.0	1.000	1261522	10.2	40506
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.344	19.342	0.002	1.000	1086082	15.7	25400
4 Perfluoroheptanoic acid	363.0 > 319.0	19.380	19.378	0.002	1.000	658044	5.12	4774
* 5 13C2-PFOA	415.0 > 370.0	20.047	20.047	0.0		1057506	10.0	27287
6 Perfluorooctanoic acid	413.0 > 369.0	20.047	20.047	0.0	1.000	1150281	10.5	429
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.619	20.619	0.0	1.000	1658139	20.7	19019
* 8 13C4 PFOS	503.0 > 80.0	20.667	20.669	-0.002		2205243	28.7	57142
9 Perfluorononanoic acid	463.0 > 419.0	20.750	20.748	0.002	1.000	1245341	10.4	13210
\$ 10 13C2 PFDA	515.0 > 470.0	21.471	21.474	-0.003	1.000	917302	9.90	28753

Reagents:

LC537-L3\_00016 Amount Added: 1.00 Units: mL

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_006.d

Injection Date: 05-Dec-2016 18:25:13

Instrument ID: A6

Lims ID: STD L3

Client ID:

Operator ID: CBW

ALS Bottle#: 3

Worklist Smp#: 4

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

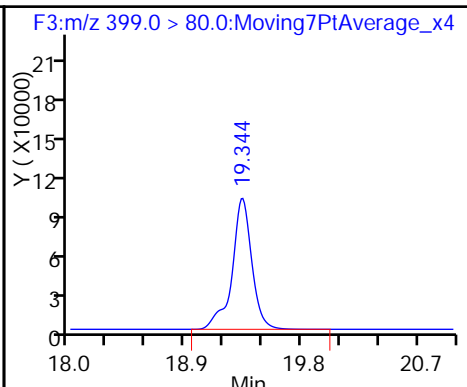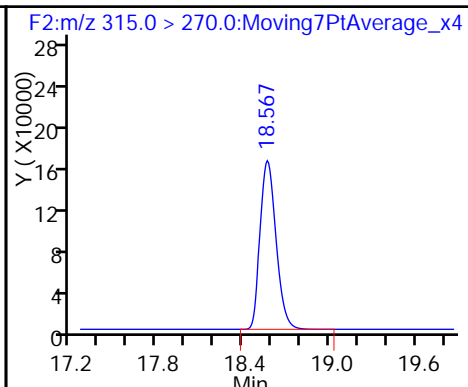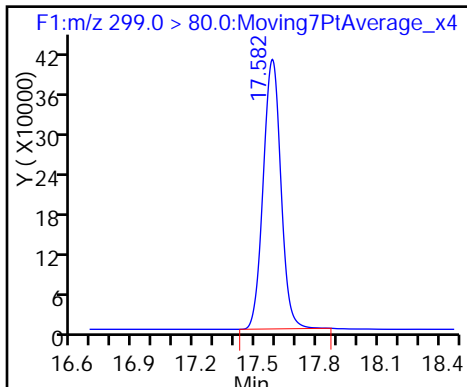

4 Perfluoroheptanoic acid

\* 5 13C2-PFOA




6 Perfluorooctanoic acid

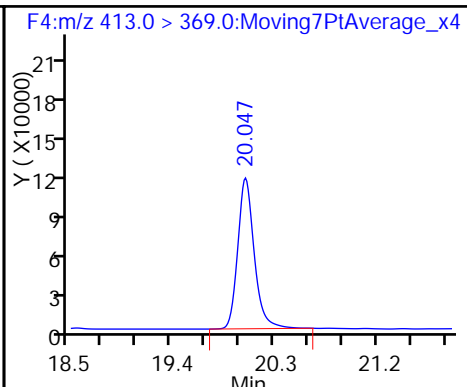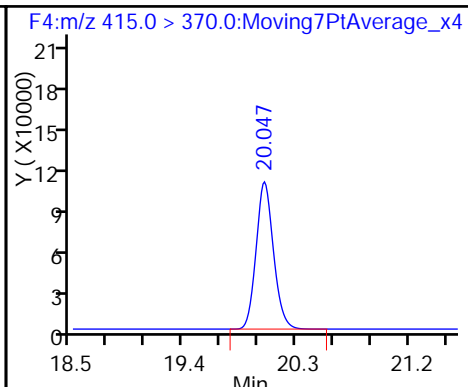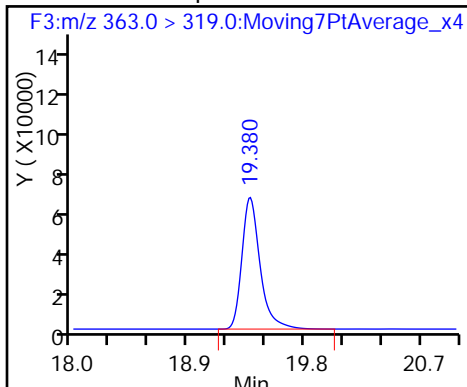


7 Perfluorooctane sulfonic acid

\* 8 13C4 PFOS

9 Perfluorononanoic acid

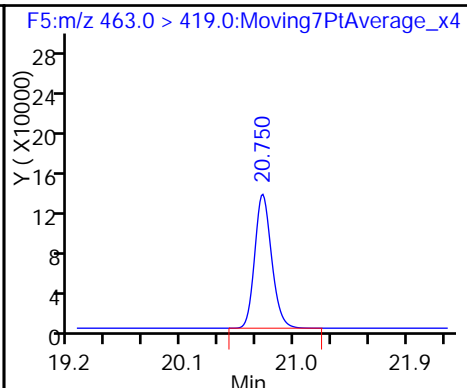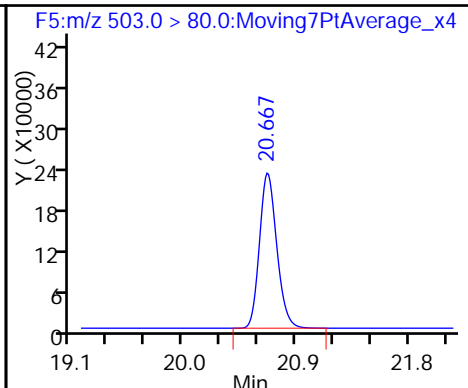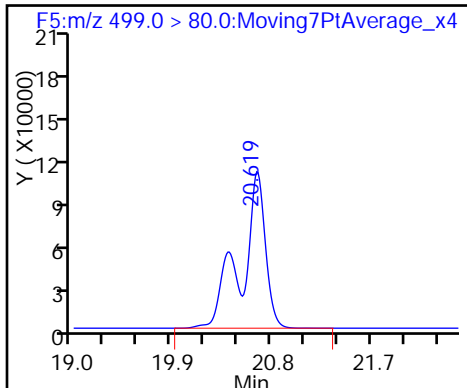

\$ 10 13C2 PFDA

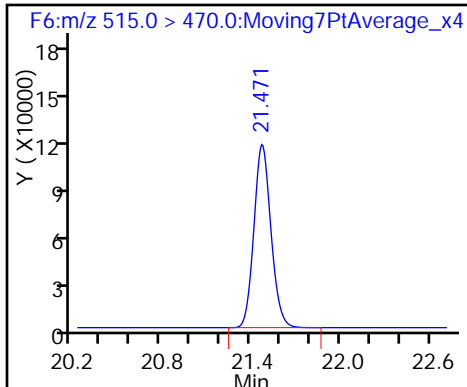

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_007.d  
 Lims ID: STD L4  
 Client ID:  
 Sample Type: ICISAV Calib Level: 4  
 Inject. Date: 05-Dec-2016 18:54:48 ALS Bottle#: 4 Worklist Smp#: 5  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: STD L4 L4  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist: chrom-537\_\_A6\*sub3  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 06-Dec-2016 16:35:37 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK024

First Level Reviewer: barnettj Date: 06-Dec-2016 13:43:03

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.579	17.581	-0.002	1.000	4401661	94.0	2768
\$ 2 13C2 PFHxA	315.0 > 270.0	18.567	18.567	0.0	1.000	1117585	10.5	28676
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.344	19.342	0.002	1.000	1938237	32.3	25196
4 Perfluoroheptanoic acid	363.0 > 319.0	19.380	19.378	0.002	1.000	1121930	10.2	12796
* 5 13C2-PFOA	415.0 > 370.0	20.047	20.047	0.0		908727	10.0	23744
6 Perfluorooctanoic acid	413.0 > 369.0	20.047	20.047	0.0	1.000	2096404	22.2	516
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.619	20.619	0.0	1.000	2969550	42.6	9704
* 8 13C4 PFOS	503.0 > 80.0	20.667	20.669	-0.002		1914415	28.7	28032
9 Perfluorononanoic acid	463.0 > 419.0	20.750	20.748	0.002	1.000	2227031	21.6	23494
\$ 10 13C2 PFDA	515.0 > 470.0	21.471	21.474	-0.003	1.000	822787	10.3	25796

Reagents:

LC537-L4\_00015 Amount Added: 1.00 Units: mL

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_007.d

Injection Date: 05-Dec-2016 18:54:48

Instrument ID: A6

Lims ID: STD L4

Client ID:

Operator ID: CBW

ALS Bottle#: 4

Worklist Smp#: 5

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

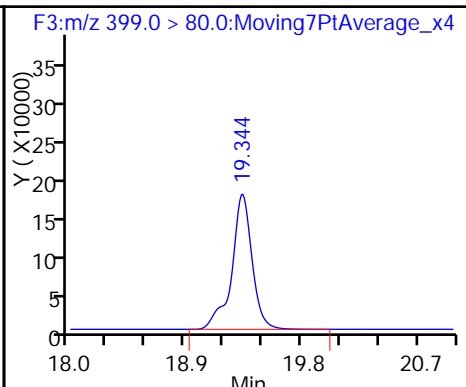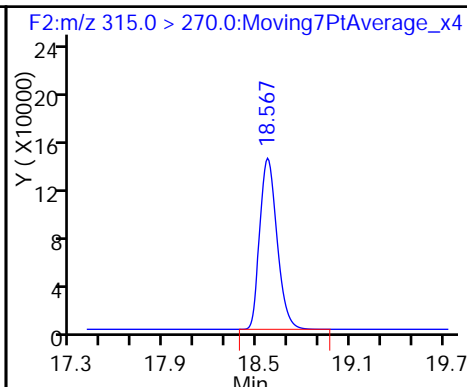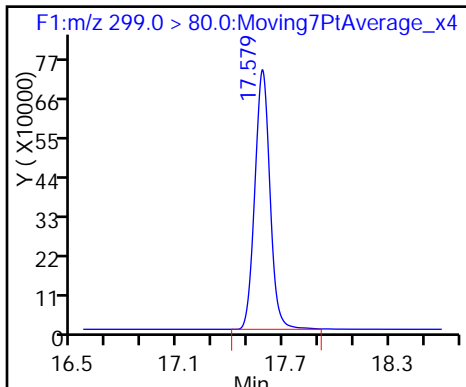

4 Perfluoroheptanoic acid

\* 5 13C2-PFOA




6 Perfluorooctanoic acid

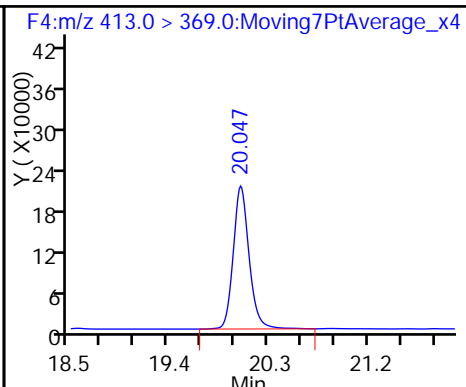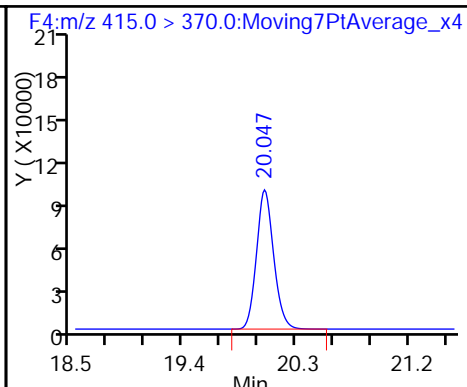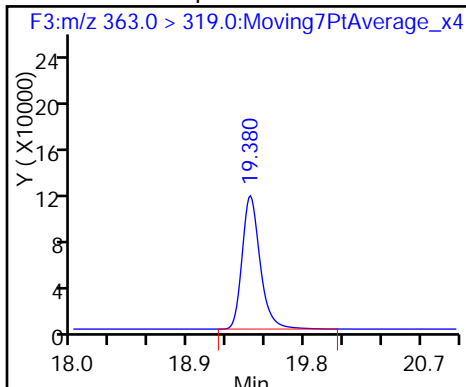


7 Perfluorooctane sulfonic acid

\* 8 13C4 PFOS

9 Perfluorononanoic acid

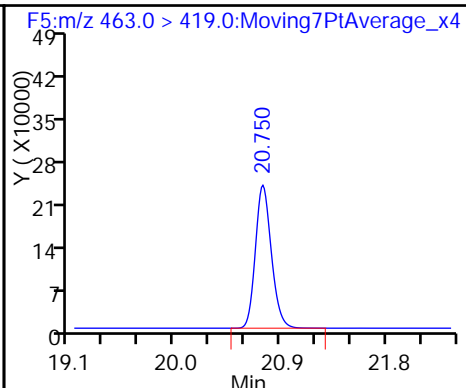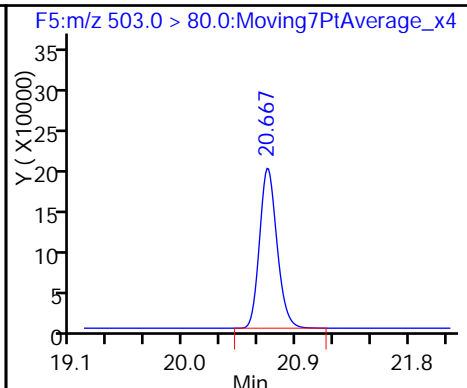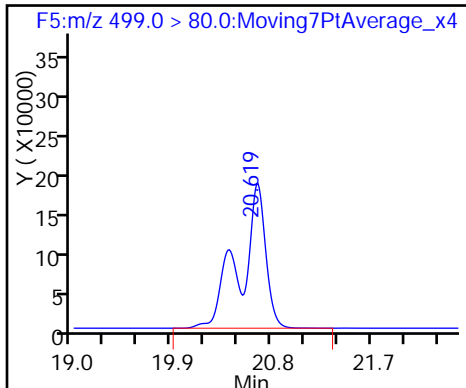

\$ 10 13C2 PFDA

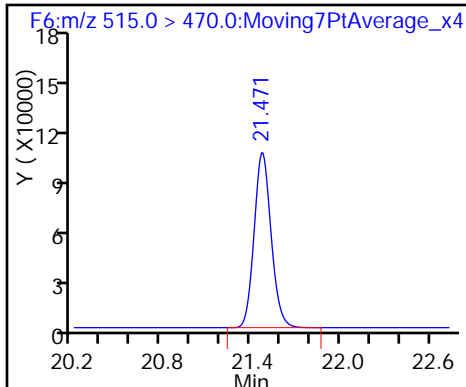

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_008.d  
 Lims ID: STD L5  
 Client ID:  
 Sample Type: IC Calib Level: 5  
 Inject. Date: 05-Dec-2016 19:24:23 ALS Bottle#: 5 Worklist Smp#: 6  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: STD L5 L5  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist: chrom-537\_\_A6\*sub3  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 06-Dec-2016 16:35:38 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK024

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.582	17.581	0.001	1.000	6630132	140.5	3208
\$ 2 13C2 PFHxA	315.0 > 270.0	18.567	18.567	0.0	1.000	1240474	11.0	39454
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.344	19.342	0.002	1.000	3077974	51.0	14553
4 Perfluoroheptanoic acid	363.0 > 319.0	19.380	19.378	0.002	1.000	1727957	14.7	6886
* 5 13C2-PFOA	415.0 > 370.0	20.047	20.047	0.0		969779	10.0	24964
6 Perfluorooctanoic acid	413.0 > 369.0	20.047	20.047	0.0	1.000	3285195	32.6	1114
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.620	20.619	0.001	1.000	4906017	69.9	10146
* 8 13C4 PFOS	503.0 > 80.0	20.679	20.669	0.010		1929192	28.7	32805
9 Perfluorononanoic acid	463.0 > 419.0	20.750	20.748	0.002	1.000	3558831	32.4	16307
\$ 10 13C2 PFDA	515.0 > 470.0	21.480	21.474	0.006	1.000	957025	11.3	30231

Reagents:

LC537-L5\_00017 Amount Added: 1.00 Units: mL

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_008.d

Injection Date: 05-Dec-2016 19:24:23

Instrument ID: A6

Lims ID: STD L5

Client ID:

Operator ID: CBW

ALS Bottle#: 5

Worklist Smp#: 6

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

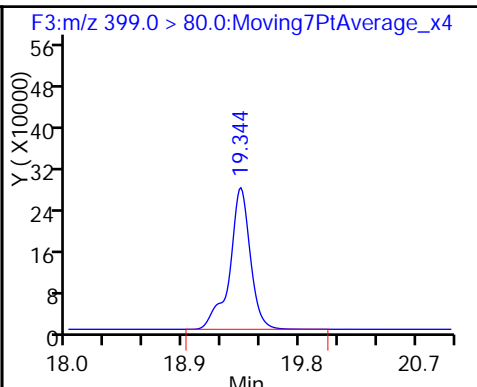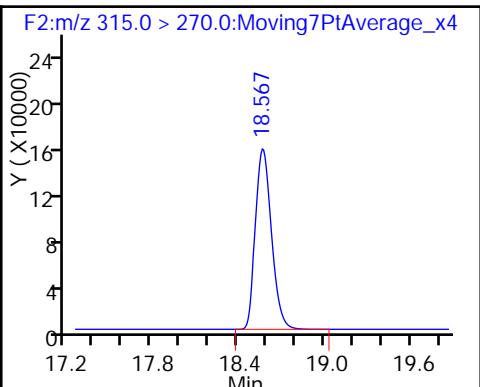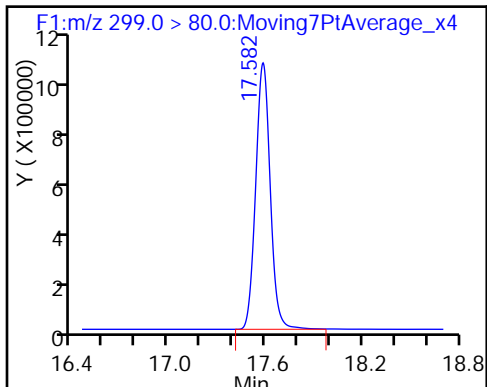

4 Perfluoroheptanoic acid

\* 5 13C2-PFOA




6 Perfluorooctanoic acid

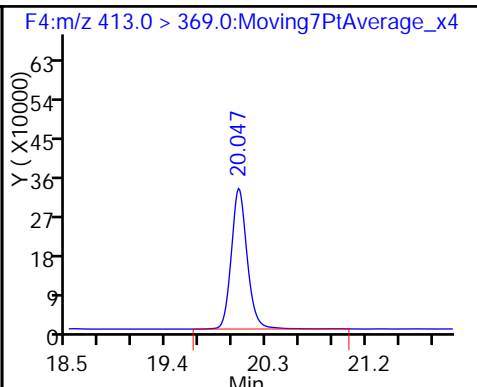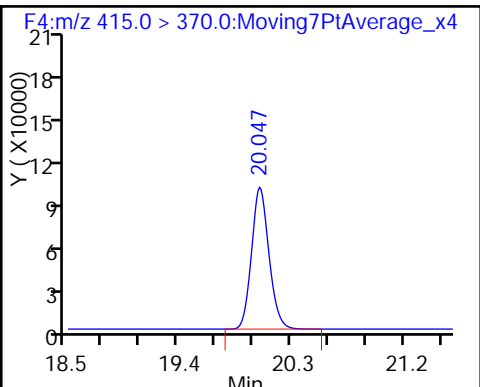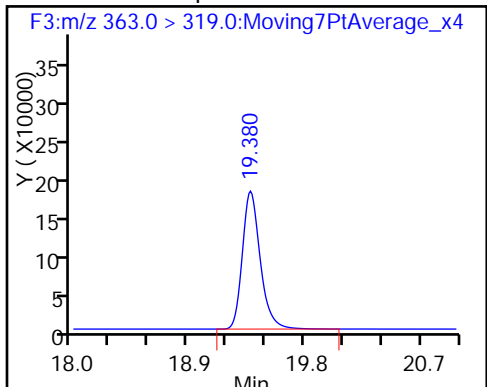


7 Perfluorooctane sulfonic acid

\* 8 13C4 PFOS

9 Perfluorononanoic acid

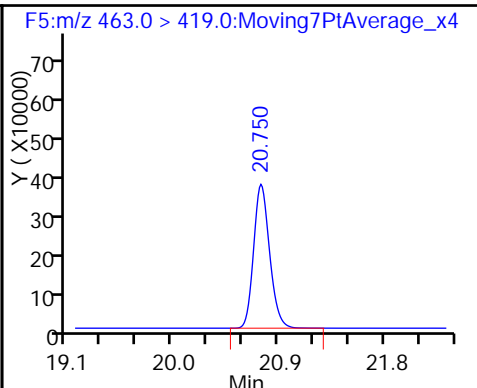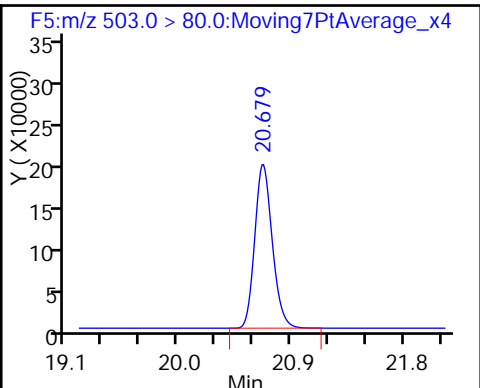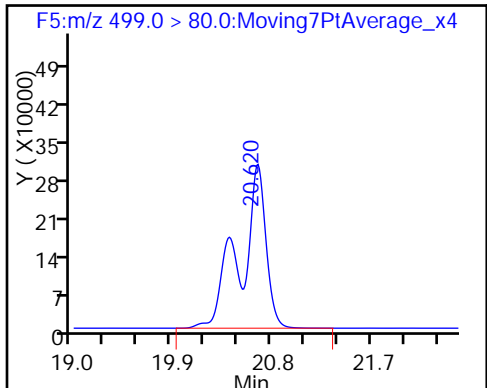

\$ 10 13C2 PFDA

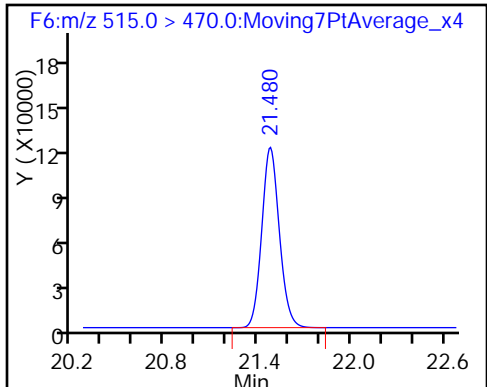

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Lims ID: STD L6  
 Client ID:  
 Sample Type: IC Calib Level: 6  
 Inject. Date: 05-Dec-2016 19:54:00 ALS Bottle#: 6 Worklist Smp#: 7  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: STD L6 L6  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist: chrom-537\_\_A6\*sub3  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 06-Dec-2016 16:35:39 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK024

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.582	17.581	0.001	1.000	7753569	166.9	8570
\$ 2 13C2 PFHxA	315.0 > 270.0	18.567	18.567	0.0	1.000	1095977	10.4	34796
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.344	19.342	0.002	1.000	3556638	59.8	31299
4 Perfluoroheptanoic acid	363.0 > 319.0	19.380	19.378	0.002	1.000	2032288	18.5	6367
* 5 13C2-PFOA	415.0 > 370.0	20.047	20.047	0.0		906416	10.0	23083
6 Perfluorooctanoic acid	413.0 > 369.0	20.047	20.047	0.0	1.000	3876381	41.1	917
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.619	20.619	0.0	1.000	5775285	83.5	12991
* 8 13C4 PFOS	503.0 > 80.0	20.667	20.669	-0.002		1899408	28.7	17628
9 Perfluorononanoic acid	463.0 > 419.0	20.750	20.748	0.002	1.000	4124664	40.1	17939
\$ 10 13C2 PFDA	515.0 > 470.0	21.471	21.474	-0.003	1.000	857144	10.8	26862

Reagents:

LC537-L6\_00014 Amount Added: 1.00 Units: mL



TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d

Injection Date: 05-Dec-2016 19:54:00

Instrument ID: A6

Lims ID: STD L6

Client ID:

Operator ID: CBW

ALS Bottle#: 6

Worklist Smp#: 7

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

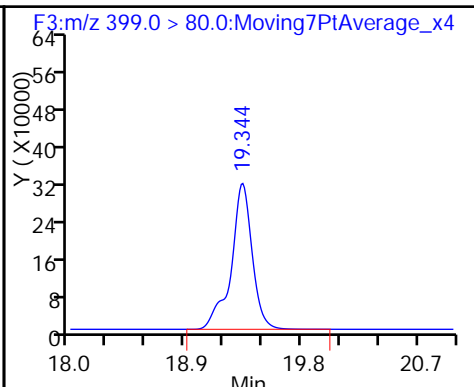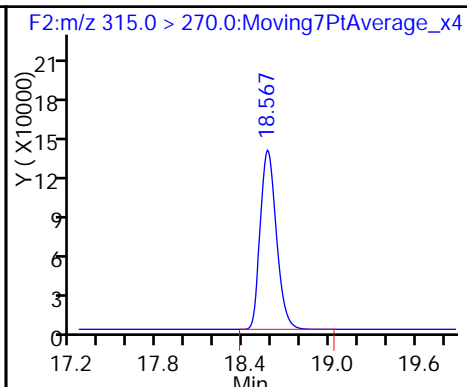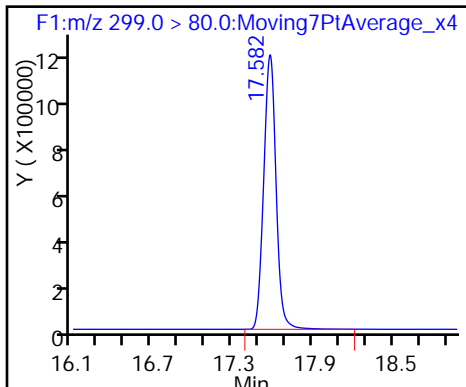

4 Perfluoroheptanoic acid

\* 5 13C2-PFOA




6 Perfluorooctanoic acid

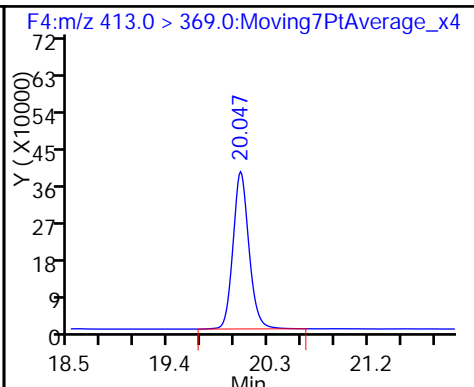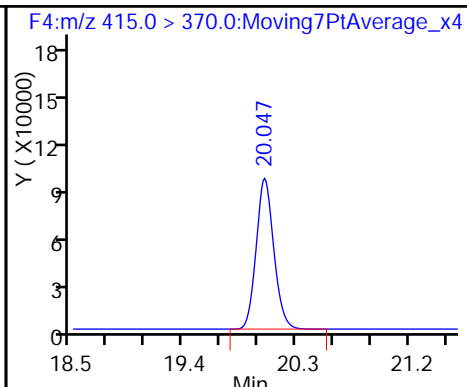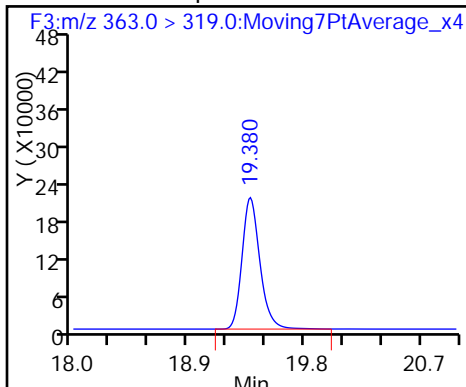


7 Perfluorooctane sulfonic acid

\* 8 13C4 PFOS

9 Perfluorononanoic acid

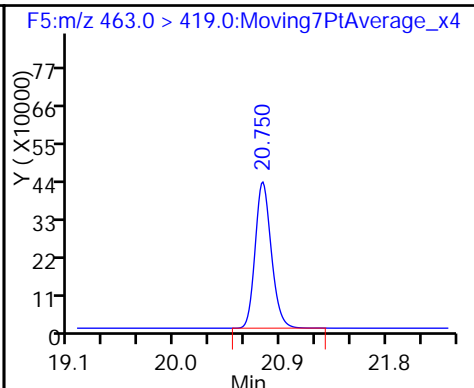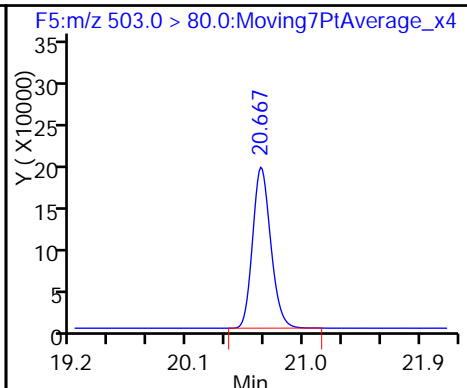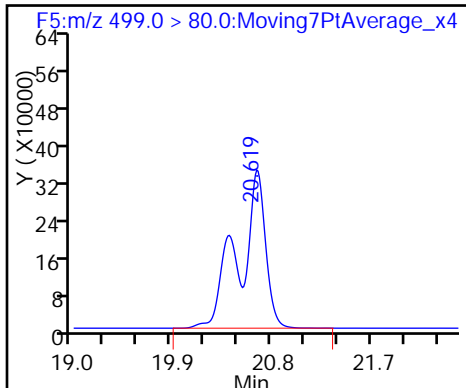

\$ 10 13C2 PFDA

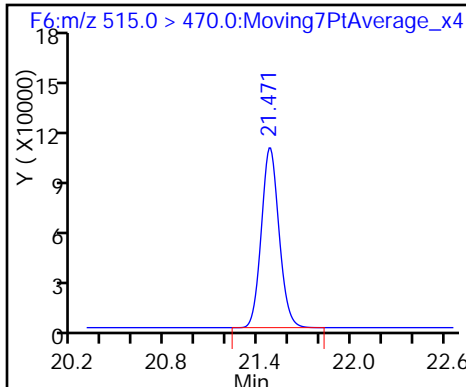

FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Lab Sample ID: CCV 320-140688/9 Calibration Date: 12/05/2016 20:53  
 Instrument ID: A6 Calib Start Date: 12/05/2016 17:26  
 GC Column: Acquity ID: 2.10 (mm) Calib End Date: 12/05/2016 19:54  
 Lab File ID: 05DEC2016A6A\_011.d Conc. Units: ng/mL

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
Perfluorobutanesulfonic acid (PFBS)	Ave	0.7015	0.6306		20.6	22.9	-10.1	50.0
Perfluorohexanesulfonic acid	Ave	0.8980	0.7822		6.72	7.72	-12.9	50.0
Perfluoroheptanoic acid	Ave	1.215	1.239		2.65	2.60	1.9	50.0
Perfluorooctanoic acid (PFOA)	Ave	1.040	0.9133		4.54	5.17	-12.2	50.0
Perfluorooctanesulfonic acid (PFOS)	Ave	1.044	0.8902		8.71	10.2	-14.7	50.0
Perfluorononanoic acid	Ave	1.134	1.093		4.83	5.01	-3.6	50.0
13C2 PFHxA	Ave	1.167	1.081		9.27	10.0	-7.3	30.0
13C2 PFDA	Ave	0.8763	0.8211		9.37	10.0	-6.3	30.0

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_011.d  
 Lims ID: CCV L2  
 Client ID:  
 Sample Type: CCVL  
 Inject. Date: 05-Dec-2016 20:53:12 ALS Bottle#: 2 Worklist Smp#: 9  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: CCV L2 CCV L2  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist: chrom-537\_\_A6\*sub3  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 06-Dec-2016 16:35:40 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK024

First Level Reviewer: barnettj Date: 06-Dec-2016 10:08:33

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.586	17.581	0.005	1.000	1186753	20.6	693
\$ 2 13C2 PFHxA	315.0 > 270.0	18.567	18.567	0.0	1.000	1108698	9.27	35970
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.344	19.342	0.002	1.000	496197	6.72	11535
4 Perfluoroheptanoic acid	363.0 > 319.0	19.380	19.378	0.002	1.000	329772	2.65	166 M
* 5 13C2-PFOA	415.0 > 370.0	20.047	20.047	0.0		1025187	10.0	21492
6 Perfluorooctanoic acid	413.0 > 369.0	20.047	20.047	0.0	1.000	484196	4.54	93.2 M
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.619	20.619	0.0	1.000	747766	8.71	8549
* 8 13C4 PFOS	503.0 > 80.0	20.667	20.669	-0.002		2358079	28.7	20478
9 Perfluorononanoic acid	463.0 > 419.0	20.750	20.748	0.002	1.000	561371	4.83	15032
\$ 10 13C2 PFDA	515.0 > 470.0	21.471	21.474	-0.003	1.000	841818	9.37	26813

### QC Flag Legend

Review Flags

M - Manually Integrated

### Reagents:

LC537-L2\_00014

Amount Added: 1.00

Units: mL

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_011.d

Injection Date: 05-Dec-2016 20:53:12

Instrument ID: A6

Lims ID: CCV L2

Client ID:

Operator ID: CBW

ALS Bottle#: 2

Worklist Smp#: 9

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

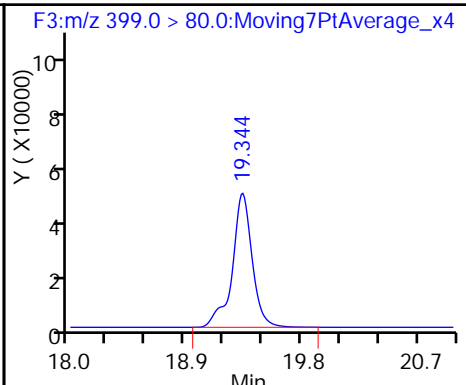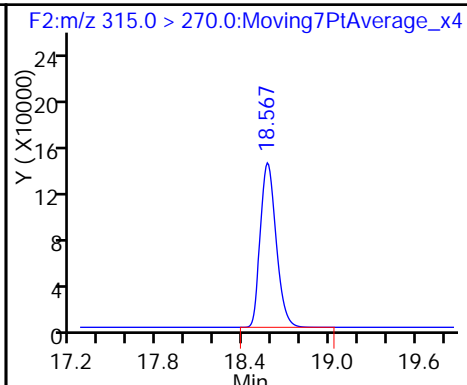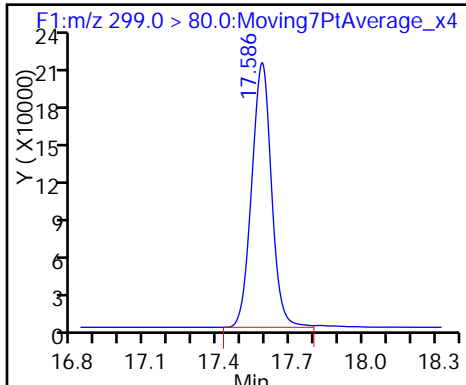

4 Perfluoroheptanoic acid (M)

\* 5 13C2-PFOA




6 Perfluorooctanoic acid (M)

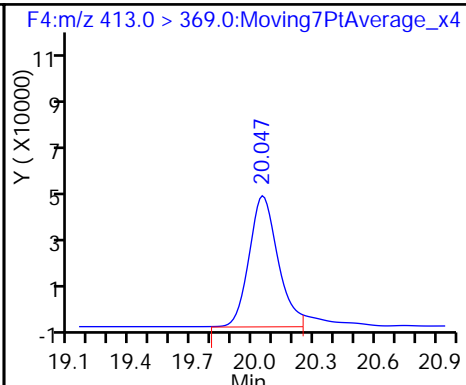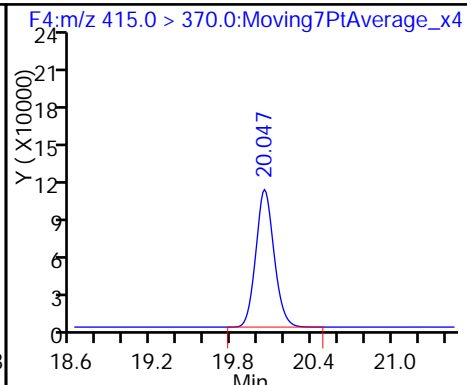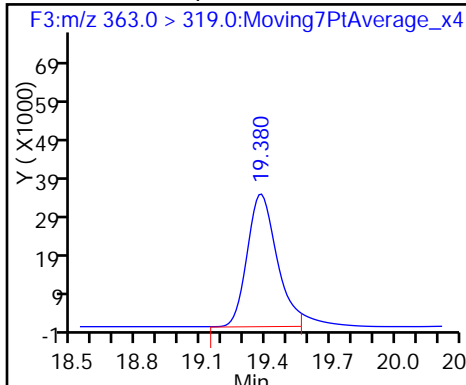


7 Perfluorooctane sulfonic acid

\* 8 13C4 PFOS

9 Perfluorononanoic acid

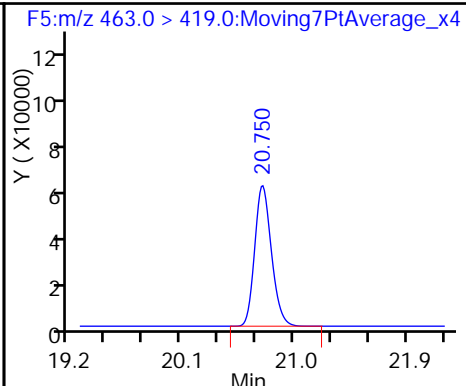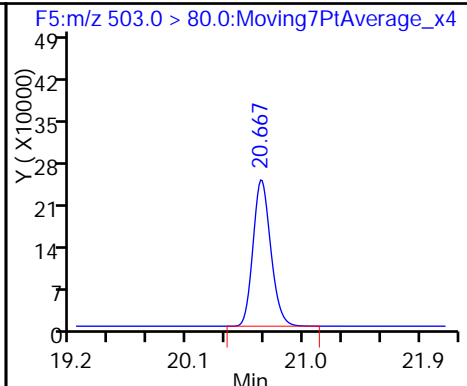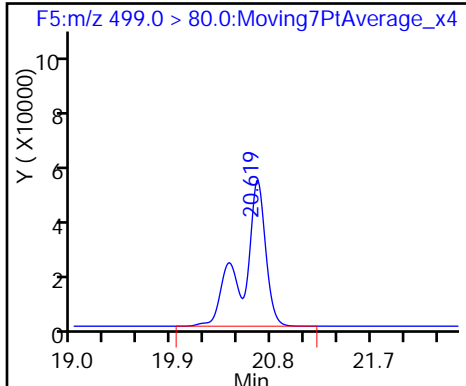

\$ 10 13C2 PFDA

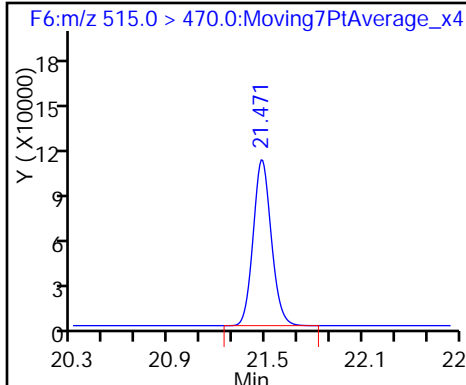

TestAmerica Sacramento


Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_011.d  
Injection Date: 05-Dec-2016 20:53:12 Instrument ID: A6  
Lims ID: CCV L2  
Client ID:  
Operator ID: CBW ALS Bottle#: 2 Worklist Smp#: 9  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Method: 537\_\_A6 Limit Group: LC 537 ICAL  
Column: Acquity BEH C18 ( 2.10 mm) Detector F3:M/RM

4 Perfluoroheptanoic acid, CAS: 375-85-9

Signal: 1


RT: 19.38  
Area: 349162  
Amount: 2.802857  
Amount Units: ng/ml

Processing Integration Results

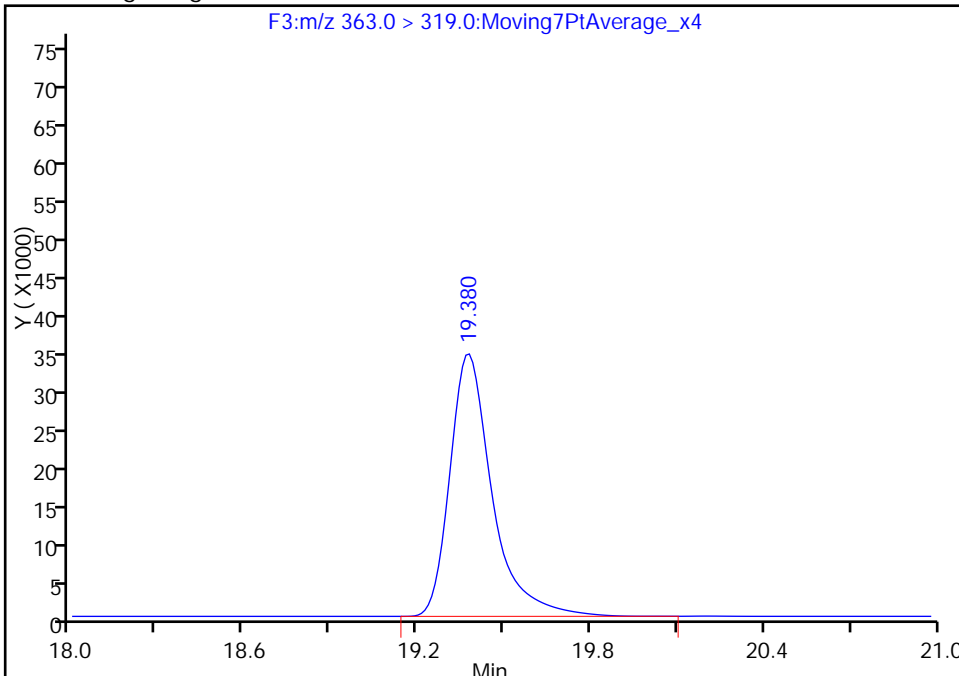

RT: 19.38  
Area: 329772  
Amount: 2.647206  
Amount Units: ng/ml

Manual Integration Results

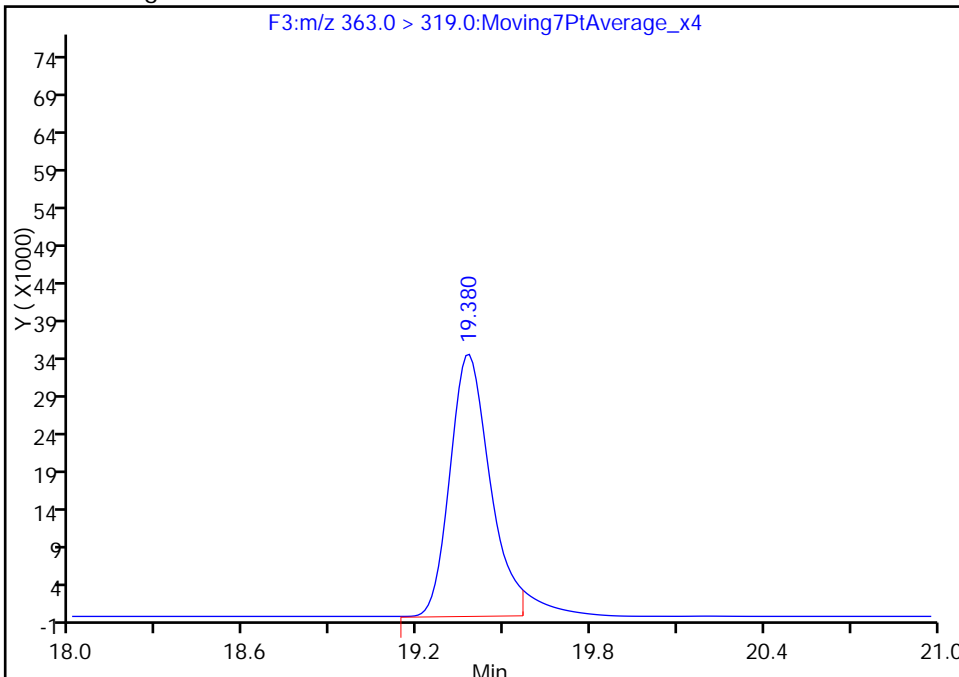

Reviewer: barnettj, 06-Dec-2016 10:08:33  
Audit Action: Manually Integrated

Audit Reason: Split Peak

TestAmerica Sacramento


Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_011.d  
Injection Date: 05-Dec-2016 20:53:12 Instrument ID: A6  
Lims ID: CCV L2  
Client ID:  
Operator ID: CBW ALS Bottle#: 2 Worklist Smp#: 9  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Method: 537\_\_A6 Limit Group: LC 537 ICAL  
Column: Acquity BEH C18 ( 2.10 mm) Detector F4:MRM

6 Perfluorooctanoic acid, CAS: 335-67-1

Signal: 1


RT: 20.05  
Area: 520603  
Amount: 4.880820  
Amount Units: ng/ml

Processing Integration Results

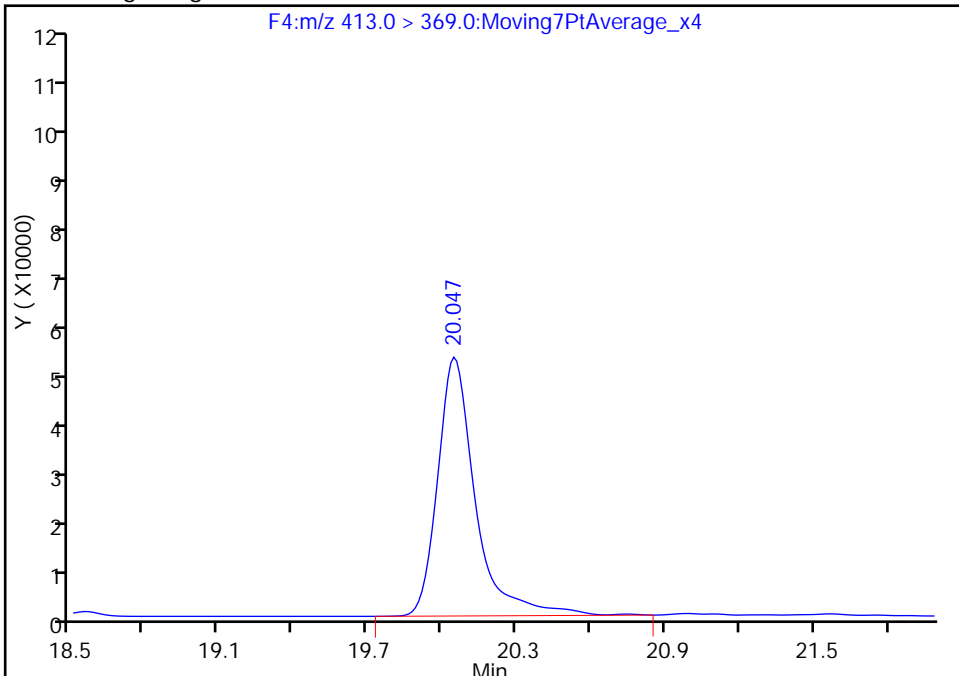

RT: 20.05  
Area: 484196  
Amount: 4.539493  
Amount Units: ng/ml

Manual Integration Results

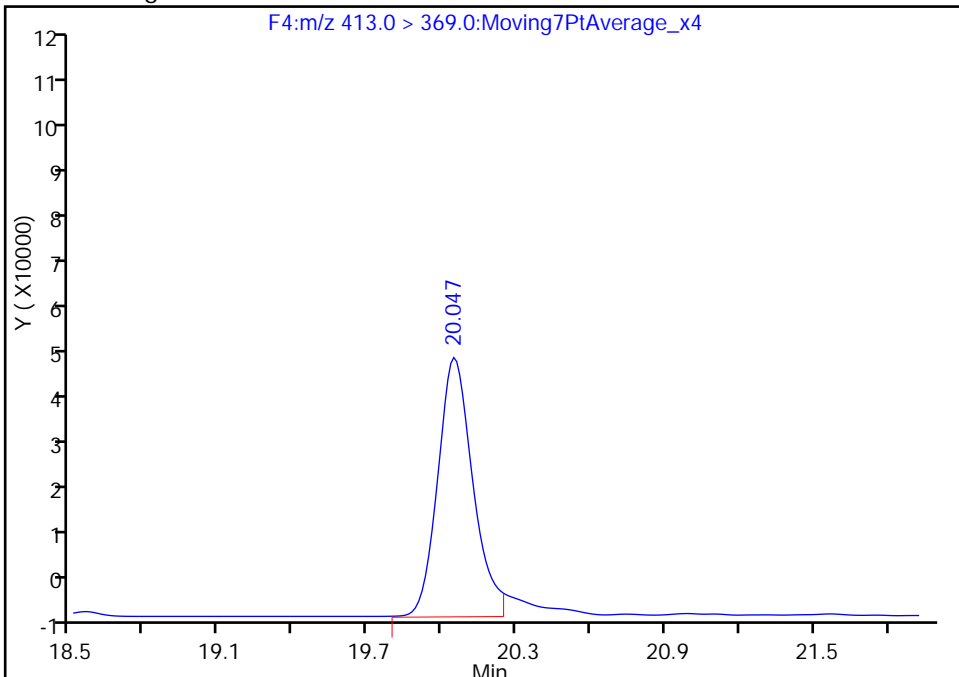

Reviewer: barnettj, 06-Dec-2016 10:08:33  
Audit Action: Manually Integrated

Audit Reason: Split Peak

FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Lab Sample ID: ICV 320-140688/11 Calibration Date: 12/05/2016 21:52  
 Instrument ID: A6 Calib Start Date: 12/05/2016 17:26  
 GC Column: Acquity ID: 2.10 (mm) Calib End Date: 12/05/2016 19:54  
 Lab File ID: 05DEC2016A6A\_013.d Conc. Units: ng/mL

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
Perfluorobutanesulfonic acid (PFBS)	Ave	0.7015	0.5756		94.2	115	-18.0	30.0
Perfluorohexanesulfonic acid	Ave	0.8980	0.6976		20.6	26.5	-22.3	30.0
Perfluoroheptanoic acid	Ave	1.215	1.155		11.9	12.5	-4.9	30.0
Perfluorooctanoic acid (PFOA)	Ave	1.040	0.9604		23.2	25.1	-7.7	30.0
Perfluorooctanesulfonic acid (PFOS)	Ave	1.044	0.8424		22.0	27.2	-19.3	30.0
Perfluorononanoic acid	Ave	1.134	0.9316		20.6	25.1	-17.9	30.0
13C2 PFHxA	Ave	1.167	1.079		9.25	10.0	-7.5	30.0
13C2 PFDA	Ave	0.8763	0.8628		9.85	10.0	-1.5	30.0



TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_013.d  
 Lims ID: ICV  
 Client ID:  
 Sample Type: ICV  
 Inject. Date: 05-Dec-2016 21:52:24 ALS Bottle#: 7 Worklist Smp#: 11  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: ICV ICV  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist:

Method: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 06-Dec-2016 16:53:23 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d

Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK024

First Level Reviewer: barnettj Date: 06-Dec-2016 16:34:53

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
--------	----	--------	--------	--------	----------	--------------	-----	-------

1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.582	17.581	0.001	1.000	4641388	94.2	8629
\$ 2 13C2 PFHxA	315.0 > 270.0	18.567	18.567	0.0	1.000	946677	9.25	29673
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.344	19.342	0.002	1.000	1298107	20.6	29738
4 Perfluoroheptanoic acid	363.0 > 319.0	19.380	19.378	0.002	1.000	1267011	11.9	9991
* 5 13C2-PFOA	415.0 > 370.0	20.047	20.047	0.0		877210	10.0	22431
6 Perfluorooctanoic acid	413.0 > 369.0	20.047	20.047	0.0	1.000	2114272	23.2	647
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.619	20.619	0.0	1.000	1612191	22.0	13496
* 8 13C4 PFOS	503.0 > 80.0	20.667	20.669	-0.002		2015178	28.7	51574
9 Perfluorononanoic acid	463.0 > 419.0	20.750	20.748	0.002	1.000	2051048	20.6	7161
\$ 10 13C2 PFDA	515.0 > 470.0	21.480	21.474	0.006	1.000	756809	9.85	23714

Reagents:

LC537-ICV\_00017 Amount Added: 1.00 Units: mL

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_013.d

Injection Date: 05-Dec-2016 21:52:24

Instrument ID: A6

Lims ID: ICV

Client ID:

Operator ID: CBW

ALS Bottle#: 7

Worklist Smp#: 11

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

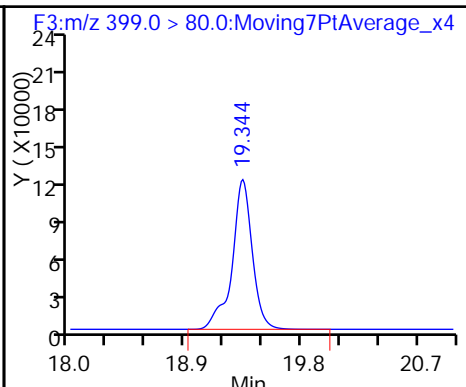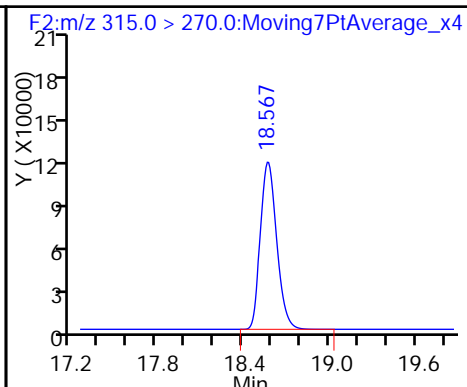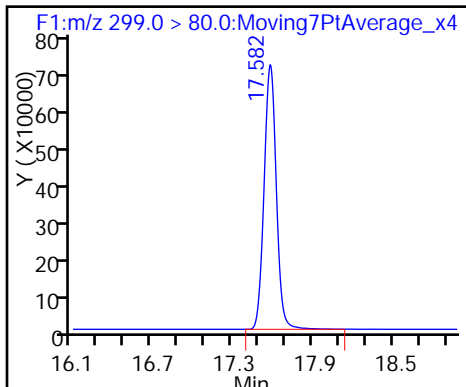

4 Perfluoroheptanoic acid

\* 5 13C2-PFOA




6 Perfluorooctanoic acid

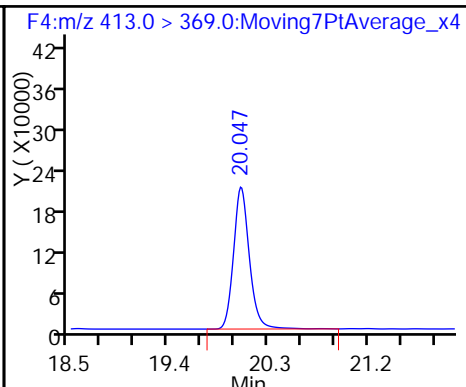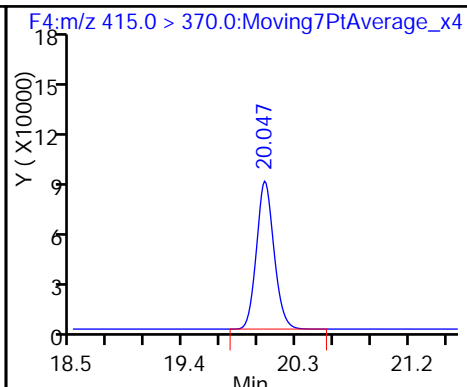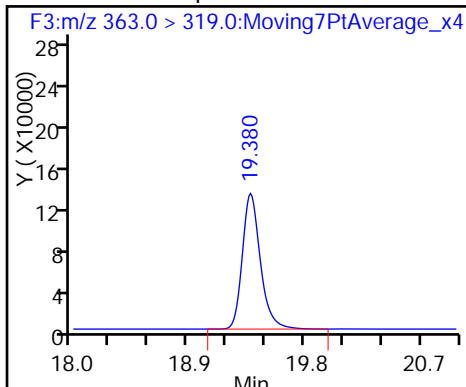


7 Perfluorooctane sulfonic acid

\* 8 13C4 PFOS

9 Perfluorononanoic acid

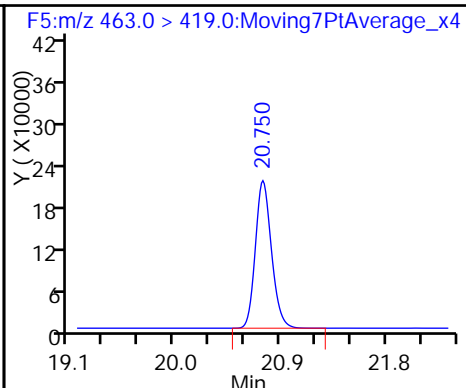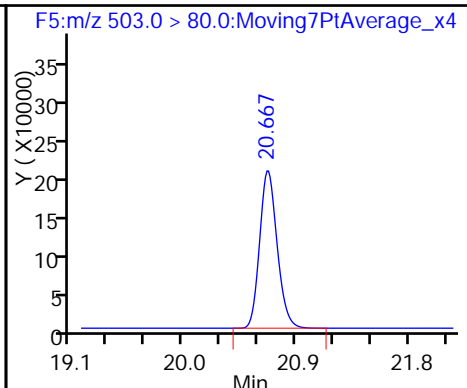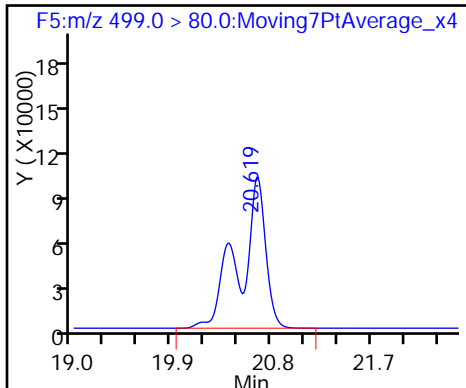

\$ 10 13C2 PFDA

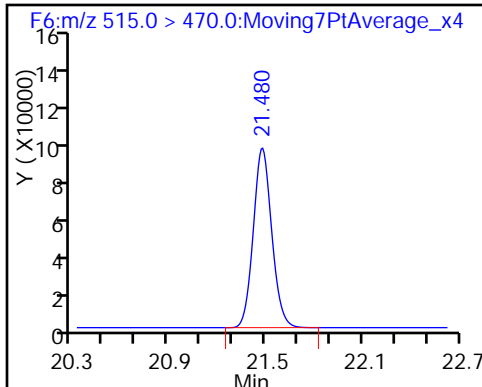

FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Lab Sample ID: CCV 320-142223/2 Calibration Date: 12/15/2016 08:03  
 Instrument ID: A6 Calib Start Date: 12/05/2016 17:26  
 GC Column: Acquity ID: 2.10 (mm) Calib End Date: 12/05/2016 19:54  
 Lab File ID: 15DEC2016A6A\_002.d Conc. Units: ng/mL

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
Perfluorobutanesulfonic acid (PFBS)	Ave	0.7015	0.6594		21.5	22.9	-6.0	50.0
Perfluorohexanesulfonic acid	Ave	0.8980	0.7732		6.64	7.72	-13.9	50.0
Perfluoroheptanoic acid	Ave	1.215	1.277		2.73	2.60	5.1	50.0
Perfluorooctanoic acid (PFOA)	Ave	1.040	0.9458		4.70	5.17	-9.1	50.0
Perfluorooctanesulfonic acid (PFOS)	Ave	1.044	0.8565		8.38	10.2	-18.0	50.0
Perfluorononanoic acid	Ave	1.134	1.091		4.82	5.01	-3.8	50.0
13C2 PFHxA	Ave	1.167	1.058		9.07	10.0	-9.3	30.0
13C2 PFDA	Ave	0.8763	0.8545		9.75	10.0	-2.5	30.0

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161214-37858.b\15DEC2016A6A\_002.d  
 Lims ID: CCV L2  
 Client ID:  
 Sample Type: CCVL  
 Inject. Date: 15-Dec-2016 08:03:38 ALS Bottle#: 2 Worklist Smp#: 2  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: CCV L2 CCV L2  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist: chrom-537\_\_A6\*sub3  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161214-37858.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 15-Dec-2016 15:21:12 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK005

First Level Reviewer: barnettj Date: 15-Dec-2016 15:05:28

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.599	17.599	0.0	1.000	845450	21.5	510
\$ 2 13C2 PFHxA	315.0 > 270.0	18.585	18.585	0.0	1.000	651205	9.07	21190
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.344	19.344	0.0	1.000	334175	6.64	7885
4 Perfluoroheptanoic acid	363.0 > 319.0	19.380	19.380	0.0	1.000	204042	2.73	115 M
* 5 13C2-PFOA	415.0 > 370.0	20.035	20.035	0.0		615516	10.0	64218
6 Perfluorooctanoic acid	413.0 > 369.0	20.047	20.047	0.0	1.000	301061	4.70	99.5 M
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.619	20.619	0.0	1.000	490144	8.38	6418
* 8 13C4 PFOS	503.0 > 80.0	20.667	20.667	0.0		1606471	28.7	33387
9 Perfluorononanoic acid	463.0 > 419.0	20.738	20.738	0.0	1.000	336444	4.82	9088
\$ 10 13C2 PFDA	515.0 > 470.0	21.471	21.471	0.0	1.000	525960	9.75	16666

### QC Flag Legend

Review Flags

M - Manually Integrated

### Reagents:

LC537-L2\_00014

Amount Added: 1.00

Units: mL

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161214-37858.b\15DEC2016A6A\_002.d

Injection Date: 15-Dec-2016 08:03:38

Instrument ID: A6

Lims ID: CCV L2

Client ID:

Operator ID: CBW

ALS Bottle#: 2

Worklist Smp#: 2

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

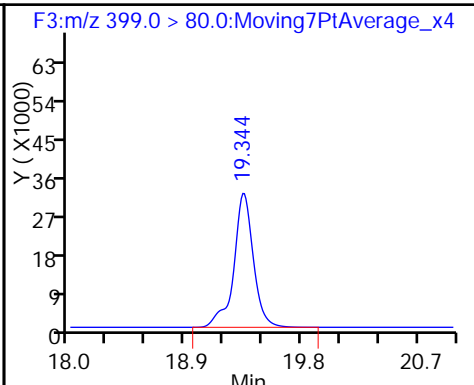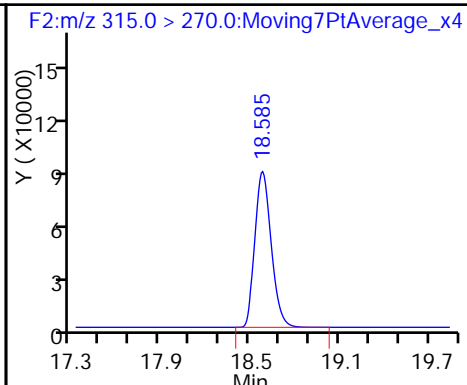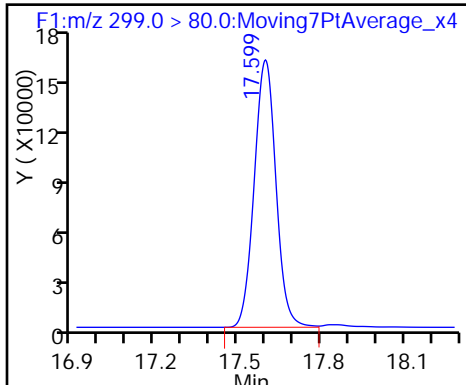

4 Perfluoroheptanoic acid (M)

\* 5 13C2-PFOA




6 Perfluorooctanoic acid (M)

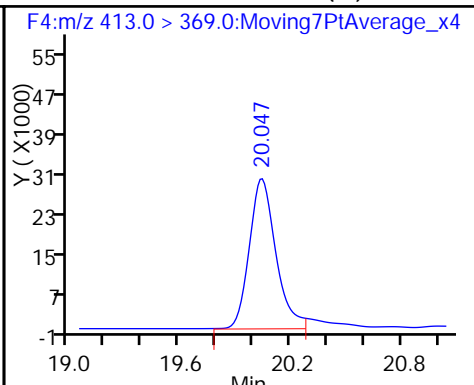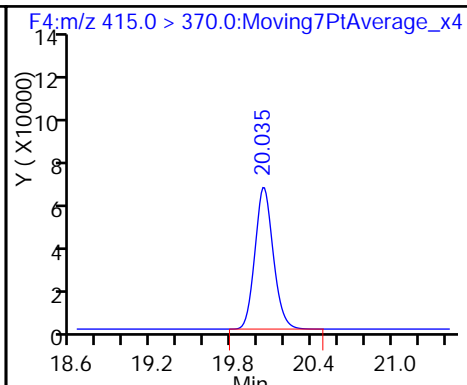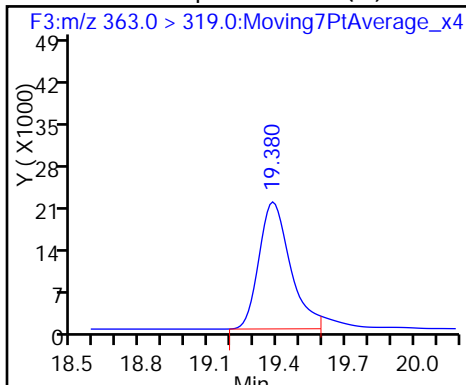


7 Perfluorooctane sulfonic acid

\* 8 13C4 PFOS

9 Perfluorononanoic acid

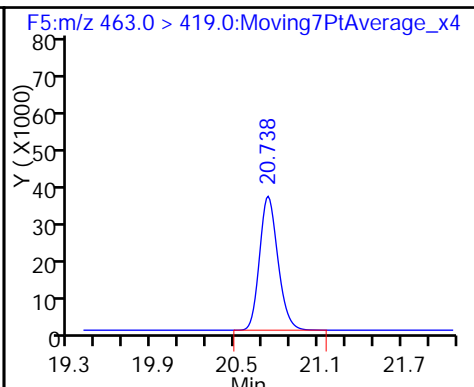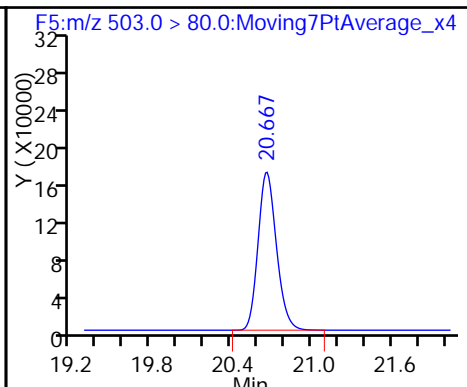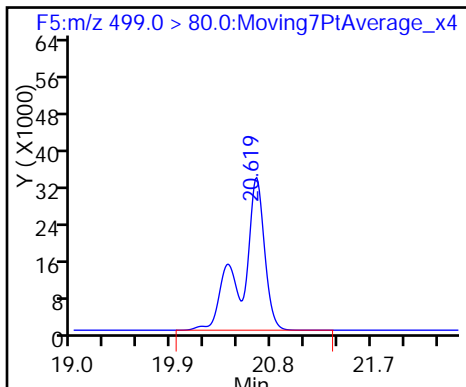

\$ 10 13C2 PFDA

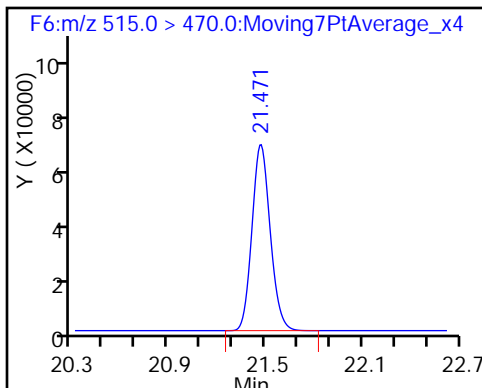

TestAmerica Sacramento


Data File: \\ChromNA\Sacramento\ChromData\A6\20161214-37858.b\15DEC2016A6A\_002.d  
Injection Date: 15-Dec-2016 08:03:38 Instrument ID: A6  
Lims ID: CCV L2  
Client ID:  
Operator ID: CBW ALS Bottle#: 2 Worklist Smp#: 2  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Method: 537\_\_A6 Limit Group: LC 537 ICAL  
Column: Acquity BEH C18 ( 2.10 mm) Detector F3:MRM

4 Perfluoroheptanoic acid, CAS: 375-85-9

Signal: 1


RT: 19.38  
Area: 213658  
Amount: 2.856649  
Amount Units: ng/ml

Processing Integration Results

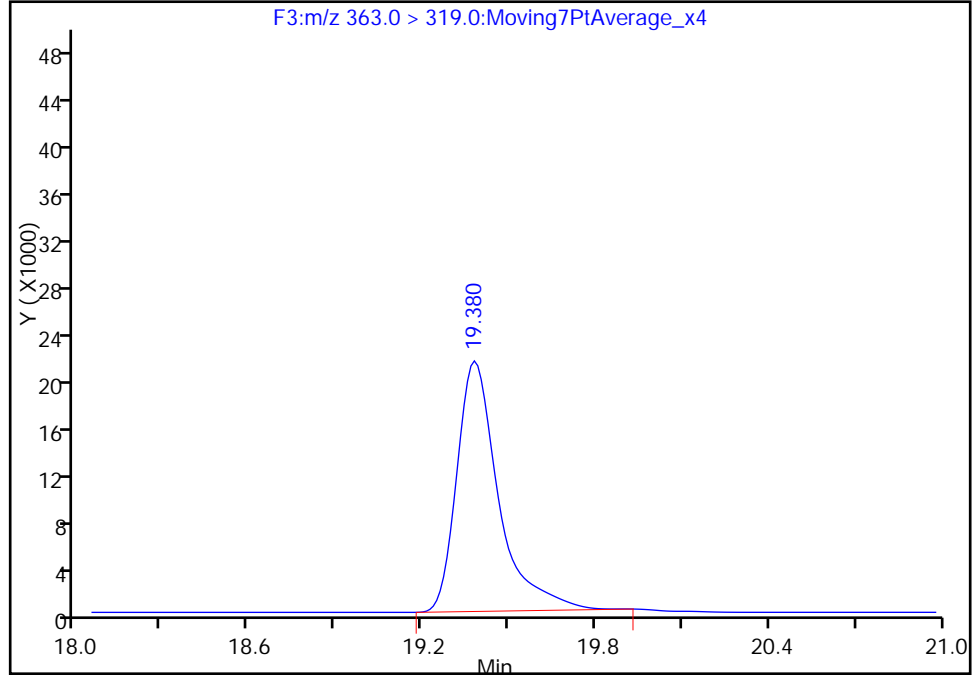

RT: 19.38  
Area: 204042  
Amount: 2.728081  
Amount Units: ng/ml

Manual Integration Results

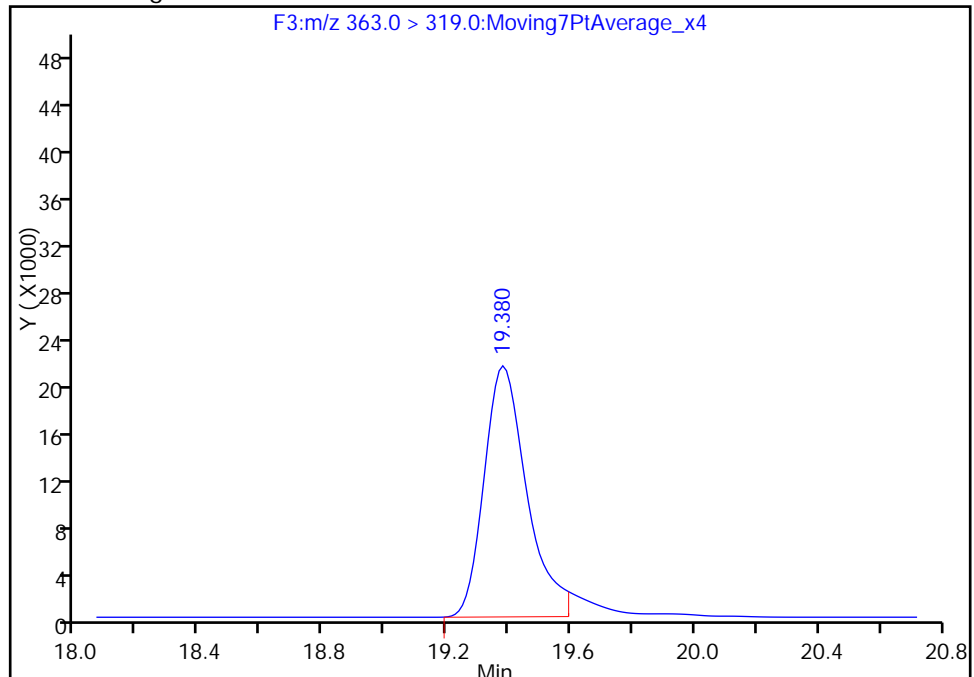

Reviewer: barnettj, 15-Dec-2016 15:05:28  
Audit Action: Manually Integrated

Audit Reason: Split Peak

TestAmerica Sacramento


Data File: \\ChromNA\Sacramento\ChromData\A6\20161214-37858.b\15DEC2016A6A\_002.d  
Injection Date: 15-Dec-2016 08:03:38 Instrument ID: A6  
Lims ID: CCV L2  
Client ID:  
Operator ID: CBW ALS Bottle#: 2 Worklist Smp#: 2  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Method: 537\_\_A6 Limit Group: LC 537 ICAL  
Column: Acquity BEH C18 ( 2.10 mm) Detector F4:M/RM

6 Perfluorooctanoic acid, CAS: 335-67-1

Signal: 1


RT: 20.05  
Area: 303357  
Amount: 4.737006  
Amount Units: ng/ml

Processing Integration Results

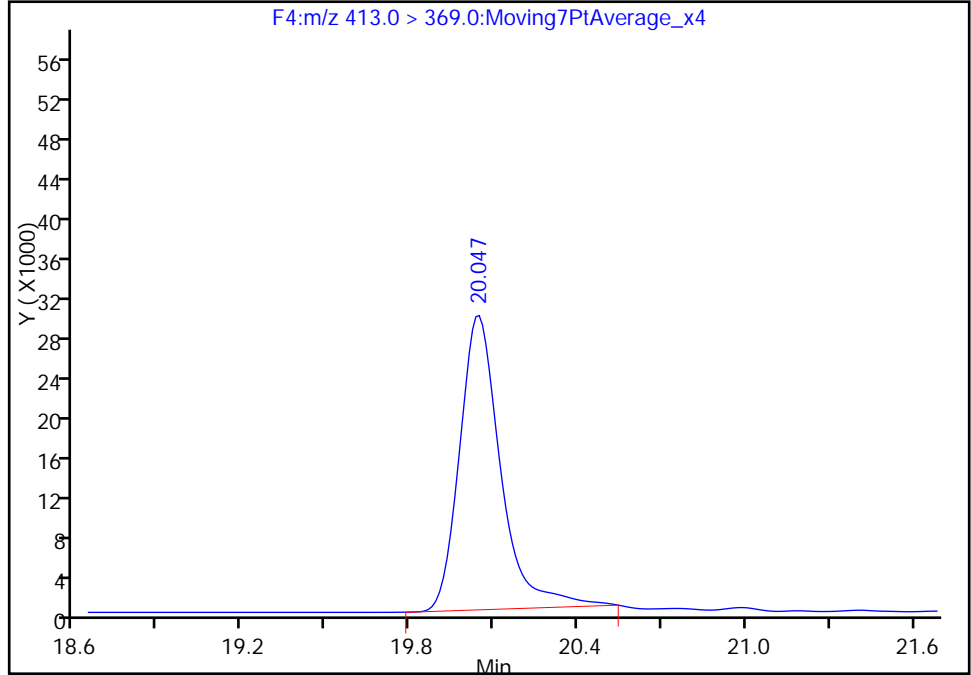

RT: 20.05  
Area: 301061  
Amount: 4.701153  
Amount Units: ng/ml

Manual Integration Results

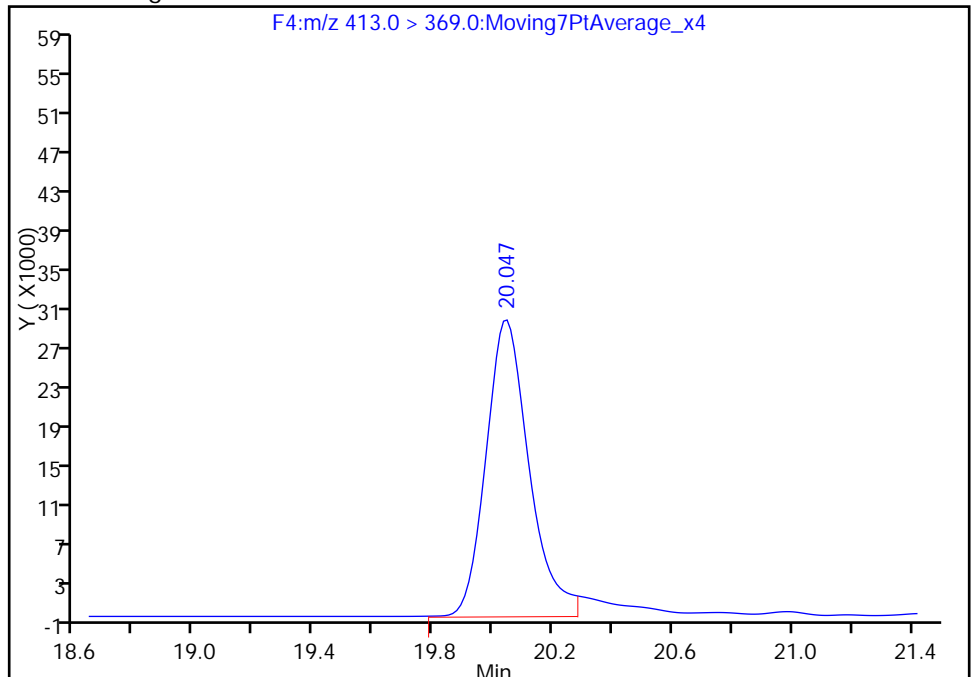

Reviewer: barnettj, 15-Dec-2016 15:05:28  
Audit Action: Manually Integrated

Audit Reason: Split Peak



FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Lab Sample ID: CCV 320-142809/35 Calibration Date: 12/18/2016 08:02  
 Instrument ID: A6 Calib Start Date: 12/05/2016 17:26  
 GC Column: Acquity ID: 2.10 (mm) Calib End Date: 12/05/2016 19:54  
 Lab File ID: 15DEC2016A6A\_147.d Conc. Units: ng/mL

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
Perfluorobutanesulfonic acid (PFBS)	Ave	0.7015	0.7192		46.2	45.1	2.5	30.0
Perfluorohexanesulfonic acid	Ave	0.8980	0.9716		16.5	15.2	8.2	30.0
Perfluoroheptanoic acid	Ave	1.215	1.346		5.67	5.12	10.8	30.0
Perfluorooctanoic acid (PFOA)	Ave	1.040	1.086		10.6	10.2	4.3	30.0
Perfluorooctanesulfonic acid (PFOS)	Ave	1.044	1.078		20.8	20.1	3.2	30.0
Perfluorononanoic acid	Ave	1.134	1.251		10.9	9.87	10.3	30.0
13C2 PFHxA	Ave	1.167	1.217		10.4	10.0	4.3	30.0
13C2 PFDA	Ave	0.8763	0.9157		10.5	10.0	4.5	30.0

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_147.d  
 Lims ID: CCV L3  
 Client ID:  
 Sample Type: CCVIS  
 Inject. Date: 18-Dec-2016 08:02:05 ALS Bottle#: 3 Worklist Smp#: 35  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: CCV L3 CCV L3  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist: chrom-537\_\_A6\*sub3  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 19-Dec-2016 09:17:18 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK034

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.570	17.570	0.0	1.000	1830394	46.2	968
\$ 2 13C2 PFHxA	315.0 > 270.0	18.548	18.548	0.0	1.000	921117	10.4	30310
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.308	19.308	0.0	1.000	833525	16.5	19547
4 Perfluoroheptanoic acid	363.0 > 319.0	19.344	19.344	0.0	1.000	521525	5.67	10784
* 5 13C2-PFOA	415.0 > 370.0	19.999	19.999	0.0		756954	10.0	19161
6 Perfluorooctanoic acid	413.0 > 369.0	19.999	19.999	0.0	1.000	837328	10.6	563
* 8 13C4 PFOS	503.0 > 80.0	20.619	20.619	0.0		1618248	28.7	41918
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.619	20.619	0.0	1.000	1224213	20.8	20146
9 Perfluorononanoic acid	463.0 > 419.0	20.691	20.691	0.0	1.000	934676	10.9	32754
\$ 10 13C2 PFDA	515.0 > 470.0	21.418	21.418	0.0	1.000	693170	10.5	21454

Reagents:

LC537-L3\_00016 Amount Added: 1.00 Units: mL

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_147.d

Injection Date: 18-Dec-2016 08:02:05

Instrument ID: A6

Lims ID: CCV L3

Client ID:

Operator ID: CBW

ALS Bottle#: 3

Worklist Smp#: 35

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

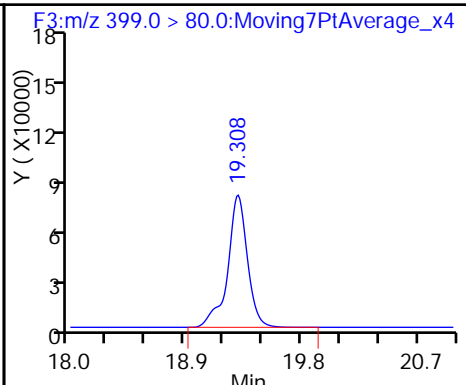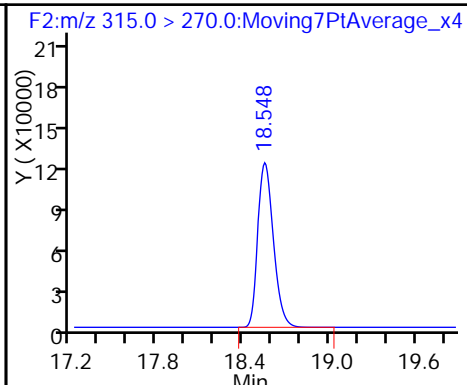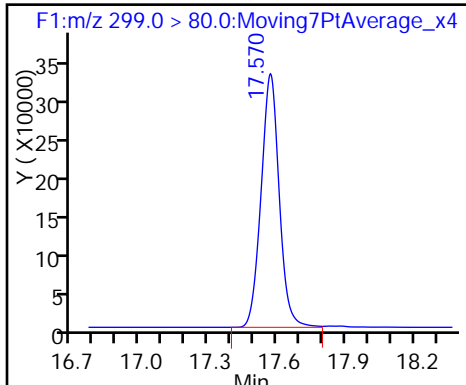

4 Perfluoroheptanoic acid

\* 5 13C2-PFOA




6 Perfluorooctanoic acid

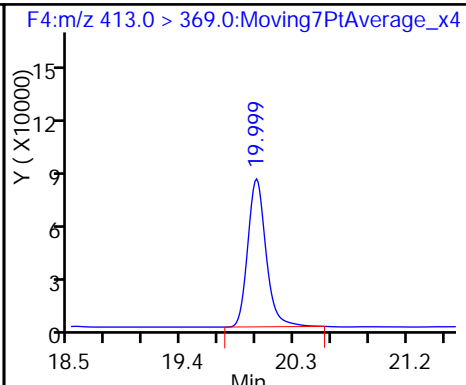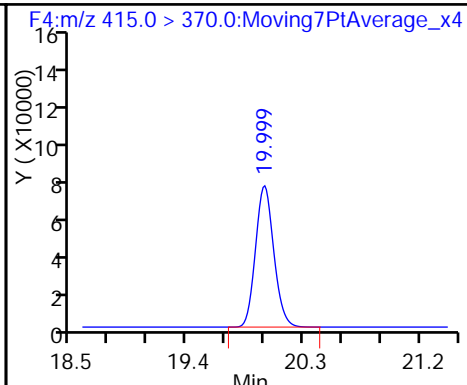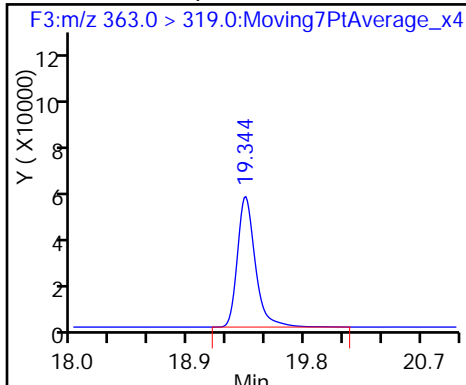


\* 8 13C4 PFOS

7 Perfluorooctane sulfonic acid

9 Perfluorononanoic acid

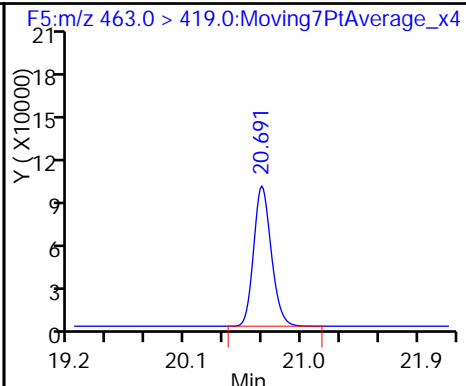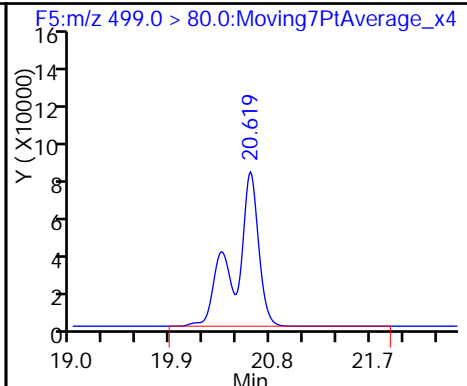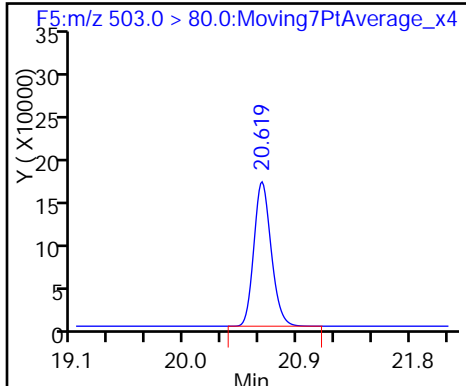

\$ 10 13C2 PFDA

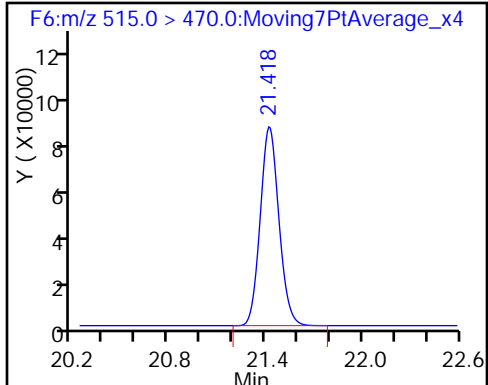

FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Lab Sample ID: CCV 320-142809/47 Calibration Date: 12/18/2016 13:57  
 Instrument ID: A6 Calib Start Date: 12/05/2016 17:26  
 GC Column: Acquity ID: 2.10 (mm) Calib End Date: 12/05/2016 19:54  
 Lab File ID: 15DEC2016A6A\_159.d Conc. Units: ng/mL

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
Perfluorobutanesulfonic acid (PFBS)	Ave	0.7015	0.7229		139	135	3.1	30.0
Perfluorohexanesulfonic acid	Ave	0.8980	1.032		52.2	45.4	14.9	30.0
Perfluoroheptanoic acid	Ave	1.215	1.224		15.4	15.3	0.7	30.0
Perfluorooctanoic acid (PFOA)	Ave	1.040	1.117		32.7	30.4	7.4	30.0
Perfluorooctanesulfonic acid (PFOS)	Ave	1.044	1.207		69.5	60.1	15.6	30.0
Perfluorononanoic acid	Ave	1.134	1.236		32.1	29.5	9.0	30.0
13C2 PFHxA	Ave	1.167	1.317		11.3	10.0	12.9	30.0
13C2 PFDA	Ave	0.8763	1.001		11.4	10.0	14.2	30.0

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_159.d  
 Lims ID: CCV L5  
 Client ID:  
 Sample Type: CCVIS  
 Inject. Date: 18-Dec-2016 13:57:17 ALS Bottle#: 5 Worklist Smp#: 47  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: CCV L5 CCV L5  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist: chrom-537\_\_A6\*sub3  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 19-Dec-2016 09:17:47 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK034

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
--------	----	--------	--------	--------	----------	--------------	-----	-------

1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.576	17.576	0.0	1.000	6747178	138.8	7798
\$ 2 13C2 PFHxA	315.0 > 270.0	18.549	18.549	0.0	1.000	1282676	11.3	42364
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.308	19.308	0.0	1.000	3247496	52.2	73666
4 Perfluoroheptanoic acid	363.0 > 319.0	19.344	19.344	0.0	1.000	1821123	15.4	21131
* 5 13C2-PFOA	415.0 > 370.0	19.999	19.999	0.0		974301	10.0	24660
6 Perfluorooctanoic acid	413.0 > 369.0	19.999	19.999	0.0	1.000	3310451	32.7	1135
* 8 13C4 PFOS	503.0 > 80.0	20.619	20.619	0.0		1988143	28.7	22664
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.619	20.619	0.0	1.000	5026998	69.5	10359
9 Perfluorononanoic acid	463.0 > 419.0	20.691	20.691	0.0	1.000	3548285	32.1	21808
\$ 10 13C2 PFDA	515.0 > 470.0	21.427	21.427	0.0	1.000	975188	11.4	30448

Reagents:

LC537-L5\_00017 Amount Added: 1.00 Units: mL

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_159.d

Injection Date: 18-Dec-2016 13:57:17

Instrument ID: A6

Lims ID: CCV L5

Client ID:

Operator ID: CBW

ALS Bottle#: 5

Worklist Smp#: 47

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

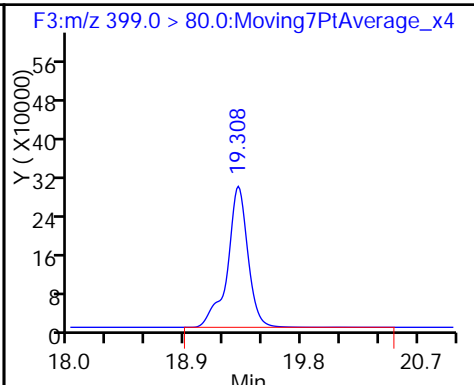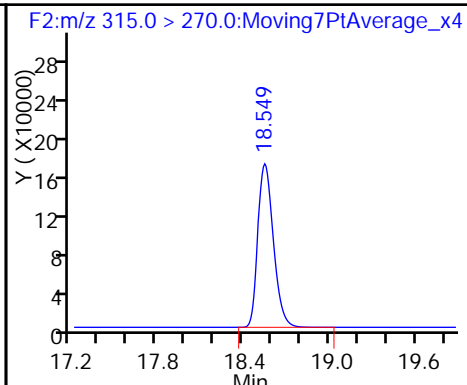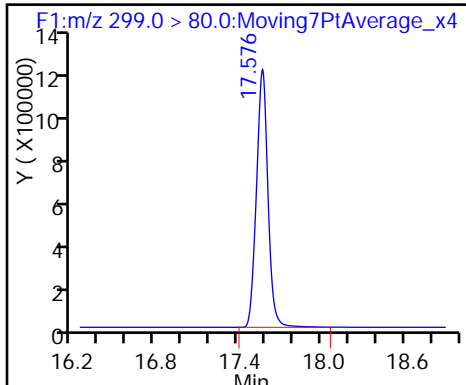

4 Perfluoroheptanoic acid

\* 5 13C2-PFOA




6 Perfluorooctanoic acid

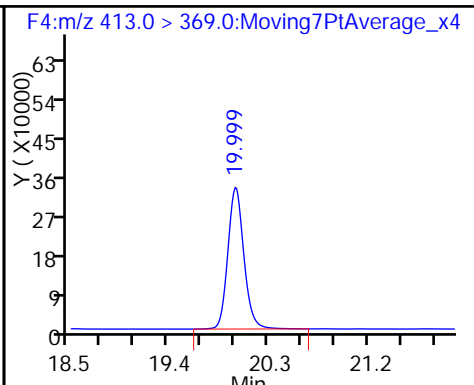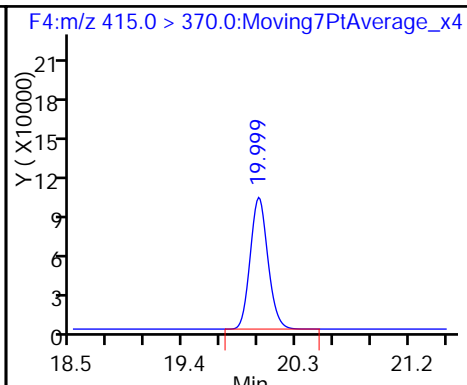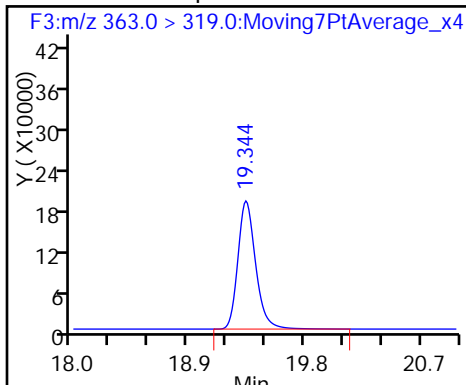


\* 8 13C4 PFOS

7 Perfluorooctane sulfonic acid

9 Perfluorononanoic acid

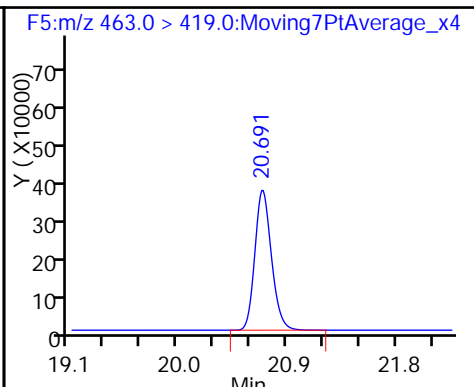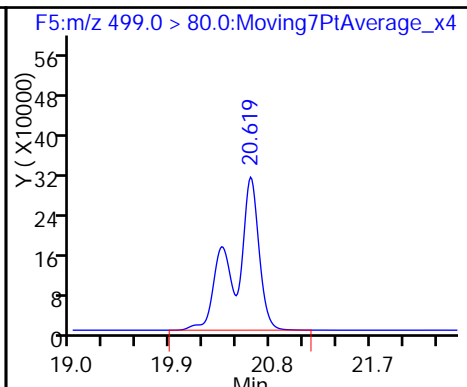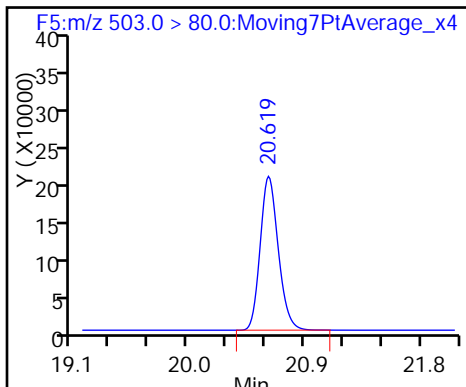

\$ 10 13C2 PFDA

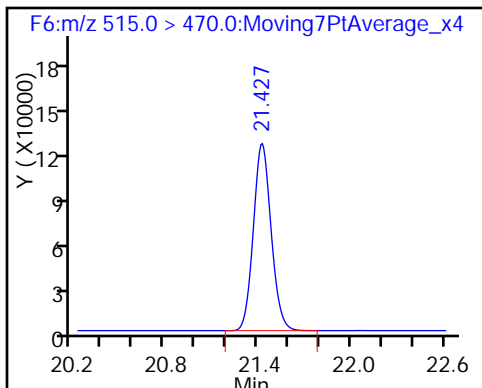

FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Lab Sample ID: CCV 320-142884/3 Calibration Date: 12/19/2016 09:45  
 Instrument ID: A6 Calib Start Date: 12/05/2016 17:26  
 GC Column: Acquity ID: 2.10 (mm) Calib End Date: 12/05/2016 19:54  
 Lab File ID: 19DEC2016A6A\_003.d Conc. Units: ng/mL

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
Perfluorobutanesulfonic acid (PFBS)	Ave	0.7015	0.6117		20.0	22.9	-12.8	50.0
Perfluorohexanesulfonic acid	Ave	0.8980	0.6980		6.00	7.72	-22.3	50.0
Perfluoroheptanoic acid	Ave	1.215	1.305		2.79	2.60	7.4	50.0
Perfluorooctanoic acid (PFOA)	Ave	1.040	0.9478		4.71	5.17	-8.9	50.0
Perfluorooctanesulfonic acid (PFOS)	Ave	1.044	0.8410		8.23	10.2	-19.4	50.0
Perfluorononanoic acid	Ave	1.134	1.022		4.51	5.01	-9.9	50.0
13C2 PFHxA	Ave	1.167	1.048		8.99	10.0	-10.1	30.0
13C2 PFDA	Ave	0.8763	0.8176		9.33	10.0	-6.7	30.0

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\19DEC2016A6A\_003.d  
 Lims ID: CCV L2  
 Client ID:  
 Sample Type: CCVL  
 Inject. Date: 19-Dec-2016 09:45:37 ALS Bottle#: 2 Worklist Smp#: 3  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: CCV L2 CCV L2  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist: chrom-537\_\_A6\*sub3  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 19-Dec-2016 13:22:17 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK025

First Level Reviewer: westendorfc Date: 19-Dec-2016 11:23:20

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.586	17.586	0.0	1.000	701934	20.0	516
\$ 2 13C2 PFHxA	315.0 > 270.0	18.549	18.549	0.0	1.000	642432	8.99	21202
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.308	19.308	0.0	1.000	269978	6.00	6558
4 Perfluoroheptanoic acid	363.0 > 319.0	19.344	19.344	0.0	1.000	207650	2.79	752
* 5 13C2-PFOA	415.0 > 370.0	19.999	19.999	0.0		612795	10.0	15691
6 Perfluorooctanoic acid	413.0 > 369.0	19.999	19.999	0.0	1.000	300363	4.71	247
* 8 13C4 PFOS	503.0 > 80.0	20.619	20.619	0.0		1437810	28.7	37350
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.619	20.619	0.0	1.000	430785	8.23	7306
9 Perfluorononanoic acid	463.0 > 419.0	20.691	20.691	0.0	1.000	313809	4.51	8317
\$ 10 13C2 PFDA	515.0 > 470.0	21.418	21.418	0.0	1.000	501037	9.33	15767

Reagents:

LC537-L2\_00014 Amount Added: 1.00 Units: mL



TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\19DEC2016A6A\_003.d

Injection Date: 19-Dec-2016 09:45:37

Instrument ID: A6

Lims ID: CCV L2

Client ID:

Operator ID: CBW

ALS Bottle#: 2

Worklist Smp#: 3

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

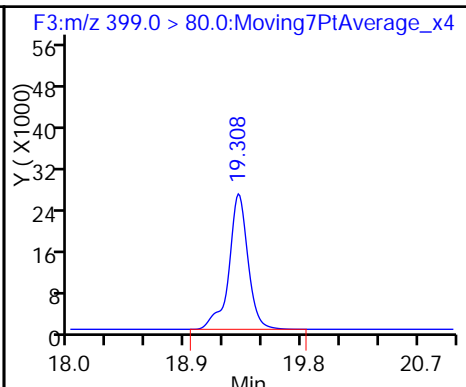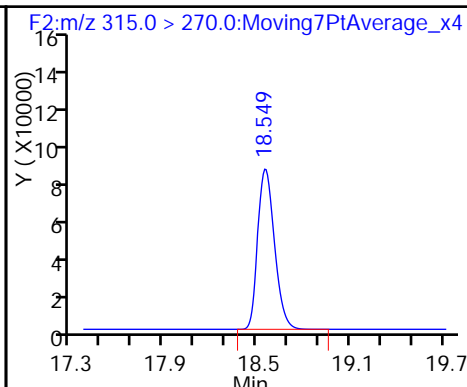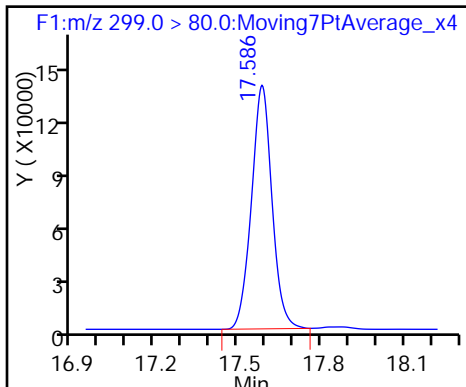

4 Perfluoroheptanoic acid

\* 5 13C2-PFOA




6 Perfluorooctanoic acid

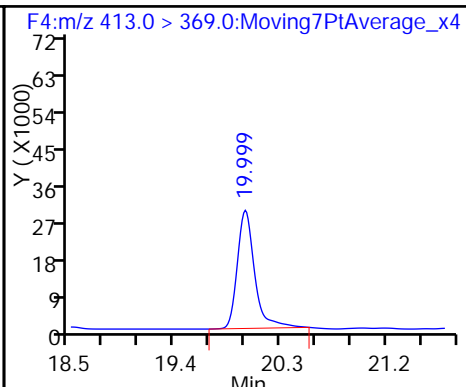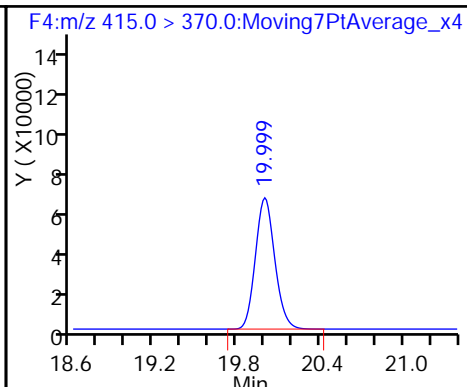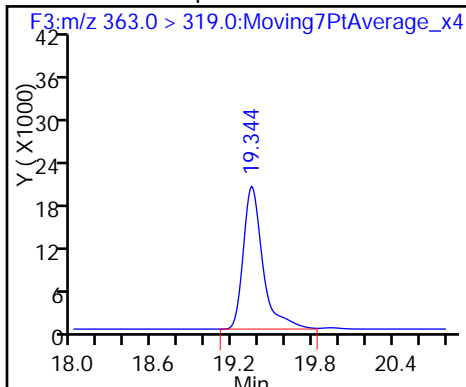


\* 8 13C4 PFOS

7 Perfluorooctane sulfonic acid

9 Perfluorononanoic acid

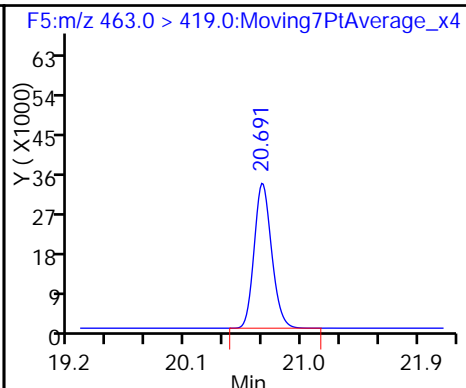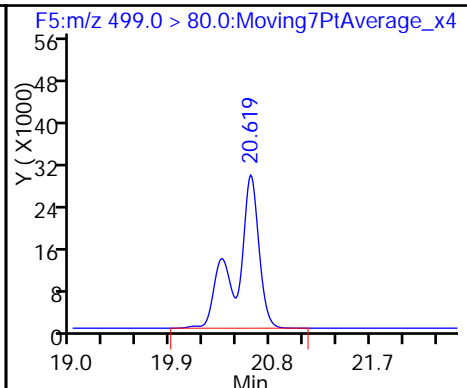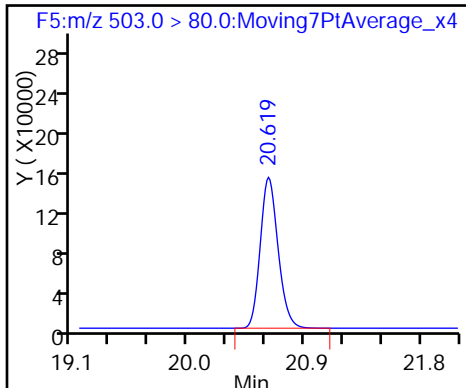

\$ 10 13C2 PFDA

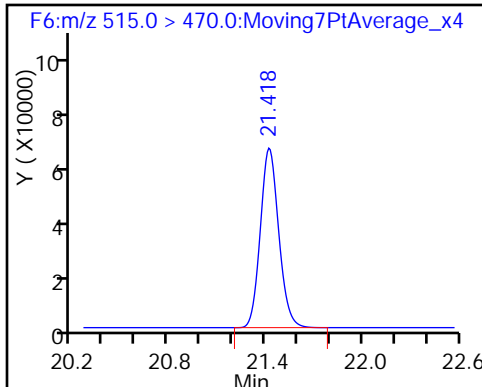

FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Lab Sample ID: CCV 320-142886/7 Calibration Date: 12/19/2016 11:51  
 Instrument ID: A6 Calib Start Date: 12/05/2016 17:26  
 GC Column: Acquity ID: 2.10 (mm) Calib End Date: 12/05/2016 19:54  
 Lab File ID: 19DEC2016A6A\_007.d Conc. Units: ng/mL

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
Perfluorobutanesulfonic acid (PFBS)	Ave	0.7015	0.7401		142	135	5.5	30.0
Perfluorohexanesulfonic acid	Ave	0.8980	1.015		51.3	45.4	13.0	30.0
Perfluoroheptanoic acid	Ave	1.215	1.230		15.5	15.3	1.2	30.0
Perfluorooctanoic acid (PFOA)	Ave	1.040	1.104		32.3	30.4	6.1	30.0
Perfluorooctanesulfonic acid (PFOS)	Ave	1.044	1.223		70.4	60.1	17.1	30.0
Perfluorononanoic acid	Ave	1.134	1.243		32.3	29.5	9.6	30.0
13C2 PFHxA	Ave	1.167	1.298		11.1	10.0	11.3	30.0
13C2 PFDA	Ave	0.8763	0.9668		11.0	10.0	10.3	30.0

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\19DEC2016A6A\_007.d  
 Lims ID: CCV L5  
 Client ID:  
 Sample Type: CCVIS  
 Inject. Date: 19-Dec-2016 11:51:43 ALS Bottle#: 5 Worklist Smp#: 7  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: CCV L5 CCV L5  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist: chrom-537\_\_A6\*sub3  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 19-Dec-2016 13:40:16 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK025

First Level Reviewer: westendorfc Date: 19-Dec-2016 12:55:25

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.589	17.589	0.0	1.000	5844162	142.0	13832
\$ 2 13C2 PFHxA	315.0 > 270.0	18.567	18.567	0.0	1.000	1017325	11.1	32771
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.320	19.320	0.0	1.000	2701403	51.3	35008
4 Perfluoroheptanoic acid	363.0 > 319.0	19.356	19.356	0.0	1.000	1472401	15.5	19025
* 5 13C2-PFOA	415.0 > 370.0	20.011	20.011	0.0		783659	10.0	19774
6 Perfluorooctanoic acid	413.0 > 369.0	20.011	20.011	0.0	1.000	2631135	32.3	1134
* 8 13C4 PFOS	503.0 > 80.0	20.619	20.619	0.0		1682080	28.7	28852
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.619	20.619	0.0	1.000	4310938	70.4	14931
9 Perfluorononanoic acid	463.0 > 419.0	20.702	20.702	0.0	1.000	2869432	32.3	49963
\$ 10 13C2 PFDA	515.0 > 470.0	21.418	21.418	0.0	1.000	757644	11.0	23668

Reagents:

LC537-L5\_00017 Amount Added: 1.00 Units: mL

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\19DEC2016A6A\_007.d

Injection Date: 19-Dec-2016 11:51:43

Instrument ID: A6

Lims ID: CCV L5

Client ID:

Operator ID: CBW

ALS Bottle#: 5

Worklist Smp#: 7

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

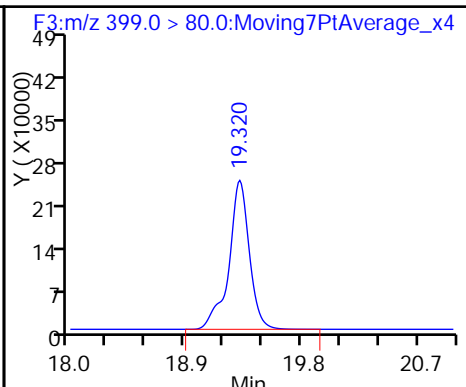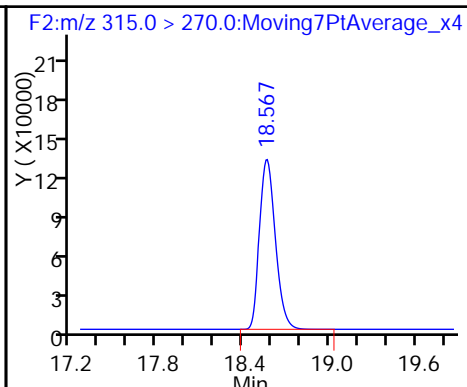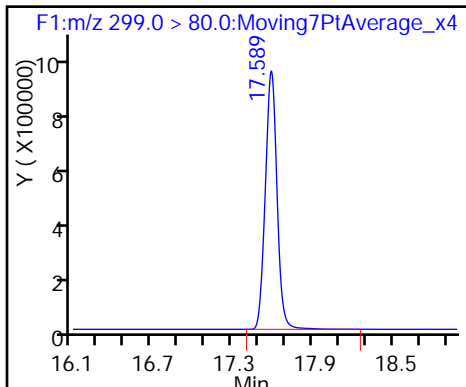

4 Perfluoroheptanoic acid

\* 5 13C2-PFOA




6 Perfluorooctanoic acid

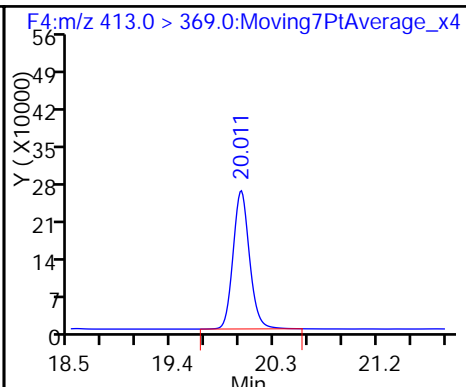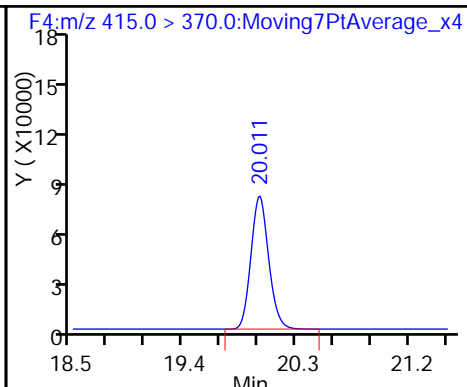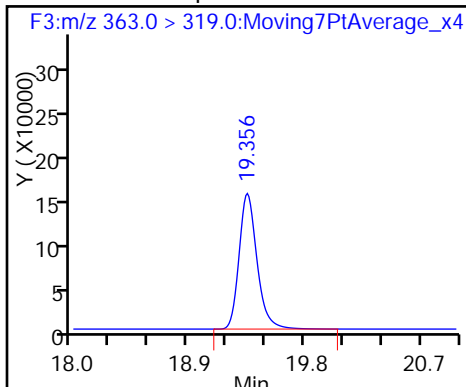


\* 8 13C4 PFOS

7 Perfluorooctane sulfonic acid

9 Perfluorononanoic acid

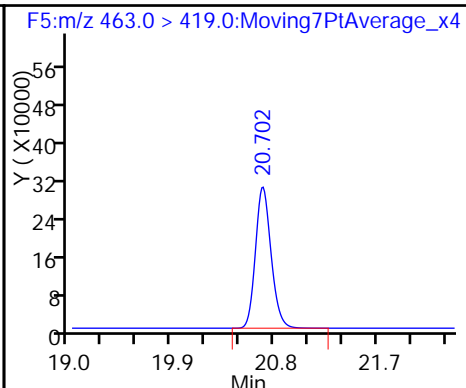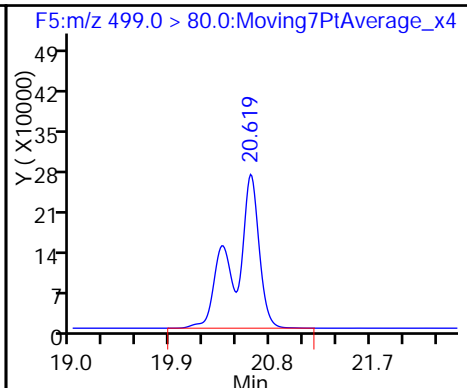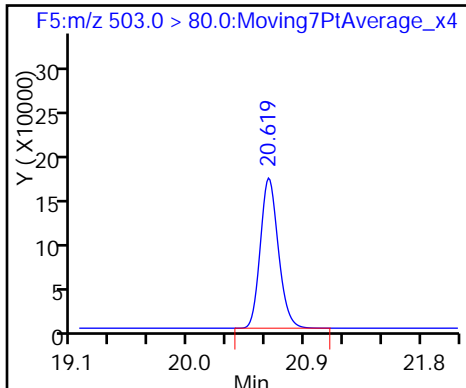

\$ 10 13C2 PFDA

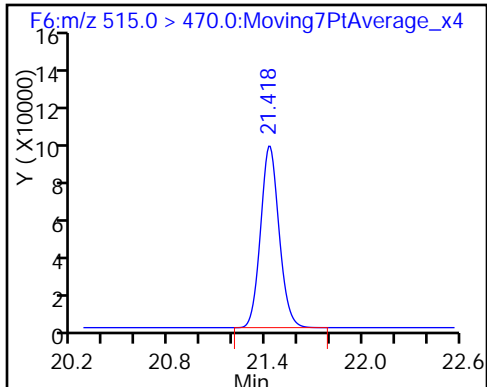

FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Lab Sample ID: CCV 320-142886/20 Calibration Date: 12/19/2016 18:32  
 Instrument ID: A6 Calib Start Date: 12/05/2016 17:26  
 GC Column: Acquity ID: 2.10 (mm) Calib End Date: 12/05/2016 19:54  
 Lab File ID: 19DEC2016A6A\_020.d Conc. Units: ng/mL

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
Perfluorobutanesulfonic acid (PFBS)	Ave	0.7015	0.7992		51.4	45.1	13.9	30.0
Perfluorohexanesulfonic acid	Ave	0.8980	1.000		16.9	15.2	11.4	30.0
Perfluoroheptanoic acid	Ave	1.215	1.276		5.37	5.12	5.0	30.0
Perfluorooctanoic acid (PFOA)	Ave	1.040	1.021		10.0	10.2	-1.9	30.0
Perfluorooctanesulfonic acid (PFOS)	Ave	1.044	1.113		21.5	20.1	6.6	30.0
Perfluorononanoic acid	Ave	1.134	1.200		10.4	9.87	5.8	30.0
13C2 PFHxA	Ave	1.167	1.256		10.8	10.0	7.6	30.0
13C2 PFDA	Ave	0.8763	0.9083		10.4	10.0	3.7	30.0

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\19DEC2016A6A\_020.d  
 Lims ID: CCV L3  
 Client ID:  
 Sample Type: CCVIS  
 Inject. Date: 19-Dec-2016 18:32:21 ALS Bottle#: 3 Worklist Smp#: 20  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: CCV L3 CCV L3  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Sublist: chrom-537\_\_A6\*sub3  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 20-Dec-2016 13:39:17 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK032

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.599	17.599	0.0	1.000	2306810	51.4	1842
\$ 2 13C2 PFHxA	315.0 > 270.0	18.567	18.567	0.0	1.000	1007576	10.8	32573
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.320	19.320	0.0	1.000	973342	16.9	22287
4 Perfluoroheptanoic acid	363.0 > 319.0	19.356	19.356	0.0	1.000	523886	5.37	13426
* 5 13C2-PFOA	415.0 > 370.0	20.011	20.011	0.0		802443	10.0	20585
6 Perfluorooctanoic acid	413.0 > 369.0	20.011	20.011	0.0	1.000	834771	10.0	375
* 8 13C4 PFOS	503.0 > 80.0	20.619	20.619	0.0		1835426	28.7	26920
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.631	20.631	0.0	1.000	1433927	21.5	22811
9 Perfluorononanoic acid	463.0 > 419.0	20.702	20.702	0.0	1.000	950270	10.4	14441
\$ 10 13C2 PFDA	515.0 > 470.0	21.427	21.427	0.0	1.000	728855	10.4	23033

Reagents:

LC537-L3\_00016 Amount Added: 1.00 Units: mL

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\19DEC2016A6A\_020.d

Injection Date: 19-Dec-2016 18:32:21

Instrument ID: A6

Lims ID: CCV L3

Client ID:

Operator ID: CBW

ALS Bottle#: 3

Worklist Smp#: 20

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

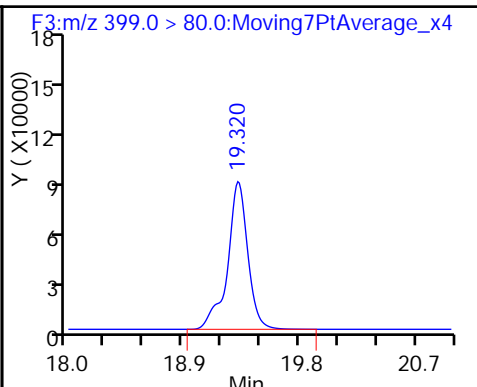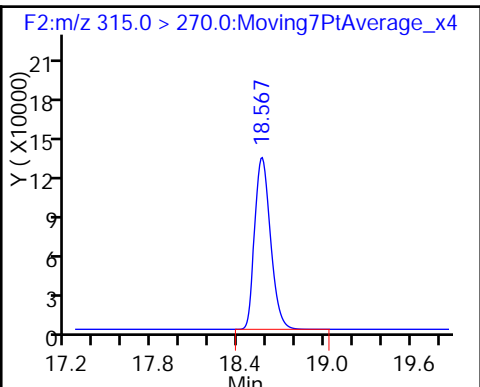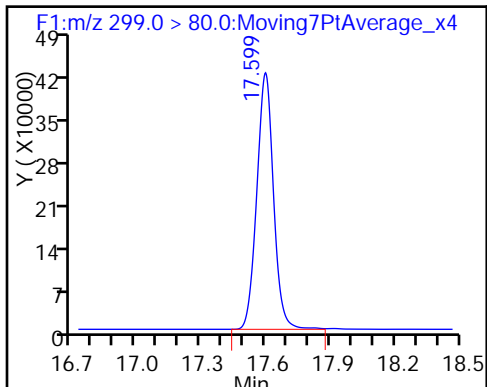

4 Perfluoroheptanoic acid

\* 5 13C2-PFOA




6 Perfluorooctanoic acid

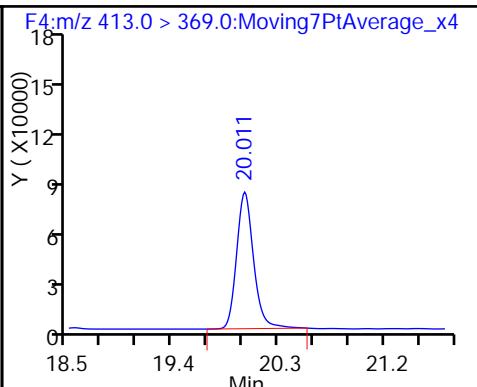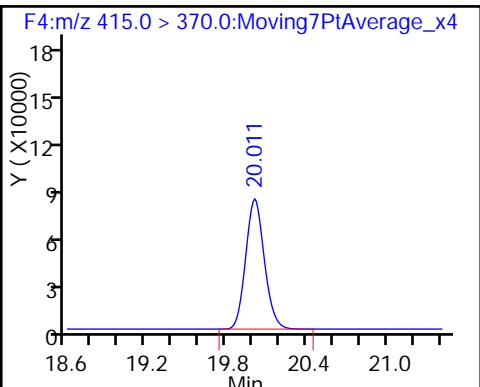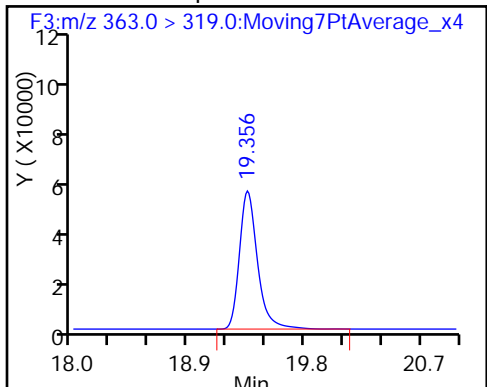


\* 8 13C4 PFOS

7 Perfluorooctane sulfonic acid

9 Perfluorononanoic acid

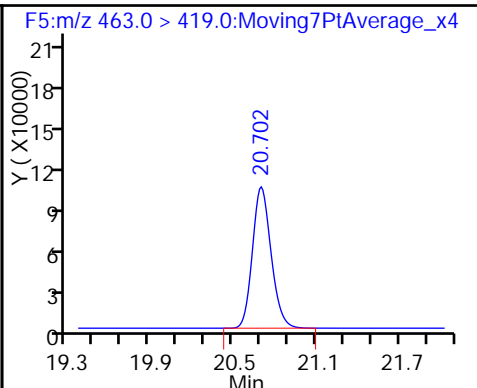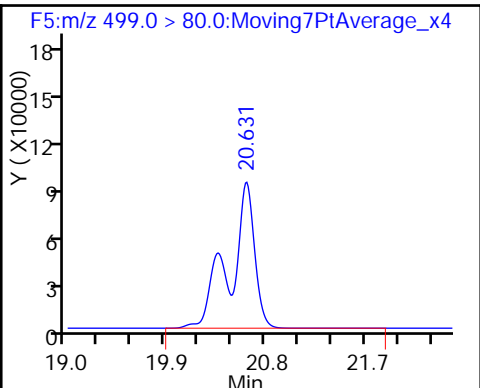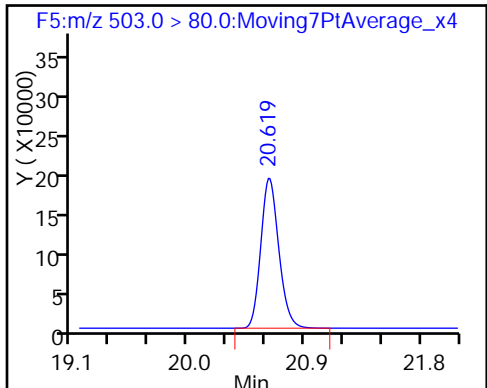

\$ 10 13C2 PFDA

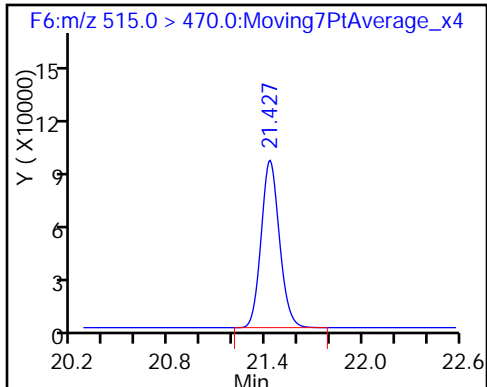

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Client Sample ID: \_\_\_\_\_ Lab Sample ID: MB 320-141642/1-A  
 Matrix: Water Lab File ID: 19DEC2016A6A\_010.d  
 Analysis Method: 537 Date Collected: \_\_\_\_\_  
 Extraction Method: 537 Date Extracted: 12/12/2016 10:03  
 Sample wt/vol: 250 (mL) Date Analyzed: 12/19/2016 13:19  
 Con. Extract Vol.: 1.0 (mL) Dilution Factor: 1  
 Injection Volume: 10 (uL) GC Column: Acquity ID: 2.1 (mm)  
 % Moisture: \_\_\_\_\_ GPC Cleanup: (Y/N) N  
 Analysis Batch No.: 142886 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	LOQ	LOD	DL
1763-23-1	Perfluorooctanesulfonic acid (PFOS)	0.048	U	0.060	0.048	0.016
335-67-1	Perfluorooctanoic acid (PFOA)	0.024	U	0.030	0.024	0.0094
375-73-5	Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.14	0.11	0.048

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL00993	13C2 PFHxA	110		70-130
STL00996	13C2 PFDA	111		70-130



TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\19DEC2016A6A\_010.d  
 Lims ID: MB 320-141642/1-A  
 Client ID:  
 Sample Type: MB  
 Inject. Date: 19-Dec-2016 13:19:16 ALS Bottle#: 34 Worklist Smp#: 10  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: mb 320-141642/1-a BOX 23  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 20-Dec-2016 13:38:59 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK032

First Level Reviewer: barnettj Date: 20-Dec-2016 09:39:34

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
--------	----	--------	--------	--------	----------	--------------	-----	-------

\$ 2 13C2 PFHxA	315.0 > 270.0	18.567	18.567	0.0	1.000	956133	11.0	30810
* 5 13C2-PFOA	415.0 > 370.0	20.011	20.011	0.0		747991	10.0	19095
* 8 13C4 PFOS	503.0 > 80.0	20.631	20.619	0.012		2054845	28.7	30485
9 Perfluorononanoic acid								M
463.0 > 419.0	20.702	20.702	0.0	1.000	1507	0.0178	40.8	M
\$ 10 13C2 PFDA	515.0 > 470.0	21.427	21.418	0.009	1.000	730640	11.1	22919

QC Flag Legend

Review Flags

M - Manually Integrated

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\19DEC2016A6A\_010.d

Injection Date: 19-Dec-2016 13:19:16

Instrument ID: A6

Lims ID: MB 320-141642/1-A

Client ID:

Operator ID: CBW

ALS Bottle#: 34

Worklist Smp#: 10

Injection Vol: 10.0 ul




Dil. Factor: 1.0000

Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid (ND) \$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid (ND)

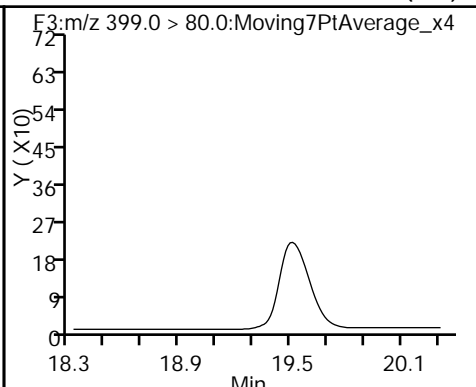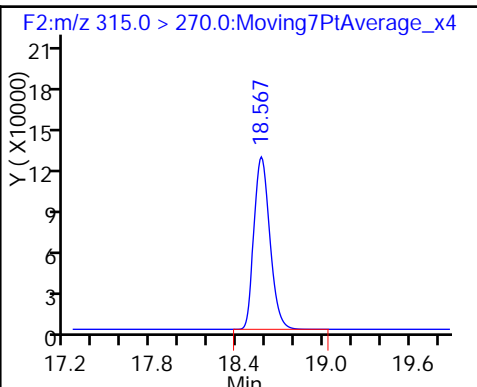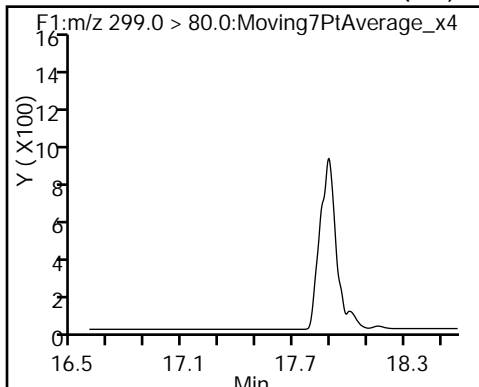

4 Perfluoroheptanoic acid (ND)

\* 5 13C2-PFOA




6 Perfluorooctanoic acid (ND)

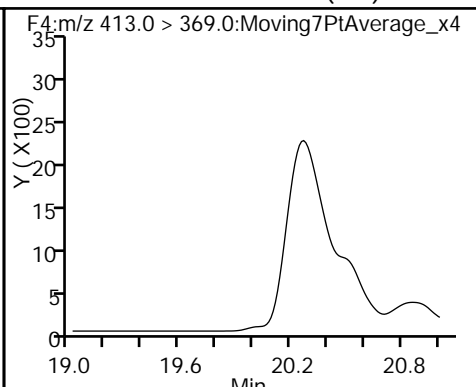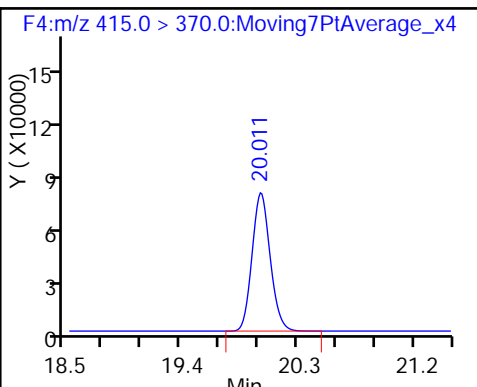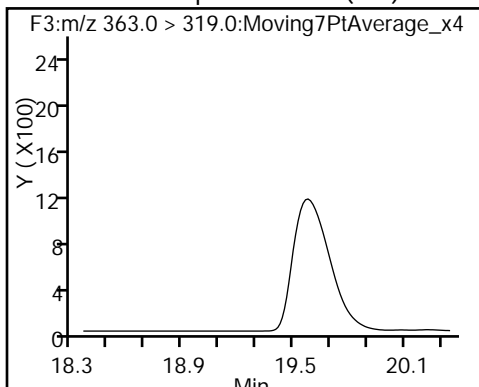


\* 8 13C4 PFOS

7 Perfluorooctane sulfonic acid (ND)

9 Perfluorononanoic acid (M)

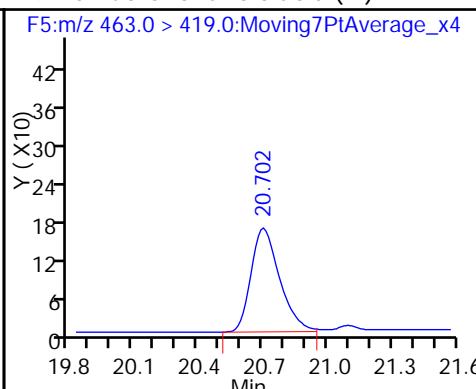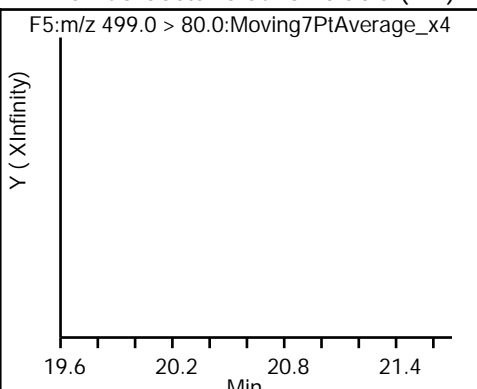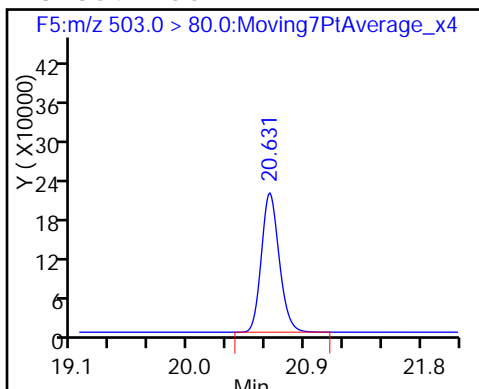

\$ 10 13C2 PFDA

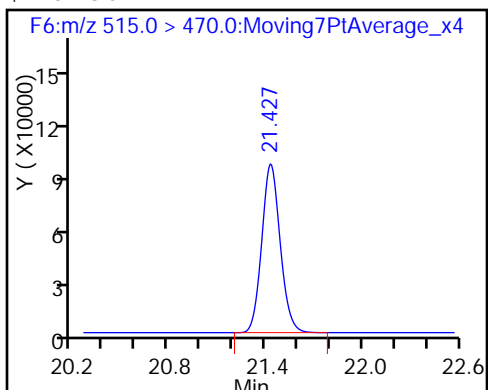

TestAmerica Sacramento  
Recovery Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\19DEC2016A6A\_010.d  
 Lims ID: MB 320-141642/1-A  
 Client ID:  
 Sample Type: MB  
 Inject. Date: 19-Dec-2016 13:19:16 ALS Bottle#: 34 Worklist Smp#: 10  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: mb 320-141642/1-a BOX 23  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35°C  
 Operator ID: CBW Instrument ID: A6  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 20-Dec-2016 13:38:59 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK032

First Level Reviewer: barnettj Date: 20-Dec-2016 09:39:34

Compound	Amount Added	Amount Recovered	% Rec.
\$ 2 13C2 PFHxA	10.0	11.0	109.58
\$ 10 13C2 PFDA	10.0	11.1	111.47

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Client Sample ID: \_\_\_\_\_ Lab Sample ID: LCS 320-141642/2-A  
 Matrix: Water Lab File ID: 19DEC2016A6A\_011.d  
 Analysis Method: 537 Date Collected: \_\_\_\_\_  
 Extraction Method: 537 Date Extracted: 12/12/2016 10:03  
 Sample wt/vol: 250 (mL) Date Analyzed: 12/19/2016 13:48  
 Con. Extract Vol.: 1.0 (mL) Dilution Factor: 1  
 Injection Volume: 10 (uL) GC Column: Acquity ID: 2.1 (mm)  
 % Moisture: \_\_\_\_\_ GPC Cleanup: (Y/N) N  
 Analysis Batch No.: 142886 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	LOQ	LOD	DL
1763-23-1	Perfluorooctanesulfonic acid (PFOS)	0.150		0.060	0.048	0.016
335-67-1	Perfluorooctanoic acid (PFOA)	0.0704		0.030	0.024	0.0094
375-73-5	Perfluorobutanesulfonic acid (PFBS)	0.334		0.14	0.11	0.048

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL00993	13C2 PFHxA	118		70-130
STL00996	13C2 PFDA	123		70-130

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\19DEC2016A6A\_011.d  
 Lims ID: LCS 320-141642/2-A  
 Client ID:  
 Sample Type: LCS  
 Inject. Date: 19-Dec-2016 13:48:50 ALS Bottle#: 35 Worklist Smp#: 11  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: lcs 320-141642/2-a  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 20-Dec-2016 13:38:59 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK032

First Level Reviewer: barnettj Date: 20-Dec-2016 09:37:03

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.602	17.589	0.013	1.000	3699394	83.4	1485
\$ 2 13C2 PFHxA	315.0 > 270.0	18.576	18.567	0.009	1.000	984523	11.8	31648
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.332	19.320	0.012	1.000	1648661	29.0	37509
4 Perfluoroheptanoic acid	363.0 > 319.0	19.368	19.356	0.012	1.000	775683	8.92	13323
* 5 13C2-PFOA	415.0 > 370.0	20.023	20.011	0.012		715286	10.0	18339
6 Perfluorooctanoic acid	413.0 > 369.0	20.023	20.011	0.012	1.000	1309461	17.6	494
* 8 13C4 PFOS	503.0 > 80.0	20.631	20.619	0.012		1813273	28.7	15426
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.631	20.619	0.012	1.000	2481036	37.6	14338
9 Perfluorononanoic acid	463.0 > 419.0	20.714	20.702	0.012	1.000	1468145	18.1	22048
\$ 10 13C2 PFDA	515.0 > 470.0	21.435	21.418	0.017	1.000	772336	12.3	24041

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\19DEC2016A6A\_011.d

Injection Date: 19-Dec-2016 13:48:50

Instrument ID: A6

Lims ID: LCS 320-141642/2-A

Client ID:

Operator ID: CBW

ALS Bottle#: 35

Worklist Smp#: 11

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

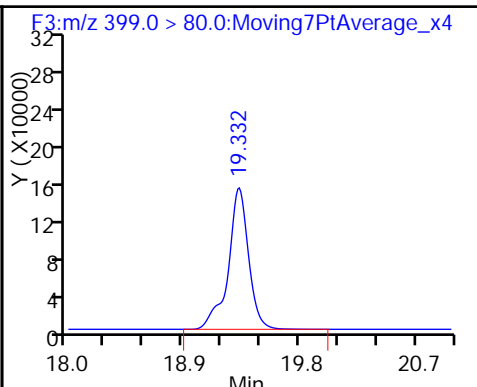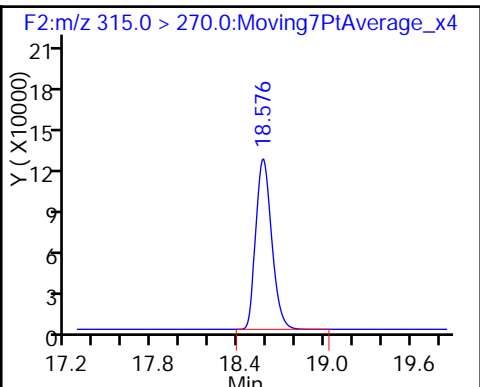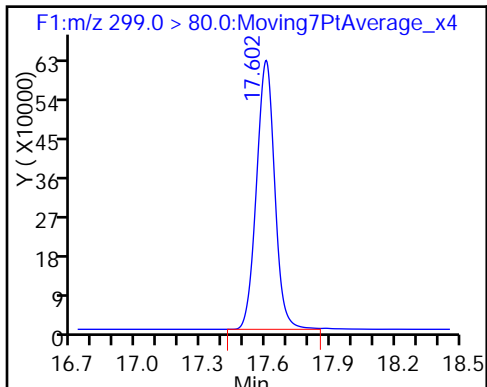

4 Perfluoroheptanoic acid

\* 5 13C2-PFOA




6 Perfluorooctanoic acid

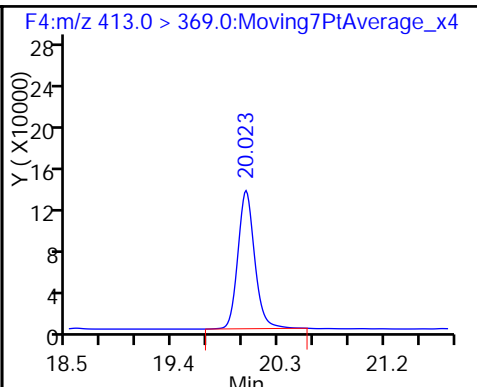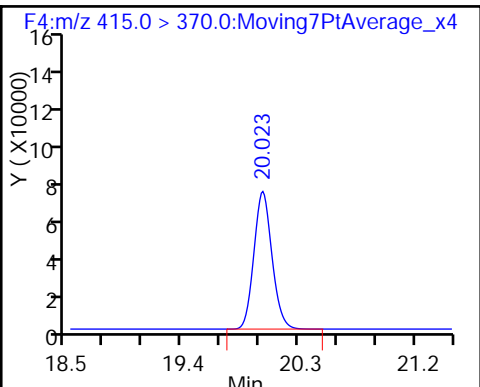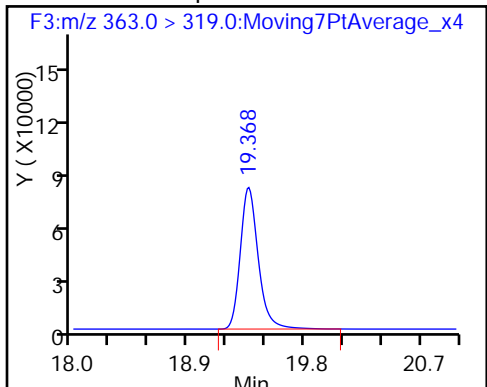


\* 8 13C4 PFOS

7 Perfluorooctane sulfonic acid

9 Perfluorononanoic acid

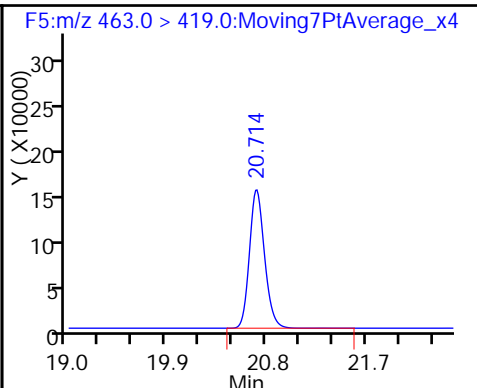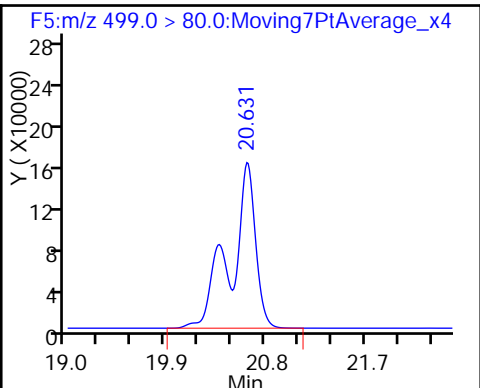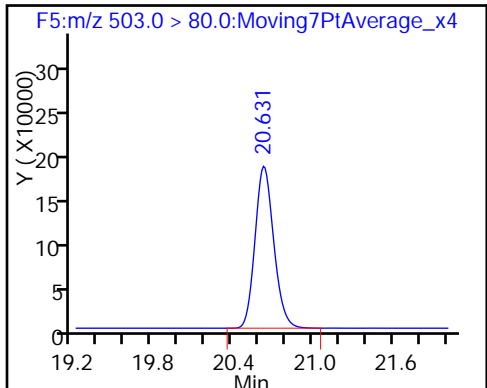

\$ 10 13C2 PFDA

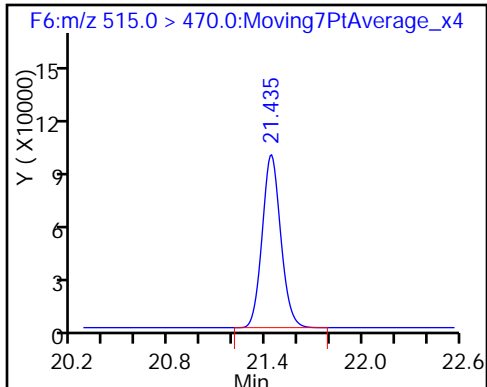

TestAmerica Sacramento  
Recovery Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\19DEC2016A6A\_011.d  
 Lims ID: LCS 320-141642/2-A  
 Client ID:  
 Sample Type: LCS  
 Inject. Date: 19-Dec-2016 13:48:50 ALS Bottle#: 35 Worklist Smp#: 11  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: lcs 320-141642/2-a  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 20-Dec-2016 13:38:59 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK032

First Level Reviewer: barnettj Date: 20-Dec-2016 09:37:03

Compound	Amount Added	Amount Recovered	% Rec.
\$ 2 13C2 PFHxA	10.0	11.8	117.99
\$ 10 13C2 PFDA	10.0	12.3	123.22

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Client Sample ID: WI-AF-3RW30-1216 MS Lab Sample ID: 320-24224-1 MS  
 Matrix: Water Lab File ID: 15DEC2016A6A\_150.d  
 Analysis Method: 537 Date Collected: 12/07/2016 14:45  
 Extraction Method: 537 Date Extracted: 12/12/2016 10:03  
 Sample wt/vol: 258 (mL) Date Analyzed: 12/18/2016 09:30  
 Con. Extract Vol.: 1.0 (mL) Dilution Factor: 1  
 Injection Volume: 10 (uL) GC Column: Acquity ID: 2.1 (mm)  
 % Moisture: \_\_\_\_\_ GPC Cleanup: (Y/N) N  
 Analysis Batch No.: 142809 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	LOQ	LOD	DL
1763-23-1	Perfluorooctanesulfonic acid (PFOS)	0.0334	J	0.058	0.047	0.015
335-67-1	Perfluorooctanoic acid (PFOA)	0.0189	J	0.029	0.023	0.0091
375-73-5	Perfluorobutanesulfonic acid (PFBS)	0.0700	J	0.14	0.11	0.046

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL00993	13C2 PFHxA	99		70-130
STL00996	13C2 PFDA	122		70-130



TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_150.d  
 Lims ID: 320-24224-A-1-B MS  
 Client ID: WI-AF-3RW30-1216  
 Sample Type: MS  
 Inject. Date: 18-Dec-2016 09:30:52 ALS Bottle#: 45 Worklist Smp#: 38  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: 320-24224-a-1-b ms  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 19-Dec-2016 15:25:13 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK025

First Level Reviewer: westendorfc Date: 19-Dec-2016 09:11:56

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
--------	----	--------	--------	--------	----------	--------------	-----	-------

1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.557	17.570	-0.013	1.000	686569	18.1	1628
\$ 2 13C2 PFHxA	315.0 > 270.0	18.530	18.548	-0.018	1.000	664119	9.94	22264
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.297	19.308	-0.011	1.000	316413	6.50	7323
4 Perfluoroheptanoic acid	363.0 > 319.0	19.320	19.344	-0.024	1.000	200013	2.87	51.5 M
* 5 13C2-PFOA	415.0 > 370.0	19.986	19.999	-0.013		572717	10.0	14467
6 Perfluorooctanoic acid	413.0 > 369.0	19.986	19.999	-0.013	1.000	290037	4.87	299
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.608	20.631	-0.023	1.000	487061	8.61	8203
* 8 13C4 PFOS	503.0 > 80.0	20.596	20.619	-0.023		1554853	28.7	40287
9 Perfluorononanoic acid	463.0 > 419.0	20.667	20.691	-0.024	1.000	311572	4.80	1373
\$ 10 13C2 PFDA	515.0 > 470.0	21.400	21.418	-0.018	1.000	610527	12.2	19209

QC Flag Legend

Review Flags

M - Manually Integrated

TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_150.d

Injection Date: 18-Dec-2016 09:30:52

Instrument ID: A6

Lims ID: 320-24224-A-1-B MS

Client ID: WI-AF-3RW30-1216

Operator ID: CBW

ALS Bottle#: 45

Worklist Smp#: 38

Injection Vol: 10.0 ul

Dil. Factor: 1.0000




Method: 537\_A6

Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

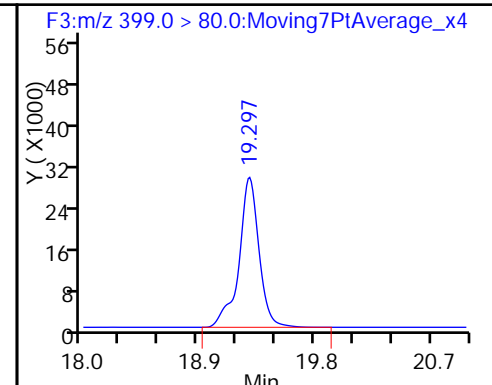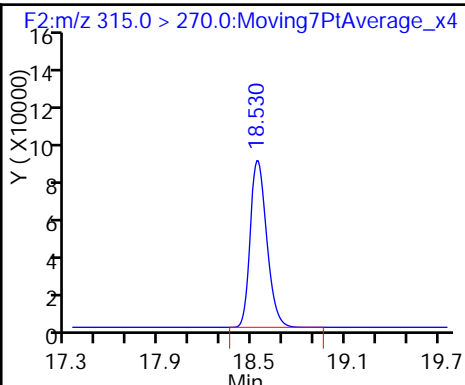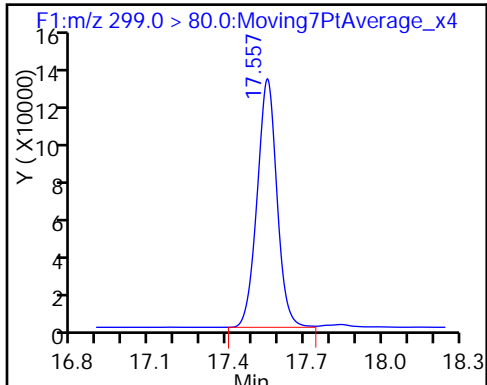

4 Perfluoroheptanoic acid (M)

\* 5 13C2-PFOA




6 Perfluorooctanoic acid

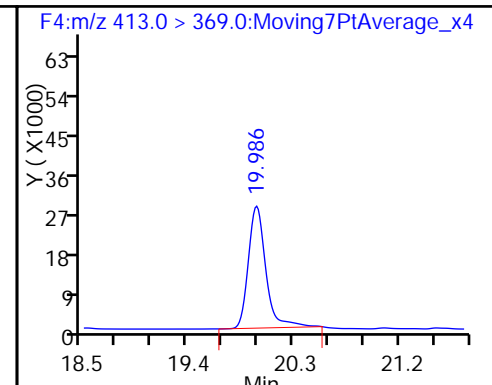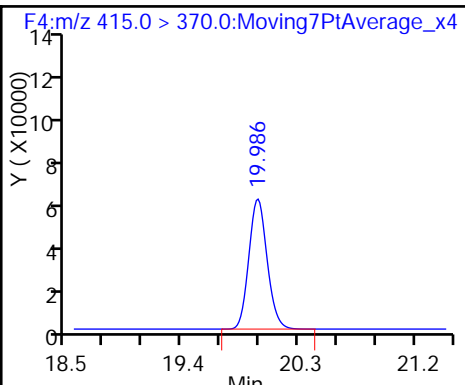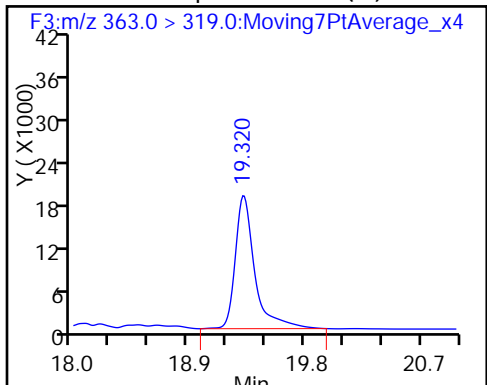


7 Perfluorooctane sulfonic acid

\* 8 13C4 PFOS

9 Perfluorononanoic acid

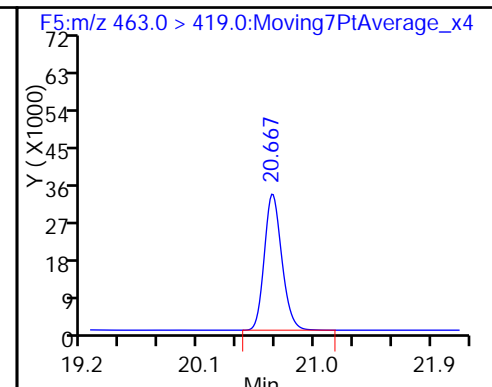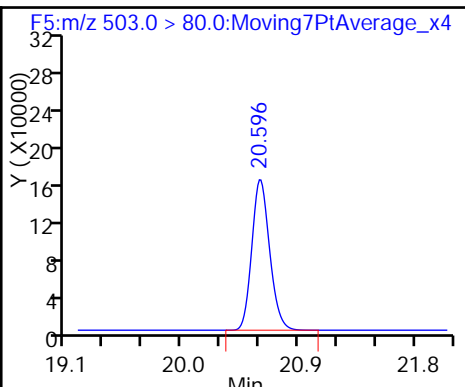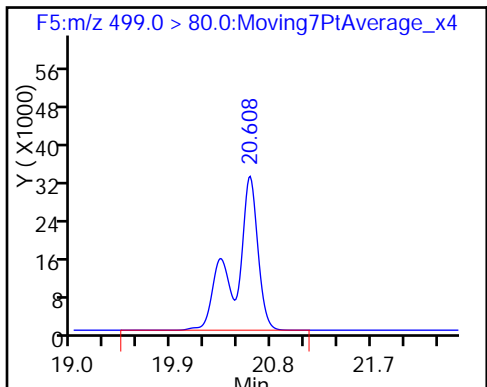

\$ 10 13C2 PFDA

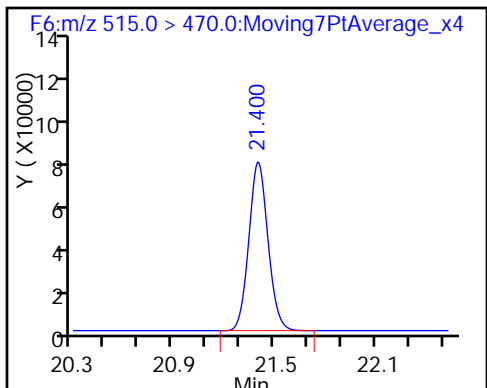

TestAmerica Sacramento  
Recovery Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_150.d  
 Lims ID: 320-24224-A-1-B MS  
 Client ID: WI-AF-3RW30-1216  
 Sample Type: MS  
 Inject. Date: 18-Dec-2016 09:30:52 ALS Bottle#: 45 Worklist Smp#: 38  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: 320-24224-a-1-b ms  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35°C  
 Operator ID: CBW Instrument ID: A6  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 19-Dec-2016 15:25:13 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK025

First Level Reviewer: westendorfc Date: 19-Dec-2016 09:11:56

Compound	Amount Added	Amount Recovered	% Rec.
\$ 2 13C2 PFHxA	10.0	9.94	99.41
\$ 10 13C2 PFDA	10.0	12.2	121.65

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Client Sample ID: WI-AF-3RW30-1216 MSD Lab Sample ID: 320-24224-1 MSD  
 Matrix: Water Lab File ID: 15DEC2016A6A\_151.d  
 Analysis Method: 537 Date Collected: 12/07/2016 14:45  
 Extraction Method: 537 Date Extracted: 12/12/2016 10:03  
 Sample wt/vol: 248.6(mL) Date Analyzed: 12/18/2016 10:00  
 Con. Extract Vol.: 1.0(mL) Dilution Factor: 1  
 Injection Volume: 10(uL) GC Column: Acquity ID: 2.1(mm)  
 % Moisture: \_\_\_\_\_ GPC Cleanup: (Y/N) N  
 Analysis Batch No.: 142809 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	LOQ	LOD	DL
1763-23-1	Perfluorooctanesulfonic acid (PFOS)	0.0337	J	0.060	0.048	0.016
335-67-1	Perfluorooctanoic acid (PFOA)	0.0165	J	0.030	0.024	0.0095
375-73-5	Perfluorobutanesulfonic acid (PFBS)	0.0706	J	0.14	0.11	0.048

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL00993	13C2 PFHxA	96		70-130
STL00996	13C2 PFDA	113		70-130

TestAmerica Sacramento  
Target Compound Quantitation Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_151.d  
 Lims ID: 320-24224-A-1-C MSD  
 Client ID: WI-AF-3RW30-1216  
 Sample Type: MSD  
 Inject. Date: 18-Dec-2016 10:00:28 ALS Bottle#: 46 Worklist Smp#: 39  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: 320-24224-a-1-c msd  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35\*C  
 Operator ID: CBW Instrument ID: A6  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 19-Dec-2016 15:25:13 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK025

First Level Reviewer: westendorfc Date: 19-Dec-2016 09:12:06

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ng/ml	S/N	Flags
1 Perfluorobutanesulfonic acid	299.0 > 80.0	17.544	17.570	-0.026	1.000	679008	17.6	1304
\$ 2 13C2 PFHxA	315.0 > 270.0	18.521	18.548	-0.027	1.000	648857	9.58	17776
3 Perfluorohexanesulfonic acid	399.0 > 80.0	19.285	19.308	-0.023	1.000	309497	6.25	1399
4 Perfluoroheptanoic acid	363.0 > 319.0	19.309	19.344	-0.035	1.000	202808	2.88	72.3
* 5 13C2-PFOA	415.0 > 370.0	19.972	19.999	-0.027		580380	10.0	14488
6 Perfluorooctanoic acid	413.0 > 369.0	19.972	19.999	-0.027	1.000	247300	4.10	171
7 Perfluorooctane sulfonic acid	499.0 > 80.0	20.596	20.631	-0.035	1.000	482751	8.39	8153
* 8 13C4 PFOS	503.0 > 80.0	20.596	20.619	-0.023		1580820	28.7	41102
9 Perfluorononanoic acid	463.0 > 419.0	20.667	20.691	-0.024	1.000	285024	4.33	1330
\$ 10 13C2 PFDA	515.0 > 470.0	21.400	21.418	-0.018	1.000	574147	11.3	18018




TestAmerica Sacramento

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_151.d  
Injection Date: 18-Dec-2016 10:00:28 Instrument ID: A6  
Lims ID: 320-24224-A-1-C MSD  
Client ID: WI-AF-3RW30-1216  
Operator ID: CBW ALS Bottle#: 46 Worklist Smp#: 39  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Method: 537\_A6 Limit Group: LC 537 ICAL

1 Perfluorobutanesulfonic acid

\$ 2 13C2 PFHxA




3 Perfluorohexanesulfonic acid

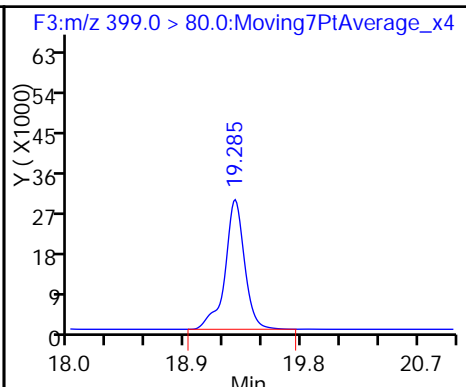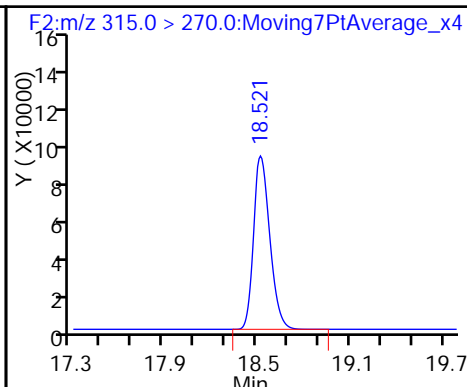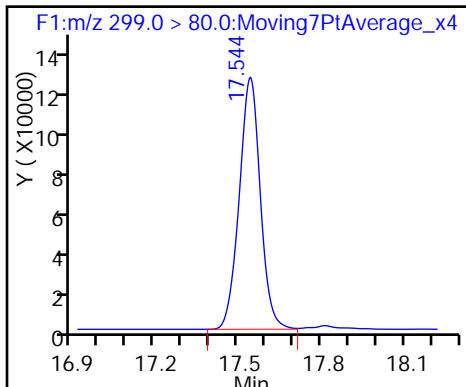

4 Perfluoroheptanoic acid

\* 5 13C2-PFOA




6 Perfluorooctanoic acid

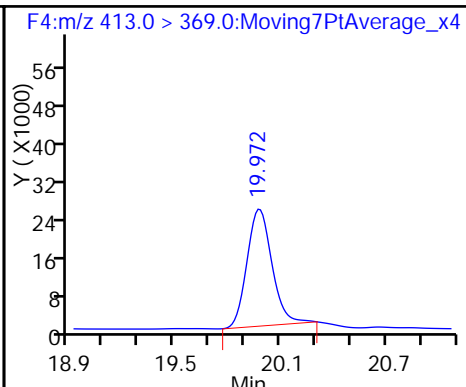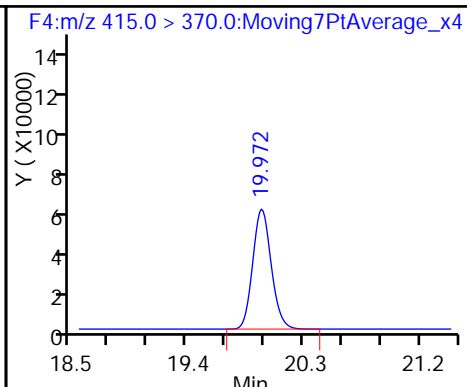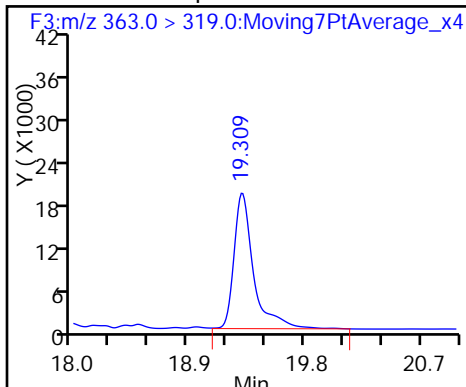


7 Perfluorooctane sulfonic acid

\* 8 13C4 PFOS

9 Perfluorononanoic acid

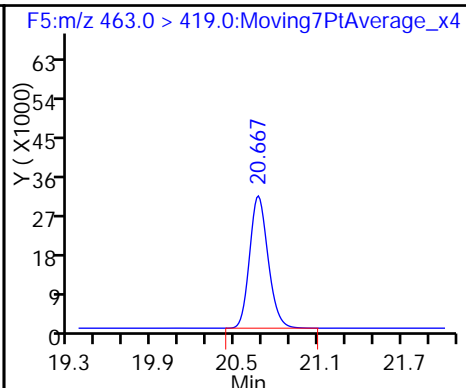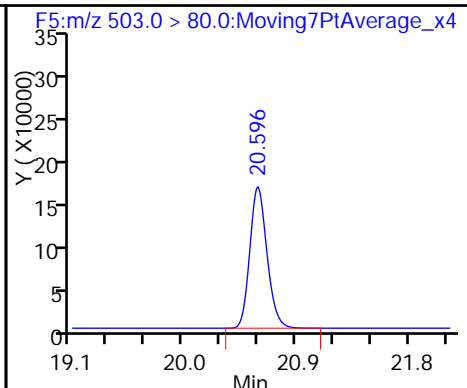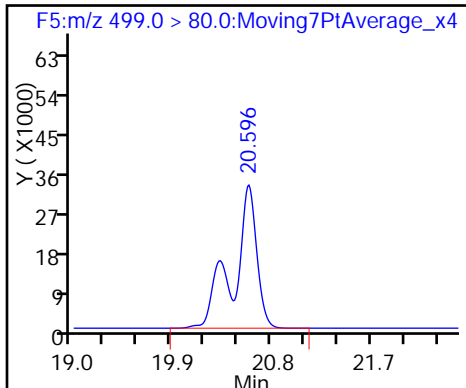

\$ 10 13C2 PFDA

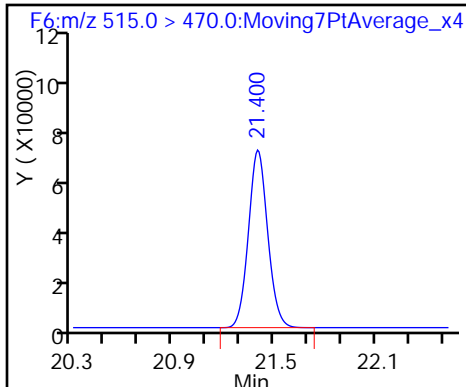

TestAmerica Sacramento  
Recovery Report

Data File: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\15DEC2016A6A\_151.d  
 Lims ID: 320-24224-A-1-C MSD  
 Client ID: WI-AF-3RW30-1216  
 Sample Type: MSD  
 Inject. Date: 18-Dec-2016 10:00:28 ALS Bottle#: 46 Worklist Smp#: 39  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: 320-24224-a-1-c msd  
 Misc. Info.: Acquity BEH 1.7um, 3X150mm T=35°C  
 Operator ID: CBW Instrument ID: A6  
 Method: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b\537\_\_A6.m  
 Limit Group: LC 537 ICAL  
 Last Update: 19-Dec-2016 15:25:13 Calib Date: 05-Dec-2016 19:54:00  
 Integrator: Picker  
 Quant Method: Internal Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Sacramento\ChromData\A6\20161205-37524.b\05DEC2016A6A\_009.d  
 Column 1 : Acquity BEH C18 ( 2.10 mm) Det: F1:MRM  
 Process Host: XAWRK025

First Level Reviewer: westendorfc Date: 19-Dec-2016 09:12:06

Compound	Amount Added	Amount Recovered	% Rec.
\$ 2 13C2 PFHxA	10.0	9.58	95.84
\$ 10 13C2 PFDA	10.0	11.3	112.89

LCMS ANALYSIS RUN LOG

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Instrument ID: A6 Start Date: 12/05/2016 17:26

Analysis Batch Number: 140688 End Date: 12/06/2016 02:48

LAB SAMPLE ID	CLIENT SAMPLE ID	DATE ANALYZED	DILUTION FACTOR	LAB FILE ID	COLUMN ID
STD 320-140688/2 IC		12/05/2016 17:26	1	05DEC2016A6A_00 4.d	Acquity 2.1 (mm)
STD 320-140688/3 IC		12/05/2016 17:55	1	05DEC2016A6A_00 5.d	Acquity 2.1 (mm)
STD 320-140688/4 IC		12/05/2016 18:25	1	05DEC2016A6A_00 6.d	Acquity 2.1 (mm)
STD 320-140688/5 ICISAV		12/05/2016 18:54	1	05DEC2016A6A_00 7.d	Acquity 2.1 (mm)
STD 320-140688/6 IC		12/05/2016 19:24	1	05DEC2016A6A_00 8.d	Acquity 2.1 (mm)
STD 320-140688/7 IC		12/05/2016 19:54	1	05DEC2016A6A_00 9.d	Acquity 2.1 (mm)
ZZZZZ		12/05/2016 20:23	1		Acquity 2.1 (mm)
CCV 320-140688/9 CCVL		12/05/2016 20:53	1	05DEC2016A6A_01 1.d	Acquity 2.1 (mm)
ZZZZZ		12/05/2016 21:22	1		Acquity 2.1 (mm)
ICV 320-140688/11		12/05/2016 21:52	1	05DEC2016A6A_01 3.d	Acquity 2.1 (mm)
ZZZZZ		12/05/2016 22:22	1		Acquity 2.1 (mm)
ZZZZZ		12/05/2016 22:51	1		Acquity 2.1 (mm)
ZZZZZ		12/05/2016 23:21	1		Acquity 2.1 (mm)
ZZZZZ		12/05/2016 23:50	1		Acquity 2.1 (mm)
ZZZZZ		12/06/2016 00:20	1		Acquity 2.1 (mm)
ZZZZZ		12/06/2016 00:49	1		Acquity 2.1 (mm)
ZZZZZ		12/06/2016 01:19	1		Acquity 2.1 (mm)
ZZZZZ		12/06/2016 01:49	1		Acquity 2.1 (mm)
ZZZZZ		12/06/2016 02:18	1		Acquity 2.1 (mm)
CCV 320-140688/21 CCVIS		12/06/2016 02:48	1		Acquity 2.1 (mm)



LCMS ANALYSIS RUN LOG

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Instrument ID: A6 Start Date: 12/15/2016 08:03

Analysis Batch Number: 142223 End Date: 12/15/2016 14:55

LAB SAMPLE ID	CLIENT SAMPLE ID	DATE ANALYZED	DILUTION FACTOR	LAB FILE ID	COLUMN ID
CCV 320-142223/2 CCVL		12/15/2016 08:03	1	15DEC2016A6A_00 2.d	Acquity 2.1(mm)
CCV 320-142223/3 CCVIS		12/15/2016 08:33	1		Acquity 2.1(mm)
ZZZZZ		12/15/2016 09:02	1		Acquity 2.1(mm)
ZZZZZ		12/15/2016 09:32	1		Acquity 2.1(mm)
ZZZZZ		12/15/2016 10:02	1		Acquity 2.1(mm)
ZZZZZ		12/15/2016 10:31	1		Acquity 2.1(mm)
ZZZZZ		12/15/2016 10:57	1		Acquity 2.1(mm)
ZZZZZ		12/15/2016 11:30	1		Acquity 2.1(mm)
ZZZZZ		12/15/2016 12:00	1		Acquity 2.1(mm)
ZZZZZ		12/15/2016 12:29	1		Acquity 2.1(mm)
ZZZZZ		12/15/2016 12:59	1		Acquity 2.1(mm)
ZZZZZ		12/15/2016 13:29	1		Acquity 2.1(mm)
ZZZZZ		12/15/2016 14:25	1		Acquity 2.1(mm)
CCV 320-142223/15 CCVIS		12/15/2016 14:55	1		Acquity 2.1(mm)

LCMS ANALYSIS RUN LOG

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Instrument ID: A6 Start Date: 12/18/2016 08:02

Analysis Batch Number: 142809 End Date: 12/18/2016 13:57

LAB SAMPLE ID	CLIENT SAMPLE ID	DATE ANALYZED	DILUTION FACTOR	LAB FILE ID	COLUMN ID
CCV 320-142809/35 CCVIS		12/18/2016 08:02	1	15DEC2016A6A_14 7.d	Acquity 2.1(mm)
ZZZZZ		12/18/2016 08:31	1		Acquity 2.1(mm)
320-24224-1		12/18/2016 09:01	1	15DEC2016A6A_14 9.d	Acquity 2.1(mm)
320-24224-1 MS		12/18/2016 09:30	1	15DEC2016A6A_15 0.d	Acquity 2.1(mm)
320-24224-1 MSD		12/18/2016 10:00	1	15DEC2016A6A_15 1.d	Acquity 2.1(mm)
320-24224-2		12/18/2016 10:30	1	15DEC2016A6A_15 2.d	Acquity 2.1(mm)
ZZZZZ		12/18/2016 10:59	1		Acquity 2.1(mm)
ZZZZZ		12/18/2016 11:29	1		Acquity 2.1(mm)
ZZZZZ		12/18/2016 11:58	1		Acquity 2.1(mm)
ZZZZZ		12/18/2016 12:28	1		Acquity 2.1(mm)
ZZZZZ		12/18/2016 12:58	1		Acquity 2.1(mm)
ZZZZZ		12/18/2016 13:27	1		Acquity 2.1(mm)
CCV 320-142809/47 CCVIS		12/18/2016 13:57	1	15DEC2016A6A_15 9.d	Acquity 2.1(mm)

LCMS ANALYSIS RUN LOG

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Instrument ID: A6 Start Date: 12/19/2016 09:45

Analysis Batch Number: 142884 End Date: 12/19/2016 11:51

LAB SAMPLE ID	CLIENT SAMPLE ID	DATE ANALYZED	DILUTION FACTOR	LAB FILE ID	COLUMN ID
CCV 320-142884/3 CCVL		12/19/2016 09:45	1	19DEC2016A6A_00 3.d	Acquity 2.1(mm)
CCV 320-142884/4 CCVIS		12/19/2016 10:15	1		Acquity 2.1(mm)
ZZZZZ		12/19/2016 10:44	1		Acquity 2.1(mm)
CCV 320-142884/7 CCVIS		12/19/2016 11:51	1		Acquity 2.1(mm)

LCMS ANALYSIS RUN LOG

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Instrument ID: A6 Start Date: 12/19/2016 11:51

Analysis Batch Number: 142886 End Date: 12/19/2016 18:32

LAB SAMPLE ID	CLIENT SAMPLE ID	DATE ANALYZED	DILUTION FACTOR	LAB FILE ID	COLUMN ID
CCV 320-142886/7 CCVIS		12/19/2016 11:51	1	19DEC2016A6A_00 7.d	Acquity 2.1(mm)
MB 320-141642/1-A		12/19/2016 13:19	1	19DEC2016A6A_01 0.d	Acquity 2.1(mm)
LCS 320-141642/2-A		12/19/2016 13:48	1	19DEC2016A6A_01 1.d	Acquity 2.1(mm)
ZZZZZ		12/19/2016 14:18	1		Acquity 2.1(mm)
ZZZZZ		12/19/2016 14:48	1		Acquity 2.1(mm)
ZZZZZ		12/19/2016 15:17	1		Acquity 2.1(mm)
ZZZZZ		12/19/2016 15:47	1		Acquity 2.1(mm)
ZZZZZ		12/19/2016 16:33	1		Acquity 2.1(mm)
ZZZZZ		12/19/2016 17:03	1		Acquity 2.1(mm)
ZZZZZ		12/19/2016 17:33	1		Acquity 2.1(mm)
ZZZZZ		12/19/2016 18:02	1		Acquity 2.1(mm)
CCV 320-142886/20 CCVIS		12/19/2016 18:32	1	19DEC2016A6A_02 0.d	Acquity 2.1(mm)

LCMS BATCH WORKSHEET

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Batch Number: 141642 Batch Start Date: 12/12/16 10:03 Batch Analyst: Sharifi, Nooshin

Batch Method: 537 Batch End Date: 12/12/16 20:52

Lab Sample ID	Client Sample ID	Method Chain	Basis	GrossWeight	TareWeight	InitialAmount	FinalAmount	ReceivedpH	LC537-IS 00026
MB 320-141642/1		537, 537				250 mL	1.0 mL	7 SU	20 uL
LCS 320-141642/2		537, 537				250 mL	1.0 mL	7 SU	20 uL
320-24224-A-1	WI-AF-3RW30-1216	537, 537	T	279.24 g	27.66 g	251.6 mL	1.0 mL	7 SU	20 uL
320-24224-A-1 MS	WI-AF-3RW30-1216	537, 537	T	285.52 g	27.56 g	258 mL	1.0 mL	7 SU	20 uL
320-24224-A-1 MSD	WI-AF-3RW30-1216	537, 537	T	276.41 g	27.79 g	248.6 mL	1.0 mL	7 SU	20 uL
320-24224-A-2	WI-AF-3FB30-1216	537, 537	T	283.94 g	26.42 g	257.5 mL	1.0 mL	7 SU	20 uL

Lab Sample ID	Client Sample ID	Method Chain	Basis	LC537-LSP 00016	LC537-MSP 00014	LC537-SU 00024	AnalysisComment		
MB 320-141642/1		537, 537				50 uL	Chlorine ND		
LCS 320-141642/2		537, 537			50 uL	50 uL	Chlorine ND		
320-24224-A-1	WI-AF-3RW30-1216	537, 537	T			50 uL	Chlorine ND		
320-24224-A-1 MS	WI-AF-3RW30-1216	537, 537	T	50 uL		50 uL	Chlorine ND		
320-24224-A-1 MSD	WI-AF-3RW30-1216	537, 537	T	50 uL		50 uL	Chlorine ND		
320-24224-A-2	WI-AF-3FB30-1216	537, 537	T			50 uL	Chlorine ND		

The pound sign (#) in the amount added field denotes that the reagent was used undiluted. All calculations are performed using the stated concentration for this reagent.

LCMS BATCH WORKSHEET

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1

SDG No.: \_\_\_\_\_

Batch Number: 141642 Batch Start Date: 12/12/16 10:03 Batch Analyst: Sharifi, Nooshin

Batch Method: 537 Batch End Date: 12/12/16 20:52

Batch Notes	
Manifold ID	2,4
Methanol ID	789820
Pipette ID	MD05306
Analyst ID - IS Reagent Drop	VPM
Analyst ID - IS Reagent Drop Witness	JER
Analyst ID - SU Reagent Drop	NSH
Analyst ID - SU Reagent Drop Witness	VPM
Analyst ID - TA Reagent Drop	NSH
Analyst ID - TA Reagent Drop Witness	VPM
SPE Cartridge ID	6332578-03
Trizma ID	SLBR4303V
Reagent Water ID	12/08/16

Basis	Basis Description
T	Total/NA

The pound sign (#) in the amount added field denotes that the reagent was used undiluted. All calculations are performed using the stated concentration for this reagent.

Job No: 24 223, 24224 Instrument ID & Date: AL 12-18-16 ICAL Batch: 140688 142993  
 Extraction Batch: 141642 Worklist #: 37983, 37999 TALS Batch: 142808, 142809, 142886

Review Items	-- Level 1 --			Level 2
	Yes	No	N/A	
<b>Initial Calibration</b>				
1. Is ICAL verified and locked in Chrom & TALS?	✓			✓
2. Is ICV properly linked in TALS?	✓			✓
<b>Continuing Calibration</b>				
1. Low-range CCV injected at start of analytical run? CCV injected after every 10 samples and at the end of the analytical run and alternated between Low-range, Mid-range and High-range?	✓			✓
2. If sequence was not after an ICAL was a low and mid range CCV injected at the start of the analytical run?	✓			✓
3. Native compounds and surrogates in control? Low-range within ±50% of true value Mid and High-range within ±30% of true value	✓			✓
4. Internal Standard areas in control? Areas ≥ 50% of average area of the ICAL and 70-140% of the most recent CCV.	✓			✓
<b>Client Samples &amp; QC Sample Results</b>				
1. Were preparation and analysis done within holding times?	✓			✓
2. Are Chromatograms reviewed and spectra verified?	✓			✓
3. Are positive results within calibration range?	✓			✓
4. Dilutions due to target cpds? _____ Dilutions due to non-targets? _____			✓	
5. All target compounds in MB < 1/3 RL ? (Requires NCM if "no.")	✓			✓
6. Are target constituents in LCS/LCSD within method control limits?	✓			✓
7. Internal Standard areas in control for all samples and QC reported? ±50% from the average area of the ICAL and 70-140% of the most recent CCV	✓			✓
8. Do results (e.g., dilutions/trip blanks) make sense?	✓			✓
9. Are MS/MSD recoveries and RPDs within method control limits?	✓			✓
10. Are all QC samples properly linked in TALS?	✓			✓
11. All manual integrations appropriate and completely documented?	✓			✓
12. Are nonconformances documented as NCMs?	✓			✓
13. Are all Chrom graphics uploaded?	✓			✓

1<sup>st</sup> Level Reviewer / Date: JRB 12-20-16 2<sup>nd</sup> Level Reviewer / Date: Murray 12/20/2016

NCM # and Comments: \_\_\_\_\_  
 \_\_\_\_\_  
 \_\_\_\_\_

Instrument ID & Date: AL 12-5-16 Worklist#: 37524

ICAL Batch: 140688, 140689 Calibration ID number: 26888, 26889

Review Items	--- Level 1 ---			Level 2
	Yes	No	N/A	
<b>Initial Calibration</b>				
1. Mass calibration, as needed, verified by full scan of PFC stock standard. All PFC ions used for quantitation are within 0.3 m/z of true mass?	✓			✓
2. Responses increase with increasing concentration?	✓			✓
3. Fit used (circle): <u>Average</u> Linear (1/x <sup>2</sup> ) Linear Quadratic (6 points minimum)				
4. Meets fit criteria? Intercept ≤ 1/2 RL RSD ≤ 30% for Average R <sup>2</sup> ≥ 0.990 for Linear R <sup>2</sup> ≥ 0.990 for Quadratic NOTE: "Force through Zero" must be used and weighted if needed	✓			✓
5. If quadratic fit used the curve does not "bend over".			✓	
6. Feed calibration points into the calculated curve. Are points ≤MRL within ±50% of true value? Are points >MRL within ±30% of true value?	✓			✓
7. Any carryover from the high calibration point must be < 1/3 RL	✓			✓
8. Asymmetry check meets criteria for the first two eluting peaks?. (0.8 - 1.5).	✓			✓
9. Is the asymmetry check scanned and linked in TALS to the calibration point?	✓			✓
10. Is ICV (2 <sup>nd</sup> source) ± 30% of true value?	✓			✓
11. Is ICV (2 <sup>nd</sup> source) internal standards ±50% of average area of the ICAL?	✓			✓
12. ICAL locked in Chrom and uploaded to TALS?	✓			✓
13. ICAL locked in TALS and scanned?				✓

1<sup>st</sup> Level Reviewer / Date: JRB 12-6-16

2<sup>nd</sup> Level Reviewer / Date: R. H. H. 12/7/16

NCM # and Comments: \_\_\_\_\_  
 \_\_\_\_\_  
 \_\_\_\_\_  
 \_\_\_\_\_



TestAmerica Laboratories  
Worklist QC Batch Report

Worklist Name: 15DEC2016B\_A6 537      Worklist Number: 37983  
 Instrument Name: A6      Chrom Method: 537\_\_A6  
 Data Directory: \\ChromNA\Sacramento\ChromData\A6\20161219-37983.b  
 QC Batching: Enabled      Limit Group Batching: Enabled

QC Batch: 1	LC 537 ICAL Raw Batch: 142797
# 1 CCV L5	# 1 CCV L5
# 2 320-24036-A-4-A	# 2 320-24036-A-4-A
# 3 320-24036-A-6-A	# 3 320-24036-A-6-A
# 4 320-24036-A-7-A	# 4 320-24036-A-7-A
# 5 LCS 320-142443/2-A	# 5 LCS 320-142443/2-A
# 6 320-24190-B-2-A	# 6 320-24190-B-2-A
# 7 320-24190-B-3-A	# 7 320-24190-B-3-A
# 8 320-24190-B-4-A	# 8 320-24190-B-4-A
# 9 320-24190-B-5-A	# 9 320-24190-B-5-A
#10 320-24190-B-6-A	#10 320-24190-B-6-A
#11 320-24190-B-7-A	#11 320-24190-B-7-A
#12 CCV L3	#12 CCV L3

QC Batch: 2	LC 537 ICAL Raw Batch: 142798
#12 CCV L3	#12 CCV L3
#13 RB	#13 RB
#14 MB Trizma	#14 MB Trizma
#15 QC 1 Trizma	#15 QC 1 Trizma
#16 QC 2 Trizma	#16 QC 2 Trizma
#17 QC 3 Trizma	#17 QC 3 Trizma
#18 QC 4 Trizma	#18 QC 4 Trizma
#19 MB 6341059-02 Cartridge QC	#19 MB 6341059-02 Cartridge QC
#20 LCS 6341059-02 Cartridge QC	#20 LCS 6341059-02 Cartridge QC
#21 MSP_00016 QC	#21 MSP_00016 QC
#22 HSP_00013 QC	#22 HSP_00013 QC
#23 CCV L5	#23 CCV L5

QC Batch: 3	LC 537 ICAL Raw Batch: 142808
#23 CCV L5	#23 CCV L5
#24 RB	#24 RB
#25 MB 320-141642/1-A	#25 MB 320-141642/1-A
#26 LCS 320-141642/2-A	#26 LCS 320-141642/2-A
#27 320-24223-A-1-A	#27 320-24223-A-1-A
#28 320-24223-A-2-A	#28 320-24223-A-2-A
#29 320-24223-A-3-A	#29 320-24223-A-3-A
#30 320-24223-A-3-B MS	#30 320-24223-A-3-B MS
#31 320-24223-A-3-C MSD	#31 320-24223-A-3-C MSD
#32 320-24223-A-4-A	#32 320-24223-A-4-A
#33 320-24223-A-5-A	#33 320-24223-A-5-A
#34 320-24223-A-6-A	#34 320-24223-A-6-A
#35 CCV L3	#35 CCV L3

QC Batch: 4	LC 537 ICAL Raw Batch: 142809
#35 CCV L3	#35 CCV L3
#36 RB	#36 RB
#37 320-24224-A-1-A	#37 320-24224-A-1-A
#38 320-24224-A-1-B MS	#38 320-24224-A-1-B MS

QC Batch: 4	LC 537 ICAL Raw Batch: 142809
#39 320-24224-A-1-C MSD	#39 320-24224-A-1-C MSD
#40 320-24224-A-2-A	#40 320-24224-A-2-A
#41 MB 320-141743/1-A	#41 MB 320-141743/1-A
#42 LCS 320-141743/2-A	#42 LCS 320-141743/2-A
#43 LCSD 320-141743/3-A	#43 LCSD 320-141743/3-A
#44 320-23932-A-1-A	#44 320-23932-A-1-A
#45 320-23932-A-2-A	#45 320-23932-A-2-A
#46 320-23932-A-3-A	#46 320-23932-A-3-A
#47 CCV L5	#47 CCV L5

QC Batch: 5	LC 537 ICAL Raw Batch: 142801	LC 537 CS ICAL Raw Batch: 142977
#47 CCV L5	#47 CCV L5	
#48 RB	#48 RB	
#49 320-23932-A-4-A	#49 320-23932-A-4-A	
#50 320-23932-A-5-A	#50 320-23932-A-5-A	
#51 320-23932-A-6-A	#51 320-23932-A-6-A	
#52 320-23932-A-7-A	#52 320-23932-A-7-A	
#53 320-23932-A-8-A	#53 320-23932-A-8-A	
#54 320-23932-A-9-A	#54 320-23932-A-9-A	
#55 320-23932-A-10-A	#55 320-23932-A-10-A	
#56 320-23932-A-11-A	#56 320-23932-A-11-A	
#57 320-23932-A-12-A	#57 320-23932-A-12-A	
#58 CCV L3	#58 CCV L3	#58 CCV L3

QC Batch: 6	LC 537 ICAL Raw Batch: 142810	LC 537 CS ICAL Raw Batch: 142811
#58 CCV L3	#58 CCV L3	#58 CCV L3
#59 RB		#59 RB
#60 MB 320-141901/1-A		#60 MB 320-141901/1-A
#61 LCS 320-141901/2-A		#61 LCS 320-141901/2-A
#62 LCSD 320-141901/3-A		#62 LCSD 320-141901/3-A
#63 320-24151-A-1-A		#63 320-24151-A-1-A
#64 320-24151-A-2-A		#64 320-24151-A-2-A
#65 320-24151-A-3-A		#65 320-24151-A-3-A
#66 320-24151-A-3-B MS		#66 320-24151-A-3-B MS
#67 320-24151-A-3-C DU		#67 320-24151-A-3-C DU
#68 320-24151-A-4-A		#68 320-24151-A-4-A
#69 320-24148-A-1-A		#69 320-24148-A-1-A
#70 CCV L5		#70 CCV L5

QC Batch: 7	LC 537 CS ICAL Raw Batch: 142813
#70 CCV L5	#70 CCV L5
#71 RB	#71 RB
#72 320-24148-A-2-A	#72 320-24148-A-2-A
#73 320-24148-A-3-A	#73 320-24148-A-3-A
#74 320-24148-A-3-B MS	#74 320-24148-A-3-B MS
#75 320-24148-A-3-C DU	#75 320-24148-A-3-C DU
#76 320-24148-A-4-A	#76 320-24148-A-4-A
#77 280-91746-A-1-A	#77 280-91746-A-1-A
#78 CCV L3	#78 CCV L3

TestAmerica Laboratories  
Worklist QC Batch Report

Worklist Name: 19DEC2016A\_A6 537      Worklist Number: 37999  
 Instrument Name: A6      Chrom Method: 537\_A6  
 Data Directory: \\ChromNA\Sacramento\ChromData\A6\20161219-37999.b  
 QC Batching: Enabled      Limit Group Batching: Enabled

QC Batch: 1	LC 537 ICAL Raw Batch: 142884	LC 537 CS ICAL Raw Batch: 142885
# 1 RB # 2 RB # 3 CCV L2 # 4 CCV L3 # 5 RB # 6 320-24284-A-5-A # 7 CCV L5	# 1 RB # 2 RB # 3 CCV L2   # 7 CCV L5	# 3 CCV L2 # 4 CCV L3 # 5 RB # 6 320-24284-A-5-A # 7 CCV L5

QC Batch: 2	LC 537 ICAL Raw Batch: 142886	LC 537 CS ICAL Raw Batch: 142941
# 7 CCV L5 # 8 RB # 9 QC IS #10 MB 320-141642/1-A #11 LCS 320-141642/2-A #12 320-24223-A-1-A #13 320-24223-A-2-A #14 320-24223-A-3-A #15 320-24223-A-3-B MS #16 320-24223-A-3-C MSD #17 320-24223-A-4-A #18 320-24223-A-5-A #19 320-24223-A-6-A #20 CCV L3	# 7 CCV L5 # 8 RB # 9 QC IS #10 MB 320-141642/1-A #11 LCS 320-141642/2-A #12 320-24223-A-1-A #13 320-24223-A-2-A #14 320-24223-A-3-A #15 320-24223-A-3-B MS #16 320-24223-A-3-C MSD #17 320-24223-A-4-A #18 320-24223-A-5-A #19 320-24223-A-6-A #20 CCV L3	# 7 CCV L5

QC Batch: 3	LC 537 ICAL Raw Batch: 142991	LC 537 CS ICAL Raw Batch: 142992
#20 CCV L3 #21 RB #22 320-23932-A-4-A #23 320-23932-A-5-A #24 320-23932-A-6-A #25 320-23932-A-7-A #26 320-23932-A-8-A #27 320-23932-A-9-A #28 320-23932-A-10-A #29 320-23932-A-11-A #30 320-23932-A-12-A #31 CCV L5	#20 CCV L3 #21 RB #22 320-23932-A-4-A #23 320-23932-A-5-A #24 320-23932-A-6-A #25 320-23932-A-7-A #26 320-23932-A-8-A #27 320-23932-A-9-A #28 320-23932-A-10-A #29 320-23932-A-11-A #30 320-23932-A-12-A #31 CCV L5	#31 CCV L5

QC Batch: 4	LC 537 ICAL Raw Batch: 142993	LC 537 CS ICAL Raw Batch: 142994
#31 CCV L5 #32 RB #33 LCS 320-141743/2-A #34 320-24148-A-2-A #35 320-24148-A-3-A #36 320-24148-A-3-B MS #37 320-24148-A-3-C DU #38 320-24148-A-4-A	#31 CCV L5 #32 RB #33 LCS 320-141743/2-A	#31 CCV L5  #34 320-24148-A-2-A #35 320-24148-A-3-A #36 320-24148-A-3-B MS #37 320-24148-A-3-C DU #38 320-24148-A-4-A

QC Batch: 4	LC 537 ICAL Raw Batch: 142993	LC 537 CS ICAL Raw Batch: 142994
#39 280-91746-A-1-A #40 320-24151-A-3-B MS #41 CCV L3	#41 CCV L3	#39 280-91746-A-1-A #40 320-24151-A-3-B MS #41 CCV L3

QC Batch: 5	LC 537 ICAL Raw Batch: 142995	LC 537 CS ICAL Raw Batch: 142996
#41 CCV L3 #42 RB #43 MB 320-137113/1-A #44 320-20510-A-4-B #45 320-20510-A-5-B #46 320-20510-A-6-B #47 MB 320-142439/1-A #48 LLCS 320-142439/2-A #49 LLCSD 320-142439/3-A #50 320-24316-A-1-A #51 320-24316-A-2-A #52 320-24316-A-3-A #53 CCV L5	#41 CCV L3 #42 RB #43 MB 320-137113/1-A #44 320-20510-A-4-B #45 320-20510-A-5-B #46 320-20510-A-6-B #47 MB 320-142439/1-A #48 LLCS 320-142439/2-A #49 LLCSD 320-142439/3-A #50 320-24316-A-1-A #51 320-24316-A-2-A #52 320-24316-A-3-A #53 CCV L5	#41 CCV L3

QC Batch: 6	LC 537 ICAL Raw Batch: 142997
#53 CCV L5 #54 RB #55 320-24316-A-4-A #56 320-24316-A-5-A #57 320-24316-A-6-A #58 320-24316-A-7-A #59 320-24316-A-8-A #60 320-24316-A-9-A #61 320-24316-A-10-A #62 320-24316-A-11-A #63 320-24316-A-12-A #64 320-24316-A-13-A #65 CCV L3 #66 RB	#53 CCV L5 #54 RB #55 320-24316-A-4-A #56 320-24316-A-5-A #57 320-24316-A-6-A #58 320-24316-A-7-A #59 320-24316-A-8-A #60 320-24316-A-9-A #61 320-24316-A-10-A #62 320-24316-A-11-A #63 320-24316-A-12-A #64 320-24316-A-13-A #65 CCV L3 #66 RB

23

# Aqueous Extraction Analysis Sheet

(To Accompany Samples to Instruments)

Batch Number: 320-141642

Method Code: 320-537\_Prep-320

Analyst: Sharifi, Nooshin

Batch Open: 12/12/2016 10:03:00AM

Batch End: 12/12/16 20:52

Screened AY 12/16/16 No DV needed

AG 12/17/16

Push

## Extraction of Perfluorinated Alkyl Acids

Input Sample Lab ID (Analytical Method)	SDG (Job #)	GrossWt TareWt	InitAmnt		PHs Adj1 Adj2	Due Date	Analytical TAT	Div Rank	Comments	Output Sample Lab ID
			FinAmnt	Rcvd						
1 MB-320-141642/1 N/A	N/A		250 mL	7		N/A	N/A	N/A	Chlorine ND	MB 320-141642/1-A
			1.0 mL							
2 LCS-320-141642/2 N/A	N/A		250 mL	7		N/A	N/A	N/A	Chlorine ND	LCS 320-141642/2-A
			1.0 mL							
3 320-24223-A-1 (537_DOD5)	N/A (320-24223-1)	272.15 g	245 mL	7		12/19/16	10_Day_Rush	4	Chlorine ND	320-24223-A-1-A
		27.15 g	1.0 mL							
4 320-24223-A-2 (537_DOD5)	N/A (320-24223-1)	276.01 g	249.8 mL	7		12/19/16	10_Day_Rush	4	Chlorine ND	320-24223-A-2-A
		26.17 g	1.0 mL							
5 320-24223-A-3 (537_DOD5)	N/A (320-24223-1)	267.09 g	241.2 mL	7		12/19/16	10_Day_Rush	4	Chlorine ND	320-24223-A-3-A
		25.89 g	1.0 mL							
6 320-24223-A-3-MS (537_DOD5)	N/A (320-24223-1)	265.63 g	239.5 mL	7		12/19/16	10_Day_Rush	4	Chlorine ND	320-24223-A-3-B
		26.13 g	1.0 mL							MS
7 320-24223-A-3-MSD (537_DOD5)	N/A (320-24223-1)	260.80 g	234.8 mL	7		12/19/16	10_Day_Rush	4	Chlorine ND	320-24223-A-3-C
		25.98 g	1.0 mL							MSD
8 320-24223-A-4 (537_DOD5)	N/A (320-24223-1)	280.83 g	254.1 mL	7		12/19/16	10_Day_Rush	4	Chlorine ND	320-24223-A-4-A
		26.75 g	1.0 mL							
9 320-24223-A-5 (537_DOD5)	N/A (320-24223-1)	281.90 g	255.5 mL	7		12/19/16	10_Day_Rush	4	Chlorine ND	320-24223-A-6-A
		26.38 g	1.0 mL							
10 320-24223-A-6 (537_DOD5)	N/A (320-24223-1)	282.00 g	255.6 mL	7		12/19/16	10_Day_Rush	4	Chlorine ND	320-24223-A-6-A
		26.41 g	1.0 mL							

# Aqueous Extraction Analysis Sheet

(To Accompany Samples to Instruments)








Batch Number: 320-141642

Analyst: Sharifi, Nooshin

Batch Open: 12/12/2016 10:03:00AM

Method Code: 320-537\_Prep-320

Batch End:

11	320-24224-A-1 (537_DOD5)	N/A (320-24224-1)	279.24 g	251.6 mL	7	12/19/16	10_Day_Rush	4	Chlorine ND	
			27.66 g	1.0 mL						
12	320-24224-A-1-MS (537_DOD5)	N/A (320-24224-1)	285.52 g	258 mL	7	12/19/16	10_Day_Rush	4	Chlorine ND	
			27.56 g	1.0 mL						
13	320-24224-A-1-MSD (537_DOD5)	N/A (320-24224-1)	276.41 g	248.6 mL	7	12/19/16	10_Day_Rush	4	Chlorine ND	
			27.79 g	1.0 mL						
14	320-24224-A-2 (537_DOD5)	N/A (320-24224-1)	283.94 g	257.5 mL	7	12/19/16	10_Day_Rush	4	Chlorine ND	
			26.42 g	1.0 mL						

# Aqueous Extraction Analysis Sheet

(To Accompany Samples to Instruments)

Batch Number: 320-141642

Method Code: 320-537\_Prep-320

Analyst: Sharifi, Nooshin

Batch Open: 12/12/2016 10:03:00AM

Batch End:

Batch Notes	
Manifold ID	2,4
Trizma ID	SLBR4303V
SPE Cartridge ID	6332578-03
Methanol ID	789820
Reagent Water ID	12/08/16
Pipette ID	MD05306
Analyst ID - TA Reagent Drop	NSH
Analyst ID - TA Reagent Drop Witness	VPM
Analyst ID - SU Reagent Drop	NSH
Analyst ID - SU Reagent Drop Witness	VPM
Analyst ID - IS Reagent Drop	VPM 799828 (20 uL)
Analyst ID - IS Reagent Drop Witness	JER
Batch Comment	

## Comments

--

# Aqueous Extraction Analysis Sheet

(To Accompany Samples to Instruments)

Analyst: Sharifi, Nooshin

Batch Number: 320-141642

Method Code: 320-537\_Prep-320

Batch Open: 12/12/2016 10:03:00AM

Batch End:

## Reagent Additions Worksheet

Lab ID	Reagent Code	Amount Added	Final Amount	By	Witness
MB 320-141642/1	LC537-SU_00024	50 uL	1.0 mL	YSH 12-12-16 	VPM 12-12-16 
LCS 320-141642/2	LC537-MSP_00014	50 uL	1.0 mL		
LCS 320-141642/2	LC537-SU_00024	50 uL	1.0 mL		
320-24223-A-1	LC537-SU_00024	50 uL	1.0 mL		
320-24223-A-2	LC537-SU_00024	50 uL	1.0 mL		
320-24223-A-3	LC537-SU_00024	50 uL	1.0 mL		
320-24223-A-3 MS	LC537-LSP_00016	50 uL	1.0 mL		
320-24223-A-3 MS	LC537-SU_00024	50 uL	1.0 mL		
320-24223-A-3 MSD	LC537-LSP_00016	50 uL	1.0 mL		
320-24223-A-3 MSD	LC537-SU_00024	50 uL	1.0 mL		
320-24223-A-4	LC537-SU_00024	50 uL	1.0 mL		
320-24223-A-5	LC537-SU_00024	50 uL	1.0 mL		
320-24223-A-6	LC537-SU_00024	50 uL	1.0 mL		
320-24224-A-1	LC537-SU_00024	50 uL	1.0 mL		
320-24224-A-1 MS	LC537-LSP_00016	50 uL	1.0 mL		
320-24224-A-1 MS	LC537-SU_00024	50 uL	1.0 mL		
320-24224-A-1 MSD	LC537-LSP_00016	50 uL	1.0 mL		
320-24224-A-1 MSD	LC537-SU_00024	50 uL	1.0 mL		



# Aqueous Extraction Analysis Sheet

(To Accompany Samples to Instruments)

Analyst: Sharifi, Nooshin

Batch Number: 320-141642

Batch Open: 12/12/2016 10:03:00AM

Method Code: 320-537\_Prep-320

Batch End:

320-24224-A-2	LC537-SU_00024	50 uL	1.0 mL	NSA 12-12-16	VPM 12/12/16
---------------	----------------	-------	--------	--------------	--------------

Reagent	Other Reagents:	Amount/Units	Lot#:

Preparation Batch Number(s): 141042

Test: 5070005 DUSH

Earliest Holding Time: 12-21-16

Sample List Tab		1 <sup>st</sup> Level Reviewer	2 <sup>nd</sup> Level Reviewer
Samples identified to the correct method			✓
All necessary NCMs filed (including holding time)		NA	NA
Method/sample/login/QAS checked and correct			✓
Worksheet Tab		1 <sup>st</sup> Level Reviewer	2 <sup>nd</sup> Level Reviewer
All samples properly preserved		✓	✓
Weights in anticipated range and not targeted		✓	✓
All additional test requirements performed, documented, and uploaded to TALS correctly (e.g. final amount, initial amount, turbidity, and CI Check)		✓	✓
The pH is transcribed correctly in TALS		✓	✓
All additional information transcribed into TALS is correct and raw data is attached		✓	✓
Comments are transcribed correctly in TALS		✓	✓
Reagents Tab		1 <sup>st</sup> Level Reviewer	2 <sup>nd</sup> Level Reviewer
All necessary reagents not expired and entered into TALS		✓	✓
All spike amounts correct and added to necessary samples and QC		✓	✓
Batch Information		1 <sup>st</sup> Level Reviewer	2 <sup>nd</sup> Level Reviewer
Date and time accurate and entered into TALS correctly		✓	✓
All necessary 'batch information' complete and entered into TALS correctly		✓	✓

1<sup>st</sup> Level Reviewer: VPM

Date: 12/12/16

2<sup>nd</sup> Level Reviewer: [Signature]

Date: 12/12/16

Comments: \_\_\_\_\_

# Shipping and Receiving Documents

Regulatory Program:  DW  NPDES  RCRA  Other: \_\_\_\_\_

Project Manager: Katie Tippin  
Tel/Fax: (757) 671-6258

Client Contact  
Tiffany Hill

Project Chemist  
1100 NE Circle Blvd Ste 300 Corvallis, OR 97330  
(541) 768-3109  
(541) 908-3794  
Project Name: CTO-08  
Site: NAS Whidbey Island  
P O #: 100067106050 - 679580 06.FIFS

Site Contact: Eric Epple  
Lab Contact: Laura Turpen

Date: 12/8/2016  
Carrier: FedEx

COC No: 7 of 1 COCs

Sampler:  
For Lab Use Only:  
Walk-in Client:  
Lab Sampling:  
Job / SDG No.:

Analysis Turnaround Time  
 CALENDAR DAYS  WORKING DAYS  
TAT if different from Below 7-Day \_\_\_\_\_  
 2 weeks  
 1 week  
 2 days  
 1 day

Sample Identification	Sample Date	Sample Time	Sample Type (C-Comp, G-Grab)	Matrix	# of Cont.	Filtered Sample (Y/N)		Perform MS/MSD (Y/N)		USEPA Method 537 (PFOA, PFOS, and PFBS)	Sample Specific Notes:
						Y	N	Y	N		
WI-AF-3RW30-1216	12/07/16	1445	G	DW	2	N	N	X	X		
WI-AF-3RW30-1216-MS	12/07/16	1445	G	DW	2	N	N	X	X		
WI-AF-3RW30-1216-SD	12/07/16	1445	G	DW	2	N	N	X	X		
* WI-AF-3FB30-1216	12/07/16	1446	G	DW	2	N	N	X	X		
* containing ID											
WI-AF-FB30-1216											
06/12/2016											



Preservation Used: 1= Ice, 2= HCl; 3= H2SO4; 4=HNO3; 5=NaOH; 6= Other Trizma  
Possible Hazard Identification: Please List any EPA Waste Codes for the sample in the Comments Section if the lab is to dispose of the sample.

Special Instructions/QC Requirements & Comments:  
 Non-Hazard  Flammable  Skin Irritant  Poison B  Unknown

Custody Seals Intact: <input type="checkbox"/> Yes <input type="checkbox"/> No	Custody Seal No.:	Cooler Temp. (°C): Obs'd: _____ Cor'd: _____	Therm ID No.:
Relinquished by: Eric Epple	Company: CH2M	Received by: Laura Turpen	Company: SA SAE
Relinquished by: _____	Company: _____	Received by: _____	Company: _____
Relinquished by: _____	Company: _____	Received in Laboratory by: _____	Company: _____

# Login Sample Receipt Checklist

Client: CH2M Hill Constructors, Inc.

Job Number: 320-24224-1

**Login Number: 24224**

**List Source: TestAmerica Sacramento**

**List Number: 1**

**Creator: Turpen, Troy**

Question	Answer	Comment
Radioactivity wasn't checked or is $\leq$ background as measured by a survey meter.	True	
The cooler's custody seal, if present, is intact.	True	(E88)
Sample custody seals, if present, are intact.	N/A	
The cooler or samples do not appear to have been compromised or tampered with.	True	
Samples were received on ice.	True	
Cooler Temperature is acceptable.	True	
Cooler Temperature is recorded.	True	
COC is present.	True	
COC is filled out in ink and legible.	True	
COC is filled out with all pertinent information.	True	
Is the Field Sampler's name present on COC?	N/A	
There are no discrepancies between the containers received and the COC.	False	Refer to Job Narrative for details.
Samples are received within Holding Time (excluding tests with immediate HTs)	True	
Sample containers have legible labels.	True	
Containers are not broken or leaking.	True	
Sample collection date/times are provided.	True	
Appropriate sample containers are used.	True	
Sample bottles are completely filled.	True	
Sample Preservation Verified.	N/A	
There is sufficient vol. for all requested analyses, incl. any requested MS/MSDs	True	
Containers requiring zero headspace have no headspace or bubble is $<6\text{mm}$ (1/4").	True	
Multiphasic samples are not present.	True	
Samples do not require splitting or compositing.	True	
Residual Chlorine Checked.	N/A	

Account ID	Product	Account Type	Account Name	Account Address	Account City	Account State	Account Zip	Account Phone	Account Email	Account Status	Account Open Date	Account Close Date	Account Balance	Account Interest	Account Fee	Account Rate	Account Term	Account Class	Account Category	Account Sub-Category	Account Description	Account Remarks	Account Note	Account Label	Account Code	Account Value	Account Unit	Account Quantity	Account Price	Account Cost	Account Profit	Account Margin	Account Commission	Account Net Income	Account Gross Income	Account Tax	Account Expense	Account Total
1000000000	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...	...

**DATA VALIDATION SUMMARY REPORT  
WHIDBEY ISLAND, WASHINGTON**

Client: CH2M HILL, Inc., Corvallis, Oregon  
 SDG: 320-24224-1  
 Laboratory: Test America, Sacramento, California  
 Site: Whidbey Island, CTO-0008, Washington  
 Date: January 5, 2017

PFCs			
EDS ID	Client Sample ID	Laboratory Sample ID	Matrix
1	WI-AF-3RW30-1216	320-24224-1	Water
1MS	WI-AF-3RW30-1216MS	320-24224-1MS	Water
1MSD	WI-AF-3RW30-1216MSD	320-24224-1MSD	Water
2	WI-AF-3FB30-1216	320-24224-2	Water

A full data validation was performed on the analytical data for one water sample and one aqueous field blank sample collected on December 7, 2016 by CH2M HILL at the Whidbey Island site in Washington. The samples were analyzed under the EPA Method “Determination of Selected Perfluorinated Alkyl Acids in Drinking Water by Solid Phase Extraction and Liquid Chromatography/Tandem Mass Spectrometry (LC/MS/MS)”.

Specific method references are as follows:

Analysis  
PFCs

Method References  
USEPA Method 537 Rev 1.1 Modified

The data have been validated according to the protocols and quality control (QC) requirements of the analytical method, and the U.S. Department of Defense (DoD) Quality Systems Manual (QSM), Version 5.0 (DoD 2013) and the USEPA National Functional Guidelines for Organic Data Review as follows:

- The USEPA “Contract Laboratories Program National Functional Guidelines for Superfund Organic Methods Data Review,” August 2014;
- and the reviewer's professional judgment.

The following data quality indicators were reviewed for this report:

***Organics***

- Date Completeness, Case Narrative & Custody Documentation
- Holding times
- Gas Chromatography/Mass Spectrometry (GC/MS) Tuning

- Initial and continuing calibration summaries
- Method blank and field QC blank contamination
- Surrogate Spike recoveries
- Matrix Spike/Matrix Spike Duplicate (MS/MSD) recoveries
- Laboratory Control Sample/Laboratory Control Sample Duplicate (LCS/LCSD) recoveries
- Internal standard area and retention time summary forms
- Target Compound Identification
- Compound Quantitation
- Field Duplicate sample precision

A full (Level IV) data validation was performed with this review including a recalculation of 10% of the detected results in the samples.

### **Data Usability Assessment**

There were no rejections of data.

Overall the data is acceptable for the intended purposes. There were no qualifications.

### **Perfluorinated Compounds (PFCs)**

### **Data Completeness, Case Narrative & Custody Documentation**

- The case narrative and chain-of-custody documentation were included in the data package as required. All criteria were met.

### **Holding Times**

- All samples were extracted within 14 days for water samples and analyzed within 28 days.

### **GC/MS Tuning**

- All criteria were met.

### **Initial Calibration**

- All relative standard deviation (%RSD) and/or correlation coefficients criteria were met.



### Continuing Calibration

- All percent difference (%D) and RRF criteria were met.

### Method Blank

- The method blanks were free of contamination.

### Field QC Blank

- The field blank samples were free of contamination.

### Surrogate Spike Recoveries

- All samples exhibited acceptable surrogate %R values.

### Matrix Spike/Matrix Spike Duplicate (MS/MSD) Recoveries

- The MS/MSD samples exhibited acceptable %R and RPD values.

### Laboratory Control Samples/Laboratory Control Sample Duplicates

- The LCS/LCSD samples exhibited acceptable percent recoveries (%R) and RPD values.

### Internal Standard (IS) Area Performance

- All internal standards met response and retention time (RT) criteria.

### Target Compound Identification

- All mass spectra and quantitation criteria were met.

### Compound Quantitation

- All criteria were met.

### Field Duplicate Sample Precision

- Field duplicate samples were not collected.

Please contact the undersigned at (757) 564-0090 if you have any questions or need further information.

Signed: Nancy Weaver  
Nancy Weaver  
Senior Chemist

Dated: 1/6/17

Data Qualifier	Definition
U	The analyte was analyzed for, but was not detected above the level of the reported sample quantitation limit.
J	The analyte is an estimated quantity. The associated numerical value is the approximate concentration of the analyte in the sample.
NJ	The analysis has been "tentatively identified" or "presumptively" as present and the associated numerical value is the estimated concentration in the samples.
UJ	The analyte was analyzed for but was not detected. The reported quantitation limit is approximate and may be inaccurate or imprecise.
R	The data are unusable. The sample results are rejected due to serious deficiencies in meeting QC criteria. The analyte may or may not be present in the samples.



»

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Client Sample ID: WI-AF-3RW30-1216 Lab Sample ID: 320-24224-1  
 Matrix: Water Lab File ID: 15DEC2016A6A\_149.d  
 Analysis Method: 537 Date Collected: 12/07/2016 14:45  
 Extraction Method: 537 Date Extracted: 12/12/2016 10:03  
 Sample wt/vol: 251.6(mL) Date Analyzed: 12/18/2016 09:01  
 Con. Extract Vol.: 1.0(mL) Dilution Factor: 1  
 Injection Volume: 10(uL) GC Column: Acquity ID: 2.1(mm)  
 % Moisture: \_\_\_\_\_ GPC Cleanup: (Y/N) N  
 Analysis Batch No.: 142809 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	LOQ	LOD	DL
1763-23-1	Perfluorooctanesulfonic acid (PFOS)	0.048	U	0.060	0.048	0.015
335-67-1	Perfluorooctanoic acid (PFOA)	0.024	U	0.030	0.024	0.0094
375-73-5	Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.14	0.11	0.047

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL00993	13C2 PFHxA	100		70-130
STL00996	13C2 PFDA	117		70-130

MW 1/5/17

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

2


Lab Name: TestAmerica Sacramento Job No.: 320-24224-1  
 SDG No.: \_\_\_\_\_  
 Client Sample ID: WI-AF-3FB30-1216 Lab Sample ID: 320-24224-2  
 Matrix: Water Lab File ID: 15DEC2016A6A\_152.d  
 Analysis Method: 537 Date Collected: 12/07/2016 14:46  
 Extraction Method: 537 Date Extracted: 12/12/2016 10:03  
 Sample wt/vol: 257.5(mL) Date Analyzed: 12/18/2016 10:30  
 Con. Extract Vol.: 1.0(mL) Dilution Factor: 1  
 Injection Volume: 10(uL) GC Column: Acquity ID: 2.1(mm)  
 % Moisture: \_\_\_\_\_ GPC Cleanup: (Y/N) N  
 Analysis Batch No.: 142809 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	LOQ	LOD	DL
1763-23-1	Perfluorooctanesulfonic acid (PFOS)	0.047	U	0.058	0.047	0.015
335-67-1	Perfluorooctanoic acid (PFOA)	0.023	U	0.029	0.023	0.0091
375-73-5	Perfluorobutanesulfonic acid (PFBS)	0.11	U	0.14	0.11	0.046

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL00993	13C2 PFHxA	121		70-130
STL00996	13C2 PFDA	126		70-130

NW 1/5/17



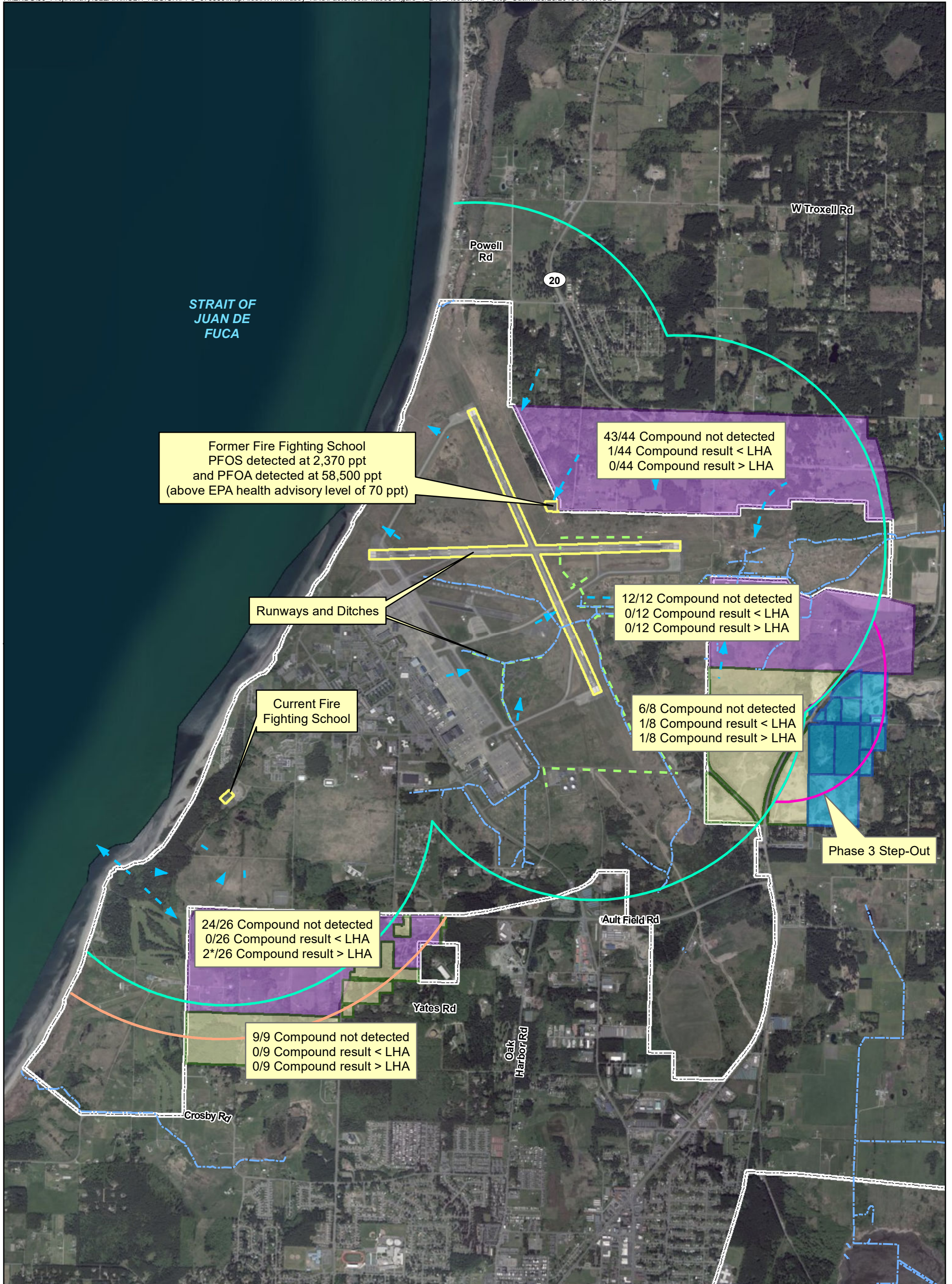

Former Fire Fighting School  
 PFOS detected at 2,370 ppt  
 and PFOA detected at 58,500 ppt  
 (above EPA health advisory level of 70 ppt)

43/44 Compound not detected  
 1/44 Compound result < LHA  
 0/44 Compound result > LHA

Runways and Ditches

12/12 Compound not detected  
 0/12 Compound result < LHA  
 0/12 Compound result > LHA

Current Fire Fighting School

6/8 Compound not detected  
 1/8 Compound result < LHA  
 1/8 Compound result > LHA

Phase 3 Step-Out

24/26 Compound not detected  
 0/26 Compound result < LHA  
 2\*/26 Compound result > LHA

9/9 Compound not detected  
 0/9 Compound result < LHA  
 0/9 Compound result > LHA

**Legend**

- 1 Mile Zone
- Half-mile Step-out Downgradient
- - - Surface Water
- - - Drainage Ditch
- Half-mile Step-out Downgradient
- Suspected Source Area
- Parcels in Phase 1 Sampling Area
- Parcels Identified in Phase 2 Sampling Area
- Parcels Identified in Phase 3 Sampling Area

- Base Boundary
- - - Inferred Groundwater Flow Direction

\* Second result above the EPA health advisory is from a duplicate sample collected from the well with the first exceedance near Ault Field.

Note:  
 PFOA and PFOS results reflected on figure,  
 PFBS results discussed in Table 2 and text.



0 0.225 0.45  
 Miles

1 inch = 0.45 mile  
 Imagery Source: Esri

Figure 2  
 Results for Drinking Water Well Sampling  
 Ault Field  
 Naval Air Station Whidbey Island  
 Oak Harbor, Washington